



Accessory catalogue

2023





Easier workflows, every day, every step of the patient pathway

Dear Customer,

Quality medical equipment, innovative digital solutions, cost-effective management, these are cornerstones that make it easier for you to provide dedicated healthcare to your patients. Now, imagine improving the quality of healthcare you deliver by simplifying your hospital's chain of acquiring the technology and services needed. Whether it's help in: planning ergonomic workstations, managing costs of care, or in improving staff satisfaction, we make it possible for you to enable high quality therapy – all from a single source. The same detailed attention that we design into our ventilators, anaesthesia workstations, incubators and monitoring products, goes into the engineering and manufacturing of our accessories and consumables, as they are an integral part of life-saving medical equipment. Improved clinical, operational and economic outcomes through the entire patient pathway is our common goal.

Freeing you of routine tasks by streamlining and facilitating your hospital processes and workflows, reducing staff stress and unnecessary costs by offering you tailor-made solutions, at the end of the day it means giving you more time to care for your patients' wellbeing and recovery. We provide you with the necessary tools to give high quality care through complete therapeutic solutions. For us, ensuring a commitment to outstanding quality means rigorously testing our products so that they fulfill all standard requirements. Together with you, we work to design intelligent, holistic workspaces that satisfy staff, manage costs, and promote patient safety and healing. Medical technology and solutions are what we do – every day, everywhere – this is what we expect of ourselves, this is why we are **Your Specialist in Acute Care.**

Helping to shape a more sustainable future

For over 130 years, we have committed ourselves to providing excellence in the quality of our products. But we don't believe in stopping here. We also take responsibility in creating sustainable products to lessen our impact on the environment.

We have agreed with nearly all of our suppliers to a code of conduct which defines specific requirements regarding workers and environmental protection. Furthermore, we have established an Air-to-Sea programme which continuously switches transport volume from air to sea freight, thus reducing CO₂ emissions. To combine patient safety and sustainability, the majority of our accessories are free of natural latex and do not contain PVC or phthalates (DEHP amongst others). As described in our Sustainability Report, we have embraced numerous measures to create the best possible conditions for working in an environmentally friendly manner. In recognition of our continuous efforts, we have, once again, been awarded with the EcoVadis Gold Medal for 2022.

For us, sustainability reaches beyond our company grounds. It means creating the best possible conditions for our customers, business partners, employees – and our planet as a whole.

As Your Specialist in Acute Care, this is our commitment.



On this way you'll find your accessories:

VENTILATION ACCESSORIES	06	WARD EQUIPMENT	258
1. Breathing circuits		1. Oxygen therapy	
2. Breathing system filters and HMEs		2. Medical aspiration	
3. CO ₂ cuvettes		3. Accessory for gas supply	
4. Catheter mounts and airway connectors			
5. Non-invasive patient interfaces and accessories		RAIL ACCESSORIES AND	
6. Oxymixer		EXAMINATION LIGHT	336
7. Device related accessories		1. Rail accessories	
8. Pack2Go®/Set2go		2. Examination light	
9. Humidifier and nebulizer		3. Hinged arms	
10. Spare parts for Resuscitators		4. Device related accessories	
SODA LIME	96	NEONATAL CARE AND THERMOREGULATIONS	344
		1. Warming therapy	
SENSORS	102	2. Ventilation accessories	
		3. Sensors	
MONITORING ACCESSORIES	108	4. Monitoring accessories	
1. Infinity® Monitoring: Gamma, Delta, Kappa and Vista series		5. Ward equipment	
2. Infinity® Acute Care System		6. Rail accessories and examination light	
3. Infinity® M300			
4. Vista 120		INDEX BY ORDER NUMBER	438
5. Vista 120 SPOTCHECK			
6. Multi Vendor Adapter		INDEX BY PRODUCT NAME	449



Ventilation Accessories

Making breathing easier

Comfortable, stable oxygen delivery is of the highest importance when administering ventilation. You need reliable lung-protective strategies that are ready for use at a moment's notice, including easily integrated ventilation accessories that work seamlessly in the background – so you can stay focused on your patients' well-being.

Our wide range of accessories enhance patient comfort whether it is for non-invasive therapies, humidification solutions, or heated breathing circuits. Patient comfort and ease of use were the focal point when designing our ventilation accessories. Moreover, our disposable accessory range of breathing circuits and filters for active protection help to prevent cross contamination and greatly improve infection prevention programmes in your clinical environment.

As a healthcare professional, we understand you value quality when delivering ventilation therapy. This is why we test our products extensively for both quality and compatibility with our devices.

We are Your Specialist in Acute Care every step of the treatment pathway: from the emergency department or surgical area to the intensive care and neonate units.

CONTENT VENTILATION ACCESSORIES

1. Breathing circuits	09
1.1 Disposable breathing circuits/sets and related accessories	09
1.2 Reusable breathing circuits/sets and related accessories	39
2. Breathing system filters and HMEs	54
3. CO₂ cuvettes	62
4. Catheter mounts and airway connectors	63
5. Non-invasive patient interfaces and accessories	65
5.1 NIV masks and accessories	65
5.2 Anesthesia masks and accessories	72
5.3 Oxygen (O ₂) and aerosol therapy masks, venturi mask kits, nasal cannulas and accessories	79
5.4 High-Flow Therapy Nasal Cannula and accessories	82
6. Oxymixer	83
7. Device related accessories	85
8. Pack2Go®/Set2go	89
9. Nebulizer	92
10. Spare parts for Resuscitators	95

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

1. BREATHING CIRCUITS

1.1 DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

With VentStar® Family Dräger offers a wide range of disposable breathing circuits for different applications as well as for different patient groups (neo, ped, adult). Customers can choose between various kinds of materials. All disposable breathing circuits are made of high quality material and all gas leading components are PVC-free (no harm through softeners like e.g. DEHP). Special focus at Dräger is the interaction between Dräger device and related breathing circuit what is documented in our declarations of system compatibility. Declarations of system compatibility underlie intensive testing in various scenarios which are based on our long-time experience concerning interaction of device with accessory. With respect to the newest device generation we enlarged our portfolio with ID breathing circuits which support the Infinity® ID functionality of the relevant Dräger devices.

A	Breathing circuit Vent Set, disposable, basic, 1.5 m, 25 pcs., not made with natural rubber latex	MP00300
B	Breathing circuit VentStar®, disposable, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP00349
	Breathing circuit VentStar® ICU 180, disposable, basic, y-piece detachable 1.8 m, 25 pcs., not made with natural rubber latex	MP17009
C	Breathing circuit VentStar®, disposable, basic, 2.5 m, 10 pcs., not made with natural rubber latex	MP00350
	Breathing circuit VentStar® ICU 250, disposable, basic, y-piece detachable, 2.5 m, 10 pcs., not made with natural rubber latex	MP17010

ADULTS
PEDIATRICS
NEONATES



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES


<p>A Breathing circuit VentStar®, disposable, P, basic, w/o Luer Lock, 1.5 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00338</p>
<p>B Breathing circuit VentStar®, disposable, P, basic, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00351</p>
<p>Breathing circuit VentStar® ICU (P) 180, disposable, P, basic, y-piece detachable, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP17011</p>
<p>C Breathing circuit VentStar®, disposable, P, basic, 2.5 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00352</p>
<p>Breathing circuit VentStar® ICU (P) 250, disposable, P, basic, y-piece detachable, 2.5 m, 25 pcs., not made with natural rubber latex</p>	<p>MP17012</p>
<p>D Breathing circuit VentStar®, disposable, N, basic, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00353</p>



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES


<p>A Breathing circuit VentStar®, disposable, water trap, w/o Luer Lock, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00334</p>
<p>Breathing circuit VentStar® WT ICU 180, disposable, water trap, y-piece detachable, w/o Luer Lock, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP17008</p>
<p>B Breathing circuit VentStar®, disposable, water trap, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00337</p>
<p>C Breathing circuit VentStar®, disposable, P, water trap, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00361</p>
<p>D Breathing circuit VentStar®, disposable, P, water trap, w/o Luer Lock, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00362</p>
<p>Breathing circuit VentStar® WT ICU (P) 180, disposable, P, water trap, y-piece detachable, w/o Luer Lock, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP17011</p>
<p>E Breathing circuit VentStar®, disposable, N, water trap, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00363</p>

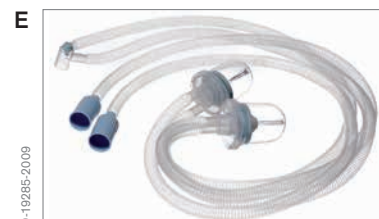

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES


<p>A Breathing circuit VentStar®, disposable, coaxial, 1.7 m, 20 pcs., not made with natural rubber latex</p>	<p>MP03373</p>
<p>B Breathing circuit VentStar®, disposable, coaxial, 1.9 m, 20 pcs., not made with natural rubber latex</p>	<p>MP03374</p>
<p>C Breathing circuit VentStar®, disposable, coaxial, 2.3 m, 20 pcs., not made with natural rubber latex</p>	<p>MP03376</p>
<p>D Breathing circuit VentStar®, disposable, coaxial, water trap, 2.1 m, 10 pcs., not made with natural rubber latex</p>	<p>MP00373</p>
<p>E Breathing circuit VentStar®, disposable, P, coaxial, 1.5 m, 20 pcs., not made with natural rubber latex</p>	<p>MP00379</p>
<p>Breathing circuit VentStar®, disposable, coaxial, GSL, 2.3 m, 20 pcs., not made with natural rubber latex</p>	<p>MP03378</p>
<p>Breathing circuit VentStar®, disposable, coaxial, water trap, GSL, 2.1 m, 10 pcs., not made with natural rubber latex</p>	<p>MP03379</p>

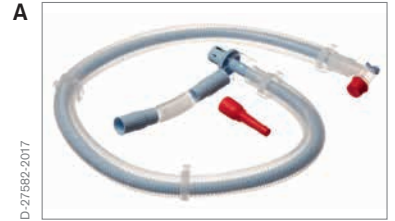

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

<p>A Breathing circuit Vent Set, disposable, extendable, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00305</p>
<p>B Breathing circuit VentStar®, disposable, extendable, 2.2 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00355</p>
<p>C Breathing circuit VentStar®, disposable, smoothbore, for MRI, 3 m, 10 pcs., not made with natural rubber latex</p>	<p>MP04800</p>
<p>D Breathing circuit VentStar®, disposable, N, MRI, 3 m, 10 pcs., not made with natural rubber latex</p>	<p>MP04802</p>

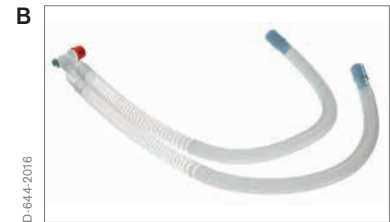

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES


<p>A Breathing bag set VentStar®, disposable, 2 L, 0.8 m, 25 pcs., not made with natural rubber latex</p>	<p><u>MP00222</u></p>
<p>Breathing bag set VentStar® (Cone 22OD), disposable, 2 L, 0.8 m, 25 pcs., not made with natural rubber latex</p>	<p><u>MP16017</u></p>
<p>B Breathing bag set VentStar®, disposable, 2 L, 1.2 m, 25 pcs., not made with natural rubber latex</p>	<p><u>MP00380</u></p>
<p>Breathing bag set VentStar® (Cone 22OD), disposable, 2 L, 1.2 m, 25 pcs., not made with natural rubber latex</p>	<p><u>MP16018</u></p>
<p>B Breathing bag set VentStar®, disposable, 2 L, 1.5 m, 25 pcs., not made with natural rubber latex</p>	<p><u>MP00381</u></p>
<p>Breathing bag set VentStar® (Cone 22OD), disposable, 2 L, 1.5 m, 25 pcs., not made with natural rubber latex</p>	<p><u>MP16019</u></p>



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

A	Breathing bag set VentStar®, disposable, P, 1 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP00384
	Breathing bag set VentStar® (Cone 22OD), disposable, P, 1 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP16020
B	Breathing bag set VentStar®, disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP00383
	Breathing bag set VentStar® (Cone 22OD), disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP16021
C	① VentStar® Breathing Bag 3 L, 25 pcs.	MP02713
	② VentStar® Breathing Bag 2 L, 25 pcs.	MP02712
	③ VentStar® Breathing Bag 1 L, 25 pcs.	MP02711
	④ VentStar® Breathing Bag 0.5 L, 25 pcs.	MP02710

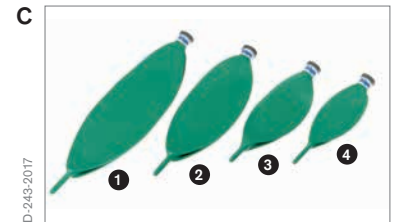

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

<p>A Anesthesia set, disposable, basic, 1.5 m/0.8 m, 25 pcs., not made with natural rubber latex</p> <hr/> <p>Anesthesia set, disposable, basic, 1.5 m/0.8 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex</p> <hr/> <p>B Anesthesia set, disposable, extendable, 1.8 m/0.8 m, 25 pcs., not made with natural rubber latex</p> <hr/> <p>Anesthesia set, disposable, extendable, 1.8 m/0.8 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex</p> <hr/> <p>C Anesthesia set VentStar®, disposable, w/o Luer Lock, 1.8 m/0.8 m, 25 pcs., not made with natural rubber latex</p> <hr/> <p>Anesthesia set VentStar®, disposable, w/o Luer Lock, 1.8 m/0.8 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex</p> <hr/>	<p>MP00301</p> <hr/> <p>MP16028</p> <hr/> <p>MP00303</p> <hr/> <p>MP16029</p> <hr/> <p>MP00330</p> <hr/> <p>MP16030</p> <hr/>
---	---



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

<p>A Anesthesia set VentStar®, disposable, P, w/o Luer Lock, 1.5 m/1.1 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00331</p>
<p>Anesthesia set VentStar®, disposable, P, w/o Luer Lock, 1.5 m/1.1 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex</p>	<p>MP16031</p>
<p>B Anesthesia set VentStar®, disposable, N, 1.8 m/1.1 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00333</p>
<p>Anesthesia set VentStar®, disposable, N, 1.8 m/1.1 m, Bag cone 22OD, 25 pcs., not made with natural rubber latex</p>	<p>MP16032</p>



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

<p>A Anesthesia set VentStar®, disposable, basic, 2 L, 1.8 m/1.5 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00365</p>
<p>Anesthesia set VentStar®, disposable, basic, 2 L, 1.8 m/1.5 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex</p>	<p>MP16033</p>
<p>B Anesthesia set VentStar®, disposable, basic, 2 L, 2.5 m/1.5 m, 10 pcs., not made with natural rubber latex</p>	<p>MP00366</p>
<p>Anesthesia set VentStar®, disposable, basic, 2 L, 2.5 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex</p>	<p>MP16034</p>
<p>C Anesthesia set VentStar®, disposable, water trap, w/o Luer Lock, 2 L, 1.8 m/1.5 m, 10 pcs., not made with natural rubber latex</p>	<p>MP00371</p>
<p>Anesthesia set VentStar®, disposable, water trap, w/o Luer Lock, 2 L, 1.8 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex</p>	<p>MP16035</p>



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

<p>A Anesthesia set VentStar®, disposable, water trap, 2 L, 1.8 m/1.5 m, 10 pcs., not made with natural rubber latex</p>	<p>MP00372</p>
<p>Anesthesia set VentStar®, disposable, water trap, 2 L, 1.8 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex</p>	<p>MP16036</p>
<p>B Anesthesia set VentStar®, disposable, water trap, w/o Luer Lock, 2 L, 2.8 m/1.5 m, 10 pcs., not made with natural rubber latex</p>	<p>MP02720</p>
<p>Anesthesia set VentStar®, disposable, water trap, w/o Luer Lock, 2 L, 2.8 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex</p>	<p>MP16039</p>
<p>C Anesthesia set VentStar®, disposable, water trap, 2 L, 2.8 m/1.5 m, 10 pcs., not made with natural rubber latex</p>	<p>MP02721</p>
<p>Anesthesia set VentStar®, disposable, water trap, 2 L, 2.8 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex</p>	<p>MP16040</p>



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES


<p>A Anesthesia set VentStar®, disposable, P, water trap, 1 L, 1.8 m/1.1 m, 10 pcs., not made with natural rubber latex</p>	<p><u>MP00374</u></p>
<p>Anesthesia set VentStar®, disposable, P, water trap, 1 L, 1.8 m/1.1 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex</p>	<p><u>MP16037</u></p>
<p>B Anesthesia set VentStar®, disposable, P, water trap, w/o Luer Lock, 1 L, 1.8 m/1.1 m, 10 pcs., not made with natural rubber latex</p>	<p><u>MP00375</u></p>
<p>Anesthesia set VentStar®, disposable, P, water trap, w/o Luer Lock, 1 L, 1.8 m/1.1 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex</p>	<p><u>MP16038</u></p>



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES


A	Anesthesia set VentStar®, disposable, coaxial, water trap, 2 L, 2.1 m/1.2 m, 10 pcs., not made with natural rubber latex	MP03381
	Anesthesia set VentStar®, disposable, coaxial, water trap, 2 L, 2.1/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16044
	Anesthesia set VentStar®, disposable, coaxial, 2 L, 2.3 m/1.2 m, 10 pcs., not made with natural rubber latex	MP03383
	Anesthesia set VentStar, disposable, coaxial, 2 L, 2.3 m/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16042
B	Anesthesia set VentStar®, disposable, coaxial, 2 L, 1.9 m/1.2 m, 10 pcs., not made with natural rubber latex	MP03384
	Anesthesia set VentStar, disposable, coaxial, 2 L, 1.9 m/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16041
C	Anesthesia set VentStar®, disposable, coaxial, GSL, 2 L, 2.3 m /1.2 m, 10 pcs., not made with natural rubber latex	MP00358
	Anesthesia set VentStar®, disposable, coaxial, GSL, 2 L, 2.3 m /1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16049
D	Anesthesia set VentStar®, coaxial, water trap, GSL, 2 L, 2.1 m /1.2 m, 10 pcs., not made with natural rubber latex	MP00359
	Anesthesia set VentStar®, coaxial, water trap, GSL, 2 L, 2.1 m /1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16050

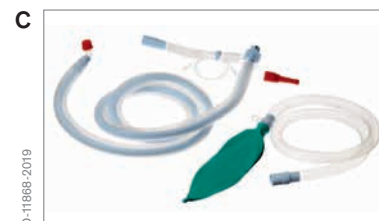

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

- A Breathing circuit VentStar®, disposable, N, heated, 1.2 m, 10 pcs., not made with natural rubber latex** **MP00308**
- B Breathing circuit VentStar®, disposable, N, heated, w/o chamber, 1.2 m, 10 pcs., not made with natural rubber latex** **MP00309**
- C Humidifier chamber VentStar®, disposable, 10 pcs.** **MP00234**

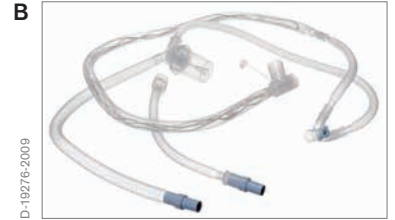

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES


A	Breathing circuit VentStar® Helix, disposable, N, dual heated, 1.7 m, 10 pcs., not made with natural rubber latex	MP02650
B	Breathing circuit VentStar® Helix, disposable, N, heated, 1.7 m, 10 pcs., not made with natural rubber latex	MP02608
C	Humidifier chamber VentStar®, disposable, N, 10 pcs.	MP02590
D	Double connector 11ID/15OD, disposable, 10 pcs.	MP02655
E	T-adapter neonatal for NO therapy, 10 pcs.	8422215
F	Aerogen T-Adapter, 11.5 mm, neonates, 10 pcs.	MP02689



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

A <u>VentStar® Helix, dual heated, 10 pcs.</u>	<u>MP02606</u>
B <u>VentStar® Helix, heated, 10 pcs.</u>	<u>MP02607</u>
C <u>Humidifier chamber, 10 pcs.</u>	<u>MP02605</u>

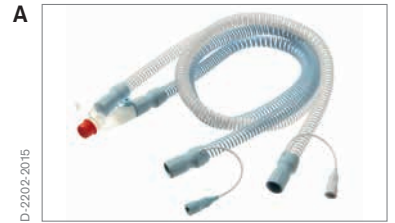

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

GAS MEASUREMENT ACCESSORIES

A	Sample line, disposable, 10 pcs.	8290286
B	Watertrap WaterLock® 2, 12 pcs.	6872130
	Water trap Infinity® ID WaterLock® 2, 12 pcs.	6872020
	Sample gas return line for Fabius®, with filter	M32692
C	Bacterial filter for Dräger patient monitors, disposable	8402868



Aquapor H300

respiratory humidifier for optimally conditioned breathing gas




DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

	<u>Aquapor H300 – Configuration</u>	<u>MP01000</u>
A	<u>Aquapor H300 Basic unit 230V</u>	<u>MP07101</u>
	<u>Aquapor H300 Basic unit 115V</u>	<u>MP07102</u>
	<u>VentStar® Aquapor dual heated, Adult and Paediatric, 10 pcs.</u>	<u>MP07130</u>
B	<u>VentStar® Aquapor, heated, Adult and Paediatric, 10 pcs.</u>	<u>MP07131</u>
C	<u>VentStar® Aquapor, dual heated, N plus, 10 pcs.</u>	<u>MP07150</u>
D	<u>VentStar® Aquapor, heated, N plus, 10 pcs.</u>	<u>MP07151</u>
	<u>HI-Flow Star Kit Aquapor A, heated inspiratory breathing circuit for the H300 humidifier, heated inspiratory hose, humidifier chamber, breathing hose when used with a ventilator, O₂ hose when used with an O₂ blender, bed clip, disposable, 10 pcs.</u>	<u>MP17090</u>
E	<u>HI-Flow Star Valve Kit w/ Pop-Up Valve and Connector 22 mm, disposable, 15 pcs.</u>	<u>MP05507</u>
F	<u>H300 Heater wire adapter, insp</u>	<u>MP07106</u>
G	<u>H300 Heater wire adapter, dual</u>	<u>MP07107</u>
	<u>H300 Temp. sensor cable*, 160 cm</u>	<u>MP07108</u>
	<u>H300 Temp. sensor cable</u>	<u>MP07109</u>

*Tradeware

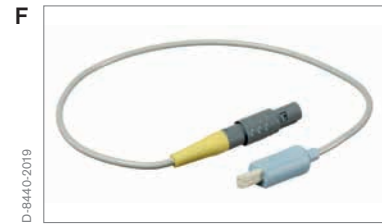

Not all articles are available worldwide.


OVERVIEW VENTILATION SODA LIME SENSORS MONITORING WARD EQUIPMENT RAIL ACCESSORIES NEONATAL CARE INDEX ORDER NO. INDEX PRODUCT

FISHER & PAYKEL HUMIDIFIER AND RELATED ACCESSORIES

FISHER & PAYKEL HUMIDIFIER AND RELATED ACCESSORIES

	F&P configuration for MR810, MR850 and F&P 950	8414144
A	F&P 950 sensor cartridge	MP16002
B	F&P 950 heater wire adaptor	MP16003
C	F&P 950 A80 breathing circuit (A)	MP16004
D	F&P 950 neonatal disposable breathing circuit	MP16005
E	Humidifier chamber* for F&P MR 290, disposable, 10 pcs.	8418282
F	Humidifier chamber* for F&P MR 340, reusable, N/P	8411047
G	Humidifier chamber* for F&P MR 370, reusable, A	8412217
	Reusable Breathing circuit (N), 900MR782*	MP17128
	Reusable Breathing circuit (A), 900MR784*	MP17129

*Tradeware

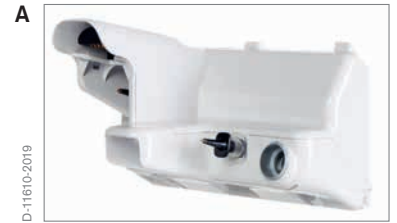

FISHER & PAYKEL HUMIDIFIER AND RELATED ACCESSORIES

HUMIDIFIER MOUNTING MATERIAL

A	Rail clamp* for F&P humidifier and Seattle PAP bubble box	8411074
B	Holder* for F&P humidifier U-shaped, for trolley, with rail clamp	MP01002
C	Water bag holder* 900 MR 290	MP01008
	Holder, lateral (F&P 900 MR 309), including 8411074	MP01009
	Universal holder with DIN rail	G93140
	Humidifier holder for side rail	8416325
	Humidifier holder chassis, swiveling	G93111

FISHER & PAYKEL HOSE HEATER

D	Hose heater* for F&P MR 700/730 and MR 850, for reusable breathing circuit, 0.7 m	8411154
D	Hose heater* for F&P MR 700/730 and MR 850, for reusable breathing circuit, 1.1 m (Low -Flow)	8411045
E	Hose heater* for F&P MR 700/730 and MR 850, for reusable breathing circuit, 1.1 m	8414993
F	Hose heater* for F&P MR 700/730 and MR 850, for reusable breathing circuit, 1.5 m	8414969

*Tradeware



FISHER & PAYKEL HUMIDIFIER AND RELATED ACCESSORIES

	Temperature probe/flow sensor* for F&P MR 850, for HFV reusable breathing circuits, 0.7 m	8418507
A	Temperature probe/flow sensor* for F&P MR 850, for all breathing circuits, 1.1 m	8414989
B	Temperature probe/flow sensor* for F&P MR 850, for all breathing circuits, 1.5 m	8414966

*Tradeware



FISHER & PAYKEL HUMIDIFIER AND RELATED ACCESSORIES

A	Hose heater adapter* for F&P MR 850, for reusable breathing circuit	8414968
B	Hose heater adapter for F&P MR 850, for single-heated disposable breathing circuit	8414992
C	Hose heater adapter for F&P MR 850, for dual-heated disposable breathing circuit	MP01006
D	Draw wire* for hose heater 1.5 m, for reusable breathing circuit	8411050
E	Temperature probe adapter* 900 MR 532, for F&P humidifiers	MP01007

*Tradeware



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

FOR OXYLOG® 2000+, 3000, 3000+ AND VE300 ADVANCED (WITH FLOW OPTION)

- A Breathing circuit VentStar®**
 Oxylog 3000F, 2000p, VE300, disposable, 1.5 m (4.92 ft), 5 pcs., for Oxylog VE300 only with Option Plus **5703041**
- B Breathing circuit VentStar®**
 Oxylog 3000F, 2000p, VE300, disposable, 3 m (9.8 ft), 5 pcs., for Oxylog VE300 only with Option Plus **MP00335**

FOR OXYLOG® 1000

- C Breathing circuit VentStar® Oxylog 1000,**
 disposable, 1.8 m, 25 pcs., not made with natural rubber latex **2M86837**



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

FOR OXYLOG® 3000: ≥ SW 1.20, 3000+ AND VE300 ADVANCED (WITH FLOW OPTION)

- | | |
|--|-----------------------|
| <p>A Breathing circuit VentStar® Oxylog 3000, 3000+, disposable, P, 1.9 m, 5 pcs., Oxylog 3000: ≥ SW 1.20, VE300 (Software ≥ 2.0, Pediatric Option), not made with natural rubber latex</p> | <p>5704964</p> |
|--|-----------------------|

FOR OXYLOG® VE300 (w/o FLOW OPTION)

- | | |
|---|-----------------------|
| <p>B Breathing circuit VentStar® Oxylog VE300, not for VE300 with Option plus, disposable, 1.5 m, 5 pcs., not made with natural rubber latex</p> | <p>MP01370</p> |
| <p>C Breathing circuit VentStar® Oxylog VE300, disposable, 3 m, 5 pcs., not made with natural rubber latex</p> | <p>MP01371</p> |

ADULTS
PEDIATRICS
NEONATES




DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

FOR CARINA®

<p>A Breathing circuit VentStar® Carina® LeakV, disposable, 1.5 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00312</p>
<p>B Breathing circuit VentStar® Carina® ExpV, disposable, 1.5 m, 10 pcs., not made with natural rubber latex</p>	<p>MP00313</p>
<p>C Breathing circuit VentStar® Carina® ExpV, disposable, 2 m, 10 pcs., not made with natural rubber latex</p>	<p>5704962</p>

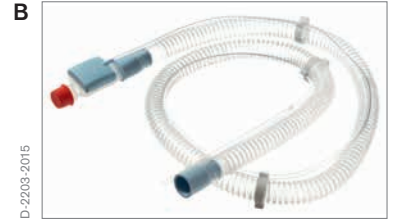

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

A	Breathing circuit Infinity® ID, disposable, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP01350
B	Breathing circuit Infinity® ID, disposable, basic, 2.5 m, 25 pcs., not made with natural rubber latex	MP01348
C	Breathing circuit Infinity® ID, disposable, smoothbore, 1.6 m, 25 pcs., not made with natural rubber latex	MP04803
D	Breathing circuit Infinity® ID, disposable, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	MP01326
E	Breathing circuit Infinity® ID, disposable, smoothbore, 3 m, 10 pcs., not made with natural rubber latex	MP04804
F	Breathing circuit Infinity® ID, disposable, extendable, 2.2 m, 25 pcs., not made with natural rubber latex	MP01347



Not all articles are available worldwide.

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

A	Breathing circuit Infinity® ID, disposable, coaxial, 1.9 m, 20 pcs., not made with natural rubber latex	MP03375
B	Breathing circuit Infinity® ID, disposable, coaxial, 2.3 m, 20 pcs., not made with natural rubber latex	MP03377
B	Breathing circuit Infinity® ID, disposable, coaxial, GSL, 2.3 m, 20 pcs., not made with natural rubber latex	MP00357
C	Breathing circuit Infinity® ID, disposable, water trap, 2.1 m, 10 pcs., not made with natural rubber latex	MP00378
C	Breathing circuit Infinity® ID, disposable, water trap, GSL, 2.1 m, 10 pcs., not made with natural rubber latex	MP00367
	Anesthesia set Infinity® ID, disposable, coaxial, water trap, 2 L, 2.1/1.2 m, 10 pcs., not made with natural rubber latex	MP03388
	Anesthesia set Infinity® ID, disposable, coaxial, water trap, 2 L, 2.1/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16045
	Anesthesia set Infinity® ID, disposable, coaxial, GSL, 2.3 m/1.2 m, 10 pcs., not made with natural rubber latex	MP00368
	Anesthesia set Infinity ID®, disposable, coaxial, GSL, 2.3 m/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16051
	Anesthesia set Infinity® ID, disposable, coaxial, water trap, GSL, 2.1 m/1.2 m, 10 pcs., not made with natural rubber latex	MP00370
	Anesthesia set Infinity ID®, disposable, coaxial, water trap, GSL, 2.1 m/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16052


ADULTS
PEDIATRICS
NEONATES



DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

<p>A Breathing circuit Infinity® ID, disposable, P, basic, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP01340</p>
<p>B Breathing circuit Infinity® ID, disposable, P, water trap, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP01341</p>
<p>C Breathing circuit Infinity® ID, disposable, N, basic, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP01327</p>
<p>Breathing circuit Infinity® ID, disposable, N, water trap, 1.8 m, 25 pcs., not made with natural rubber latex</p>	<p>MP01328</p>

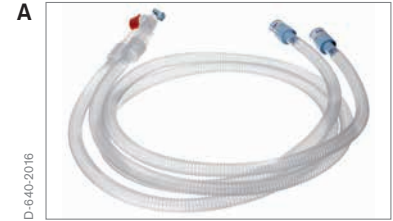

DISPOSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

ADULTS
PEDIATRICS
NEONATES

<p>A Breathing bag set Infinity® ID, disposable, 2 L, 1.2 m, 25 pcs., not made with natural rubber latex</p>	<p>MP00341</p>
<p>B Breathing bag set Infinity® ID, disposable, 2 L, 1.5 m, 25 pcs., not made with natural rubber latex</p>	<p>MP01330</p>
<p>C Breathing bag set Infinity® ID, disposable, P, 1 L, 1.1 m, 25 pcs., not made with natural rubber latex</p>	<p>MP01343</p>
<p>D Breathing bag set Infinity® ID, disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex</p>	<p>MP01329</p>



REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

1.2 REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES


REUSABLE BREATHING SYSTEMS: BLUESET VENTILATION

A	BlueSet Ventilation (N)	8412344
Consists of:		
1	Y-piece, reusable, N, 90°, with suction opening	8403075
2	Corrugated hose, reusable, N/P, 14 cm, connector 10 mm/10 mm, Ø 10 mm	8403333
3	Temperature sensor mount, reusable, N	8403056
4	Water trap with outlet, reusable N/P	8404760
5	Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/10 mm, Ø 10 mm	2165848
6	Water trap, reusable, N/P	8409627
7	Red cap, 7 mm	8401645
8	Breathing hose, silicone, reusable, N/P, 40 cm, connector 10 mm/22 mm, Ø 10 mm	2165856
9	Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/22 mm, Ø 10 mm	2165821
	Dummy Plug	8410049




REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

A	BlueSet Ventilation (P/A)	8413146
	Consists of:	
1	Adapter, reusable, 22 mm/22 mm	M25647
2	Adapter, reusable, 22 mm/15 mm, 90°	8290274
3	Hose clamp, 10 pcs.	8403566
4	Y-piece with sensor mount, reusable, A	8405435
5	Breathing hose, silicone, reusable, A, 35 cm, connector 22 mm/22 mm, Ø 22 mm	2165619
6	Breathing hose, silicone, reusable, A, 60 cm, connector 22 mm/22 mm, Ø 22 mm	2165627
7	Water trap, reusable, A	8404985
	Dummy Plug	8410049

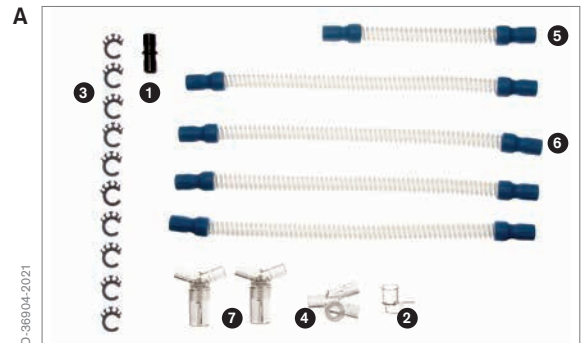

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES


A	BlueSet Ventilation HME (P/A)	8412860
	Consists of:	
1	Breathing hose, silicone, reusable, A, 1.5 m, connector 22 mm/22 mm, Ø 22 mm	2166046
2	Y-piece, polyamide, reusable, straight, A	M25650

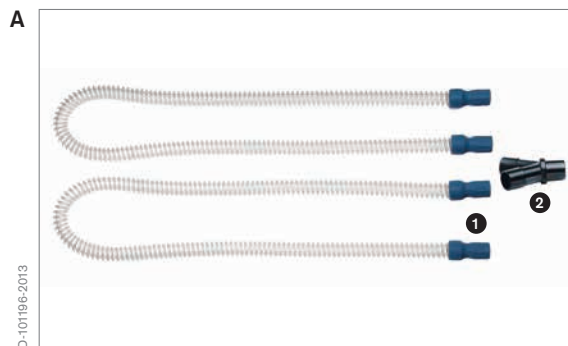

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES


REUSABLE ANESTHESIA SYSTEMS: BLUESET ANESTHESIA

A	BlueSet Anesthesia (P/A)	M30146
	Consists of:	
1	Breathing hose, silicone, reusable, A, 1.5 m, connector 22 mm/22 mm, Ø 22 mm	2166046
2	Breathing hose, silicone, reusable, A, 1.1 m, connector 22 mm/22 mm, Ø 22 mm	2166038
3	Breathing bag, silicone, reusable, 2.3 L, not made with natural rubber latex	2166062
4	Y-piece 90° with Luer Lock, reusable, A	M33278
5	Adapter, reusable, 22 mm/22 mm	M25647

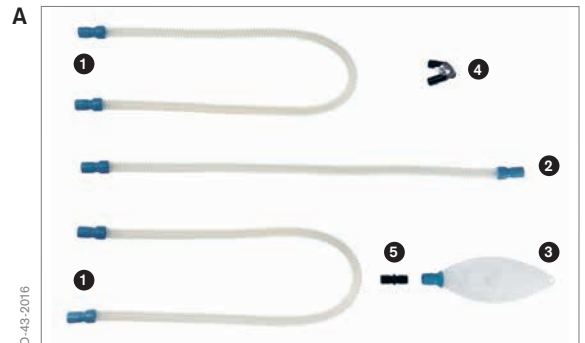


REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

A	BlueSet Anesthesia water trap (P/A)	MK04248
	Consists of:	
1	Breathing hose, silicone, reusable, A, 60 cm, connector 22 mm/22 mm, Ø 22 mm	2165627
2	Breathing hose, silicone, reusable, A, 1.1 m, connector 22 mm/22 mm, Ø 22 mm	2166038
3	Adapter, reusable, 22 mm/22 mm	M25647
4	Water trap, reusable, A	8404985
5	Y-piece 90° with Luer Lock, reusable, A	M33278
6	Breathing bag, silicone, reusable, 2.3 L, not made with natural rubber latex	2166062

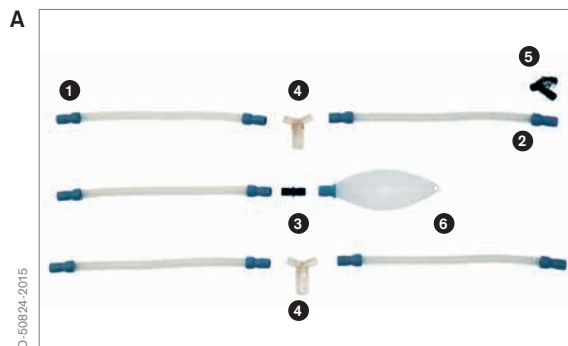

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES


A	BlueSet Anesthesia (N)	M27542
	Consists of:	
①	Breathing hose, silicone, reusable, N/P, 1.5 m, connector 10 mm/22 mm, Ø 10 mm	2165961
②	Breathing hose, silicone, reusable, N/P, 1.1 m, connector 22 mm/22 mm, Ø 10 mm	2165678
③	Breathing bag, silicone, reusable, 0.5 L, not made with natural rubber latex	2166054
④	Adapter, reusable, 22 mm/22 mm	M25647
⑤	Y-piece, reusable, N, 90°, with suction opening	8403075
⑥	Elbow, reusable, N/P, metal, Luer Lock	M34568
⑦	Set of 5 red caps, 7 mm	8402953

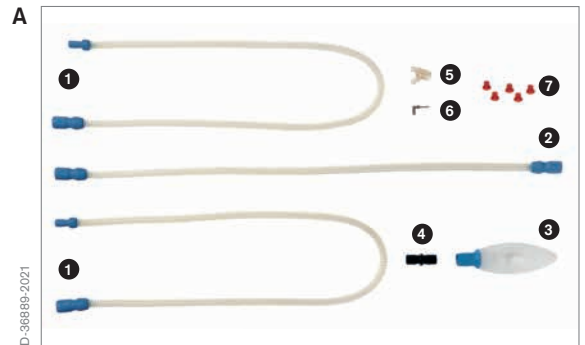

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

A	Breathing hose, silicone, reusable, A, 1.5 m, connector 22 mm/22 mm, Ø 22 mm	2166046
B	Ventilation valve for Oxylog® 1000, reusable	2M86800
C	PEEP valve for Oxylog® 1000, reusable, adjustable from 0 – 10 mbar	8407475
D	PEEP valve for Oxylog® 1000*, disposable, adjustable from 5 – 20 mbar	2M86832
	Breathing hose for Oxylog® 1000, reusable, silicone, 3 m	2112760

*Tradeware

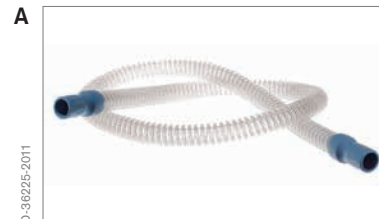

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

FOR OXYLOG® 2000, 2000+, 3000, 3000+ AND VE300 ADVANCED (WITH FLOW OPTION)

A Breathing hose set for Oxylog® 2000/3000 family, reusable, w/ flow measurement line, 1.5 m 8412068

Consists of:

Breathing hose, silicone, reusable, A, 1.5 m, connector 22 mm/22 mm, Ø 22 mm 2166046

Flow measuring line for Oxylog® 2000/3000 family, 1.5 m 8410929

B Breathing hose set for Oxylog® 2000/3000 family, reusable, w/ flow measurement line, 3 m 8412913

Consists of:

Breathing hose for Oxylog® 1000, reusable, silicone, 3 m 2112760

Flow measuring line for Oxylog® 2000/3000 family, 3 m 8413293

C Ventilation valve for Oxylog® 2000/3000, reusable 8412001

D Pressure difference flow sensor, PPSU, sterilizable 8412034

E Elbow, reusable 8412235



FOR OXYLOG® VE300

Hose Oxylog® VE300, 150, reusable MP01390

Hose Oxylog® VE300, 300, reusable MP01391

Connector pressure Oxylog® VE300 MP01372


E Elbow, reusable 8412235

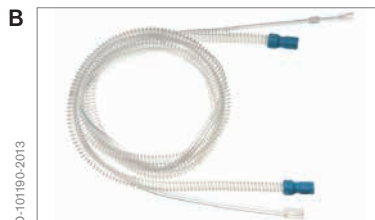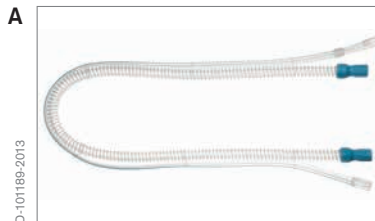

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

Silicone breathing hose (adult)
with standard sleeves 22 mm/22 mm (0.9 in/0.9 in)

A	Breathing hose, silicone, reusable, A, 35 cm, connector 22 mm/22 mm, Ø 22 mm	2165619
B	Breathing hose, silicone, reusable, A, 60 cm, connector 22 mm/22 mm, Ø 22 mm	2165627
C	Breathing hose, silicone, reusable, A, 1.1 m, connector 22 mm/22 mm, Ø 22 mm	2166038
D	Breathing hose, silicone, reusable, A, 1.5 m, connector 22 mm/22 mm, Ø 22 mm	2166046

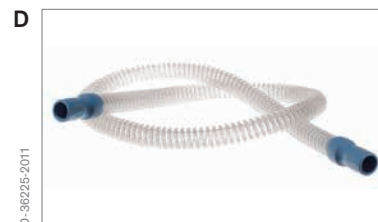


REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

Silicone breathing hose (pediatric)
with standard sleeves 10 mm (0.4 in) on both sides
and 10 mm (0.4 in) inner diameter

**A Breathing hose, silicone, reusable, N/P, 60 cm,
connector 10 mm/10 mm, Ø 10 mm** **2165848**

Silicone breathing hose (pediatric)
with standard sleeves 22 mm (0.9 in) on both sides
and 10 mm (0.4 in) inner diameter



**B Breathing hose, silicone, reusable, N/P, 1.1 m,
connector 22 mm/22 mm, Ø 10 mm** **2165678**

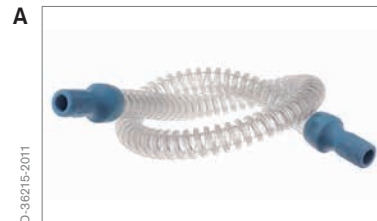

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

Silicone breathing hose (pediatric) with standard sleeves with one 10 mm (0.4 in) and one 22 mm (0.9 in) sleeve and 10 mm (0.4 in) inner diameter

A	Breathing hose, silicone, reusable, N/P, 40 cm, connector 10 mm/22 mm, Ø 10 mm	2165856
B	Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/22 mm, Ø 10 mm	2165821
C	Breathing hose, silicone, reusable, N/P, 1.1 m, connector 10 mm/22 mm, Ø 10 mm	2165651
D	Breathing hose, silicone, reusable, N/P, 1.5 m, connector 10 mm/22 mm, Ø 10 mm	2165961

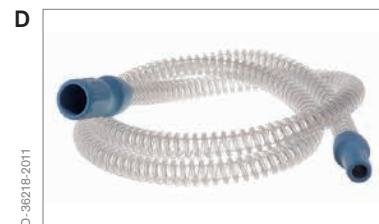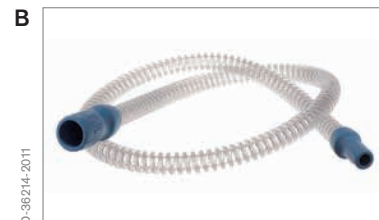

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

A	Y-piece, polyamide, reusable, straight, A	M25650
B	Y-piece 90° with Luer Lock, reusable, A	M33278
C	Y-piece for F&P, reusable, A*	8412888
D	Y-piece with sensor mount, reusable, A	8405435
E	Y-piece straight, reusable, N/P	M27077
F	Y-piece, reusable, N, 90°, with suction opening	8403075
G	Temperature probe adapter 900 MR 532, for F&P humidifiers*	MP01007


*Tradeware



REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

A	Elbow, reusable	8412235
B	Elbow for F&P, reusable*	8418289
C	Adapter, reusable, 22 mm/22 mm	M25647
	Adapter straight, reusable, A, 18 mm/11 mm	8403076
	Adapter, reusable, 22 mm/22 mm, 90°	M25649
D	Adapter, reusable, 22 mm/11 mm, N/P	8409897
E	Elbow, reusable, N/P, metal, Luer Lock	M34568
	Temperature sensor mount, reusable, N	8403056
	T-piece connector, reusable	8411044
F	Water trap, reusable, A	8404985
G	Water trap, reusable, N/P	8409627
	Red cap, 7 mm	8401645


*Tradeware

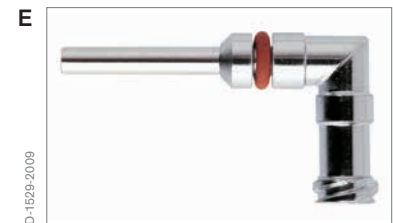

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

Transparent breathing bags of silicone elastomer

A	① Breathing bag, silicone, reusable, 0.5 L, not made with natural rubber latex	<u>2166054</u>
	② Breathing bag, silicone, reusable, 1.5 L, not made with natural rubber latex	<u>2166127</u>
	③ Breathing bag, silicone, reusable, 2.3 L, not made with natural rubber latex	<u>2166062</u>
	④ Breathing bag, silicone, reusable, 3 L, not made with natural rubber latex	<u>2166070</u>

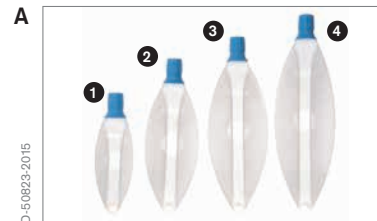

REUSABLE BREATHING CIRCUITS/SETS AND RELATED ACCESSORIES

A	Test lung SelfTestLung, latex-free, reusable	MP20020
B	Test lung SelfTestLung, reusable, max. 1 L, not made with natural rubber latex	MP02400
C	Demo lung SmartLung, max.1 L, reusable, not made with natural rubber latex	MP02405
D	Test Lung, reusable, max. 1.5 L, not made with natural rubber latex	8403201
E	Neonatal Test Lung, reusable, 50 ml, not made with natural rubber latex	MP20050



Filter/HME

avoiding nosocomial infections and support of lung protective ventilation

Our new filter portfolio

High performance filtration rates
to reduce risk of nosocomial infections

Writable surface
to facilitate clinical workflows

Performance tested
in accordance with Dräger standards and accredited laboratories

Clear color coding
to quickly identify right filter

Excellent moistening performance
to support lung protective ventilation

D-4013-2021

To be able to rely on the quality of medical equipment is a needed prerequisite to fully concentrate on the application of therapies to take care of patient recovery and wellbeing. We provide you with high quality filters for a wide range of applications and patient groups.

BREATHING SYSTEM FILTERS AND HMEs

2. BREATHING SYSTEM FILTERS AND HMEs

ELECTROSTATIC FILTER AND HME

The TwinStar® Plus filter/HME from Dräger is a combination of an efficient HME and a high-performance breathing system filter, thereby significantly contributing to infection prophylaxis in ventilation treatment.

- Combination of efficient HME and high-performance breathing system filter
- High retention rates for bacteria and viruses
- High return of moisture
- Minimal workload and low cost
- Sampling connector with tethered cap
- Transparent housing for visual control
- Standardized connectors for safe connection to other components
- Clear labeling and blue color coding for quick identification

A Filter + HME TwinStar® 90 Plus, disposable, ~90 ml deadspace, 100 pcs. **MP05800**

B Filter + HME TwinStar® 60A Plus, disposable, ~60 ml deadspace, angled, 100 pcs. **MP05810**

C Filter + HME TwinStar® 55 Plus, disposable, ~55 ml deadspace, 100 pcs. **MP05805**

D Filter + HME TwinStar® 25 Plus, disposable, ~25 ml deadspace, 100 pcs. **MP05815**

E Filter + HME TwinStar® 9 Plus, disposable, ~9 ml deadspace, 100 pcs. **MP05820**

MECHANICAL FILTER AND HME

F Mechanical filter + HME TwinStar® HEPA Plus, disposable, ~90 ml deadspace, 100 pcs. **MP05801**



BREATHING SYSTEM FILTERS

ELECTROSTATIC FILTER

The CareStar® Plus breathing system filter significantly contributes to infection prophylaxis in anesthesia and ventilation treatment through high retention rates.

- High-performance electrostatic filter
- Patient-side or device-side use
- Minimal workload and low cost
- Sampling connector with tethered cap
- Transparent housing for visual control
- Standardized connectors for safe connection to other components
- Clear labeling and yellow color coding for quick identification

A Electrostatic filter CareStar® 35 Plus, disposable, ~35 ml deadspace, 100 pcs.

MP05755

B Electrostatic filter CareStar® 20 Plus, disposable, ~20 ml deadspace, 100 pcs.

MP05770



BREATHING SYSTEM FILTERS

MECHANICAL FILTER

The SafeStar® Plus mechanical high performance breathing system filter shows excellent retention rates, thereby significantly contributing to infection prophylaxis in anesthesia and ventilation treatment.

- High-performance mechanical HEPA filter
- High hydrophobicity
- Patient-side or device-side use
- Sampling connector with tethered cap
- Transparent housing for visual control
- Standardized connectors for safe connection to other components
- Clear labeling and red color coding for quick identification

A	Mechanical filter SafeStar® 90 Plus, disposable, ~90 ml deadspace, 100 pcs.	MP05785
B	Mechanical filter SafeStar® 60A Plus, disposable, ~60 ml deadspace, angled, 100 pcs.	MP05795
C	Mechanical filter SafeStar® 55 Plus, disposable, ~55 ml deadspace, 100 pcs.	MP05790



HMEs

HME

HumidStar® Plus moisture exchangers from Dräger partially take over the function of the upper airways and efficiently humidify the inspired air.

- High return of moisture
- Minimal workload and low cost
- Sampling connector with tethered cap
- Transparent housing for visual control
- Standardized connectors for safe connection to other components
- Clear labeling and green color coding for quick identification

A	HME HumidStar® 55 Plus*, disposable, ~55 ml deadspace, 100 pcs.	MP05730
B	HME HumidStar® 25 Plus*, disposable, ~25 ml deadspace, 100 pcs.	MP05735
C	HME HumidStar® 2 Plus, disposable, 2 ml deadspace, 100 pcs.	MP05845
D	HME HumidStar® Trach Plus, disposable, 6 ml deadspace, 100 pcs.	MP05750



BREATHING SYSTEM FILTERS AND HMEs

SPECIFICATIONS

	TwinStar® 90 Plus	TwinStar® 60A Plus	TwinStar® 55 Plus
Part number	MP05800	MP05810	MP05805
Deadspace	90 ml	60 ml	55 ml
Recommended tidal volume	300 to 1500 ml	300 to 1500 ml	300 to 1500 ml
Bacterial retention**	≥ 99.99%	≥ 99.99%	≥ 99.99%
Viral retention**	≥ 99.9%	≥ 99.9%	≥ 99.9%
Filtration method	electrostatic	electrostatic	electrostatic
Moisture loss*	≤ 5.9 mg/l at VT = 500 ml	≤ 6.3 mg/l at VT = 500 ml	≤ 9.4 mg/l at VT = 500 ml
Moisture return	≥ 38.1 mg/l at VT = 500 ml	≥ 37.7 mg/l at VT = 500 ml	≥ 34.6 mg/l at VT = 500 ml
Resistance	≤ 1.0 mbar at 30 l/min	≤ 1.3 mbar at 30 l/min	≤ 1.3 mbar at 30 l/min
Maximum duration of use	24 hours	24 hours	24 hours
Product	PVC-free/ Latex free	PVC-free/ Latex free	PVC-free/ Latex free

	TwinStar® 25 Plus	TwinStar® 9 Plus	TwinStar® HEPA Plus
Part number	MP05815	MP05820	MP05801
Deadspace	25 ml	9 ml	90 ml
Recommended tidal volume	100 to 500 ml	30 to 150 ml	300 to 1500 ml
Bacterial retention**	≥ 99,98%	≥ 99.99%	≥ 99.9999%
Viral retention**	≥ 99.9%	≥ 99.9%	≥ 99.9999%
Filtration method	electrostatic	electrostatic	mechanical
Moisture loss*	≤ 11.8 mg/l at VT = 250 ml	≤ 10.30 mg/l at VT = 50 ml	≤ 10.9 mg/l at VT = 500 ml
Moisture return	≥ 32.2 mg/l at VT = 250 ml	≥ 33.7 mg/l at VT = 50 ml	≥ 33.1 mg/l at VT = 500 ml
Resistance	≤ 1.8 mbar at 30 l/min	≤ 1.5 mbar at 15 l/min	≤ 1.6 mbar at 30 l/min
Maximum duration of use	24 hours	24 hours	24 hours
Product	PVC-free/ Latex free	PVC-free/ Latex free	PVC-free/ Latex free

Cleanroom manufactured according to EN ISO 14644-1:2014.

*According to DIN EN ISO 9360-1 2009.

**Test has been performed by Nelson Lab, USA, in 2021.

BREATHING SYSTEM FILTERS

SPECIFICATIONS

	CareStar® 35 Plus	CareStar® 20 Plus
Part number	MP05755	MP05770
Deadspace	35 ml	20 ml
Recommended tidal volume	300 to 1500 ml	100 to 500 ml
Bacterial retention**	≥ 99.99%	≥ 99.99%
Viral retention**	≥ 99.9%	≥ 99.9%
Filtration method	electrostatic	electrostatic
Moisture loss*	-	-
Moisture return	-	-
Resistance	≤ 0.9 mbar at 30 l/min	≤ 1.3 mbar at 30 l/min
Maximum duration of use	24 hours	24 hours
Product	PVC-free/ Latex free	PVC-free/ Latex free

	SafeStar® 90 Plus	Safe Star® 60A Plus	SafeStar® 55 Plus
Part number	MP05785	MP05795	MP05790
Deadspace	90 ml	60 ml	55 ml
Recommended tidal volume	300 to 1500 ml	300 to 1500 ml	300 to 1500 ml
Bacterial retention**	≥ 99.9999%	≥ 99.999%	≥ 99.999%
Viral retention**	≥ 99.999%	≥ 99.9999%	≥ 99.999%
Filtration method	mechanical	mechanical	mechanical
Moisture loss*	-	-	-
Moisture return	-	-	-
Resistance	≤ 1.3 mbar at 30 l/min	≤ 2.0 mbar at 30 l/min	≤ 2.0 mbar at 30 l/min
Maximum duration of use	24 hours	24 hours	24 hours
Product	PVC-free/ Latex free	PVC-free/ Latex free	PVC-free/ Latex free

Cleanroom manufactured according to EN ISO 14644-1:2014.

*According to DIN EN ISO 9360-1 2009.

**Test has been performed by Nelson Lab, USA, in 2021.

HMEs

SPECIFICATIONS

	HumidStar® 55 Plus	HumidStar® 25 Plus	HME HumidStar® 2 Plus
Part number	MP05730	MP05735	MP05845
Deadspace	55 ml	25 ml	2 ml
Recommended tidal volume	300 to 1,500 ml	100 to 500 ml	10 to 50 ml
Bacterial retention	-	-	-
Viral retention	-	-	-
Filtration method	-	-	-
Moisture loss*	≤ 7.8 mg/l at VT = 500 ml	≤ 9.3 mg/l at VT = 250 ml	≤ 11.5 mg/l at VT = 45 ml
Moisture return	≥ 36.2 mg/l at VT = 500 ml	≥ 34.7 mg/l at VT = 250 ml	≥ 32.5 mg/l at VT = 45 ml
Resistance	≤ 0.6 mbar at 30 l/min	≤ 0.3 mbar at 30 l/min	≤ 1.2 mbar at 15 l/min
Maximum duration of use	24 hours	24 hours	24 hours
Product	PVC-free/ Latex free	PVC-free/ Latex free	PVC-free/ Latex free

	HME HumidStar® Trach Plus
Part number	MP05750
Deadspace	6 ml
Recommended tidal volume	100 to 1,500 ml
Bacterial retention	-
Viral retention	-
Filtration method	-
Moisture loss*	≤ 10.8 mg/l at Vt = 250 ml
Moisture return	≥ 33.2 mg/l at Vt = 250 ml
Resistance	≤ 0.3 mbar at 60 l/min
Maximum duration of use	24 hours
Product	PVC-free/ Latex free

Cleanroom manufactured according to EN ISO 14644-1:2014.

*According to DIN EN ISO 9360-1 2009.

CO₂ CUVETTES

3. CO₂ CUVETTES

CO₂ cuvettes are designed for the monitoring of CO₂ concentration during patient ventilation.

CO₂ CUVETTES, DISPOSABLE

A	CO₂ cuvette disposable, adult, 10 pcs.	MP01062
B	CO₂ cuvette disposable, pediatric, 10 pcs.	MP01063

CO₂ CUVETTES, REUSABLE

C	CO₂ cuvette reusable, adult	6870279
D	CO₂ cuvette reusable, pediatric	6870280

CO₂ CABLES

E	CO₂ mainstream sensor, for Evita® 2 dura, Evita® 4 Edition and Evita® XL	6871500
F	MCable Mainstream CO₂, for Evita V-family and Babylog VN-family, Infinity® M540, Vista120, Vista120 S, and Delta VF9.1 and higher	6871950



CATHETER MOUNTS AND AIRWAY CONNECTORS

4. CATHETER MOUNTS AND AIRWAY CONNECTORS

- Smoothbore products provide low resistance to flow while external coil reduces kinks and blockages
- Extendable product provides a lightweight and flexible solution
- Double swivel connector increases ease of use
- Double cap to maintain PEEP during suctioning or bronchoscopy
- Each catheter mount and airway connector is equipped with a red safety cap so that the system is completely closed and protected until directly prior to use on the patient
- 22F/15M connector on device side as a two-in-one solution
- Standardized connectors for safe connection to other components
- Luer Lock port for gas sampling
- Disposable convenience enhances workflow
- PVC-free design of gas leading components eliminates risk of softeners such as DEHP

A	Catheter mount ErgoStar® CM 40, disposable, 10 cm, double swivel elbow, 50 pcs.	MP01840
B	Catheter mount ErgoStar® CM 45, disposable, 15 cm, double swivel elbow, 50 pcs.	MP01845
C	Catheter mount ErgoStar® CM 50, disposable, 10 – 16 cm, extendable, 50 pcs.	MP01850
D	Catheter mount ErgoStar® CM 55, disposable, 15 cm, straight, 50 pcs.	MP01855
E	Catheter mount ErgoStar® CM 60, disposable, 10 cm, small, double swivel elbow, 50 pcs.	MP01860
F	Airway connector ErgoStar® AC 90, disposable, double-swivel elbow w/ double cap, 100 pcs.	MP01890
G	Airway connector ErgoStar® AC 95, disposable, w/ Luer Lock sampling port, 100 pcs.	MP01895



Not all articles are available worldwide.

COMBINATIONS OF FILTER AND TUBE EXTENTION

COMBINATIONS OF FILTER AND TUBE EXTENTION

CombiStar A Filter HME, combination of airway connector MP01845 and filter TwinStar® 55 Plus (MP05805), 25 pcs.	MP12060
CombiStar Filter HME flex, combination of flexible tube extention MP01850 and filter TwinStar® 55 Plus (MP05805), 25 pcs.	MP12061
CombiStar F-HME HEPA flex, combination of flexible tube extention MP01850 and filter TwinStar® HEPA Plus (MP05801), 25 pcs.	MP12062
CombiStar mechanical Filter flex, combination of flexible tube extention MP01850 and filter SafeStar® 55 Plus (MP05790), 25 pcs.	MP12063

NIV MASKS AND ACCESSORIES

5. NON-INVASIVE PATIENT INTERFACES AND ACCESSORIES

5.1 NIV MASKS AND ACCESSORIES

ClassicStar® NIV oronasal masks cover the nose and mouth of the patient. Therefore, they provide effective therapy even in case of mouth breathing.


The mask design ensures high patient comfort combined with anatomical fit and effective seal due to individually adjustable mask cushion. The 6-point headgear provides a stable and secure fixation. Adjustable forehead support and forehead pad allow for distributing the pressure from the bridge of nose towards the forehead for more patient comfort.

ClassicStar® NIV masks with standard elbow (SE) have been successfully tested for compatibility with Dräger ventilators which have NIV-option. Variants with Anit-Asphyxia Valve (AAV) are intended for Multi-Vendor-Accessory business.

All ClassicStar® NIV masks are for single patient use.

Sizing gauges for easy mask size identification as well as a pumpball are delivered together with the masks.



A	NIV oronasal mask ClassicStar® SE, disposable, size S	MP01573
	NIV oronasal mask ClassicStar® SE, disposable, size M	MP01574
	NIV oronasal mask ClassicStar® SE, disposable, size L	MP01575
B	NIV oronasal mask ClassicStar® AAV, disposable, size S	MP01570
	NIV oronasal mask ClassicStar® AAV, disposable, size M	MP01571
	NIV oronasal mask ClassicStar® AAV, disposable, size L	MP01572

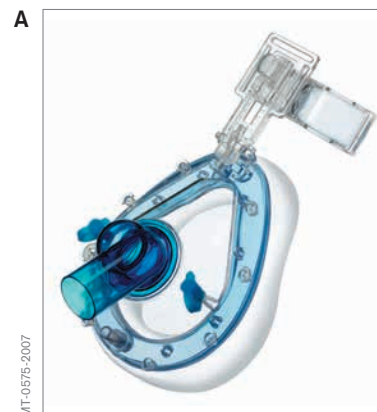

NIV MASKS AND ACCESSORIES

Oronasal masks are often the first-line strategy in the initial management of acute respiratory failure as they cover nose and mouth. Nasal masks are often chosen if NIV therapy has to be prolonged due to improved comfort and reduced complications associated with oronasal masks. Similar to ClassicStar® oronasal masks, the nasal mask cushion can be individually adjusted to the patient's face providing high patient comfort and effective seal. The 4-point-headgear ensures a stable and secure fixation reducing mask displacement. All ClassicStar® NIV nasal masks are for single patient use. Sizing gauges for easy mask size identification as well as a pumpball are delivered together with the masks.

A	<u>NIV nasal mask ClassicStar® non-vented, disposable, size S</u>	<u>MP01623</u>
	<u>NIV nasal mask ClassicStar® non-vented, disposable, size M</u>	<u>MP01624</u>
	<u>NIV nasal mask ClassicStar® non-vented, disposable, size L</u>	<u>MP01625</u>
B	<u>NIV nasal mask ClassicStar® vented, disposable, size S</u>	<u>MP01620</u>
	<u>NIV nasal mask ClassicStar® vented, disposable, size M</u>	<u>MP01621</u>
	<u>NIV nasal mask ClassicStar® vented, disposable, size L</u>	<u>MP01622</u>

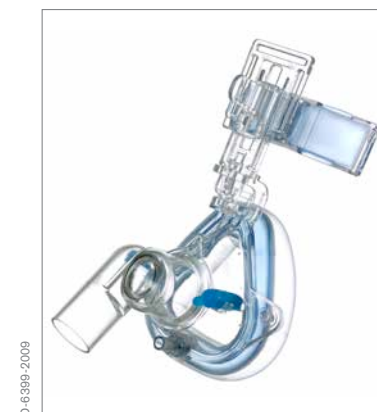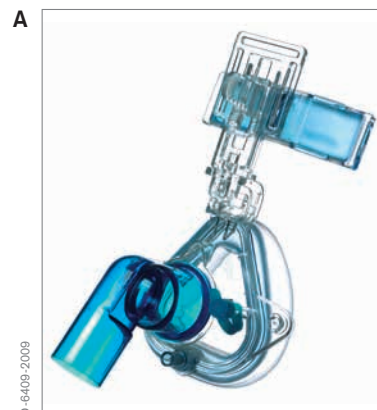

NIV MASKS AND ACCESSORIES

The ClassicStar® NIV total-face mask seals around the perimeter of the face, covering mouth, nose and eyes. These areas are less pressure sensitive and have smoother facial contours compared to the areas a oronasal masks covers (nasal bridge). The mask is suitable for patients with nasal bridge breakdown or facial irregularities. The air filled cushion can be adjusted individually to the patients face, providing high patient comfort and effective seal.

The 5-point headgear which includes a headstrap at the top of the mask provides a stable and secure attachment, thus preventing mask displacement and reducing pressure. Quick release clips are integrated for quick removal of the mask.

Sizing gauges for easy mask size identification as well as a pumpball are delivered together with the masks.



A	NIV total-face mask ClassicStar® SE, disposable, size S	MP02910
	NIV total-face mask ClassicStar® SE, disposable, size M	MP02911
	NIV total-face mask ClassicStar® SE, disposable, size L	MP02912
B	NIV total-face mask ClassicStar® AAV, disposable, size S	MP02913
	NIV total-face mask ClassicStar® AAV, disposable, size M	MP02914
	NIV total-face mask ClassicStar® AAV, disposable, size L	MP02915
C	NIV total-face mask ClassicStar® AAV without exhalation ports, disposable, size S	MP02916
	NIV total-face mask ClassicStar® AAV without exhalation ports, disposable, size M	MP02917
	NIV total-face mask ClassicStar® AAV without exhalation ports, disposable, size L	MP02918



NIV MASKS AND ACCESSORIES

The ClassicStar® plus mask series features a soft and anatomically shaped silicone cushion as sealing interface to the patients face. The mask has an easy to adjust forehead pad with silicone cushion which allows optimal adaption to the sensitive nasal bridge. The mask is available in four sizes (S, M, L, XL) as oronasal and nasal mask. ClassicStar® plus NIV masks with standard elbow (SE) have been successfully tested for compatibility with Dräger ventilators which have NIV-option. Variants with Anti-Asphyxia-Valve (AAV) are intended for Multi-Vendor-Accessory business.

A	NIV oronasal mask ClassicStar® plus SE, disposable, size S, 5 pcs.	MP04711
	NIV oronasal mask ClassicStar® plus SE, disposable, size M, 5 pcs.	MP04712
	NIV oronasal mask ClassicStar® plus SE, disposable, size L, 5 pcs.	MP04713
	NIV oronasal mask ClassicStar® plus SE, disposable, size XL, 5 pcs.	MP04714
B	NIV oronasal mask ClassicStar® plus AAV, disposable, size S, 5 pcs.	MP04721
	NIV oronasal mask ClassicStar® plus AAV, disposable, size M, 5 pcs.	MP04722
	NIV oronasal mask ClassicStar® plus AAV, disposable, size L, 5 pcs.	MP04723
	NIV oronasal mask ClassicStar® plus AAV, disposable, size XL, 5 pcs.	MP04724
C	NIV oronasal mask ClassicStar® plus AAV without exhalation ports, disposable, size S, 5 pcs.	MP04731
	NIV oronasal mask ClassicStar® plus AAV without exhalation ports, disposable, size M, 5 pcs.	MP04732
	NIV oronasal mask ClassicStar® plus AAV without exhalation ports, disposable, size L, 5 pcs.	MP04733
	NIV oronasal mask ClassicStar® plus AAV without exhalation ports, disposable, size XL, 5 pcs.	MP04734

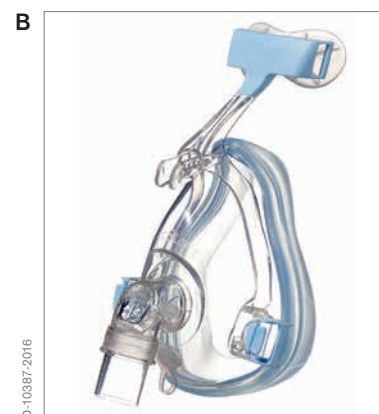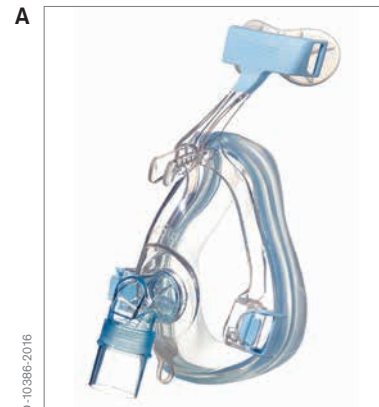

NIV MASKS AND ACCESSORIES

A	NIV nasal mask ClassicStar® plus non-vented, disposable, size S, 5 pcs.	MP04741
	NIV nasal mask ClassicStar® plus non-vented, disposable, size M, 5 pcs.	MP04742
	NIV nasal mask ClassicStar® plus non-vented, disposable, size L, 5 pcs.	MP04743
	NIV nasal mask ClassicStar® plus non-vented, disposable, size XL, 5 pcs.	MP04744
B	NIV nasal mask ClassicStar® plus vented, disposable, size S, 5 pcs.	MP04751
	NIV nasal mask ClassicStar® plus vented, disposable, size M, 5 pcs.	MP04752
	NIV nasal mask ClassicStar® plus vented, disposable, size L, 5 pcs.	MP04753
	NIV nasal mask ClassicStar® plus vented, disposable, size XL, 5 pcs.	MP04754



NIV MASKS AND ACCESSORIES

The NovaStar® plus NiV mask series cover nose and mouth of the patient. Therefore, they provide effective therapy even in case of mouth breathing. This range of products provides ultimate comfort and an effective seal due to the individual adaption of the silicone-gel cushion with a pliable embedded ring. NovaStar® plus NiV masks with standard elbow (SE) have been successfully tested for compatibility with Dräger ventilators which have NiV-option. Variants with Anti-Asphyxia-Valve (AAV) are intended for Multi-Vendor-Accessory business. All NovaStar® plus NiV masks are for multi patient use. Replacement headgears with magnetic clips are available individually (cf. accessories for NiV masks).

A	<u>NIV oronasal mask NovaStar® plus, SE, reusable, size S, 3 pcs.</u>	<u>MP05051</u>
	<u>NIV oronasal mask NovaStar® plus, SE, reusable, size M, 3 pcs.</u>	<u>MP05052</u>
	<u>NIV oronasal mask NovaStar® plus, SE, reusable, size L, 3 pcs.</u>	<u>MP05053</u>
	<u>NIV oronasal mask NovaStar® plus, SE, reusable, size XL, 3 pcs.</u>	<u>MP05054</u>
B	<u>NIV oronasal mask NovaStar® plus, AAV, reusable, size S, 3 pcs.</u>	<u>MP05061</u>
	<u>NIV oronasal mask NovaStar® plus, AAV, reusable, size M, 3 pcs.</u>	<u>MP05062</u>
	<u>NIV oronasal mask NovaStar® plus, AAV, reusable, size L, 3 pcs.</u>	<u>MP05063</u>
	<u>NIV oronasal mask NovaStar® plus, AAV, reusable, size XL, 3 pcs.</u>	<u>MP05064</u>
C	<u>NIV oronasal mask NovaStar® plus, AAV, without exhalation ports, reusable, size S, 3 pcs.</u>	<u>MP05071</u>
	<u>NIV oronasal mask NovaStar® plus, AAV, without exhalation ports, reusable, size M, 3 pcs.</u>	<u>MP05072</u>
	<u>NIV oronasal mask NovaStar® plus, AAV, without exhalation ports, reusable, size L, 3 pcs.</u>	<u>MP05073</u>
	<u>NIV oronasal mask NovaStar® plus, AAV, without exhalation ports, reusable, size XL, 3 pcs.</u>	<u>MP05074</u>



NIV MASKS AND ACCESSORIES

The **MiniMe2*** nasal mask is designed for children ages 2 – 12. It features an ergonomic silicone-gel cushion and a shapeable shell mask to fit individual needs and thus ensuring an effective seal and a comfortable fit.

A	Pediatric nasal mask MiniMe2 non-vented, disposable, size S	MP01555
	Pediatric nasal mask MiniMe2 non-vented, disposable, size L	MP01556
B	Pediatric nasal mask MiniMe2 vented, disposable, size S	MP01553
	Pediatric nasal mask MiniMe2 vented, disposable, size L	MP01554



*Tradeware



NIV MASKS AND ACCESSORIES

ACCESSORIES FOR NIV MASKS

- | | | |
|----------|---|----------------|
| A | Two-sided pump ball for air management | MP01590 |
| B | NovaStar® plus Headgear inclusive magnetic clips, disposable, 3 pcs. | MP05080 |

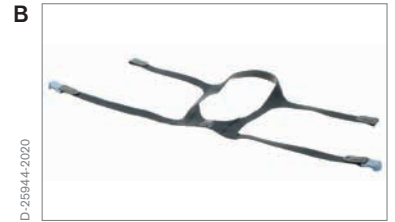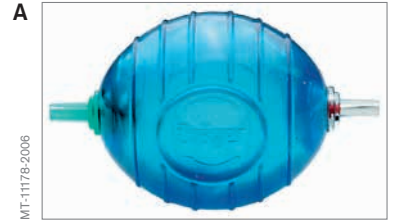

ANESTHESIA MASKS AND ACCESSORIES

5.2 ANESTHESIA MASKS AND ACCESSORIES

ANESTHESIA MASKS SEDISTAR®

The disposable anesthesia face mask SediStar® features a soft and anatomically shaped lip as the sealing interface to the patient's face. The mask is BPA and PVC free.

Optimal fit


Due to the ergonomic design and sealing lip the mask offers an optimized fit for different facial anatomies.

Color-coding

Every mask size is color coded by a different label, which makes it easy to identify the correct mask size right of the shelf.

Additional mask accessories

Easy application is supported through additional mask accessories: hook ring and headgears for fixation. Two different headgears in small and large designs offer patients and physicians a soft and elastic support and holding position.

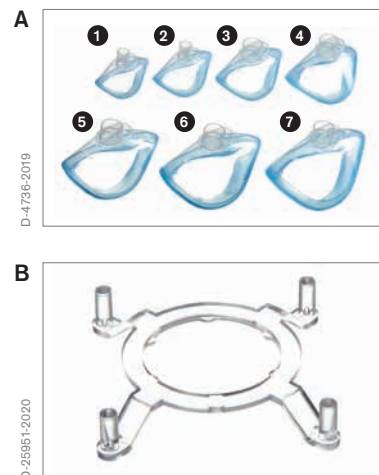

A	1 Anesthesia mask SediStar®, disposable, w/o hook ring, size 1, infant S, 30 pcs.	MP01701
	2 Anesthesia mask SediStar®, disposable, w/o hook ring, size 2, infant L, 30 pcs.	MP01702
	3 Anesthesia mask SediStar®, disposable, w/o hook ring, size 3, child, 30 pcs.	MP01703
	4 Anesthesia mask SediStar®, disposable, w/o hook ring, size 4, adult S, 30 pcs.	MP01704
	5 Anesthesia mask SediStar®, disposable, w/o hook ring, size 5, adult M, 30 pcs.	MP01705
	6 Anesthesia mask SediStar®, disposable, w/o hook ring, size 6, adult L, 30 pcs.	MP01706
	7 Anesthesia mask SediStar®, disposable, w/o hook ring, size 7, adult XL, 30 pcs.	MP01707
B	Hook ring for SediStar® anesthesia masks, disposable, 22 mm ID, for sizes 3 – 7, 30 pcs.	MP01708

ANESTHESIA MASKS AND ACCESSORIES

ANESTHESIA MASKS LITESTAR®

The re-designed recessed grip as well as the teardrop shape of the of our disposable anesthesia mask LiteStar® support an anatomical fit. The transparent mask body allows for continuous observation of the patient's condition.


Color-coding

Every mask size is color coded by a different label, which makes it easy to identify the correct mask size.

Additional mask accessories

Easy application is supported through additional mask accessories: hook ring and headgears for fixation. Two different headgears in small and large designs offer patients and physicians a soft and elastic support and holding position.

A	① Anesthesia mask LiteStar®, disposable, w/o hook ring, size 1, infant S, 30 pcs.	MP11001
	② Anesthesia mask LiteStar®, disposable, w/o hook ring, size 2, infant L, 30 pcs.	MP11002
	③ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 3, child, 30 pcs.	MP11003
	④ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 4, adult S, 30 pcs.	MP11004
	⑤ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 5, adult M, 30 pcs.	MP11005
	⑥ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 6, adult L, 30 pcs.	MP11006
	⑦ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 7, adult XL, 30 pcs.	MP11007
B	Hook ring for LiteStar® anesthesia masks, disposable, 22 mm ID, for sizes 3 – 7, 30 pcs.	MP11008

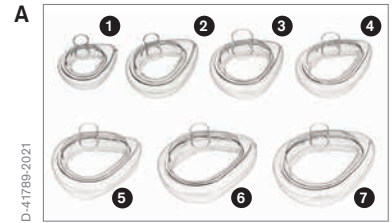

ANESTHESIA MASKS AND ACCESSORIES

ANESTHESIA MASKS COMFORTSTAR®

Patient comfort and ease of use were the top priorities when developing the Dräger disposable ComfortStar® anesthesia mask portfolio. Patient comfort is supported by teardrop or round-shaped masks and adjustable (inflatable) mask cushion which provide anatomical fit with effective seal. Built-in ridges (sizes 4 – 6) and a flexible mask body provide healthcare professionals an easy and comfortable grip. ComfortStar® scented anesthesia masks have the same excellent benefits of the ComfortStar® standard masks, but in addition, they have added aromas. The appropriate color scheme provides optical stimulation and is intended to reduce pre-operative anxiety in patients.

All masks are pre-equipped with ‚Traffic Light‘ colored hook rings. They are color coded for easier size identification. Smaller sizes (0 – 1) have reddish color, mid-sizes (2 – 3) yellowish color and large sizes (4 – 6) greenish color. All ComfortStar® masks are compatible with all Dräger anesthesia devices.

Pumpball as well as different headgears for fixation are available (cf. accessories for anesthesia masks).



A	① Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 0, premature baby, 20 pcs.	MP01510
	② Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 1, neonatal, 20 pcs.	MP01511
	③ Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 2, baby, 20 pcs.	MP01512
	④ Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 3, infant, 20 pcs.	MP01513
	⑤ Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 4, adult S, 20 pcs.	MP01514
	⑥ Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 5, adult M, 20 pcs.	MP01515
	⑦ Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 6, adult L, 20 pcs.	MP01516

ANESTHESIA MASKS AND ACCESSORIES

A	<p>① Anesthesia mask ComfortStar®, disposable, hook ring, strawberry, size 1, neonatal, 20 pcs. <u>MP01521</u></p> <p>② Anesthesia mask ComfortStar®, disposable, hook ring, strawberry, size 2, baby, 20 pcs. <u>MP01522</u></p> <p>③ Anesthesia mask ComfortStar®, disposable, hook ring, strawberry, size 3, infant, 20 pcs. <u>MP01523</u></p> <p>④ Anesthesia mask ComfortStar®, disposable, hook ring, strawberry, size 4, adult S, 20 pcs. <u>MP01524</u></p>
B	<p>① Anesthesia mask ComfortStar®, disposable, hook ring, bubble gum, size 1, neonatal, 20 pcs. <u>MP01531</u></p> <p>② Anesthesia mask ComfortStar®, disposable, hook ring, bubble gum, size 2, baby, 20 pcs. <u>MP01532</u></p> <p>③ Anesthesia mask ComfortStar®, disposable, hook ring, bubble gum, size 3, infant, 20 pcs. <u>MP01533</u></p> <p>④ Anesthesia mask ComfortStar®, disposable, hook ring, bubble gum, size 4, adult S, 20 pcs. <u>MP01534</u></p>
C	<p>① Anesthesia mask ComfortStar®, disposable, hook ring, mint, size 3, infant, 20 pcs. <u>MP01543</u></p> <p>② Anesthesia mask ComfortStar®, disposable, hook ring, mint, size 4, adult S, 20 pcs. <u>MP01544</u></p> <p>③ Anesthesia mask ComfortStar®, disposable, hook ring, mint, size 5, adult M, 20 pcs. <u>MP01545</u></p> <p>④ Anesthesia mask ComfortStar®, disposable, hook ring, mint, size 6, adult L, 20 pcs. <u>MP01546</u></p>



ANESTHESIA MASKS AND ACCESSORIES

SILICONE MASKS

Silicone masks are made of medical grade silicone because this is generally considered the material of choice for its superior ability to conform to the patient's facial anatomy. The field of application includes anesthesia and mask ventilation for pediatric and adult patients. The masks are available in various shapes and sizes.


The smaller sizes have a rather round shape in order to provide comfortable fit and effective seal for the faces of infants.

The bigger sizes have an anatomically formed teardrop shape to ensure convenient fit and seal for facial contours of adult patients.

Silicone masks have a flexible, soft and transparent mask body.

The mask body design provides a comfortable grip. Shape and material allow for safe handling by doctors and nurses.

The silicone masks are reusable.

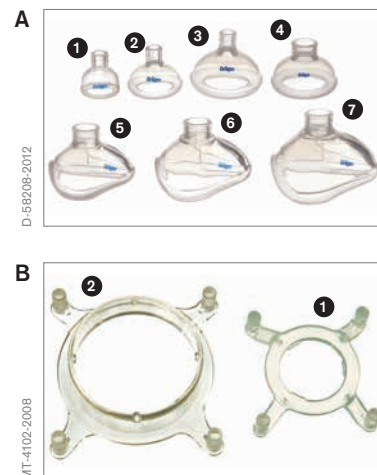

A	① Silicone mask, reusable, size 00, premature baby	MP02900
	② Silicone mask, reusable, size 0, small neonate	MP01910
	③ Silicone mask, reusable, size 1, large neonate	MP01911
	④ Silicone mask, reusable, size 2, small child	MP01912
	⑤ Silicone mask, reusable, size 3, large child	MP01913
	⑥ Silicone mask, reusable, size 4, small adult	MP01914
	⑦ Silicone mask, reusable, size 5, large adult	MP01915
B	① Hook ring for silicone masks, reusable, 15 mm OD, for sizes 00 – 1, 10 pcs.	MP01916
	② Hook ring for silicone masks, reusable, 22 mm ID, for sizes 2 – 5, 10 pcs.	MP01917



Additional mask accessories, e. g. reusable hook rings as well as headgears are available for fixation (cf. accessories for anesthesia masks).

ANESTHESIA MASKS AND ACCESSORIES

ACCESSORIES FOR ANESTHESIA MASKS

There are three different headgears in small and large designs offering patients and physicians a soft, elastic, lightweight fixation of disposable anesthesia/ventilation masks.

A	① Headgear, disposable, soft fabric, small	MP01559
	② Headgear, disposable, soft fabric, large	MP01560
B	① Headgear, disposable, soft rubber, small	MP01561
	② Headgear, disposable, soft rubber, large	MP01562
C	Two-sided pump ball for air management	MP01590

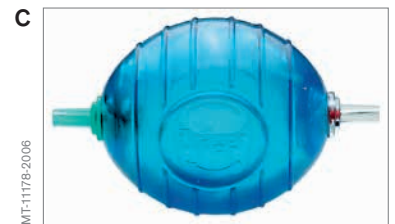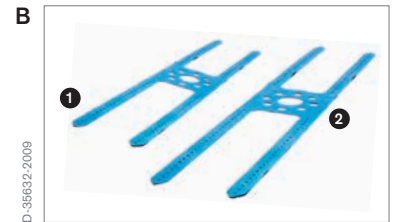

OXYGEN (O₂) AND AEROSOL THERAPY MASKS, VENTURI MASK KITS, NASAL CANNULAS AND ACCESSORIES

5.3 OXYGEN (O₂) AND AEROSOL THERAPY MASKS, VENTURI MASK KITS, NASAL CANNULAS AND ACCESSORIES

Oxygen therapy is the administration of oxygen at a concentration higher than in ambient air. The purpose is to prevent and/or treat the symptoms and manifestation of hypoxia. Aerosol therapy is used to administer medicine in the patient airways/lungs in order to reduce inflammation caused by different respiratory conditions. The aerosol and oxygen therapy masks are for disposable use.


A	O2Star™ Non-rebreather oxygen mask with connection tube 2.1 m, disposable, size S, 50 pcs.	MP01930
	O2Star™ Non-rebreather oxygen mask with connection tube 2.1 m, disposable, size M/L, 50 pcs.	MP01920
B	O2Star™ Medium concentration oxygen mask with connection tube 2.1 m, disposable, size S, 50 pcs.	MP01931
	O2Star™ Medium concentration oxygen mask with connection tube 2.1 m, disposable, size M/L, 50 pcs.	MP01921
C	O2Star™ Oxygen and aerosol mask with 22 mm OD, disposable, size S, 50 pcs.	MP01934
	O2Star™ Oxygen and aerosol mask with 22 mm OD, disposable, size M/L, 50 pcs.	MP01926
D	O2Star™ Tracheostomy mask, disposable, size S, 50 pcs.	MP01932
	O2Star™ Tracheostomy mask, disposable, size M/L, 50 pcs.	MP01924
E	O2Star™ Face tent, disposable, size S, 50 pcs.	MP01933
	O2Star™ Face tent, disposable, size M/L, 50 pcs.	MP01925
F	O2Star™ Venturi mask kit, adjustable venturi valve, disposable, size S, 25 pcs.	MP01935
G	O2Star™ Venturi mask kit, adjustable venturi valve, disposable, size M/L, 25 pcs.	MP01922



Not all articles are available worldwide.

OXYGEN (O₂) AND AEROSOL THERAPY MASKS, VENTURI MASK KITS, NASAL CANNULAS AND ACCESSORIES


A	O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size S, 50 pcs.	<u>MP02920</u>
	O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size M, 50 pcs.	<u>MP02921</u>
	O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size L, 50 pcs.	<u>MP02922</u>
	O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size XL, 50 pcs.	<u>MP02923</u>
B	O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size S, 50 pcs.	<u>MP02924</u>
	O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size M, 50 pcs.	<u>MP02925</u>
	O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size L, 50 pcs.	<u>MP02926</u>
	O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size XL, 50 pcs.	<u>MP02927</u>

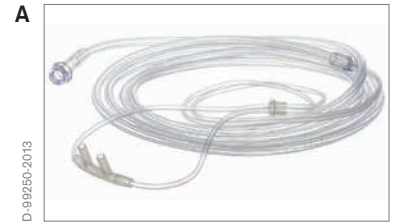

**OXYGEN (O₂) AND AEROSOL THERAPY MASKS,
VENTURI MASK KITS, NASAL CANNULAS AND ACCESSORIES**

**ACCESSORIES FOR OXYGEN (O₂)
AND AEROSOL THERAPY MASKS,
VENTURI MASK KITS, NASAL CANNULAS**

A	O2Star™ Connector kit, 6 mm OD/22 mm OD + 6 mm OD/22 mm ID with connection tube, disposable, 5 pcs.	MP01940
B	O2Star™ Connector 22 mm ID / 22 mm ID, disposable, 5 pcs.	MP01942
C	O2Star™ Oxygen connection tube, 2.1 m, disposable, 25 pcs.	MP01927
	O2Star™ Oxygen connection tube, 2.1 m, with 9/16" connector, disposable, 25 pcs.	MP01944
D	O2Star™ Oxygen connection tube, 4 m, disposable, 25 pcs.	MP01943
	O2Star™ Oxygen connection tube, 4 m, with 9/16" connector, disposable, 25 pcs.	MP01945

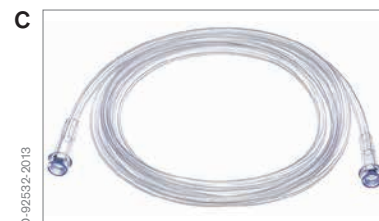


HIGH-FLOW THERAPY NASAL CANNULA AND ACCESSORIES

5.4 HIGH-FLOW THERAPY NASAL CANNULA AND ACCESSORIES

High Flow Therapy allows to provide warm humidified air with enriched oxygen concentration at flow rates up to 60 l per minute to patients via nasal cannulas. Patient comfort is one of the advantages of High Flow therapy over other treatment forms for patient with compromised respiration and can in many cases successfully avoid using less comfortable or even invasive forms of respiratory support. The HI-Flow Star Nasal Cannula and the inspiratory breathing circuit offer advanced features designed for enhanced patient safety and comfort such as integrated clips in order to provide effective breathing circuit tension relief, a flexible, breathable hose minimizing internal condensation in the high flow nasal cannula as well as soft materials for the prongs and the headgear.

Thanks to its symmetrical design, the nasal cannula can be connected from either the left or the right side of the bed. This helps caregivers save time and effort. The HI-Flow Star products are disposable. Oxygen and compressed air are required which are provided by an O₂ mixer such as the Oxymixer from Dräger, as well as a humidification source for the hose system.

A	① HI-Flow Star Nasal Cannula, disposable, adult S, Flow Range < 60 l/min, 25 pcs.	MP05511
A	② HI-Flow Star Nasal Cannula, disposable, adult M, Flow Range < 60 l/min, 25 pcs.	MP05512
A	③ HI-Flow Star Nasal Cannula, disposable, adult L, Flow Range < 60 l/min, 25 pcs.	MP05513
B	HI-Flow Star System, heated inspiratory breathing circuit, consists of: heated inspiratory hose, humidifier chamber, breathing hose when used with a ventilator, O ₂ hose when used with an O ₂ blender, bed clip, disposable, 25 pcs.	MP05601
C	HI-Flow Star Kit Aquapor A, heated inspiratory breathing circuit for the H300 humidifier, heated inspiratory hose, humidifier chamber, breathing hose when used with a ventilator, O ₂ hose when used with an O ₂ blender, bed clip, disposable, 10 pcs.	MP17090
D	HI-Flow Star Valve Kit w/ Pop-Up Valve and Connector 22 mm, disposable, 15 pcs.	MP05507

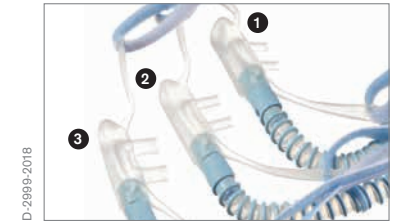

OXYMIXER*

6. OXYMIXER

Oxymixer enables a continuous and individually adjusted mixture of medical air and medical oxygen to be administered to infants, children and adults. The O₂ concentration is easily adjusted to a highly accurate flow – even with flow rates of less than 3 l/min.

Oxymixer configuration

MP04200

Free choice of Low- or High Flow Blender with:
1 or 2 flowmeter, NIST or DISS gas inlet,
Monitor, Bleed swith, Nozzle, Sound damper.
The rail fixation is included.

A Oxymixer High Flow

MP04201

Including: Rail clamp, Flowmeter 0 – 32 l/min right,
Flowmeter 0 – 85 l/min left.

B Oxymixer High Flow with monitor

MP04202

Including: Rail clamp, Flowmeter 0 – 32 l/min right,
Flowmeter 0 – 85 l/min left, Monitor, Oxygen sensor,
Oxygen sensor cable, Sensor receptacle.

C Multi-P trolley* equipped

MP13020

Including base with 5 castors, 38 mm Ø pole, 3 rail clamps, basket.

D Multi-P trolley*

MP13021

Including base with 5 castors and 38 mm Ø pole for accessories.

Several hoses to connect Oxymixer with the central supply.
For CS hoses see chapter 4, "Accessory for gas supply".

*Tradeware



OXYMIXER*

A	<u>Oxymixer sound damper 9/16"</u>	<u>MP04218</u>
B	<u>Oxymixer sound damper 22 and 15 mm</u>	<u>MP04219</u>
C	<u>Adapter 9/16" to 15/22mm</u>	<u>MP05305</u>
D	<u>Adapter, with 6 mm nozzle, e.g., for disposable breathing masks, 9/16" connection</u>	<u>MP00416</u>
E	<u>Oxymixer oxygen sensor (spare part)</u>	<u>MP04221</u>

*Tradeware



DEVICE RELATED ACCESSORIES

7. DEVICE RELATED ACCESSORIES

EXPIRATORY VALVES FOR EVITA® 2 DURA, 4 EDITION, XL

A	Expiratory valve for Evita® 2 dura, 4 edition, XL, disposable, 10 pcs.	8414776
B	Expiratory valve for Evita® 2 dura, 4 edition, XL, disposable, 10 pcs.	MP02600
	Expiratory valve for Evita®, reusable	8410580
	Membrane for Evita® 2 dura, Evita® 4, Evita® XL, for expiratory valve	8412015
C	Water trap for Evita® 2 dura, Evita® 4, Evita® XL, reusable	8413125

EXPIRATORY VALVES FOR BABYLOG®

	Expiratory valve, for neonatal/pediatric patients (Non-RFID), reusable	8415296
D	Expiratory valve for Babylog® VN 500, reusable (RFID)	8415270
	Expiratory valve for Babylog® 8000	8408950



DEVICE RELATED ACCESSORIES

**DISPOSABLE EXPIRATORY VALVES
FOR CARINA®, SAVINA® AND V-FAMILY**

A	Infinity® ID, expiratory valve, disposable, 10 pcs.	MP01060/ MP20108
B	Expiratory valve, disposable, 10 pcs.	MP01061/ MP20109
C	Expiratory valve for Carina®, disposable, 25 pcs.	MP00220
D	Leakage valve for Carina® SyncVent, disposable, 5 pcs.	MP00224

**REUSABLE EXPIRATORY VALVES
FOR SAVINA® AND V-FAMILY**

E	Expiratory valve for Evita® Infinity® V500, reusable	8416750
	Expiratory valve for Savina®, reusable	8413660
F	Expiratory valve, reusable	8417050



DEVICE RELATED ACCESSORIES

A Expiratory filter, disposable, 20 pcs.

MP01781

A

D-43830-2012



DEVICE RELATED ACCESSORIES


A Fluxmed™, Pressure Probe with esophageal balloon, disposable **MP20110**

ACCESSORIES FOR MX 300-i

B O₂ sensor* for MX 300 and MX 300-i, reusable **MP01442**

C T-piece adapter* for MX 300 and MX 300-i, reusable **MP01445**

*Tradeware

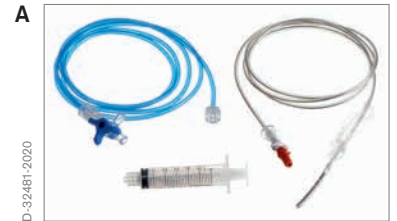

PACK2GO®/SET2GO

8. PACK2GO®

Dräger has developed a comprehensive range of disposable accessories. These accessories can now also be ordered together in a Pack2Go® package, in order to benefit the hospital's costs and logistical/clinical processes. Pack2Go® package solutions include dedicated Dräger and general accessories, which can be used as a system solution with Dräger devices.

- | | | |
|----------|--|----------------|
| A | Pack2Go® Evita® invasive ventilation pack with breathing circuit VentStar®, basic, + cuvette, 10 pcs. | MP20370 |
| B | Pack2Go® Evita® invasive ventilation pack with breathing circuit VentStar®, basic, 1.8 m, 10 pcs. | MP20371 |
| C | Pack2Go® Evita® invasive ventilation pack with breathing circuit VentStar®, coaxial, + cuvette, 10 pcs. | MP20377 |
| D | Pack2Go® Evita® invasive ventilation pack with breathing circuit VentStar®, coaxial, 1.8 m, 10 pcs. | MP20378 |



PACK2GO®/SET2GO

- | | | |
|----------|---|----------------|
| A | Pack2Go® Carina® NIV ventilation pack
with oronasal mask ClassicStar® SE, size S, 10 pcs. | MP20372 |
| | Pack2Go® Carina® NIV ventilation pack
with oronasal mask ClassicStar® SE, size M, 10 pcs. | MP20373 |
| | Pack2Go® Carina® NIV ventilation pack
with oronasal mask ClassicStar® SE, size L, 10 pcs. | MP20374 |
| B | Pack2Go® invasive ventilation pack
with breathing circuit VentStar®, coaxial, 1.8 m, 10 pcs. | MP20375 |
| C | Pack2Go® invasive ventilation pack
with breathing circuit VentStar®, basic, 1.8 m, 10 pcs. | MP20376 |
| D | Pack2Go® anesthesia pack, A,
with anesthesia set + watertrap, 10 pcs. | MP20379 |



PACK2GO®/SET2GO

SET2GO

- A** With Dräger's Set2go solutions we offer you totally customised sets consisting of an individual selection of single-patient-use accessories. Dräger Set2go makes using and ordering consumables and accessories much simpler for you. You just need one article number to get all the items in a single package and accomplish your specific task.

Please contact your local Dräger sales organisation to create your individual accessory solution.



NEBULIZER

9. NEBULIZER

A	Pneum. Medication Nebulizer N	8411030
B	Pneum. Medication Nebulizer	8412935
C	Nebulizer hose	8412985



NEBULIZER

AEROGEN PRO*

Aerogen Pro works with an electronic micropump. Aerogen Pro is ideally suited to nebulizer therapy requirements for inline mechanical ventilation and it accommodates all patient age ranges from neonates to adults. It operates without changing patient ventilator parameters, without heating or degrading medications and it can be refilled without interrupting ventilation.

A	Medication nebulizer Aerogen Pro	MP01010
B	Nebulizer unit with filler cap for Aerogen Pro	MP01030
C	T-piece with plug for Aerogen Pro, A, ISO 22 mm, 5 pcs.	MP01031
D	T-piece set with plug for Aerogen Pro, N, 15 mm (0.6 in), w/ adapter 10 mm (0.4 in), 5 pcs.	MP01032
E	T-connector for Aerogen Pro, 12 mm, with plug, 5 pcs.	MP01037
F	Controller cable for Aerogen Pro	MP01036
G	Filler cap for Aerogen Pro, silicone, 5 pcs.	MP01035

*Tradeware



NEBULIZER

AEROGEN SOLO*

Aerogen Solo Starter Kit

Consists of: Control module, 2 Aerogen Solo Nebulizers with filler caps, adult T-adapter with silicone plug, Universal mounting bracket, AC/DC adapter and Instruction manual.

A	Aerogen Solo Starter Kit, Int.	MP01071
	Aerogen Solo Starter Kit, FR, NL, IT, GE	MP01072
	Aerogen Solo Starter Kit, ES, IT, PT, GR, TK	MP01073
	Aerogen Solo Starter Kit, SW, NO, FI, DK	MP01074
	Aerogen Solo Starter Kit, UK	MP01075
	Aerogen Solo Starter Kit, US	MP01076
B	Nebulizer head, 5 pcs.	MP01078
C	Nebulizer head, 10 pcs.	MP01079
	T-Adapter, adult, 10 pcs.	MP01085
	T-Adapter, pediatric, 10 pcs.	MP01086
	Aerogen T-Adapter**, 11.5 mm, neonates, 10 pcs.	MP02689
	Kit T-Adapter, Pediatric with Neo connectors, 10 pcs.	MP01087
	Aerogen Solo USB Controller Starter Kit NE	MP01089
	Aerogen USB Controller Replacement w/o AC/DC adapter	MP01095

*Tradeware

**Made by Drägerwerk AG & Co. KGaA




SPARE PARTS FOR RESUSCITATORS

10. SPARE PARTS FOR RESUSCITATORS

A	Oxidem 3000 demand valve incl. adapter (2120496, 2120488)	2111764
	1.5 m spiral hose and M12×1 screw connector	2M86291
	1.5 m spiral hose and probe for central supply system	2M86292
	3.0 m spiral hose and probe for central supply system	2M86293
	Adapter for Oxidem for M12×1 (for existing hose systems)	8410547
B	Adapter for Ambu Mark III*	2120488
C	Adapter for Laerdal**	2120496
	Adapter for Laerdal	5790124
D	Adapter for Combibag	2M86913
E	Adapter for Resutator 2000/MR-100**	2120461
F	Sealing rings for spiral hoses, 20 pcs.	M07152

*Made by Spiracle Technology

**Made by Allied Healthcare, USA

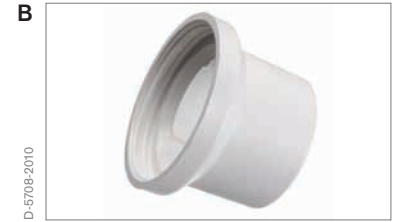



Dräger

Dräger AG, 26149 Luebeck, Germany
Dräger AG, 26149 Luebeck, Germany
電話 49 43068 0
FAX 49 43068 1201
E-MAIL: dräger@dräger.com
www.dräger.com
製造販売
ダイカルクマバ
〒100-0001 東京都千代田区千代田2-15-17 日興東急ビル
〒100-0001 東京都千代田区千代田2-15-17 日興東急ビル
アフロ
R3P033-09
Made in Germany

Soda lime

Contactless and dustless—protect your patient, protect yourself

Protecting your patients, protecting your staff – we take safety seriously. In the operating room every breath the patient takes is paramount to their safety. Safeguarding both patient and clinician during this critical time is of the highest importance. With our Drägersorb CLIC disposable absorber, you can protect your patients from breathing exhaled carbon dioxide thanks to the soda lime absorption, which filters and absorbs the patient's exhaled carbon dioxide. In addition, in just one click, you can easily and quickly change absorbers during a surgical procedure, giving you more time to care for your patient.

Improve infection prevention and maximise economic efficiency: Drägersorb CLIC is a disposal solution in which staff have no direct contact with dust and saves on costs due to its colour indicator that tells you exactly when it is used up.

CONTENT SODA LIME

Soda lime	99
Absorber with CLIC System	100

SODA LIME AND RELATED ACCESSORIES

SODA LIME

A Soda lime Drägersorb® 800+, for CO ₂ absorption, canister, 5 l, 2 pcs.	<u>MX00001</u>
B Soda lime Drägersorb® Free, for CO ₂ absorption, canister, 5 l, 2 pcs.	<u>MX50050</u>
C Canister opener, for canisters MX00001 and MX50050, 5 pcs.	<u>MX10209</u>



SODA LIME AND RELATED ACCESSORIES

ABSORBER WITH CLIC SYSTEM

A	<u>CLIC Absorber 800+, disposable absorber, 1.2 l, 6 pcs.</u>	<u>MX00004</u>
B	<u>CLIC Absorber Free, disposable absorber, 1.2 l, 6 pcs.</u>	<u>MX50100</u>
	<u>Infinity® ID CLIC Absorber 800+, disposable absorber, 1.2 l, 6 pcs.</u>	<u>MX50004</u>
	<u>Infinity® ID CLIC Absorber Free, disposable absorber, 1.2 l, 6 pcs.</u>	<u>MX50120</u>



SODA LIME AND RELATED ACCESSORIES

A	Absorber adapter CLIC, for Dräger anesthesia workstations	MX50090
B	Cleaning plate for CLIC adapter	MX50081
C	Bacteria filter IBF absorber, only use with Zeus®/Zeus® IE, 5 pcs.	MK02588
D	Dust filter Perseus®, 5 pcs.	MX50115





Sensors

Improved diagnostics means improved therapy solutions

Patient data is valuable for both diagnosing and delivering the correct therapies for your patient's pathway to recovery. This is why the reliability and quality of collected data is imperative when diagnosing and monitoring. It must provide you with accurate volume, flow, and pressure values so as to assess your patient's lung condition and reduce ventilator-associated pneumonia (VAP).

Our sensors are rigorously tested for both quality and compatibility with our devices because we know this allows you to reliably monitor every patient's vitals so as to improve your patient outcomes.

CONTENT SENSORS

Flow sensors by Dräger	105
Oxygen sensors by Dräger	107
Oxygen sensor for MX 300-i	107

SENSORS

FLOW SENSORS BY DRÄGER

Hot-wire anemometry:

Accurate and intelligent flow measurement

Spirolog[®], neonatal flow sensors and SpiroLife[®] are based on proven hot-wire anemometry techniques. Due to extremely fast and accurate measurement results, as well as extremely low inspiratory and expiratory resistances, it is clearly superior to other methods. You only need one flow sensor for all patients.


Neonatal flow sensors offer customized solutions for the sole exception of the specific needs of neonates in intensive care. The sterilizable neonatal flow sensors and SpiroLife[®] combine the advantages of hot-wire technology with the possibility of sterilizability. Therefore, they integrate seamlessly into your hospital's sterilization processes.

A	Flow sensor Spirolog[®], PC/ABS, 5 pcs.	8403735
B	Flow sensor SpiroLife[®], PPSU, sterilizable	MK01900
C	Flow sensor Infinity[®] ID, PC/ABS, 5 pcs.	6871980
D	Neonatal flow sensor ISO 15, PSU, sterilizable	8411130
E	Insert for neonatal flow sensor, PSU, sterilizable, 5 pcs.	8410179

Differential pressure measurement: proven technology for emergencies


The sterilizable reusable flow sensor for the emergency room is based on the proven principle of differential pressure measurement. Because of its robustness, it really is the ideal solution for the specialties of the emergency room. The pressure differential sensor is used on Oxylog[®] 3000 plus, Oxylog[®] 3000 and Oxylog[®] 2000 plus. Whichever Dräger flow sensor you choose – you can be sure to have made the right choice with our extensive family of flow sensors.

F	Pressure difference flow sensor, PPSU, sterilizable	8412034
----------	--	----------------

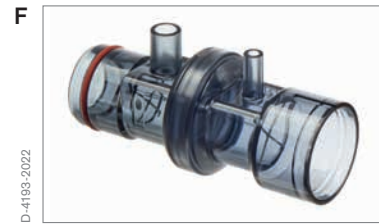

SENSORS

A	<u>Flow sensor cable, for neonatal flow sensor, e.g., Babylog®</u>	<u>8409626</u>
B	<u>MRI flow sensor extension cable, for MRI patient monitor</u>	<u>8306370</u>
C	<u>Flow sensor cable for V-series with Neo option and VN-series</u>	<u>8416600</u>

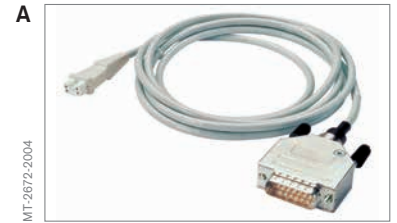

SENSORS

OXYGEN SENSORS BY DRÄGER

Oxygen for Life

Oxygen sensors have the important task of measuring and monitoring the inspiratory and/or expiratory O₂ concentration in anesthesia and respiratory devices. Oxygen sensors from Dräger satisfy the highest national and international requirements. Furthermore, they have been extensively tested and fine-tuned for proper functional performance with Dräger devices and systems.

Dräger's years of experience in the field of sensors, and their high quality has even convinced NASA. Dräger O₂ sensors are electrochemical transducers used to measure O₂ partial pressure.

They operate on the principle of a galvanic cell.

Examples of locations where the sensors are used include the Fabius® family, Evita® family and the Primus®.

To take into account the high value of oxygen measurement, we only recommend using O₂ sensors from Dräger.

A	<u>O₂ sensor capsule, for mainstream measurement</u>	<u>6850645</u>
B	<u>O₂ sensor OxyTrace VE, for Savina® ventilator</u>	<u>MX01049</u>
C	<u>O₂ sensor OxyTrace INCU, for Caleo®</u>	<u>MX01050</u>
D	<u>O₂ sensor Oxycell, for Babyleo® TN500</u>	<u>6873065</u>
E	<u>O₂ sensor housing Oxydig, for O₂ sensor capsule, with spiral cable</u>	<u>6850250</u>
F	<u>O₂ sensor "Oxycell A" for Atlan</u>	<u>6872666</u>

OXYGEN SENSOR FOR MX 300-i

G	<u>O₂ sensor* for MX 300 and MX 300-i, reusable</u>	<u>MP01442</u>
---	--	----------------

*Tradeware





Monitoring accessories

Patient-centred care: Integrated solutions that fit

To support the delivery of quality care, you need quick access to patient data with all relevant information at a glance to reliably assess your patient's status. So as to enable you to get the relevant information and to accurately take care of the patient's needs, we provide you with a broad portfolio of monitoring accessories – all from a single source.

Equally important, to help you break the chain of infections while offering you accuracy and reliability, we provide you with hygienic, single-patient use solutions such as our:

- Electrocardiography leads
- Non-invasive blood pressure (NIBP) cuffs
- Temperature probes, e.g. Tcore sensor

Assisting you in delivering a patient-centred approach, as Your Specialist in Acute Care, we have paid detailed attention in designing our accessories using quality materials that promote comfort and patient satisfaction.

CONTENT MONITORING ACCESSORIES

1. Infinity® Monitoring: Gamma, Delta, Kappa and Vista series	111	3. Infinity® M300	210
1.1 MultiMed® Pods	111	3.1 Electrocardiography (ECG)	210
1.2 Electrocardiography (ECG)	113	3.2 Pulse oximetry (SpO ₂)	213
1.3 Pulse oximetry (SpO ₂)	121	3.3 Pouch	219
1.4 Temperature	139	3.4 Codensation	219
1.5 Tcore® Temperature Monitoring System	141	4. Vista 120	220
1.6 Non-invasive blood pressure (NIBP)	142	4.1 Electrocardiography (ECG)	220
1.7 CNAP continuous non-invasive arterial pressure	146	4.2 Pulse oximetry (SpO ₂)	225
1.8 Invasive blood pressure (IBP) and cardiac output (CO)	148	4.3 Non-invasive blood pressure (NIBP)	227
1.9 Infinity® PiCCO	161	4.4 Temperature	229
1.10 Bispectral index (BIS)	164	4.5 Invasive blood pressure (IBP) and cardiac output (CO)	232
1.11 Electroencephalography (EEG)	165	4.6 CO ₂ Sidestream	234
1.12 Neuromuscular Transmission (NMT)	166	4.7 CO ₂ Mainstream	235
1.13 Capnography etCO ₂	167	5. Vista 120 SPOTCHECK	238
1.14 Scio	172	5.1 Pulse oximetry (SpO ₂)	238
1.15 Recorder paper	173	5.2 Non-invasive blood pressure (NIBP)	248
2. Infinity® Acute Care System	174	5.3 Temperature	250
2.1 Electrocardiography (ECG)	174	5.4 Capnography etCO ₂	252
2.2 Pulse oximetry (SpO ₂)	179	6. Multi Vendor Adapter	254
2.3 Temperature	191	6.1 ECG Adapter	254
2.4 Tcore® Temperature Monitoring System	193	6.2 ECG Adapter Telemetry	256
2.5 Non-invasive blood pressure (NIBP)	194	6.3 TEMP Adapter	256
2.6 Invasive blood pressure (IBP) and cardiac output (CO)	198		
2.7 etCO ₂	205		
2.8 ToFscan® Neuromuscular Transmission Monitor	208		

MULTIMED® PODS

**1. INFINITY® MONITORING:
GAMMA, DELTA, KAPPA AND VISTA SERIES**

1.1 MULTIMED® PODS

A	Multiparameter cable MultiMed® 5 Pod, for ECG, resp., SpO ₂ , and temp., 1.5 m	5950196
	Multiparameter cable MultiMed® 5 Pod, for ECG, resp., SpO ₂ , and temp., 2.5 m	3368391
B	Multiparameter cable MultiMed® 6 Pod, for ECG, resp., SpO ₂ , and temp., 2.5 m	5191221
C	Multiparameter cable MultiMed® plus, for ECG, resp., SpO ₂ , and temp., 2.5 m	MS20093
D	Multiparameter cable MultiMed® plus OR, for ECG, SpO ₂ , and temp., 2.5 m	MS20094
E	Multiparameter cable NeoMed Pod, for ECG, resp., SpO ₂ , temp., and FiO ₂ , 2.5 m	5590539

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL




Not all articles are available worldwide.

MULTIMED® PODS

A	Holder for MultiMed® 12 Pod	5202747
B	Pod connection cable, 30 cm	7257988
	Pod connection cable, 3 m	3368425
	Pod connection cable, 5 m	5195198
C	Holder for MultiMed® 5/6, MultiMed® Plus/OR, NeoMed Pods and Tcore®, swiveling 360°	MP00721
D	Holder for MultiMed® 5/6, MultiMed® Plus/OR, and NeoMed Pods	5598128
E	Universal pole mount	7485621

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

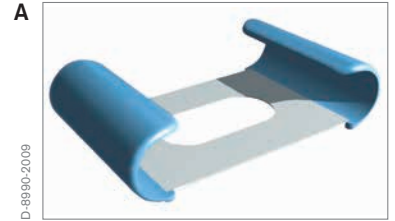

ELECTROCARDIOGRAPHY (ECG)

1.2 ELECTROCARDIOGRAPHY (ECG)

MONOLEAD® DUAL-PIN CONNECTOR (1ST GENERATION)

ECG dual-pin connectors can be connected directly to MultiMed® 5,6 and MultiMed® 12.


A	ECG cable, 3-lead, dual-pin connector MonoLead® 3, IEC1 (European color code), 2 m	MS16160	■	■	■	■	■	■
	ECG cable, 3-lead, dual-pin connector MonoLead® 3, IEC2 (AHA/US color code), 2 m	MS16233	■	■	■	■	■	■
B	ECG cable, 5-lead, dual-pin connector MonoLead® 5, IEC1 (European color code), 2.5 m	MS16161	■	■	■	■	■	■
	ECG cable, 5-lead, dual-pin connector MonoLead® 5, IEC2 (AHA/US color code), 2.5 m	MS16229	■	■	■	■	■	■

CHESTLEADS FOR USE WITH MULTIMED® 12 (1ST GENERATION)

ECG dual-pin connectors can be connected directly to MultiMed® 12.

C	ECG cable for C lead MonoLead® 5, IEC1 (European color code), 1.5 m	MS16232	■	■	■	■	■	■
	ECG cable for V lead MonoLead® 5, IEC2 (AHA/US color code), 1.5 m	MS16230	■	■	■	■	■	■

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL

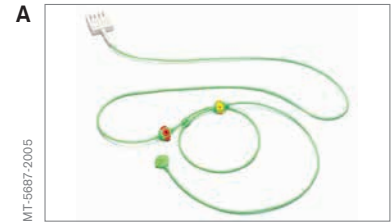

ELECTROCARDIOGRAPHY (ECG)

MONOLEAD® SINGLE-PIN CONNECTOR (2ND GENERATION)

ECG single-pin connectors can be connected directly to MultiMed® plus/MultiMed® plus OR as well as to Infinity® M300 und Infinity® M540.
Note: Not for use with Infinity® TruST™ Telemetrie.

A	ECG MonoLead®, 3-lead, EUR, 2 m	MP05120
B	ECG MonoLead®, 3-lead, AHA/US, 2 m	MP05121
C	ECG MonoLead®, 5-lead, EUR, 2.5 m	MP05122
D	ECG MonoLead®, 5-lead, AHA/US, 2.5 m	MP05123
E	ECG MonoLead®, 6-lead, EUR, 3 m	MP05124
F	ECG MonoLead®, 6-lead, AHA/US, 3 m	MP05125

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL




ELECTROCARDIOGRAPHY (ECG)

MONOLEAD® LONG, SINGLE-PIN CONNECTOR (2ND GENERATION)

A	ECG MonoLead®, 3-lead, EUR, 4.1 m	MP05130
B	ECG MonoLead®, 3-lead, AHA/US, 4.1 m	MP05131
C	ECG MonoLead®, 5-lead, EUR, 5.3 m	MP05132
D	ECG MonoLead®, 5-lead, AHA/US, 5.3 m	MP05133
E	ECG MonoLead®, 6-lead, EUR, 5.7 m	MP05134
F	ECG MonoLead®, 6-lead, AHA/US, 5.7 m	MP05135

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

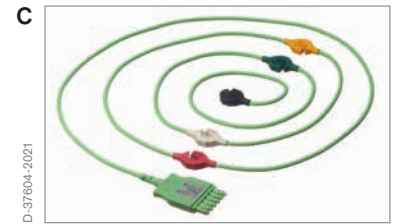

Not all articles are available worldwide.

ELECTROCARDIOGRAPHY (ECG)

SINGLE-PATIENT USE ECG LEADS SINGLE-PIN CONNECTOR

ECG single-pin connectors can be connected directly to MultiMed® plus/MultiMed® plus OR as well as to Infinity® M300 and Infinity® M540.

A ECG cable, 3-lead, single-pin connector, single-patient use, IEC1 (European color code), 1 m, 20 pcs.	MP00875
ECG cable, 3-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1 m, 20 pcs.	MP00877
B ECG cable, 5-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs.	MP00879
ECG cable, 5-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs.	MP00881
C ECG cable, 6-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs.	MP03122
ECG cable, 6-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs.	MP03123

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL



ELECTROCARDIOGRAPHY (ECG)

REUSABLE ECG LEADS SINGLE-PIN CONNECTOR

ECG single-pin connectors can be connected directly to MultiMed® plus/MultiMed® plus OR as well as to Infinity® M300 and Infinity® M540.

A	ECG cable, 3-lead, single-pin connector, IEC1 (European color code), 1 m	MP03401							
	ECG cable, 3-lead, single-pin connector, IEC2 (AHA/US color code), 1 m	MP03402							
B	ECG cable, 5-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03403							
	ECG cable, 5-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03404							
C	ECG cable, 6-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03405							
	ECG cable, 6-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03406							

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL



ELECTROCARDIOGRAPHY (ECG)

REUSABLE ECG LEADS DUAL-PIN CONNECTOR

ECG dual-pin connectors can be connected directly to MultiMed® 5,6 and MultiMed® 12.

A	ECG cable, 3-lead, dual-pin connector, IEC1 (European color code), 1 m	MP03411
	ECG cable, 3-lead, dual-pin connector, IEC2 (AHA/US color code), 1 m	MP03412
B	ECG cable, 5-lead, dual-pin connector, IEC1 (European color code), 1/1.5 m	MP03413
	ECG cable, 5-lead, dual-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03414
C	ECG cable, 6-lead, dual-pin connector, IEC1 (European color code), 1/1.5 m	MP03415
	ECG cable, 6-lead, dual-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03416

CHESTLEADS FOR USE WITH MULTIMED® 12

ECG dual-pin connectors can be connected directly to MultiMed® 12.

D	ECG cable for C lead, IEC1 (European color code), 1 m	MP03417
	ECG cable for V lead, IEC2 (AHA/US color code), 1 m	MP03418

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL



ELECTROCARDIOGRAPHY (ECG)

ELECTRODES

ECG electrodes*,
single-patient use, neonate, 300 pcs. **5195024**

NEONATAL ECG


ECG electrodes*,
single-patient use, neonate, 300 pcs. **5195024**

A Electrode adapter, to connect neonatal electrodes
to MultiMed® pod, 10 pcs. **5194779**

B ECG intermediate cable, 3-lead,
for NeoMed Pod, 1.5 m **5592162**

*Tradeware

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL


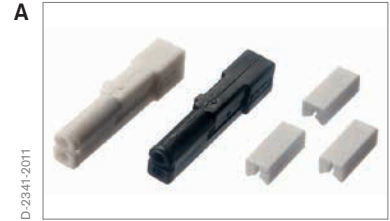


ELECTROCARDIOGRAPHY (ECG)

ECG ACCESSORIES

A	ECG adapter for single-pin ECG cable to MultiMed® 5	MS14679
	ECG adapter for single-pin ECG cable to MultiMed® 6	MS14680
B	ECG extension cable, 6-lead, 2 m	MP05115
C	ECG extension cable, dual-pin connector, 2 m	MS16492
D	ECG ESU block, for OR use of 3- or 5-lead ECG cables with MultiMed® 5	5947226
E	ECG ESU block, for OR use of 3-, 5-, or 6-lead ECG cables with MultiMed® 6	7486140

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

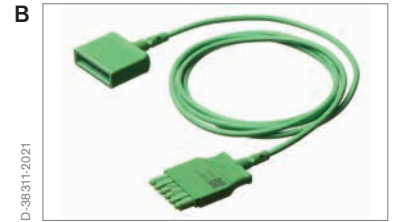

PULSE OXIMETRY (SPO₂)

1.3 PULSE OXIMETRY (SPO₂)

DRÄGER

DRÄGER/NELLCOR® CABLES

SpO₂ intermediate cables for connection to MultiMed® 5, 6, 12 and NeoMed Pods.

A SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 1 m **3368433**

SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 2 m **3375834**

SpO₂ intermediate cables for connection to MultiMed® plus/plus OR and Infinity® M300.

B SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 1.2 m **MS18683**

C SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 3 m **MS17330**

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

GAMMA / GAMMA XL							
GAMMA X XL							
DELTA / DELTA XL							
KAPPA							
KAPPA XLT							
VISTA							
VISTA XL							



Not all articles are available worldwide.

PULSE OXIMETRY (SPO₂)

DRÄGER REUSABLE SENSOR

A SpO₂ finger sensor Dräger, reusable MS13235

DRÄGER DISPOSABLE SENSORS

B SpO₂ adhesive sensor Dräger, vinyl, disposable, adult, 24 pcs. MS16449

C SpO₂ adhesive sensor Dräger, foam, disposable, adult, 24 pcs. MS16445

D SpO₂ adhesive sensor Dräger, vinyl, disposable, child, 24 pcs. MS16448

E SpO₂ adhesive sensor Dräger, foam, disposable, child, 24 pcs. MS16444

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL





PULSE OXIMETRY (SPO₂)

A SpO ₂ adhesive sensor Dräger, cloth, disposable, baby, 24 pcs.	MS16447
B SpO ₂ adhesive sensor Dräger, cloth, disposable, neonate, 24 pcs.	MS16446

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

■	■	■	■	■
■	■	■	■	■

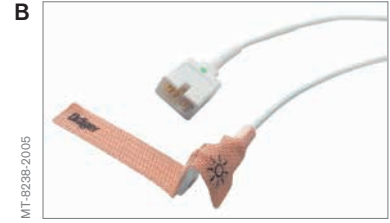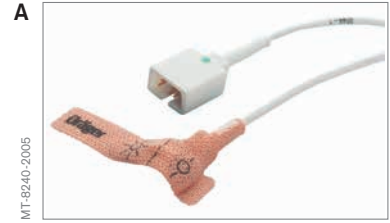

Not all articles are available worldwide.

PULSE OXIMETRY (SpO₂)

DRÄGER/NELLCOR® CABLES

SpO₂ intermediate cables for connection to MultiMed® 5, 6, 12 and NeoMed Pods.

A SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 1 m* **3368433**

SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 2 m* **3375834**

SpO₂ intermediate cables for connection to MultiMed® plus/plus OR.

B SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 1.2 m* **MS18683**

C SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 3 m* **MS17330**

MASIMO®

MASIMO® LNCS CABLES**

SpO₂ intermediate cables for connection to Infinity® M300 and Masimo® SET Pod.

D SpO₂ intermediate cable Masimo® LNCS, 3 m **MS17522**

D SpO₂ intermediate cable Masimo® LNCS, 1.2 m **MS24303**

MASIMO® RD SET CABLES**

Connects the Masimo® SET approved sensors to the Masimo® SET pod or to Infinity M300. This cable cannot be used with the MultiMed.

E SpO₂ intermediate cable Masimo® RD SET, 3 m (9.8 ft) **MS34437**

E SpO₂ intermediate cable Masimo® RD SET, 1.2 m (3.9 ft) **MS34436**

*Note: Requires Nellcor sensor software setting.

**Tradeware

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL



PULSE OXIMETRY (SPO₂)

MASIMO® LNCS REUSABLE SENSORS*

A	SpO ₂ finger sensor Masimo® LNCS DCI, reusable, adult	MP00796
B	SpO ₂ finger sensor Masimo® LNCS DCIP, reusable, child	MP00795
	SpO ₂ Masimo® sensor LNCS-TC-I Ear, reusable, pediatric/adult	MP00788
	SpO ₂ Masimo® sensor LNCS-YI Multisite, reusable, neonate/pediatric/adult	MP00789
	SpO ₂ Masimo® sensor LNCS-TF-I, reusable, pediatric/adult	MP00799
C	SpO ₂ Masimo® sensor LNCS-DBI, reusable, adult	MP02993

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



Not all articles are available worldwide.

PULSE OXIMETRY (SPO₂)

MASIMO® LNCS DISPOSABLE SENSORS*

A	SpO₂ adhesive sensor Masimo® LNCS Amtx, disposable, adult, 20 pcs.	MP00790
B	SpO₂ adhesive sensor Masimo® LNCS Inf, disposable, baby, 20 pcs.	MP00791
C	SpO₂ adhesive sensor Masimo® LNCS Neo, disposable, neonate, 20 pcs.	MP00792
D	SpO₂ adhesive sensor Masimo® LNCS Pmtx, disposable, child, 20 pcs.	MP00793
E	SpO₂ adhesive sensor Masimo® LNCS NeoPt, disposable, premature baby, 20 pcs.	MP00794

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



PULSE OXIMETRY (SPO₂)



MASIMO® LNCS TRAUMA SENSORS*

Masimo® LNCS Trauma sensors are for use with Masimo® SET technology only.

A	Masimo® LNCS Trauma sensor, adult, 20 pcs.	MP03191
B	Masimo® LNCS Trauma sensor, infant/pedi, 20 pcs.	MP03193
C	Masimo® LNCS Trauma sensor, neonate, 20 pcs.	MP03192

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

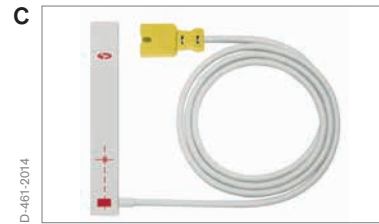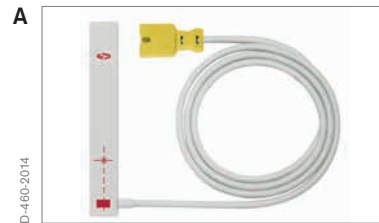

PULSE OXIMETRY (SPO₂)

MASIMO® RD SET ADAPTER CABLES*

<p>A Adapter cable RD SET to LNOP, for connection of RD SET sensors to Infinity® M300 and Masimo® SET Pod via the LNOP intermediate cables MS18680 and MS17041, as well as to MultiMed® 5/6/12, NeoMed Pods and MultiMed® Plus/Plus OR via the according intermediate cables MS13926/7492601 or MS20163**</p>	<p>MS33731</p>
<p>B Adapter cable RD SET to LNCS, for connection of RD SET sensors to Infinity® M300 and Masimo® SET Pod via the LNCS intermediate cables MS24303 or MS17522</p>	<p>MS33732</p>
<p>B Adapter cable RD SET to MultiMed®, for connection of RD SET sensors to MultiMed® 5/6/12, NeoMed Pods as well as MultiMed® Plus/Plus OR via the according intermediate cables 3368433/3375834 or MS18683/MS17330</p>	<p>MS33733</p>

*Tradeware

**Products are not available any more.
Adapter is to be used to switch from sensors to the installed base.

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



PULSE OXIMETRY (SPO₂)

MASIMO® RD SET REUSABLE SENSORS*

A	Masimo® RD SET DCI sensor, adult, reusable	MS33740	■	■	■	■
B	Masimo® RD SET DC-IP sensor, pediatric, reusable	MS33741	■	■	■	■
C	Masimo® RD SET sensor TC-I Ear, pediatric/adult, reusable	MS40067	■	■	■	■
D	Masimo® RD SET sensor YI Multisite, neonate/pediatric/adult, reusable	MS40070	■	■	■	■
E	Masimo® RD SET sensor TF-I, pediatric/adult, reusable	MS40073	■	■	■	■
F	Masimo® RD SET sensor DB-I, adult, reusable	3718728	■	■	■	■

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL





PULSE OXIMETRY (SPO₂)

MASIMO® RD SET DISPOSABLE SENSORS*

A	Masimo® RD SET sensor, adult, disposable	MS33726
B	Masimo® RD SET sensor, pediatric, disposable	MS33727
C	Masimo® RD SET sensor, infant, disposable	MS33728
D	Masimo® RD SET sensor, neonatal, disposable	MS33729
E	Masimo® RD SET sensor, neonatal preterm, disposable	MS33730

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

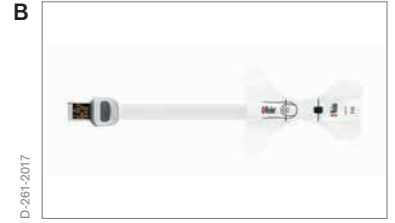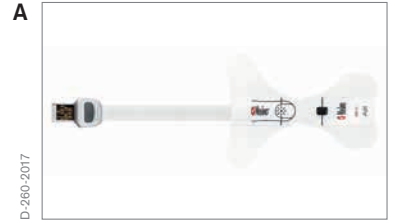

PULSE OXIMETRY (SPO₂)




**MASIMO® RD SET TRAUMA
DISPOSABLE SENSORS***

For use with Masimo® SET technology only.

A	Masimo® Trauma sensor RD SET, adult, disposable	MS40088				
B	Masimo® Trauma sensor RD SET, infant/pedi, disposable	MS40091				
C	Masimo® Trauma sensor RD SET, neonate, disposable	MS40094				

*Tradeware

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL

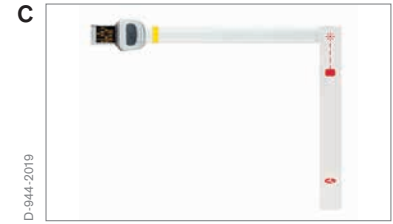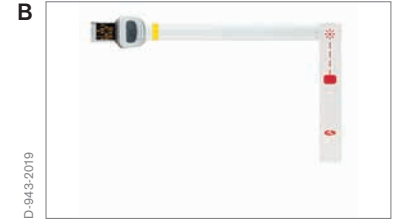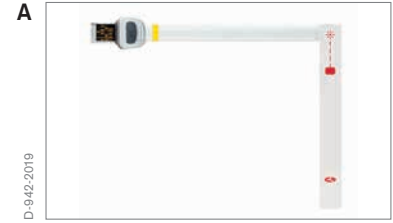



PULSE OXIMETRY (SpO₂)

MASIMO® ACCESSORIES*

A	SpO₂ adhesive wrap Masimo® Multisite-Wrap, for LNOP YI SpO₂ sensor, 100 pcs.	7496917
B	SpO₂ adhesive wrap Masimo® Petite-Wrap, for LNOP YI SpO₂ sensor, 100 pcs.	7496909
C	SpO₂ adhesive square Masimo®, 12 pcs.	7496875
D	SpO₂ foam wrap Masimo® Foam-Wrap, neonate, LNOP YI SpO₂ sensor, 12 pcs.	7496891

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

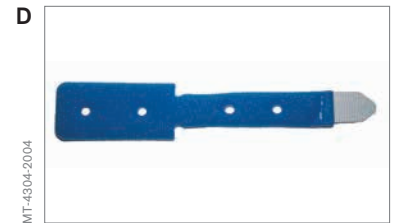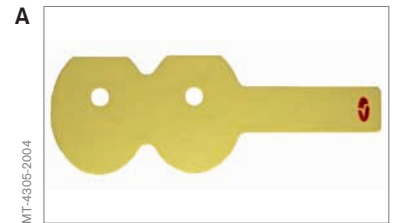

PULSE OXIMETRY (SPO₂)

NELLCOR®

DRÄGER/NELLCOR® CABLES

SpO₂ intermediate cables for connection to MultiMed® 5, 6, 12 and NeoMed Pods.

A SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 1 m **3368433**

SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 2 m **3375834**

SpO₂ intermediate cables for connection to, MultiMed® plus/plus OR, Infinity® M300 and Nellcor® Oximax.

B SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 1.2 m **MS18683**

C SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 3 m **MS17330**

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



PULSE OXIMETRY (SPO₂)

NELLCOR® REUSABLE SENSORS*

A	SpO₂ finger sensor Nellcor® DURA DS 100A	7262764
B	SpO₂ multisite sensor Nellcor® DURA-Y™, reusable	MX01004
C	Ear clip for Nellcor® DURA-Y™	MX01005
D	Adhesive wrap for Nellcor® Oxiband™ A/N, 100 pcs.	8201012
E	Adhesive wrap for Nellcor® Oxiband™ P/I, 100 pcs.	2170221

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



PULSE OXIMETRY (SPO₂)

NELLCOR® DISPOSABLE SENSORS*

A	SpO₂ adhesive sensor Nellcor® Oximax® Max-A, disposable, adult, cable 46 cm, 24 pcs.	MX50065
B	SpO₂ adhesive sensor Nellcor® Oximax® Max-AL, disposable, adult, cable 91 cm, 24 pcs.	MX50071
C	SpO₂ adhesive sensor Nellcor® Oximax® Max-I, disposable, child, 24 pcs.	MX50067
D	SpO₂ adhesive sensor Nellcor® Oximax® Max-N, disposable, neonate/adult, 24 pcs.	MX50068
E	SpO₂ adhesive sensor Nellcor® Oximax® Max-P, disposable, child, 24 pcs.	MX50066

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



TEMPERATURE

1.4 TEMPERATURE


ADAPTER CABLES

A	Adapter cable, 7-pin connector for temperature probe with ¼" connector	5198333
B	Y-cable, reusable, 7-pin connector, for 2 temperature probes	5592154

SINGLE-PATIENT USE TEMPERATURE PROBES

C	General-purpose temperature probe, single-patient use, adult, 7-pin connector, 1.6 m	MP00991
	General-purpose temperature probe, single-patient use, adult, 7-pin connector, 3 m	MP00992
D	General-purpose temperature probe, single-patient use, child, 7-pin connector, 1.6 m	MP00993
	General-purpose temperature probe, single-patient use, child, 7-pin connector, 3 m	MP00994
E	Skin temperature probe, single-patient use, adult, 7-pin connector, 1.6 m	MP00995
	Skin temperature probe, single-patient use, adult, 7-pin connector, 3 m	MP00996
F	Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 1.6 m	MP00997
	Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 3 m	MP00998

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

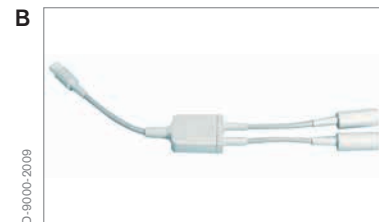

Not all articles are available worldwide.

TEMPERATURE

REUSABLE TEMPERATURE PROBES

A	General-purpose temperature probe, reusable, adult, 7-pin connector, 1.5 m	4329889
A	General-purpose temperature probe, reusable, adult, 7-pin connector, 3 m	5204644
A	General-purpose temperature probe, reusable, child, 7-pin connector, 1.5 m	4329848
A	General-purpose temperature probe, reusable, child, 7-pin connector, 3 m	5204651
B	Skin temperature probe, reusable, with 7-pin connector, 1.5 m	4329822
B	Skin temperature probe, reusable, with 7-pin connector, 3 m	5204669
C	Protective cover for general-purpose temperature probe, contains latex, 10 pcs.	7014616
	Protective cover* for general-purpose temperature probe, latexfree, 100 pcs.	MS40243

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



TCORE® TEMPERATURE MONITORING SYSTEM

1.5 TCORE® TEMPERATURE MONITORING SYSTEM

Tcore® – a new non-invasive technology for accurate core temperature monitoring Tcore® employs a unique dual-sensor heat flux technology, which, following a short ramp-up time, calculates core body temperature continuously and accurately. A simple, self-adhesive sensor – for children from 5 years of age and adults – placed on the patient's forehead is all that's required. This single-use sensor can be connected through a battery-powered adapter to all current Dräger monitors*, eliminating the need for a dedicated display.

*Exception: Vista 120 product family, patient monitors.

		GAMMA / GAMMA XL	GAMMA X XL	DELTA / DELTA XL	KAPPA	KAPPA XLT	VISTA	VISTA XL
A	Tcore® adapter, reusable	MP00999	■	■	■	■	■	■
B	Tcore® sensor, single patient use	MP00989	■	■	■	■	■	■
C	Tcore® extension cable	MP02201	■	■	■	■	■	■
D	Holder for MultiMed® 5/6, MultiMed® Plus/OR, NeoMed Pods and Tcore®, swiveling 360°	MP00721	■	■	■	■	■	■



NON-INVASIVE BLOOD PRESSURE (NIBP)

1.6 NON-INVASIVE BLOOD PRESSURE (NIBP)

NIBP HOSES

A NIBP connection tube, adult, 3.7 m, for Infinity® monitors	MP00953
B NIBP hose, neonate, 2.4 m	2870298

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL




NON-INVASIVE BLOOD PRESSURE (NIBP)

NIBP CUFFS, REUSABLE

A	① NIBP cuff XXS, reusable, infant, 8 – 13 cm/13 cm	MP00911
	② NIBP cuff XS, reusable, child, 12 – 19 cm/19 cm	MP00912
	③ NIBP cuff S, reusable, small adult, 17 – 25 cm/29 cm	MP00913
	④ NIBP cuff M, reusable, adult, 23 – 33 cm/33 cm	MP00915
	⑤ NIBP cuff M+, reusable, adult, long, 23 – 33 cm/43 cm	MP00916
	⑥ NIBP cuff M++, reusable, adult, extra long, 23 – 33 cm/53 cm	MP00917
	⑦ NIBP cuff L, reusable, large adult, 31 – 40 cm/40 cm	MP00918
	⑧ NIBP cuff L+, reusable, large adult, long, 31 – 40 cm/55 cm	MP00919
	⑨ NIBP cuff XL, reusable, thigh, 38 – 50 cm/50 cm	MP00921

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

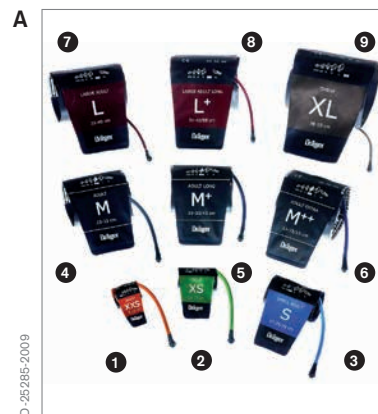


Not all articles are available worldwide.

NON-INVASIVE BLOOD PRESSURE (NIBP)

NIBP CUFFS, SINGLE-PATIENT USE

A	1	NIBP cuff XXS, single-patient use, infant, 8 – 13 cm/13 cm, 10 pcs.	MP00924						
	2	NIBP cuff XS, single-patient use, child, 12 – 19 cm/19 cm, 10 pcs.	MP00925						
	3	NIBP cuff S, single-patient use, small adult, 17 – 25 cm/29 cm, 10 pcs.	MP00926						
	4	NIBP cuff M, single-patient use, adult, 23 – 33 cm/33 cm, 10 pcs.	MP00928						
	5	NIBP cuff M+, single-patient use, adult long, 23 – 33 cm/43 cm, 10 pcs.	MP00929						
	6	NIBP cuff M++, single-patient use, adult extra long, 23 – 33 cm/53 cm, 10 pcs.	MP00930						
	7	NIBP cuff L, single-patient use, large adult, 31 – 40 cm/40 cm, 10 pcs.	MP00931						
	8	NIBP cuff L+, single-patient use, large adult long, 31 – 40 cm/55 cm, 10 pcs.	MP00932						
	9	NIBP cuff XL, single-patient use, thigh, 38 – 50 cm/50 cm, 10 pcs.	MP00934						

GAMMA / GAMMA XL
 GAMMA XXL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

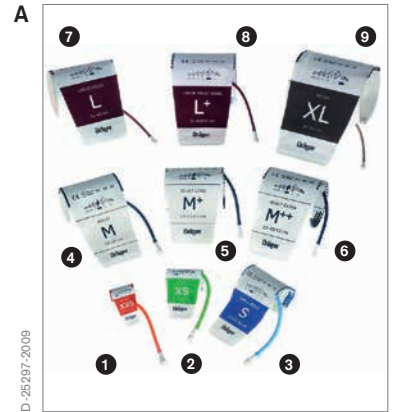

D-25297-2009

NON-INVASIVE BLOOD PRESSURE (NIBP)

NIBP CUFFS NEONATAL, DISPOSABLE

A	NIBP cuff, disposable, neonate, size 1, 3.1 – 5.7 cm, 10 pcs.	2870181
B	NIBP cuff, disposable, neonate, size 2, 4.3 – 8 cm, 10 pcs.	2870199
C	NIBP cuff, disposable, neonate, size 3, 5.8 – 10.9 cm, 10 pcs.	2870207
D	NIBP cuff, disposable, neonate, size 4, 7.1 – 13.1 cm, 10 pcs.	2870215
E	NIBP cuff, disposable, neonate, size 5, 8.3 – 15 cm, 10 pcs.	2870173

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



Not all articles are available worldwide.

CNAP CONTINUOUS NON-INVASIVE ARTERIAL PRESSURE

1.7 CNAP CONTINUOUS NON-INVASIVE ARTERIAL PRESSURE

	CNAP™ cuff controller and cradle	MS26125
A	CNAP™ intermediate cable for Infinity® CNAP™ SmartPod®, for CNAP™ controller	MS15893
B	CNAP™ cuff controller forearm strap	MS26122
C	Cuff controller cradle allows CNAP™ cuff controller to “slide” with movement of patient’s arm	MS26123


GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



CNAP CONTINUOUS NON-INVASIVE ARTERIAL PRESSURE

A	CNAP™ sensor cuff, size L	MS15896
B	CNAP™ sensor cuff, size M	MS15895
C	CNAP™ sensor cuff, size S	MS15894
D	Pod connection cable, 3 m	3368425
E	Universal pole mount	7485621

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL




INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

1.8 INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

IBP AND CO HEMO PODS

A	Infinity® MPod-QuadHemo, cardiac output, 2 temperatures, up to 4 invasive pressures	MS20725
B	Patient parameter unit Infinity® HemoMed Pod	5588822
C	Connection cable, connects HemoMed Pod to monitor, 3 m	5591925
	Connection cable, connects HemoMed Pod to monitor, 5 m	5591933
D	Universal pole mount	7485621

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

IBP QUADHEMO ADAPTERS

A	IBP transducer adapter Abbott/Medex/ICU Medical	MS22535
	IBP transducer adapter Argon/Merit Medical	MS22148
	IBP transducer adapter Baxter/Edwards	MS22147
B	IBP transducer adapter Utah Medical	MS22534
C	IBP transducer adapter cable, Dräger 7-pin	MS22533
D	IBP transducer adapter cable, Dräger 10-pin	MS22532

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL





INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

For use with Infinity® MPod-QuadHemo.

A	IBP transducer adapter cable, Codan/DPT, 2 m	2601235
B	IBP transducer adapter cable, Codan/DPT, 0.3 m	2601234
C	IBP transducer adapter cable, B.Braun, 2 m	2601233
D	IBP transducer adapter cable, B.Braun, 0.3 m	2601232




GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

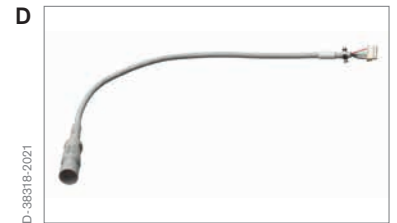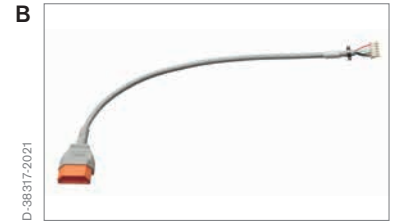

INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

A	IBP transducer adapter cable, Abbott/Medex, Smiths/TranStar, 2 m	2601231
B	IBP transducer adapter cable, Abbott/Medex, Smiths/TranStar, 0.3 m	2601230
C	IBP transducer adapter cable, Utah Medical, 2 m	2601229
D	IBP transducer adapter cable, Utah Medical, 0.3 m	2601228

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

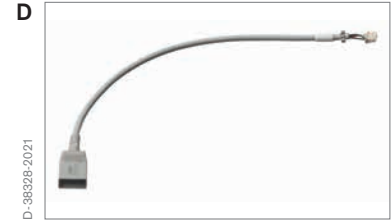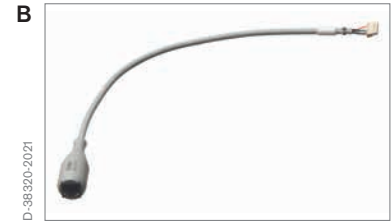



Not all articles are available worldwide.

INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

A	IBP transducer adapter cable, Dräger 7-pin, 2 m	2601227
B	IBP transducer adapter cable, Dräger 7-pin, 0.3 m	2601226
C	IBP transducer adapter cable, Dräger 10-pin, 2 m	2601225
D	IBP transducer adapter cable, Dräger 10-pin, 0.3 m	2601224

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

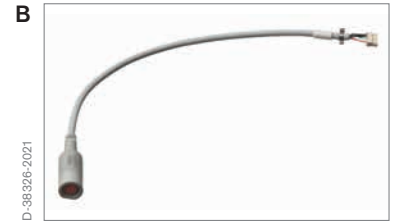



INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

A	IBP transducer adapter cable, Argon/BD/Ohmeda, 2 m	2601223
B	IBP transducer adapter cable, Argon/BD/Ohmeda, 0.3 m	2601222
C	IBP transducer adapter cable, Baxter/Edwards, 2 m	2601221
D	IBP transducer adapter cable, Baxter/Edwards, 0.3 m	2601220
E	IBP transducer adapter cable, Codan/Xtrans, 0.3 m	2601236

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

■	■	■	■	■	■
■	■	■	■	■	■
■	■	■	■	■	■
■	■	■	■	■	■
■	■	■	■	■	■

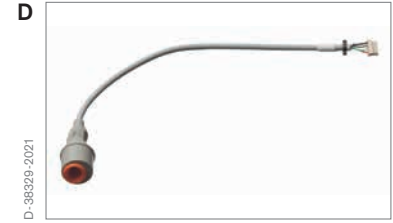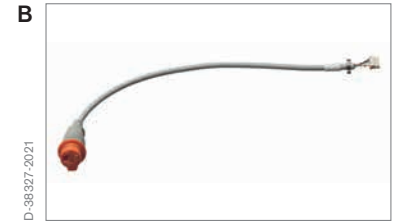


Not all articles are available worldwide.

INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

IBP HEMO POD ADAPTERS AND PLATES

A	Pod adapter, 10-pin, for pressure transducer	3375958
B	Pod adapter, for pressure transducer Abbott/Medex/ICU Medical	5196998
C	Pod adapter, for Edwards pressure transducer	5196980
D	Pod adapter, for pressure transducer Argon/Merit Medical (former Becton Dickinson/Ohmeda)	3375941
	Pod adapter SensoNor, for connecting pressure transducer cables, 7-pin	MS23598

GAMMA / GAMMA XL
 GAMMA XXL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLIT
 VISTA
 VISTA XL


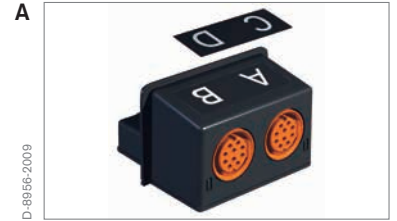


INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

A IBP pressure transducer plate, for Hemo2/4 Pods, for Abbott Transpac IV, 5 pcs.	7270460
B IBP pressure transducer plate, for Hemo2/4 Pods, for Abbott TranspacII/III/Braun, 5 pcs.	5192112
C IBP pressure transducer plate, for Hemo4 Pod, for Memscap (Capto/SensoNor), 5 pcs.	4721416
D IBP pressure transducer plate, for Hemo4 Pod, universal, 5 pcs.	4721424
IBP pressure transducer plate, for Hemo2 Pod, for Memscap (Capto/SensoNor), 5 pcs.	4721614

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

■	■	■	■	■
■	■	■	■	■
■	■	■	■	■
■	■	■	■	■
■	■	■	■	■

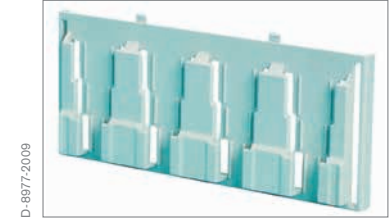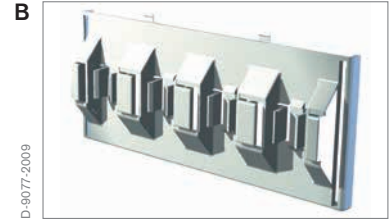

Not all articles are available worldwide.

INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

IBP Y-ADAPTERS AND CABLES

A	IBP Y-adapter, 7-pin	5592147
B	IBP Y-adapter, 10-pin	5731281
C	IBP adapter, 10-pin to 7-pin	3368383
D	IBP Y-adapter, 10-pin to 7-pin	5588095
E	IBP Y-cable, Baxter/Edwards, 3.7 m	5206607

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL




INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

IBP MEMSCAP CABLES AND TRANSDUCERS

A	IBP pressure transducer cable for Memscap (Capto/SensoNor), 3.7 m	4321563
	IBP pressure transducer for Memscap (Capto/SensoNor) 844, disposable, sterile, 50 pcs.	7864908
B	IBP Y-cable, for Memscap (Capto/SensoNor), 3.7 m	5195180
C	IBP pressure transducer set* for Memscap (Capto/SensoNor) 844, disposable, 10 pcs.	7489433
D	IBP holder* for pressure transducer Memscap (Capto/SensoNor) 844, 4 pcs.	7864155
E	IBP pressure transducer* for Memscap (Capto/SensoNor) 844, reusable	7489417
	IBP pressure transducer set* for Memscap (Capto/SensoNor) 844	7498947

*Tradeware

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL

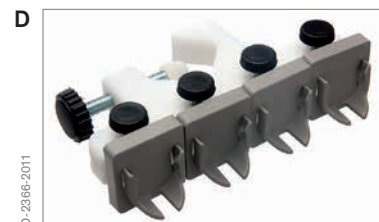


INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

IBP OHMEDA CABLES AND TRANSDUCERS

A	IBP connection cable Argon/Merit Medical, 3.7 m (12.1 ft)	3375933
B	IBP Y-cable Argon/Merit Medical 3.7 m (12.1 ft)	5206581
C	IBP pressure transducer set* Argon/Merit Medical, disposable, sterile, 5 pcs.	4528741

*Tradeware

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL

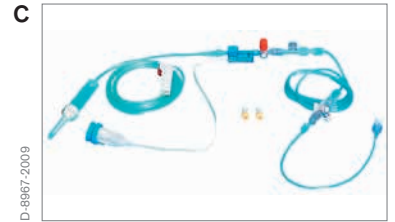

INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

CO CABLES, CATHETERS AND THERMISTORS

A	CO thermistor cable* Argon/Merit Medical	8420077
B	CO thermistor T-piece* Argon/Merit Medical	5741975
C	CO thermistor cable* Baxter/Edwards	8539983
D	CO connection cable, 1 m	3368458
E	CO catheter cable	8419160

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)



Hemodynamic Connectivity

Infinity Acute Care System Infinity M540 MS3987 2.0 M MS3988 2.9 M 	+ HemoMed 7257985 0.3 M 5591826 0.9 M 5591833 5.0 M or HemoMed 7257985 0.3 M 5591826 0.9 M 5591833 5.0 M Monitor IBP Aux/Hemo
Infinity Bedside Monitors Infinity Gamma X XL Vista XL 	+ HemoMed Aux/Hemo 2, 3
Infinity Delta/XL 	+ HemoMed Aux/Hemo 2, 3
Infinity Kappa 	+ HemoMed Aux/Hemo 2, 3
Infinity Kappa XLT *PICCO not supported 	+ HemoMed Aux/Hemo 2, 3

Dual Hemo MS20783 A B	+ MS20725 Quad Hemo A B C D	+ Hemo2 4319435 Sold out A B	+ Hemo4 4315861 Sold out A B C D	+ PICCO MS16734 Kappa XLT does not support PICCO Sold out C D	+ PICCO Gen II MS23133 Kappa XLT does not support PICCO Gen II A B C D	+ HemoMed 5589822 A B C D
--	--	--	--	--	--	--

Dräger 7 Pin 2601227 2M 2601226 2M 2601224 2M Dräger 10 Pin 2601225 2M	+ B. Braun Medical 2601222 2M 2601220 2M Argon 2601221 2M 2601223 2M Baster 2601230 2M 2601229 2M Abbott SM Translar 2601231 2M 2601228 2M Coflex DPT 2601234 2M 2601235 2M Coflex Xtrans 2601236 2M 2601235 2M	+ B. Braun 8 Pin MS23598 7 Pin 3375968 10 Pin MS23598 7 Pin 3375968 10 Pin Abbott-Medex 5198988 5198980 Biotek Enox 5198980 Argon/IBD 3375941	+ Reusable Adult 4329827 1.5m 5204669 3m Pediatric 5204644 1.5m 5204651 3m Single Patient Use Adult MP00991 1.6m MP00992 3m Pediatric MP00993 1.6m MP00994 3m	+ Argon/IBD 2 C.O. 1 2 3386458 849180 8420077 5741975 Biotek Edwards 2 8539963	+ C.O. 1 2 MS161616 MS161616 MS161619 PICCO
Vendor specific Hemo Pod adapters [AB-CD]			Temperature		
Vendor specific adapter cables [A-B-C-D]			Cardiac Output (C.O.)		



INFINITY® PiCCO

1.9 INFINITY® PiCCO


A	Infinity® PiCCO Gen II SmartPod®	MS23133	■	■		
	Infinity® PiCCO Gen II SmartPod® kit	MS32282	■	■		
B	CO thermistor cable* for Infinity® PiCCO Pod	MS16919	■	■		
	CO catheter cable for Infinity® PiCCO Pod	MS16918	■	■		
	CO cable for Infinity® PiCCO Pod	MS16916	■	■		
C	IBP transducer adapter cable, Dräger 10-pin	MS22532	■	■		
D	Pod adapter, 10-pin, for pressure transducer	3375958	■	■		
	Adapter cable* for Infinity® PiCCO SmartPod®, for 10-pin connector	MS16920	■	■		

*Tradeware


GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

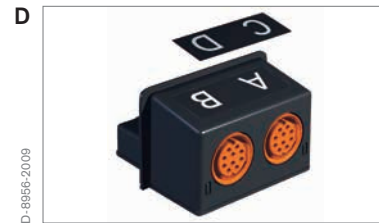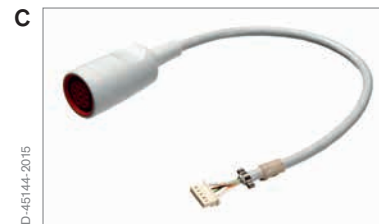

INFINITY® PiCCO

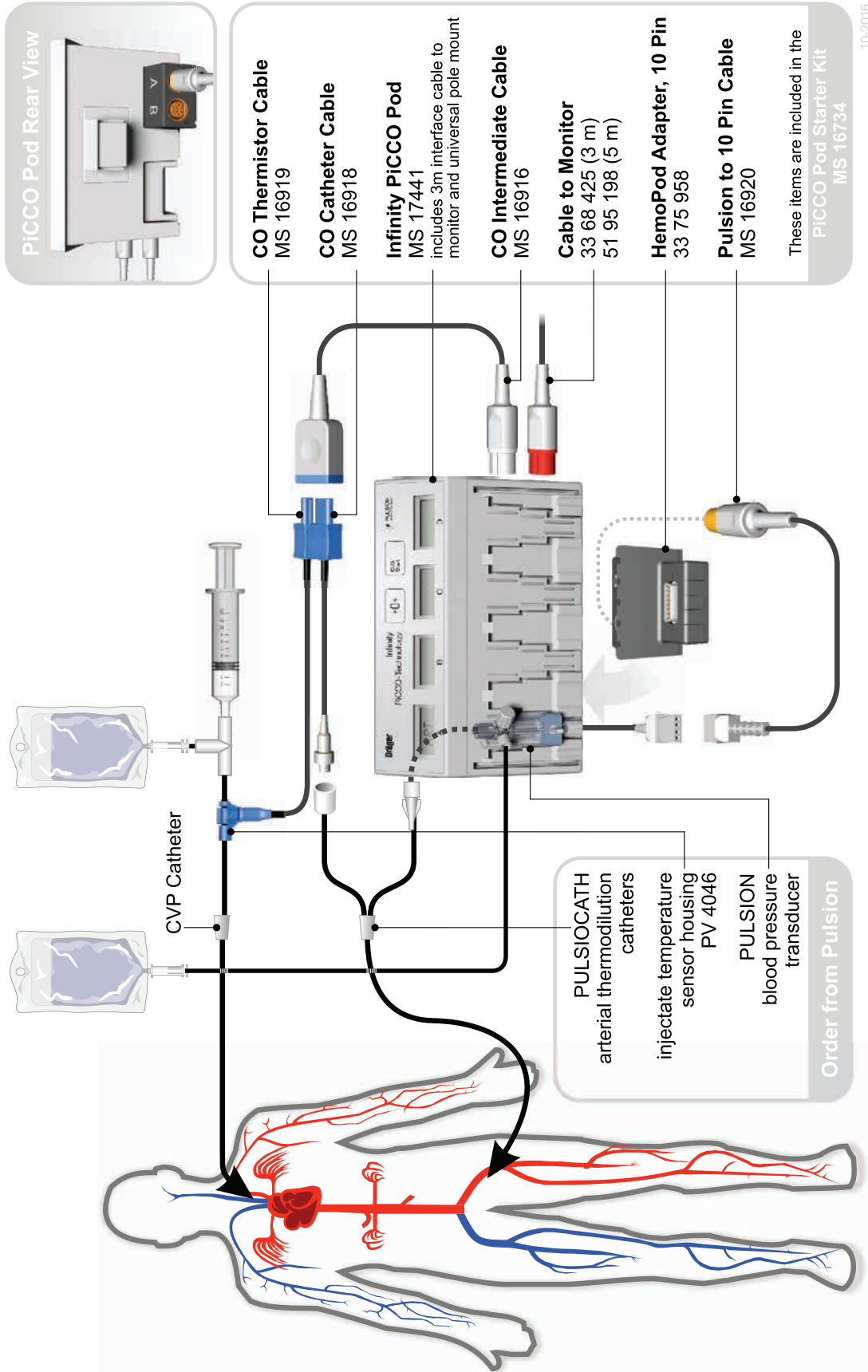

INFINITY® PiCCO GEN II SMARTPOD®

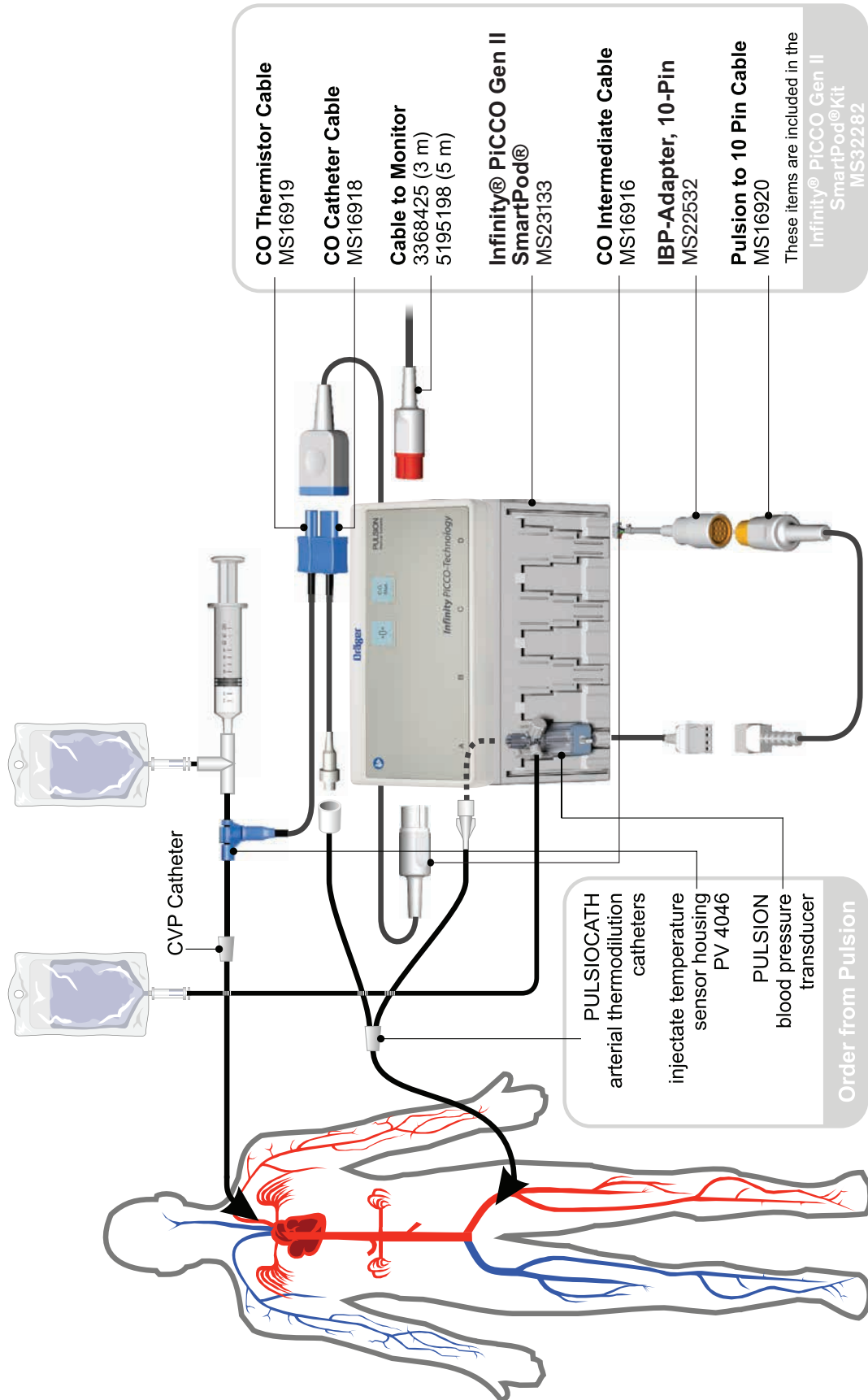

Not all articles are available worldwide.

BISPECTRAL INDEX (BIS)

1.10 BISPECTRAL INDEX (BIS)

INFINITY® BISX SMARTPOD®

A	Patient parameter unit Infinity® BISx SmartPod®	MS14796
---	--	---------

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



ELECTROENCEPHALOGRAPHY (EEG)


1.11 ELECTROENCEPHALOGRAPHY (EEG)

INFINITY® EEG POD AND ACCESSORIES

A	EEG electrode cable, color-coded, 9 pcs.	5947804
B	Block for EEG referential lead conversion	5954859
C	EEG electrodes*, disposable, 25 pcs.	MS29468

*Tradeware

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

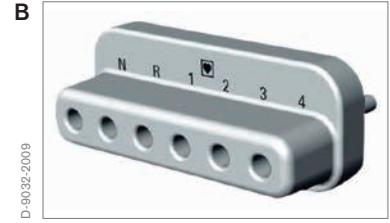


NEUROMUSCULAR TRANSMISSION (NMT)

1.12 NEUROMUSCULAR TRANSMISSION (NMT)

INFINITY® TRIDENT (NMT) SMARTPOD® AND ACCESSORIES

A	NMT accelerometer Trident®	MS15085
B	NMT thermistor Trident®	MS15053
C	NMT hand adapter Trident®	MS15086
D	NMT intermediate cable for Infinity® Trident (NMT) SmartPod®, 1.9 m	7872174
	NMT electrode cable Trident®, 2-lead, 0.3 m	MS13218
E	Bed clip for cables	MS15087

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

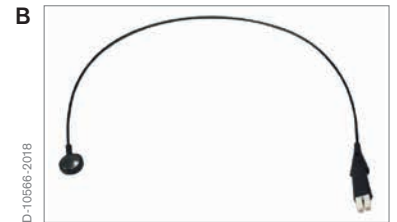

CAPNOGRAPHY ETCO₂

CO₂ CABLES

A MCable Mainstream CO₂,
for Evita V-family and Babylog VN-family,
Infinity® M540, Vista120, Vista120 S,
and Delta VF9.1 and higher **6871950**

B Y-cable CO₂ and analog/sync **MS22259**

CO₂ CUVETTES, DISPOSABLE

C CO₂ cuvette disposable, adult, 10 pcs. **MP01062**

D CO₂ cuvette disposable, pediatric, 10 pcs. **MP01063**

CO₂ CUVETTES, REUSABLE

E CO₂ cuvette reusable, adult **6870279**

F CO₂ cuvette reusable, pediatric **6870280**

GAMMA / GAMMA XL
GAMMA X XL
DELTA / DELTA XL
KAPPA
KAPPA XLT
VISTA
VISTA XL




CAPNOGRAPHY ETCO₂

ACCESSORIES

Calibration set, with test connection and test gas cylinder	8412710
Test gas cylinder 5 Vol.% CO ₂ , 95 Vol. % N ₂	6850435
A Parking holder for CO ₂ sensor	8416352
Test filter	6870281

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

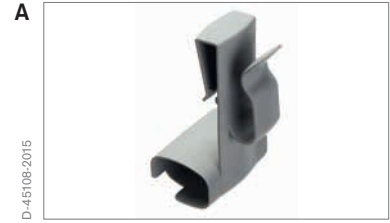

CAPNOGRAPHY ETCO₂

A	etCO₂ Microstr. Intubated Adult/Pediatric - Extended Duration, 25 pcs.	MP15001
	etCO₂ Microstr. Intubated Adult/Pediatric - Short-term use, 25 pcs.	MP15002
B	etCO₂ Microstr. Intubated Infant/Neonatal - Extended Duration, 25 pcs.	MP15003
C	etCO₂ Microstr. Nasal Adult - Extended Duration, 25 pcs.	MP15004
	etCO₂ Microstr. Nasal Adult - Short-term use, 25 pcs.	MP15005
D	etCO₂ Microstr. Nasal Pediatric - Short-term use, 25 pcs.	MP15006
E	etCO₂ Microstr. Nasal Infant/Neonatal - Extended Duration, 25 pcs.	MP15007
F	etCO₂ Microstr. Nasal/Oral Adult - Short-term use, 25 pcs.	MP15008
G	etCO₂ Microstr. Nasal/Oral Pediatric - Short-term use, 25 pcs.	MP15009


GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL



CAPNOGRAPHY ETCO₂

A	CO₂ Microstr. O₂ Nasal Adult - Extended Duration, 25 pcs.	MP15010
	CO₂ Microstr. O₂ Nasal Pediatric - Extended Duration, 25 pcs.	MP15011
	CO₂ Microstr. O₂ Nasal/Oral Adult - Short-term use, 25 pcs.	MP15012
B	CO₂ Microstr. O₂ Nasal/Oral Pediatric - Short-term use, 25 pcs.	MP15013

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

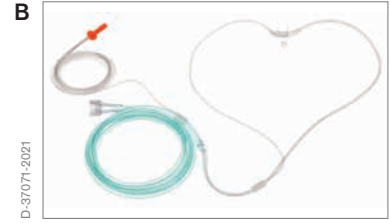

Not all articles are available worldwide.

SCIO

1.14 SCIO

SCIO ACCESSORIES

A	Sample line, disposable, 10 pcs.	8290286
B	Watertrap WaterLock® 2, 12 pcs.	6872130
	Water trap Infinity® ID WaterLock® 2, 12 pcs.	6872020
	Sample gas return line for Fabius®, with filter	M32692
C	Bacterial filter for Dräger patient monitors, disposable	8402868




RECORDER PAPER

1.15 RECORDER PAPER

A Recording paper for Infinity® R 50, 50 mm, approx. 30 m, 10 rolls	4711201
---	----------------

GAMMA / GAMMA XL
 GAMMA X XL
 DELTA / DELTA XL
 KAPPA
 KAPPA XLT
 VISTA
 VISTA XL

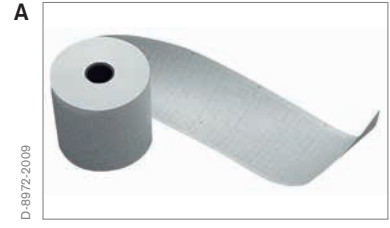

INFINITY® ACUTE CARE SYSTEM

2. INFINITY® ACUTE CARE SYSTEM

2.1 ELECTROCARDIOGRAPHY (ECG)



MONOLEAD® SINGLE-PIN CONNECTOR (2ND GENERATION)

A	ECG MonoLead®, 3-lead, EUR, 2 m	MP05120
B	ECG MonoLead®, 3-lead, AHA/US, 2 m	MP05121
C	ECG MonoLead®, 5-lead, EUR, 2.5 m	MP05122
D	ECG MonoLead®, 5-lead, AHA/US, 2.5 m	MP05123
E	ECG MonoLead®, 6-lead, EUR, 3 m	MP05124
F	ECG MonoLead®, 6-lead, AHA/US, 3 m	MP05125

MONOLEAD® LONG, SINGLE-PIN CONNECTOR (2ND GENERATION)*

A	ECG MonoLead®, 3-lead, EUR, 4.1 m	MP05130
B	ECG MonoLead®, 3-lead, AHA/US, 4.1 m	MP05131
C	ECG MonoLead®, 5-lead, EUR, 5.3 m	MP05132
D	ECG MonoLead®, 5-lead, AHA/US, 5.3 m	MP05133
E	ECG MonoLead®, 6-lead, EUR, 5.7 m	MP05134
F	ECG MonoLead®, 6-lead, AHA/US, 5.7 m	MP05135

*Not for use with extension cables MS16256 and MP05115.

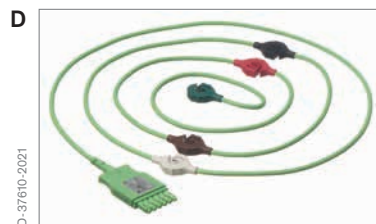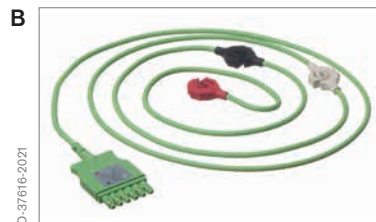

INFINITY® ACUTE CARE SYSTEM

CHESTLEADS MONOLEAD® (2ND GENERATION)

A	ECG MonoLead®, 4-lead, EUR, 1.7 m	MP05126
B	ECG MonoLead®, 4-lead, AHA/US, 1.7 m	MP05127
A	ECG MonoLead®, 4-lead, EUR, 3.7 m*	MP05128
B	ECG MonoLead®, 4-lead, AHA/US, 3.7 m*	MP05129

*Not for use with extension cables MS31435 and MP05114.



INFINITY® ACUTE CARE SYSTEM

REUSABLE ECG LEADS SINGLE-PIN CONNECTOR

A	ECG cable, 3-lead, single-pin connector, IEC1 (European color code), 1 m	MP03401
	ECG cable, 3-lead, single-pin connector, IEC2 (AHA/US color code), 1 m	MP03402
B	ECG cable, 5-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03403
	ECG cable, 5-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03404
C	ECG cable, 6-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03405
	ECG cable, 6-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03406

CHESTLEADS

D	ECG cable, 4-lead, single-pin connector, IEC1 (European color code), 1 m	MP03407
	ECG cable, 4-lead, single-pin connector, IEC2 (AHA/US color code), 1 m	MP03408



INFINITY® ACUTE CARE SYSTEM

**SINGLE-PATIENT USE ECG LEADS
SINGLE-PIN CONNECTOR**

A ECG cable, 3-lead, single-pin connector, single-patient use, IEC1 (European color code), 1 m, 20 pcs. **MP00875**

ECG cable, 3-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1 m, 20 pcs. **MP00877**

B ECG cable, 5-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs. **MP00879**

ECG cable, 5-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs. **MP00881**

C ECG cable, 6-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs. **MP03122**

ECG cable, 6-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs. **MP03123**

**SINGLE-PATIENT USE ECG LEADS
SINGLE-PIN CONNECTOR FOR CHESTLEADS**

D ECG cable, 4-lead, single-pin connector, single-patient use, IEC1 (European color code), 1 m, 20 pcs. **MP03120**

ECG cable, 4-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1 m, 20 pcs. **MP03121**



Not all articles are available worldwide.

INFINITY® ACUTE CARE SYSTEM

ELECTRODES

ECG electrodes*, single-patient use, neonate, 300 pcs. 5195024

NEONATAL ECG

ECG electrodes*, single-patient use, neonate, 300 pcs. 5195024

ECG neonatal 3-lead adapter cable, 2 m MS25951

ECG ACCESSORIES

A ECG adapter for single-pin ECG cable to MultiMed® 5 MS14679


ECG adapter for single-pin ECG cable to MultiMed® 6 MS14680

B ECG "4-lead" spacer, 10 pcs. MS22832

C ECG extension cable, 6-lead, 2 m MP05115

D ECG extension cable, 4-lead, 2 m MP05114

*Tradeware

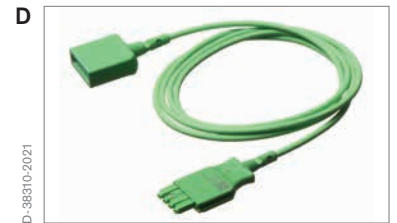

INFINITY® ACUTE CARE SYSTEM

2.2 PULSE OXIMETRY (SpO₂)

MASIMO®

For Masimo® accessories, specific part numbers for USA and Japan apply. Please contact your local Dräger Sales representative for further information.

MASIMO® LNCS CABLES*

For connection of LNCS sensors to Masimo® SET Mcable.

A	<u>SpO₂ intermediate cable Masimo® LNCS, 1.2 m</u>	<u>MP00800</u>
A	<u>SpO₂ intermediate cable Masimo® LNCS, 3 m</u>	<u>MP00797</u>
A	<u>SpO₂ intermediate cable Masimo® LNCS, 4.3 m</u>	<u>MP03194</u>

For connection of LNCS sensors to Masimo® rainbow SET Mcable.

B	<u>SpO₂ intermediate cable Masimo® Red LNCS, 1.2 m</u>	<u>MP02991</u>
B	<u>SpO₂ intermediate cable Masimo® Red LNCS, 3 m</u>	<u>MP02992</u>
B	<u>SpO₂ intermediate cable Masimo® Red LNCS, 4.3 m</u>	<u>MP03196</u>

*Tradeware



INFINITY® ACUTE CARE SYSTEM

MASIMO® LNCS REUSABLE SENSORS*

A	SpO₂ finger sensor Masimo® LNCS DCI, reusable, adult	MP00796
B	SpO₂ finger sensor Masimo® LNCS DCIP, reusable, child	MP00795
	SpO₂ Masimo® sensor LNCS-TC-I Ear, reusable, pediatric/adult	MP00788
	SpO₂ Masimo® sensor LNCS-TF-I, reusable, pediatric/adult	MP00799
	SpO₂ Masimo® sensor LNCS-YI Multisite, reusable, neonate/pediatric/adult	MP00789
C	SpO₂ Masimo® sensor LNCS-DBI, reusable, adult	MP02993

*Tradeware



INFINITY® ACUTE CARE SYSTEM

MASIMO® DISPOSABLE SENSORS*

A SpO ₂ adhesive sensor Masimo® LNCS Adtx, disposable, adult, 20 pcs.	MP00790
B SpO ₂ adhesive sensor Masimo® LNCS Inf, disposable, baby, 20 pcs.	MP00791
C SpO ₂ adhesive sensor Masimo® LNCS Neo, disposable, neonate, 20 pcs.	MP00792
D SpO ₂ adhesive sensor Masimo® LNCS Pdtx, disposable, child, 20 pcs.	MP00793
E SpO ₂ adhesive sensor Masimo® LNCS NeoPt, disposable, premature baby, 20 pcs.	MP00794

*Tradeware




INFINITY® ACUTE CARE SYSTEM

MASIMO® LNCS TRAUMA SENSORS*

Masimo® LNCS Trauma sensors are for use with Masimo® SET technology only.

A Masimo® LNCS Trauma sensor, adult, 20 pcs.	MP03191
B Masimo® LNCS Trauma sensor, infant/pedi, 20 pcs.	MP03193
C Masimo® LNCS Trauma sensor, neonate, 20 pcs.	MP03192

*Tradeware

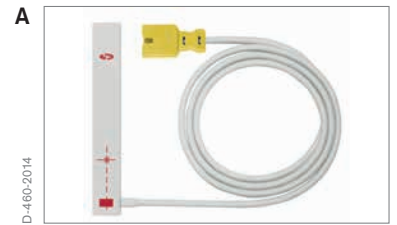

INFINITY® ACUTE CARE SYSTEM

MASIMO® RD SET CABLES*

For connection of RD SET sensors and rainbow sensors to Masimo® rainbow SET Mcable.

A RD SET intermediate cable for Masimo® SET Mcable, 1.5 m MS33736

A RD SET intermediate cable for Masimo® SET Mcable, 3.6 m MS33737

For connection of RD SET sensors to Masimo® rainbow SET Mcable.

B RD SET intermediate cable for Masimo® rainbow SET Mcable, 1.5 m MS33738

B RD SET intermediate cable for Masimo® rainbow SET Mcable, 3.6 m MS33739

MASIMO® RD SET REUSABLE SENSORS*

C Masimo® RD SET DCI sensor, adult, reusable MS33740

D Masimo® RD SET DC-IP sensor, pediatric, reusable MS33741

E Masimo® RD SET sensor TC-I Ear, pediatric/adult, reusable MS40067

F Masimo® RD SET sensor YI Multisite, neonate/pediatric/adult, reusable MS40070

G Masimo® RD SET sensor TF-I, pediatric/adult, reusable MS40073

Masimo® RD SET sensor DB-I, adult, reusable 3718728

*Tradeware



Not all articles are available worldwide.

INFINITY® ACUTE CARE SYSTEM

MASIMO® RD SET DISPOSABLE SENSORS*


A	Masimo® RD SET sensor, adult, disposable	MS33726
B	Masimo® RD SET sensor, pediatric, disposable	MS33727
C	Masimo® RD SET sensor, infant, disposable	MS33728
D	Masimo® RD SET sensor, neonatal, disposable	MS33729
E	Masimo® RD SET sensor, neonatal preterm, disposable	MS33730

MASIMO® RD SET TRAUMA DISPOSABLE SENSORS*

For use with Masimo® SET technology only.

F	Masimo® Trauma sensor RD SET, adult, disposable	MS40088
G	Masimo® Trauma sensor RD SET, infant/pedi, disposable	MS40091
	Masimo® Trauma sensor RD SET, neonate, disposable	MS40094

*Tradeware

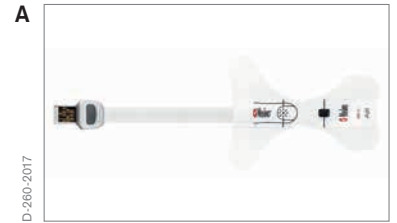

INFINITY® ACUTE CARE SYSTEM

MASIMO® RAINBOW CABLES*

For connection of RD SET sensors and rainbow sensors to Masimo® rainbow SET Mcable.

A RD SET intermediate cable for Masimo® rainbow SET Mcable, 1.5 m	MS33738
A RD SET intermediate cable for Masimo® rainbow SET Mcable, 3.6 m	MS33739





MASIMO® RAINBOW DISPOSABLE SENSORS (SPCO, SPMET, SPOC, PVI)*

B Masimo rainbow® sensor 8 Lambda (SpO ₂ , SpCO, SpMet), disposable, adult, 10 pcs.	3702841
C Masimo rainbow® sensor 8 Lambda (SpO ₂ , SpCO, SpMet), disposable, child, 10 pcs.	3702837
D Masimo rainbow® sensor 8 Lambda (SpO ₂ , SpCO, SpMet), disposable, infant, 10 pcs.	3702834
E Masimo rainbow® sensor 8 Lambda (SpO ₂ , SpCO, SpMet), disposable, adult/neonate, 10 pcs.	3702831

MASIMO® RAINBOW REUSABLE SENSORS (SPCO, SPMET, SPOC, PVI)*

F Masimo rainbow® sensor DCI (SpO ₂ , SpCO, SpMet), reusable, adult	3722471
G Masimo rainbow® sensor DCIP (SpO ₂ , SpCO, SpMet), reusable, pediatric	3722467

*Tradeware

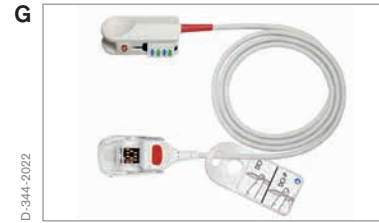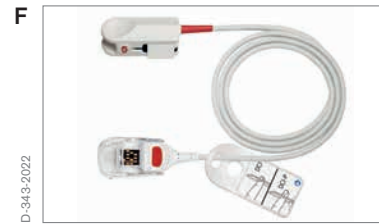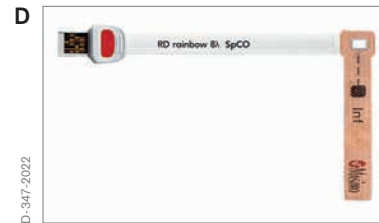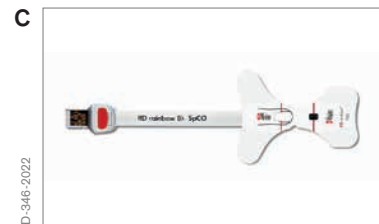

Not all articles are available worldwide.

INFINITY® ACUTE CARE SYSTEM

MASIMO® RD RAINBOW SET CABLE*

For connection of RD SET sensors and rainbow sensors to Masimo® rainbow SET Mcable.

A RD SET intermediate cable for Masimo® rainbow SET Mcable, 1.5 m MS33738

A RD SET intermediate cable for Masimo® rainbow SET Mcable, 3.6 m MS33739

MASIMO® RD RAINBOW DISPOSABLE SENSORS (SPHB, SPMET, SPOC, PVI)*

B Masimo® RD rainbow SET sensor, disposable, adult MS40076

C Masimo® RD rainbow SET sensor, disposable, pediatric MS40079

D Masimo® RD rainbow SET sensor, disposable, infant MS40082

E Masimo® RD rainbow SET sensor, disposable, neonate MS40085

*Tradeware




INFINITY® ACUTE CARE SYSTEM

MASIMO® ACCESSORIES*

A SpO ₂ adhesive wrap Masimo® Multisite-Wrap, for LNOP YI SpO ₂ sensor, 100 pcs.	7496917
B SpO ₂ adhesive wrap Masimo® Petite-Wrap, for LNOP YI SpO ₂ sensor, 100 pcs.	7496909
C SpO ₂ adhesive square Masimo®, 12 pcs.	7496875
D SpO ₂ foam wrap Masimo® Foam-Wrap, neonate, LNOP YI SpO ₂ sensor, 12 pcs.	7496891

*Tradeware

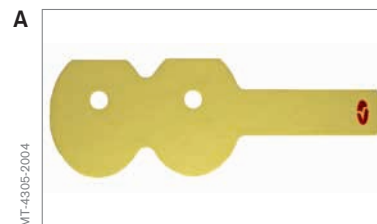

INFINITY® ACUTE CARE SYSTEM

A	Headband* LNCS TF-I, 5 pcs.	MP03197
B	Adhesive Pad* LNCS TF-I, 20 pcs.	MP03198
C	Masimo® rainbow sensor cover*, 10 pcs.	MP03203

*Tradeware



INFINITY® ACUTE CARE SYSTEM



Masimo SpO₂ and rainbow Sensor Connectivity

INFINITY PATIENT MONITORS

M540

Masimo SET MCable MS20867

Masimo rainbow SET MCable MS27003 (order via MS27000)

rainbow Sensors for SpO₂, PLS, PI, SpMet, SpOC, PVI and SpCO

Reusable

RD rainbow	MS33786 1.5 m
RD rainbow	MS33757 3.0 m
RD rainbow	MS33738 1.5 m
RD rainbow	MS33739 3.0 m
RD rainbow SET 2 Pdt	MS40076
RD rainbow SET 2 Neo	MS40085

Disposable

RD rainbow	MP09800 1.2 m
RD rainbow	MP09801 3.0 m
RD rainbow	MP09814 4.5 m
RD rainbow SET 2 Pdt	MS40079
RD rainbow SET 2 Neo	MS40085

Gamma XL (VF6, VF7)

Masimo SET Pod MS16358

Multimed Plus MS20993 1.2 m
MS20993 3.0 m
Multimed Plus OR MS20994 1.2 m
MS20994 3.0 m

Multimed 5 MS53723

Multimed 12 5589663 3.0 m

Vista XL (VF6, VF7, VF8)

Masimo SET Pod MS16359

Multimed Plus MS20993 1.2 m
MS20993 3.0 m
Multimed Plus OR MS20994 1.2 m
MS20994 3.0 m

Multimed 5 MS53723

Multimed 12 5589663 3.0 m

Delta/Delta XL/Kappa (VF6, VF7, VF8)

Masimo SET Pod MS16356

Multimed Plus MS20993 1.2 m
MS20993 3.0 m
Multimed Plus OR MS20994 1.2 m
MS20994 3.0 m

Multimed 5 MS53723

Multimed 12 5589663 3.0 m

rainbow SET for SpO₂, PLS, PI, SpMet, SpOC, PVI and SpCO

Reusable

RD rainbow	rainbow DCI 8L 372471
RD rainbow	rainbow DCIP BL 372467
RD rainbow	rainbow DCIP BL 372467
RD rainbow	Neo BL 3702831
RD rainbow	Neo BL 3702837
RD rainbow	rainbow DCI 8L 372471
RD rainbow	rainbow DCIP BL 372467
RD rainbow	rainbow DCIP BL 372467
RD rainbow	Neo BL 3702831
RD rainbow	Neo BL 3702837

Disposable

RD rainbow SET 2 Adt	MS40076
RD rainbow SET 2 Pdt	MS40079
RD rainbow SET 2 Inf	MS40082
RD rainbow SET 2 Neo	MS40085

RD SET Sensors for SpO₂, PLS, PI, PVI**

Disposable

RD SET Adt	MS33726
RD SET Inf	MS33728
RD SET Neo	MS33729
RD SET NeoP	MS33730
RD SET Pdt	MS33727
RD SET SpHb	MS33722 3.0 m
RD SET SpHb	MS33722 1.2 m
RD SET SpHb	MS33722 3.0 m
RD SET SpHb	MS33722 3.0 m

Reusable

RD SET DCI A41	MS33740
RD SET DCIP Pdt	MS33741
RD SET Inf	MS40067
RD SET Inf	MS40067
RD SET Neo	MS40070
RD SET NeoP	MS40070
RD SET Pdt	MS40070
RD SET Pdt	MS40070
RD SET SpHb	MS40070
RD SET SpHb	MS40070

LNCS Sensors for SpO₂, PLS, PI, PVI**

Disposable

LNCS Adt	MP00790
LNCS Inf	MP00791
LNCS Neo	MP00792
LNCS NeoP	MP00794
LNCS DCI	MP00796
LNCS DCI	MP00795
LNCS TC-I	MP00788
LNCS TC-I	MP00788
LNCS YI	MP00789
LNCS YI	MP00789
LNCS DB-I	MP00788
LNCS DB-I	MP00788

Reusable

LNCS Adt	MP00790
LNCS Inf	MP00791
LNCS Neo	MP00792
LNCS NeoP	MP00794
LNCS DCI	MP00796
LNCS DCI	MP00795
LNCS TC-I	MP00788
LNCS TC-I	MP00788
LNCS YI	MP00789
LNCS YI	MP00789
LNCS DB-I	MP00788
LNCS DB-I	MP00788

LNCS Trauma Sensors
Sensors are only for use with Masimo SET (not via Multimed Plus)
Trauma LNCS Adt/Pdt
Trauma LNCS Neo
Trauma LNCS NeoP

LNCS Adt	MP03191
LNCS Inf	MP03192
LNCS Neo	MP03193
LNCS NeoP	MP03192

* Requires Nellcor sensor software setting.
** PVI works only with Masimo rainbow SET MCable.

Masimo, LNCS, RD SET, SET, SpHb, SpMet, SpCO, and SpOC are registered trademarks of Masimo Corporation. Nellcor is a registered trademark of Nellcor Puritan Bennett LLC.

Mentioned trademarks are only registered in certain countries and not necessarily in the country in which this material is released.

INFINITY® ACUTE CARE SYSTEM

NELLCOR®

NELLCOR® CABLES*

A	Nellcor® intermediate cable, 1.2 m	MP00748
	Nellcor® intermediate cable, 3 m	MS20979

NELLCOR® REUSABLE SENSOR*

B	SpO₂ finger sensor Nellcor® DURA DS 100A	7262764
	SpO₂ Nellcor® FLEXMAX Adt	MS34677
	SpO₂ Nellcor® FLEXMAX Pdt	MS34678

NELLCOR® DISPOSABLE SENSORS*

C	SpO₂ adhesive sensor Nellcor® Oximax® Max-A, disposable, adult, cable 46 cm, 24 pcs.	MX50065
D	SpO₂ adhesive sensor Nellcor® Oximax® Max-AL, disposable, adult, cable 91 cm, 24 pcs.	MX50071
E	SpO₂ adhesive sensor Nellcor® Oximax® Max-I, disposable, child, 24 pcs.	MX50067
F	SpO₂ adhesive sensor Nellcor® Oximax® Max-N, disposable, neonate/adult, 24 pcs.	MX50068
G	SpO₂ adhesive sensor Nellcor® Oximax® Max-P, disposable, child, 24 pcs.	MX50066

*Tradeware




INFINITY® ACUTE CARE SYSTEM

2.3 TEMPERATURE

ADAPTER CABLE

A	Y-cable, reusable, 7-pin connector, for 2 temperature probes	5592154
----------	--	----------------

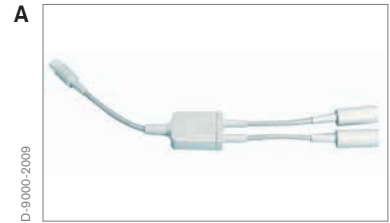

SINGLE-PATIENT USE TEMPERATURE PROBES

B	General-purpose temperature probe, single-patient use, adult, 7-pin connector, 1.6 m	MP00991
----------	--	----------------



	General-purpose temperature probe, single-patient use, adult, 7-pin connector, 3 m	MP00992
--	--	----------------



C	General-purpose temperature probe, single-patient use, child, 7-pin connector, 1.6 m	MP00993
----------	--	----------------

	General-purpose temperature probe, single-patient use, child, 7-pin connector, 3 m	MP00994
--	--	----------------

D	Skin temperature probe, single-patient use, adult, 7-pin connector, 1.6 m	MP00995
----------	---	----------------

	Skin temperature probe, single-patient use, adult, 7-pin connector, 3 m	MP00996
--	---	----------------

E	Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 1.6 m	MP00997
----------	--	----------------



	Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 3 m	MP00998
--	--	----------------



Not all articles are available worldwide.

INFINITY® ACUTE CARE SYSTEM

REUSABLE TEMPERATURE PROBES

A	General-purpose temperature probe, reusable, adult, 7-pin connector, 1.5 m	4329889
A	General-purpose temperature probe, reusable, adult, 7-pin connector, 3 m	5204644
A	General-purpose temperature probe, reusable, child, 7-pin connector, 1.5 m	4329848
A	General-purpose temperature probe, reusable, child, 7-pin connector, 3 m	5204651
B	Skin temperature probe, reusable, with 7-pin connector, 1.5 m	4329822
B	Skin temperature probe, reusable, with 7-pin connector, 3 m	5204669
C	Protective cover for general-purpose temperature probe, contains latex, 10 pcs.	7014616
	Protective cover* for general-purpose temperature probe, latexfree, 100 pcs.	MS40243

*Tradeware



INFINITY® ACUTE CARE SYSTEM

2.4 Tcore® TEMPERATURE MONITORING SYSTEM

Tcore® – a new non-invasive technology for accurate core temperature monitoring Tcore® employs a unique dual-sensor heat flux technology, which, following a short ramp-up time, calculates core body temperature continuously and accurately. A simple, self-adhesive sensor – for children from 5 years of age and adults – placed on the patient's forehead is all that's required. This single-use sensor can be connected through a battery-powered adapter to all current Dräger monitors*, eliminating the need for a dedicated display.

*Exception: Vista 120 product family, patient monitors.

A Tcore® adapter, reusable	MP00999
B Tcore® sensor, single patient use	MP00989
C Tcore® extension cable	MP02201
Tcore® mount for M540	MS33239
D Holder for MultiMed® 5/6, MultiMed® Plus/OR, NeoMed Pods and Tcore®, swiveling 360°	MP00721



INFINITY® ACUTE CARE SYSTEM

2.5 NON-INVASIVE BLOOD PRESSURE (NIBP)

NIBP HOSES


A NIBP connection tube, adult, 3.7 m, for Infinity® monitors	MP00953
B NIBP hose, neonate, 2.4 m	2870298



INFINITY® ACUTE CARE SYSTEM

NIBP CUFFS, REUSABLE


A	① NIBP cuff XXS, reusable, infant, 8 – 13 cm/13 cm	MP00911
	② NIBP cuff XS, reusable, child, 12 – 19 cm/19 cm	MP00912
	③ NIBP cuff S, reusable, small adult, 17 – 25 cm/29 cm	MP00913
	④ NIBP cuff M, reusable, adult, 23 – 33 cm/33 cm	MP00915
	⑤ NIBP cuff M+, reusable, adult, long, 23 – 33 cm/43 cm	MP00916
	⑥ NIBP cuff M++, reusable, adult, extra long, 23 – 33 cm/53 cm	MP00917
	⑦ NIBP cuff L, reusable, large adult, 31 – 40 cm/40 cm	MP00918
	⑧ NIBP cuff L+, reusable, large adult, long, 31 – 40 cm/55 cm	MP00919
	⑨ NIBP cuff XL, reusable, thigh, 38 – 50 cm/50 cm	MP00921

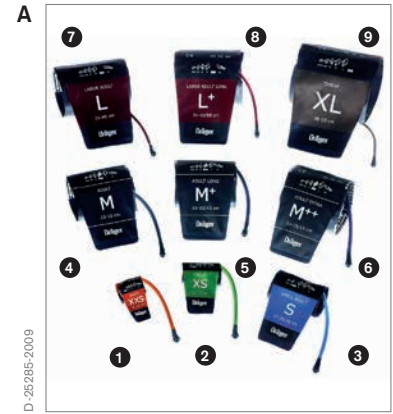

INFINITY® ACUTE CARE SYSTEM

NIBP CUFFS, SINGLE-PATIENT USE

A	① NIBP cuff XXS, single-patient use, infant, 8 – 13 cm/13 cm, 10 pcs.	MP00924
	② NIBP cuff XS, single-patient use, child, 12 – 19 cm/19 cm, 10 pcs.	MP00925
	③ NIBP cuff S, single-patient use, small adult, 17 – 25 cm/29 cm, 10 pcs.	MP00926
	④ NIBP cuff M, single-patient use, adult, 23 – 33 cm/33 cm, 10 pcs.	MP00928
	⑤ NIBP cuff M+, single-patient use, adult long, 23 – 33 cm/43 cm, 10 pcs.	MP00929
	⑥ NIBP cuff M++, single-patient use, adult extra long, 23 – 33 cm/53 cm, 10 pcs.	MP00930
	⑦ NIBP cuff L, single-patient use, large adult, 31 – 40 cm/40 cm, 10 pcs.	MP00931
	⑧ NIBP cuff L+, single-patient use, large adult long, 31 – 40 cm/55 cm, 10 pcs.	MP00932
	⑨ NIBP cuff XL, single-patient use, thigh, 38 – 50 cm/50 cm, 10 pcs.	MP00934

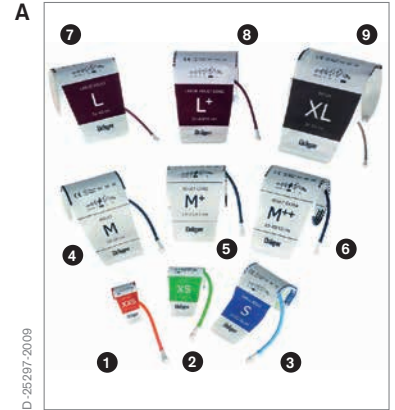

INFINITY® ACUTE CARE SYSTEM

NIBP CUFFS NEONATAL, DISPOSABLE

A NIBP cuff, disposable, neonate, size 1, 3.1 – 5.7 cm, 10 pcs.	2870181
B NIBP cuff, disposable, neonate, size 2, 4.3 – 8 cm, 10 pcs.	2870199
C NIBP cuff, disposable, neonate, size 3, 5.8 – 10.9 cm, 10 pcs.	2870207
D NIBP cuff, disposable, neonate, size 4, 7.1 – 13.1 cm, 10 pcs.	2870215
E NIBP cuff, disposable, neonate, size 5, 8.3 – 15 cm, 10 pcs.	2870173



INFINITY® ACUTE CARE SYSTEM

2.6 INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT (CO)

A	MPod-Dual Hemo, 2 invasive blood pressures for M540*	MS20783
B	Infinity® MPod-QuadHemo, cardiac output, 2 temperatures, up to 4 invasive pressures*	MS20725
Required to connect MPod-QuadHemo, Hemo2 or Hemo4 to IACS:		
C	M540 Hemo connection cable, 1.5 m	MS23897
	M540 Hemo connection cable, 3 m	MS23898
Required to connect two MPod-QuadHemo:		
	Pod connection cable, 30 cm	7257988
	Pod connection cable, 3 m	3368425
	Pod connection cable, 5 m	5195198



*Adapter cable not included.



INFINITY® ACUTE CARE SYSTEM

IBP TRANSDUCER ADAPTER CABLES


A	IBP transducer adapter Abbott/Medex/ICU Medical	MS22535
	IBP transducer adapter Argon/Merit Medical	MS22148
	IBP transducer adapter Baxter/Edwards	MS22147
B	IBP transducer adapter Utah Medical	MS22534
C	IBP transducer adapter cable, Dräger 7-pin	MS22533
D	IBP transducer adapter cable, Dräger 10-pin	MS22532

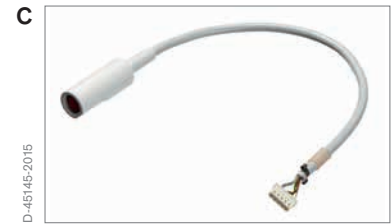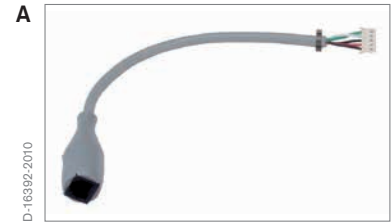

INFINITY® ACUTE CARE SYSTEM


For use with Infinity® MPod-QuadHemo.

A	<u>IBP transducer adapter cable, Codan/DPT, 2 m</u>	<u>2601235</u>
B	<u>IBP transducer adapter cable, Codan/DPT, 0.3 m</u>	<u>2601234</u>
C	<u>IBP transducer adapter cable, B.Braun, 2 m</u>	<u>2601233</u>
D	<u>IBP transducer adapter cable, B.Braun, 0.3 m</u>	<u>2601232</u>

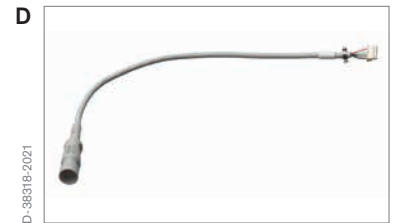



INFINITY® ACUTE CARE SYSTEM

A IBP transducer adapter cable, Abbott/Medex, Smiths/TranStar, 2 m	2601231
B IBP transducer adapter cable, Abbott/Medex, Smiths/TranStar, 0.3 m	2601230
C IBP transducer adapter cable, Utah Medical, 2 m	2601229
D IBP transducer adapter cable, Utah Medical, 0.3 m	2601228

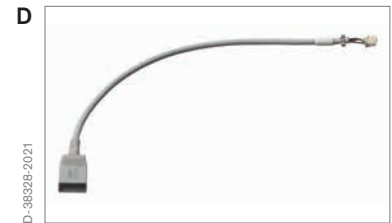

INFINITY® ACUTE CARE SYSTEM

A	IBP transducer adapter cable, Dräger 7-pin, 2 m	2601227
B	IBP transducer adapter cable, Dräger 7-pin, 0.3 m	2601226
C	IBP transducer adapter cable, Dräger 10-pin, 2 m	2601225
D	IBP transducer adapter cable, Dräger 10-pin, 0.3 m	2601224

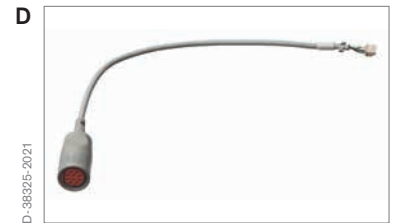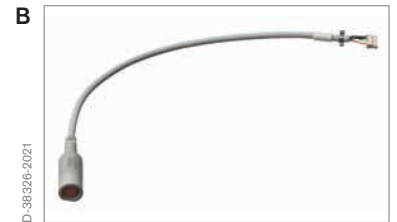

INFINITY® ACUTE CARE SYSTEM

A	IBP transducer adapter cable, Argon/BD/Ohmeda, 2 m	2601223
B	IBP transducer adapter cable, Argon/BD/Ohmeda, 0.3 m	2601222
C	IBP transducer adapter cable, Baxter/Edwards, 2 m	2601221
D	IBP transducer adapter cable, Baxter/Edwards, 0.3 m	2601220
E	IBP transducer adapter cable, Codan/Xtrans, 0.3 m	2601236



INFINITY® ACUTE CARE SYSTEM

CO CABLES AND THERMISTORS

A	CO thermistor cable* Argon/Merit Medical	8420077
B	CO thermistor T-piece* Argon/Merit Medical	5741975
C	CO thermistor cable Baxter/Edwards	8539983
D	CO connection cable, 1 m	3368458
E	CO catheter cable	8419160

CARDIAC OUTPUT KITS

	Cardiac output kit, Baxter	7265155
	Cardiac output kit, Argon/Merit Medical	7265163

*Tradeware



INFINITY® ACUTE CARE SYSTEM

2.7 ETCO₂

CO₂ CABLES

A	MCable Mainstream CO₂, for Evita V-family and Babylog VN-family, Infinity® M540, Vista120, Vista120 S, and Delta VF9.1 and higher	6871950
B	Y-cable CO₂ and analog/sync	MS22259

CO₂ CUVETTES, DISPOSABLE


C	CO₂ cuvette disposable, adult, 10 pcs.	MP01062
D	CO₂ cuvette disposable, pediatric, 10 pcs.	MP01063

CO₂ CUVETTES, REUSABLE

E	CO₂ cuvette reusable, adult	6870279
F	CO₂ cuvette reusable, pediatric	6870280

ACCESSORIES



	Calibration set, with test connection and test gas cylinder	8412710
	Test gas cylinder 5 Vol.% CO₂, 95 Vol. % N₂	6850435
G	Parking holder for CO₂ sensor	8416352
	Test filter	6870281

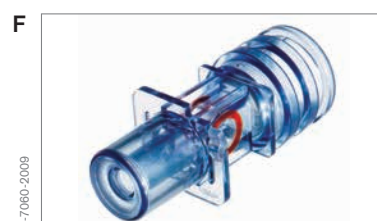

Not all articles are available worldwide.


INFINITY® ACUTE CARE SYSTEM

A	etCO₂ Microstr. Intubated Adult/Pediatric - Extended Duration, 25 pcs.	MP15001
	etCO₂ Microstr. Intubated Adult/Pediatric - Short-term use, 25 pcs.	MP15002
B	etCO₂ Microstr. Intubated Infant/Neonatal - Extended Duration, 25 pcs.	MP15003
C	etCO₂ Microstr. Nasal Adult - Extended Duration, 25 pcs.	MP15004
	etCO₂ Microstr. Nasal Adult - Short-term use, 25 pcs.	MP15005
D	etCO₂ Microstr. Nasal Pediatric - Short-term use, 25 pcs.	MP15006
E	etCO₂ Microstr. Nasal Infant/Neonatal - Extended Duration, 25 pcs.	MP15007
F	etCO₂ Microstr. Nasal/Oral Adult - Short-term use, 25 pcs.	MP15008
G	etCO₂ Microstr. Nasal/Oral Pediatric - Short-term use, 25 pcs.	MP15009

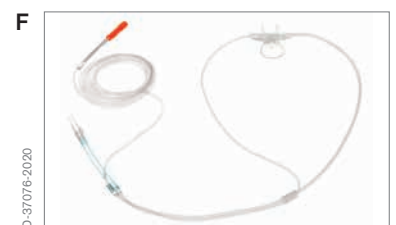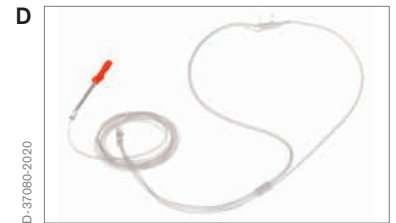

INFINITY® ACUTE CARE SYSTEM

A CO ₂ Microstr. O ₂ Nasal Adult - Extended Duration, 25 pcs.	MP15010
CO ₂ Microstr. O ₂ Nasal Pediatric - Extended Duration, 25 pcs.	MP15011
CO ₂ Microstr. O ₂ Nasal/Oral Adult - Short-term use, 25 pcs.	MP15012
B CO ₂ Microstr. O ₂ Nasal/Oral Pediatric - Short-term use, 25 pcs.	MP15013

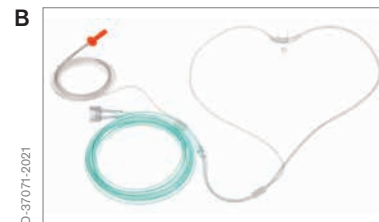

INFINITY® ACUTE CARE SYSTEM

2.8 TOFSCAN® NEUROMUSCULAR TRANSMISSION MONITOR*

ToFscan® is a nerve stimulator module for the measurement of neuromuscular transmission, via accelerometry. It provides an easy, accurate way to measure of patient muscle relaxation. ToFscan® supports evidence-based judgment about the state of a patient's neuromuscular blockade.

A	ToFscan® Kit configurator	MS32318
	Hand sensor, reusable, adult	MS32154
B	Hand sensor, reusable, adult, as of Software 1.8	MS40317
	Hand sensor, reusable, pediatric	MS32289
C	Hand sensor, reusable, pediatric, as of Software 1.8	MS40318
D	Hand sensor, reusable, small pediatric, as of Software 1.8	MS40319
	Eyebrow sensor, reusable	MS32291
E	Eyebrow sensor, reusable, as of Software 1.8	MS40321
	Tape for Eyebrow Sensor, 96 pcs.	MS34209
	Foot sensor, reusable	MS32290
F	Foot sensor, reusable, as of Software 1.8	MS40322
G	Disposable hand sensor, adult, as of Software 1.8.5	MS40320
	Intermediate cable for disposable sensor 3 m, reusable	MS40323

*Tradeware

Not all articles are available worldwide.

TOFSCAN® MONITOR
TOFSCAN® MONITOR,
SOFTWARE 1.8



INFINITY® ACUTE CARE SYSTEM

TOFSCAN® NEUROMUSCULAR TRANSMISSION MONITOR*

A	Extension cable	MS32153
B	Mount clamp 10 – 40 mm (0.4 – 1.57 in)	MS32150
	Mount clamp 20 – 60 mm (0.79 – 2.36 in)	MS32151
C	Optic serial RS 232 cable (1 meter) to connect the device to a Dräger monitorAC	MS32156
	Optic serial RS 232 cable (2.5 meter) to connect the device to a Dräger monitorAC	MS40021

*Tradeware

TOFSCAN® MONITOR
TOFSCAN® MONITOR,
SOFTWARE 1.8



INFINITY® M300

3. INFINITY® M300

3.1 ELECTROCARDIOGRAPHY (ECG)

MONOLEAD® SINGLE-PIN CONNECTOR (2ND GENERATION)

A	ECG MonoLead®, 3-lead, EUR, 2 m	MP05120
B	ECG MonoLead®, 3-lead, AHA/US, 2 m	MP05121
C	ECG MonoLead®, 5-lead, EUR, 2.5 m	MP05122
D	ECG MonoLead®, 5-lead, AHA/US, 2.5 m	MP05123
E	ECG MonoLead®, 6-lead, EUR, 3 m	MP05124
F	ECG MonoLead®, 6-lead, AHA/US, 3 m	MP05125



INFINITY® M300

ECG-LEAD SETS FOR TELEMETRY

Infinity® TruST™ Telemetry lead sets can also be used in combination with MultiMed® plus. Shielded, not for use in the OR.

A	ECG cable, 3-lead Infinity® TruST™, IEC1 (European color code), 70 cm	MP03421
	ECG cable, 3-lead Infinity® TruST™, IEC2 (AHA/US color code), 70 cm	MP03422
B	ECG cable, 5-lead Infinity® TruST™, IEC1 (European color code), 70 cm	MP03423
	ECG cable, 5-lead Infinity® TruST™, IEC2 (AHA/US color code), 70 cm	MP03424
C	ECG cable, 6-lead Infinity® TruST™, IEC1 (European color code), 70 cm	MP03425
	ECG cable, 6-lead Infinity® TruST™, IEC2 (AHA/US color code), 70 cm	MP03426

REUSABLE ECG LEADS SINGLE-PIN CONNECTOR

D	ECG cable, 3-lead, single-pin connector, IEC1 (European color code), 1 m	MP03401
	ECG cable, 3-lead, single-pin connector, IEC2 (AHA/US color code), 1 m	MP03402
E	ECG cable, 5-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03403
	ECG cable, 5-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03404
F	ECG cable, 6-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03405
	ECG cable, 6-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03406



INFINITY® M300

SINGLE-PATIENT USE ECG LEADS SINGLE-PIN CONNECTOR

A	ECG cable, 3-lead, single-pin connector, single-patient use, IEC1 (European color code), 1 m, 20 pcs.	MP00875
	ECG cable, 3-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1 m, 20 pcs.	MP00877
B	ECG cable, 5-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs.	MP00879
	ECG cable, 5-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs.	MP00881
C	ECG cable, 6-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs.	MP03122
	ECG cable, 6-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs.	MP03123



INFINITY® M300

3.2 PULSE OXIMETRY (SPO₂)

DRÄGER

DRÄGER CABLES

A	SpO ₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 1.2 m	MS18683
B	SpO ₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 3 m	MS17330



Not all articles are available worldwide.

INFINITY® M300

DRÄGER REUSABLE SENSOR

A	SpO₂ finger sensor Dräger, reusable	MS13235
----------	---	----------------

DRÄGER DISPOSABLE SENSORS

B	SpO₂ adhesive sensor Dräger, vinyl, disposable, adult, 24 pcs.	MS16449
----------	--	----------------

C	SpO₂ adhesive sensor Dräger, foam, disposable, adult, 24 pcs.	MS16445
----------	---	----------------

D	SpO₂ adhesive sensor Dräger, vinyl, disposable, child, 24 pcs.	MS16448
----------	--	----------------

E	SpO₂ adhesive sensor Dräger, foam, disposable, child, 24 pcs.	MS16444
----------	---	----------------



INFINITY® M300

MASIMO®

For Masimo® accessories, specific part numbers for USA and Japan apply. Please contact your local Dräger Sales representative for further information.

MASIMO® LNCS CABLES*

A	SpO ₂ intermediate cable Masimo® LNCS, 3 m	MS17522
A	SpO ₂ intermediate cable Masimo® LNCS, 1.2 m	MS24303

MASIMO® LNCS REUSABLE SENSORS*

B	SpO ₂ finger sensor Masimo® LNCS DCI, reusable, adult	MP00796
C	SpO ₂ finger sensor Masimo® LNCS DCIP, reusable, child	MP00795
	SpO ₂ Masimo® sensor LNCS-TC-I Ear, reusable, pediatric/adult	MP00788
	SpO ₂ Masimo® sensor LNCS-YI Multisite, reusable, neonate/pediatric/adult	MP00789
	SpO ₂ Masimo® sensor LNCS-TF-I, reusable, pediatric/adult	MP00799
D	SpO ₂ Masimo® sensor LNCS-DBI, reusable, adult	MP02993

MASIMO® LNCS DISPOSABLE SENSORS*

E	SpO ₂ adhesive sensor Masimo® LNCS Amtx, disposable, adult, 20 pcs.	MP00790
F	SpO ₂ adhesive sensor Masimo® LNCS Inf, disposable, baby, 20 pcs.	MP00791
G	SpO ₂ adhesive sensor Masimo® LNCS Pmtx, disposable, child, 20 pcs.	MP00793

*Tradeware




INFINITY® M300

MASIMO® RD SET ADAPTER CABLES*

<p>A Adapter cable RD SET to LNOP, for connection of RD SET sensors to Infinity® M300 and Masimo® SET Pod via the LNOP intermediate cables MS18680 and MS17041, as well as to MultiMed® 5/6/12, NeoMed Pods and MultiMed® Plus/Plus OR via the according intermediate cables MS13926/7492601 or MS20163**</p>	<p>MS33731</p>
<p>B Adapter cable RD SET to LNCS, for connection of RD SET sensors to Infinity® M300 and Masimo® SET Pod via the LNCS intermediate cables MS24303 or MS17522</p>	<p>MS33732</p>
<p>SpO₂ intermediate cable Masimo® RD SET, 3 m (9.8 ft)</p>	<p>MS34437</p>
<p>SpO₂ intermediate cable Masimo® RD SET, 1.2 m (3.9 ft)</p>	<p>MS34436</p>
<h3>MASIMO® RD SET DISPOSABLE SENSORS*</h3>	
<p>C Masimo® RD SET sensor, adult, disposable</p>	<p>MS33726</p>
<p>D Masimo® RD SET sensor, pediatric, disposable</p>	<p>MS33727</p>
<h3>MASIMO® RD SET REUSABLE SENSORS*</h3>	
<p>E Masimo® RD SET DCI sensor, adult, reusable</p>	<p>MS33740</p>
<p>F Masimo® RD SET DC-IP sensor, pediatric, reusable</p>	<p>MS33741</p>
<p>G Masimo® RD SET sensor DB-I, adult, reusable</p>	<p>3718728</p>

*Tradeware

**Products are not available any more. Adapter is to be used to switch from sensors to the installed base.

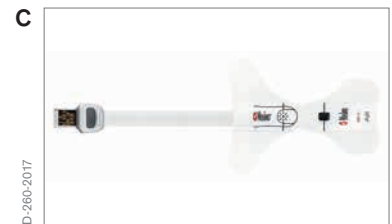


INFINITY® M300

MASIMO® ACCESSORIES*

A SpO ₂ adhesive wrap Masimo® Multisite-Wrap, for LNOP YI SpO ₂ sensor, 100 pcs.	7496917
B SpO ₂ adhesive wrap Masimo® Petite-Wrap, for LNOP YI SpO ₂ sensor, 100 pcs.	7496909
C SpO ₂ adhesive square Masimo®, 12 pcs.	7496875
D SpO ₂ foam wrap Masimo® Foam-Wrap, neonate, LNOP YI SpO ₂ sensor, 12 pcs.	7496891

*Tradeware

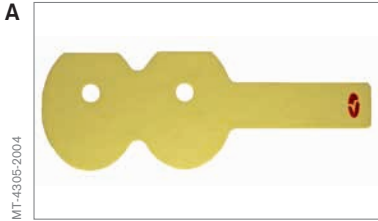

INFINITY® M300

NELLCOR®

NELLCOR® CABLES*

A SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 1.2 m **MS18683**

B SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 3 m **MS17330**

NELLCOR® REUSABLE SENSOR*

C SpO₂ finger sensor Nellcor® DURA DS 100A **7262764**

NELLCOR® DISPOSABLE SENSORS*

D SpO₂ adhesive sensor Nellcor® Oximax® Max-A, disposable, adult, cable 46 cm, 24 pcs. **MX50065**

E SpO₂ adhesive sensor Nellcor® Oximax® Max-AL, disposable, adult, cable 91 cm, 24 pcs. **MX50071**

F SpO₂ adhesive sensor Nellcor® Oximax® Max-P, disposable, child, 24 pcs. **MX50066**

G SpO₂ adhesive sensor Nellcor® Oximax® Max-N, disposable, neonate/adult, 24 pcs. **MX50068**

*Tradeware



INFINITY® M300

3.3 POUCH

A Telemetry transmitter pouch, disposable, 100 pcs. MS22905



3.4 CODENSATION

B Holder for Infinity® M300 bedside charger, swiveling 360° MP00990



Not all articles are available worldwide.

VISTA 120

4. VISTA 120

4.1 ELECTROCARDIOGRAPHY (ECG)

Not all articles are available worldwide.

VISTA 120 ECG SINGLE-PIN

A	<u>ECG defib trunk cable, 3/5-lead single-pin, IEC/AHA, 1.5 m</u>	<u>2612020</u>
	<u>ECG defib trunk cable, 3/5-lead single-pin, IEC/AHA, 2.5 m</u>	<u>2612021</u>
B	<u>ECG ESU trunk cable, 3/5-lead single-pin, IEC/AHA, 1.5 m</u>	<u>2612024</u>
	<u>ECG ESU trunk cable, 3/5-lead single-pin, IEC/AHA, 2.5 m</u>	<u>2612025</u>
C	<u>ECG defib trunk cable, Neonate, single-pin, IEC/AHA, 1.5 m</u>	<u>2612026</u>
	<u>ECG defib trunk cable, Neonate, single-pin, IEC/AHA, 2.5 m</u>	<u>2612027</u>
D	<u>ECG ESU trunk cable, Neonate, single-pin, IEC/AHA, 1.5 m</u>	<u>2612028</u>
	<u>ECG ESU trunk cable, Neonate, single-pin, IEC/AHA, 2.5 m</u>	<u>2612029</u>



VISTA 120

A	<u>ECG leadwire, 3-lead single-pin, IEC, 1 m</u>	<u>2612000</u>
	<u>ECG leadwire, 3-lead single-pin, AHA, 1 m</u>	<u>2612001</u>
	<u>ECG leadwire, 5-lead single-pin, IEC, 1.5 m</u>	<u>2612002</u>
B	<u>ECG leadwire, 5-lead single-pin, AHA, 1.5 m</u>	<u>2612003</u>



VISTA 120

INFINITY ECG CABLE

The following lead cables from the Infinity ECG portfolio can be used together with the Vista 120 single-pin connection cables.



A	ECG cable, 3-lead, single-pin connector, IEC1 (European color code), 1 m	MP03401
	ECG cable, 3-lead, single-pin connector, IEC2 (AHA/US color code), 1 m	MP03402
B	ECG cable, 5-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03403
	ECG cable, 5-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03404
C	ECG cable, 3-lead, single-pin connector, single-patient use, IEC1 (European color code), 1 m, 20 pcs.	MP00875
	ECG cable, 3-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1 m, 20 pcs.	MP00877
D	ECG cable, 5-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs.	MP00879
	ECG cable, 5-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs.	MP00881



VISTA 120

The following lead cables from the Infinity ECG portfolio can be used together with the Vista 120 single-pin connection cables.

A	ECG MonoLead[®], 3-lead, EUR, 2 m	MP05120
B	ECG MonoLead[®], 3-lead, AHA/US, 2 m	MP05121
C	ECG MonoLead[®], 5-lead, EUR, 2.5 m	MP05122
D	ECG MonoLead[®], 5-lead, AHA/US, 2.5 m	MP05123
E	ECG MonoLead[®], 3-lead, EUR, 4.1 m	MP05130
	ECG MonoLead[®], 3-lead, AHA/US, 4.1 m	MP05131
F	ECG MonoLead[®], 5-lead, EUR, 5.3 m	MP05132
	ECG MonoLead[®], 5-lead, AHA/US, 5.3 m	MP05133

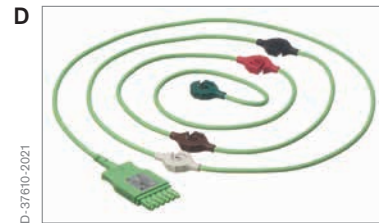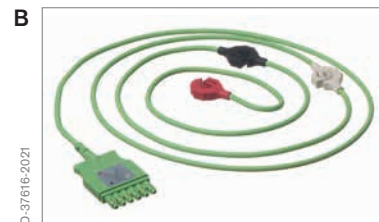

VISTA 120

VISTA 120 DUAL-PIN ECG CABLE

A	ECG trunk cable, 5-lead, IEC/AHA colour code	2606496
	ECG trunk cable, 3-lead, IEC/AHA colour code	2606495
B	ECG trunk cable, 3-lead, for neonates, IEC & AHA	2606492
C	ECG limb wires, 5-lead, Grabber, IEC colour code	2606494
	ECG limb wires, 3-lead, Grabber, IEC colour code	2606493
D	ECG limb wires, 5-lead, Grabber, AHA colour code	2606498
	ECG limb wires, 3-lead, Grabber, AHA colour code	2606497
	ECG electrodes*, neonate, single-patient use, 300 pcs.	5195024

*Tradeware



VISTA 120

4.2 PULSE OXIMETRY (SPO₂)

DRÄGER SPO₂

Not for use with Vista 120 Nellcor models.

A	SpO ₂ finger sensor, adult, 2.5 m, reusable	2606483
B	SpO ₂ finger sensor, adult, 1 m, reusable	2606484
C	SpO ₂ extension cable, 2 m, reusable	2606487
D	SpO ₂ silicone Soft-tip sensor, adult, 1 m, reusable	2606485
	SpO ₂ silicone Soft-tip sensor, pediatric, 1 m, reusable	2606486
E	SpO ₂ sensor*, adult, 0.5 m, disposable, 25 pcs.	2606208
F	SpO ₂ sensor*, pediatric, 0.5 m, disposable, 25 pcs.	2606209
G	SpO ₂ sensor*, infant, 0.5 m, disposable, 25 pcs.	2606210
	SpO ₂ sensor*, neonate, 0.5 m, disposable, 25 pcs.	2606211

Can be used with Vista 120 SC Model A and Model B.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120

NELLCOR OXIMAX SPO₂*

Only for use with Vista 120 Nellcor models.

A	Nellcor® intermediate cable, 1.2 m	MP00748
	Nellcor® intermediate cable, 3 m	MS20979
B	SpO₂ finger sensor Nellcor® DURA DS 100A	7262764
C	SpO₂ multisite sensor Nellcor® DURA-Y™, reusable	MX01004
D	Ear clip for Nellcor® DURA-Y™	MX01005
E	SpO₂ adhesive sensor Nellcor® Oximax® Max-Fast, disposable, application on forehead, 24 pcs.	MX50070

Can be used with Vista 120 SC Model C.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120

4.3 NON-INVASIVE BLOOD PRESSURE (NIBP)*

A	NIBP tube, 3 m	2601195
B	NIBP cuff, infant, 10 cm – 15 cm, reusable, Vista 120	2606151
C	NIBP cuff, small child, 13 cm – 17 cm, reusable, Vista 120	2606152
D	NIBP cuff, child, 16 cm - 21 cm, reusable, Vista 120	2606153
E	NIBP cuff, small adult, 20.5 cm – 28 cm, reusable, Vista 120	2606154
F	NIBP cuff, adult, 27 cm – 35 cm, reusable, Vista 120	2606155
G	NIBP cuff, larger adult, 34 cm – 43 cm, reusable, Vista 120	2606156

Can be used with Vista 120 SC Model A and Model B.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120

A	NIBP cuff, neonatal #1, 3 – 6 cm, disposable, 20 pcs., Vista 120	2607034
B	NIBP cuff, neonatal #2, 4 – 8 cm, disposable, 20 pcs., Vista 120	2607035
C	NIBP cuff, neonatal #3, 6 – 11 cm, disposable, 20 pcs., Vista 120	2607036
D	NIBP cuff, neonatal #4, 7 – 13 cm, disposable, 20 pcs., Vista 120	2607037
E	NIBP cuff, neonatal #5, 8 – 15 cm, disposable, 20 pcs., Vista 120	2607038

VISTA 120
VISTA 120 SC



VISTA 120

4.4 TEMPERATURE

	Temperature probe skin*, adult, 3 m, reusable	2601197
	Temperature probe*, rectal/oral, adult, 3 m, reusable	2601198
A	Temperature probe skin*, neonate/infant, 3 m, reusable	2601200
B	Temperature probe*, rectal/oral, neonatal/Infant, 3 m, reusable	2601199
C	Protective cover for general-purpose temperature probe, contains latex, 10 pcs.	7014616
	Protective cover* for general-purpose temperature probe, latexfree, 100 pcs.	MS40243

*Tradeware



VISTA 120

The following temperature probes from the Infinity® temperature portfolio can be used from software 3.11 with the Vista 120 in combination with the MP00943.

SINGLE-PATIENT USE TEMPERATURE PROBES

A	General-purpose temperature probe, single-patient use, adult, 7-pin connector, 1.6 m	MP00991
	General-purpose temperature probe, single-patient use, adult, 7-pin connector, 3 m	MP00992
B	General-purpose temperature probe, single-patient use, child, 7-pin connector, 1.6 m	MP00993
	General-purpose temperature probe, single-patient use, child, 7-pin connector, 3 m	MP00994
C	Skin temperature probe, single-patient use, adult, 7-pin connector, 1.6 m	MP00995
	Skin temperature probe, single-patient use, adult, 7-pin connector, 3 m	MP00996
D	Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 1.6 m	MP00997
	Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 3 m	MP00998



VISTA 120

REUSABLE TEMPERATURE PROBES

A	General-purpose temperature probe, reusable, adult, 7-pin connector, 1.5 m	4329889
A	General-purpose temperature probe, reusable, adult, 7-pin connector, 3 m	5204644
A	General-purpose temperature probe, reusable, child, 7-pin connector, 1.5 m	4329848
A	General-purpose temperature probe, reusable, child, 7-pin connector, 3 m	5204651
B	Skin temperature probe, reusable, with 7-pin connector, 1.5 m	4329822
B	Skin temperature probe, reusable, with 7-pin connector, 3 m	5204669
C	Protective cover for general-purpose temperature probe, contains latex, 10 pcs.	7014616
	Protective cover* for general-purpose temperature probe, latexfree, 100 pcs.	MS40243

*Tradeware



VISTA 120

4.5 INVASIVE BLOOD PRESSURE (IBP) AND CARDIAC OUTPUT* (CO)

*For supported Vista 120 models only.

A	IBP pressure transducer kit** , BD/Argon, disposable, 5 kits	2606225
B	IBP cable** Becton Dickinson/Argon	2606488
C	IBP cable** Edwards	2606489
D	IBP cable** Abbott, Medex	2606490
E	IBP cable** Utah	2606491
F	IBP cable** 7pin	2607558
G	IBP cable** 10pin	2607559

**Tradeware



VISTA 120

A	Cardiac Output cable*, 3 m	2601201
B	CO In-line Injection temperature probe* Becton Dickinson	2601190
	CO In-line Injection temperature probe housing*	2601191
C	CO Control Syringe*	2601192

*Tradeware






VISTA 120

4.6 CO₂ SIDESTREAM

DRÄGER G2 SIDESTREAM*

A	CO₂ sampling line, disposable, for Vista 120, 100 pcs.	2607550
B	CO₂ nasal cannula, adult, disposable, for Vista 120, 100 pcs.	2607551
C	CO₂ nasal cannula, pediatric, disposable, for Vista 120, 100 pcs.	2607552
D	CO₂ nasal cannula, infant, disposable, for Vista 120, 100 pcs.	2607553
E	CO₂ nasal cannula + O₂, adult, disposable, for Vista 120, 50 pcs.	2607554
F	CO₂ nasal cannula + O₂, child, disposable, for Vista 120, 50 pcs.	2607555
G	CO₂ sidestream water trap, for Vista 120, 20 pcs.	2605124

*Tradeware

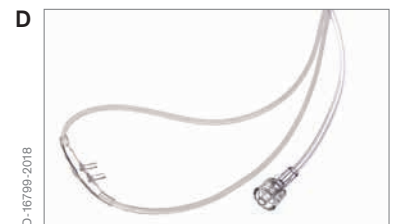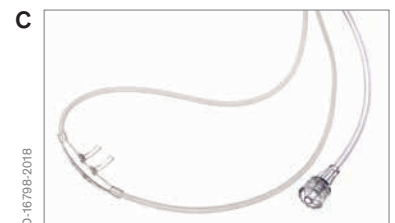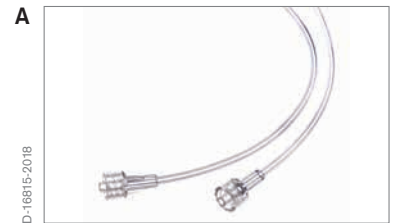

VISTA 120

4.7 CO₂ MAINSTREAM

DRÄGER CO₂ MAINSTREAM

Compatible with Vista 120 from software 3.11.

<p>A MCable Mainstream CO₂, for Evita V-family and Babylog VN-family, Infinity® M540, Vista120, Vista120 S, and Delta VF9.1 and higher</p>	<p>6871950</p>
<p>etCO₂ adapter cable M11.1</p>	<p>2607032</p>
<p>B CO₂ cuvette disposable, adult, 10 pcs.</p>	<p>MP01062</p>
<p>C CO₂ cuvette disposable, pediatric, 10 pcs.</p>	<p>MP01063</p>
<p>D CO₂ cuvette reusable, adult</p>	<p>6870279</p>
<p>E CO₂ cuvette reusable, pediatric</p>	<p>6870280</p>




VISTA 120

RESPIRONICS CO₂ MAINSTREAM*

A	etCO₂ sensor Capnostat 5, reusable	2606226
B	CO₂ airway adapter, adult, disposable, 10 pcs.	2606227
C	CO₂ airway adapter, neonatal (infant/pediatric), disposable, 10 pcs.	2606228

*Tradeware

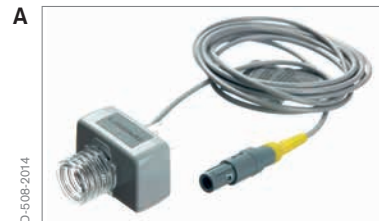

VISTA 120

	Vista 120 rechargeable Lithium-Ion battery*	2606229
A	Vista 120 printing paper*, 4 rolls	2606231
	BISx adapter cable Vista 120	2607560
	BISx pod Vista 120	2607561

*Tradeware



VISTA 120 SPOTCHECK

5. VISTA 120 SPOTCHECK
5.1 PULSE OXIMETRY (SPO₂)

DRÄGER SPO₂

A	SpO₂ extension cable, 2 m, reusable	2606487
B	SpO₂ silicone Soft-tip sensor, adult, 1 m, reusable	2606485
	SpO₂ silicone Soft-tip sensor, pediatric, 1 m, reusable	2606486
C	SpO₂ sensor*, adult, 0.5 m, disposable, 25 pcs.	2606208
D	SpO₂ sensor*, pediatric, 0.5 m, disposable, 25 pcs.	2606209
E	SpO₂ sensor*, infant, 0.5 m, disposable, 25 pcs.	2606210
F	SpO₂ sensor*, neonate, 0.5 m, disposable, 25 pcs.	2606211
G	SpO₂ finger sensor, adult, 2.5 m, reusable	2606483
	SpO₂ finger sensor, adult, 1 m, reusable	2606484

Can be used with Vista 120 SC Model A and Model B.

*Tradeware

VISTA 120
 VISTA 120 SC



VISTA 120 SPOTCHECK

NELLCOR® SPO₂*

A	Nellcor® intermediate cable, 1.2 m	MP00748
	Nellcor® intermediate cable, 3 m	MS20979
B	SpO₂ finger sensor Nellcor® DURA DS 100A	7262764
C	SpO₂ multisite sensor Nellcor® DURA-Y™, reusable	MX01004
D	Ear clip for Nellcor® DURA-Y™	MX01005
E	SpO₂ adhesive sensor Nellcor® Oximax® Max-Fast, disposable, application on forehead, 24 pcs.	MX50070

Can be used with Vista 120 SC Model C.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120 SPOTCHECK

MASIMO® LNCS CABLES*

A	<u>SpO₂ intermediate cable Masimo® LNCS, 1.2 m</u>	<u>MP00800</u>
A	<u>SpO₂ intermediate cable Masimo® LNCS, 3 m</u>	<u>MP00797</u>
A	<u>SpO₂ intermediate cable Masimo® LNCS, 4.3 m</u>	<u>MP03194</u>

Can be used with Vista 120 SC Model D and E.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120 SPOTCHECK

MASIMO® LNCS REUSABLE SENSORS*

A	SpO ₂ finger sensor Masimo® LNCS DCI, reusable, adult	MP00796
B	SpO ₂ finger sensor Masimo® LNCS DCIP, reusable, child	MP00795
	SpO ₂ Masimo® sensor LNCS-TC-I Ear, reusable, pediatric/adult	MP00788
	SpO ₂ Masimo® sensor LNCS-YI Multisite, reusable, neonate/pediatric/adult	MP00789
	SpO ₂ Masimo® sensor LNCS-TF-I, reusable, pediatric/adult	MP00799
C	SpO ₂ Masimo® sensor LNCS-DBI, reusable, adult	MP02993

Can be used with Vista 120 SC Model D and E.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120 SPOTCHECK

MASIMO® LNCS DISPOSABLE SENSORS*

A	SpO₂ adhesive sensor Masimo® LNCS Amtx, disposable, adult, 20 pcs.	MP00790
B	SpO₂ adhesive sensor Masimo® LNCS Inf, disposable, baby, 20 pcs.	MP00791
C	SpO₂ adhesive sensor Masimo® LNCS Neo, disposable, neonate, 20 pcs.	MP00792
D	SpO₂ adhesive sensor Masimo® LNCS Pmtx, disposable, child, 20 pcs.	MP00793
E	SpO₂ adhesive sensor Masimo® LNCS NeoPt, disposable, premature baby, 20 pcs.	MP00794

Can be used with Vista 120 SC Model D and E.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120 SPOTCHECK




MASIMO® LNCS TRAUMA SENSORS*

A Masimo® LNCS Trauma sensor, adult, 20 pcs.	MP03191
B Masimo® LNCS Trauma sensor, infant/pedi, 20 pcs.	MP03193
C Masimo® LNCS Trauma sensor, neonate, 20 pcs.	MP03192

Can be used with Vista 120 SC Model D and E.

*Tradeware

VISTA 120
VISTA 120 SC

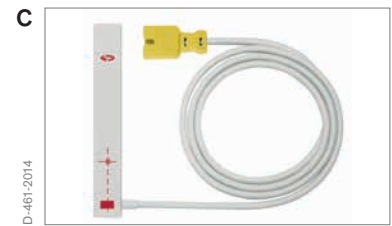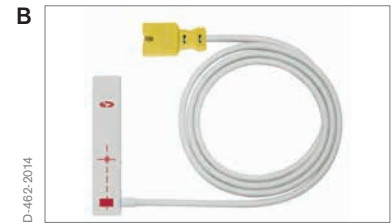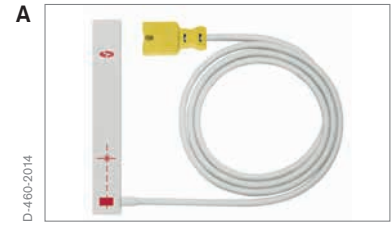

VISTA 120 SPOTCHECK

MASIMO® RD SET CABLES*

A	RD SET intermediate cable for Masimo® SET Mcable, 1.5 m	MS33736
A	RD SET intermediate cable for Masimo® SET Mcable, 3.6 m	MS33737

Can be used with Vista 120 SC Model D and E.

*Tradeware

VISTA 120
VISTA 120 SC



D-267-2017

VISTA 120 SPOTCHECK


MASIMO® RD SET DISPOSABLE SENSORS*

A	Masimo® RD SET sensor, adult, disposable	MS33726
B	Masimo® RD SET sensor, pediatric, disposable	MS33727
C	Masimo® RD SET sensor, infant, disposable	MS33728
D	Masimo® RD SET sensor, neonatal, disposable	MS33729
E	Masimo® RD SET sensor, neonatal preterm, disposable	MS33730

Can be used with Vista 120 SC Model D and E.

*Tradeware

VISTA 120
VISTA 120 SC

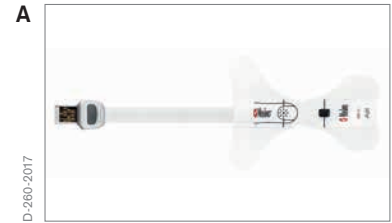

VISTA 120 SPOTCHECK

MASIMO® RD SET REUSABLE SENSORS*

A	Masimo® RD SET DCI sensor, adult, reusable	MS33740
B	Masimo® RD SET DC-IP sensor, pediatric, reusable	MS33741
C	Masimo® RD SET sensor TC-I Ear, pediatric/adult, reusable	MS40067
D	Masimo® RD SET sensor YI Multisite, neonate/pediatric/adult, reusable	MS40070
E	Masimo® RD SET sensor TF-I, pediatric/adult, reusable	MS40073

Can be used with Vista 120 SC Model D and E.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120 SPOTCHECK



MASIMO® RD SET TRAUMA DISPOSABLE SENSORS*

A	Masimo® Trauma sensor RD SET, adult, disposable	MS40088
B	Masimo® Trauma sensor RD SET, infant/pedi, disposable	MS40091
C	Masimo® Trauma sensor RD SET, neonate, disposable	MS40094

Can be used with Vista 120 SC Model D and E.

*Tradeware

VISTA 120
VISTA 120 SC

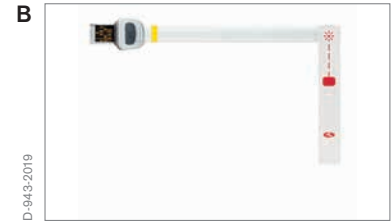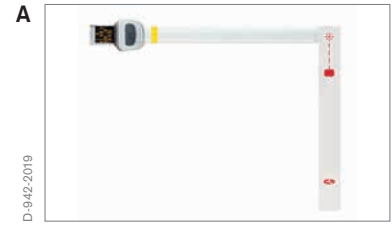

VISTA 120 SPOTCHECK

5.2 NON-INVASIVE BLOOD PRESSURE (NIBP)*

A	NIBP tube, 3 m	2601195
B	NIBP cuff, infant, 10 cm – 15 cm, reusable, Vista 120	2606151
C	NIBP cuff, small child, 13 cm – 17 cm, reusable, Vista 120	2606152
D	NIBP cuff, child, 16 cm - 21 cm, reusable, Vista 120	2606153
E	NIBP cuff, small adult, 20.5 cm – 28 cm, reusable, Vista 120	2606154
F	NIBP cuff, adult, 27 cm – 35 cm, reusable, Vista 120	2606155
G	NIBP cuff, larger adult, 34 cm – 43 cm, reusable, Vista 120	2606156

Can be used with Vista 120 SC Model A and Model B.

*Tradeware

VISTA 120
VISTA 120 SC



VISTA 120 SPOTCHECK

SUNTECH NIBP*

	APC NIBP Cuff, child, 12 cm – 19 cm, reusable, Vista 120 SC	<u>2600930</u>
A	APC NIBP Cuff, small adult, 17 cm – 25 cm, reusable, Vista 120 SC	<u>2600932</u>
	APC NIBP Cuff, adult, 23 cm – 33 cm, reusable, Vista 120 SC	<u>2600933</u>
B	APC NIBP Cuff, large adult, 31 cm – 40 cm, reusable, Vista 120	<u>2600934</u>
	OPC NIBP Cuff, child, 12 cm – 19 cm, reusable, 5 pcs., Vista 120 SC	<u>2600935</u>
	OPC NIBP Cuff, small adult, 17 cm – 25 cm, reusable, 5 pcs., Vista 120 SC	<u>2600936</u>
	OPC NIBP Cuff, adult, 23 cm – 33 cm, reusable, 5 pcs., Vista 120 SC	<u>2600937</u>
	OPC NIBP Cuff, large adult, 31 cm – 40 cm, reusable, 5 pcs., Vista 120 SC	<u>2600938</u>

Can be used with Vista 120 SpotCheck Model C, D and E.
APC = All Purpose Cuff, OPC= One Piece Cuff

*Tradeware



VISTA 120 SPOTCHECK

5.3. TEMPERATURE

EXERGEN THERMOMETER*

A	TAT-5000S Temporal Scanner, Vista 120 SC	2600946
B	TAT-5000S Probe Caps, single-use, 1,000 pcs., Vista 120 SC	2600948
	TAT-5000S Sheats, single-use, latexfree, 250 pcs., Vista 120 SC	2600949

Can be used with Vista 120 SpotCheck Model B, C, D and E.

*Tradeware



VISTA 120 SPOTCHECK

FILAC 3000*

A	Covidien Quick Temp Pod Oral	2600146
B	Filac 3000 Probe, oral, reusable, 1.2 m, Vista 120 SC	2600939
	Filac 3000 Probe, oral, reusable, 2.7 m, Vista 120 SC	2600940
C	Covidien Quick Temp Pod Rectal	2600147
D	Filac 3000 Probe, rectal, reusable, 1.2 m, Vista 120 SC	2600942
	Filac 3000 Probe, rectal, reusable, 2.7 m, Vista 120 SC	2600943
E	Filac 3000 Probe Covers, 20 pcs. Vista 120 SC	2600945

Can be used with Vista 120 SpotCheck Model B, C, D and E.

*Tradeware




VISTA 120 SPOTCHECK

5.4 CAPNOGRAPHY ETCO₂

A	etCO ₂ Microstr. Intubated Adult/Pediatric - Extended Duration, 25 pcs.	MP15001
	etCO ₂ Microstr. Intubated Adult/Pediatric - Short-term use, 25 pcs.	MP15002
B	etCO ₂ Microstr. Intubated Infant/Neonatal - Extended Duration, 25 pcs.	MP15003
C	etCO ₂ Microstr. Nasal Adult - Extended Duration, 25 pcs.	MP15004
	etCO ₂ Microstr. Nasal Adult - Short-term use, 25 pcs.	MP15005
D	etCO ₂ Microstr. Nasal Pediatric - Short-term use, 25 pcs.	MP15006
E	etCO ₂ Microstr. Nasal Infant/Neonatal - Extended Duration, 25 pcs.	MP15007
F	etCO ₂ Microstr. Nasal/Oral Adult - Short-term use, 25 pcs.	MP15008
G	etCO ₂ Microstr. Nasal/Oral Pediatric - Short-term use, 25 pcs.	MP15009


Can be used with Vista 120 SC Model E.

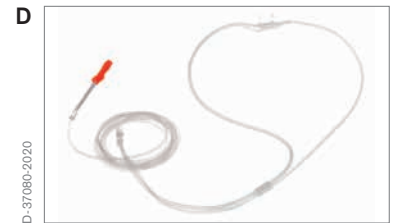

VISTA 120 SPOTCHECK

A CO ₂ Microstr. O ₂ Nasal Adult - Extended Duration, 25 pcs.	MP15010
CO ₂ Microstr. O ₂ Nasal Pediatric - Extended Duration, 25 pcs.	MP15011
CO ₂ Microstr. O ₂ Nasal/Oral Adult - Short-term use, 25 pcs.	MP15012
B CO ₂ Microstr. O ₂ Nasal/Oral Pediatric - Short-term use, 25 pcs.	MP15013

Can be used with Vista 120 SC Model E.

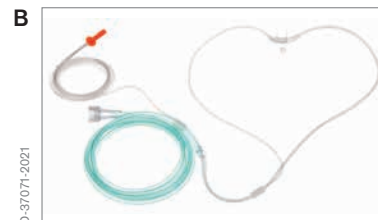

MULTI VENDOR ADAPTER

6. MULTI VENDOR ADAPTER

Multi Vendor Accessories (MVA)


6.1 ECG ADAPTER*

A	ECG Adapter, 3 Lead GE	MP04322
	ECG Adapter, 5 Lead GE	MP04320
B	ECG Adapter, 3 Lead Mindray	MP04312
	ECG Adapter, 5 Lead Mindray	MP04310
C	ECG Adapter, 3 Lead NK	MP04316
	ECG Adapter, 6 Lead NK	MP04315
D	ECG Adapter, 3 Lead Philips	MP04325
	ECG Adapter, 5 Lead Philips	MP04326
	ECG Adapter, 6 Lead Philips	MP04327
E	ECG Adapter, 3 Lead Spacelabs	MP04300
	ECG Adapter, 5 Lead Spacelabs	MP04301
F	ECG Adapter, 3 Lead Welch Allyn	MP04305
	ECG Adapter, 5 Lead Welch Allyn	MP04306

*Tradeware



MULTI VENDOR ADAPTER

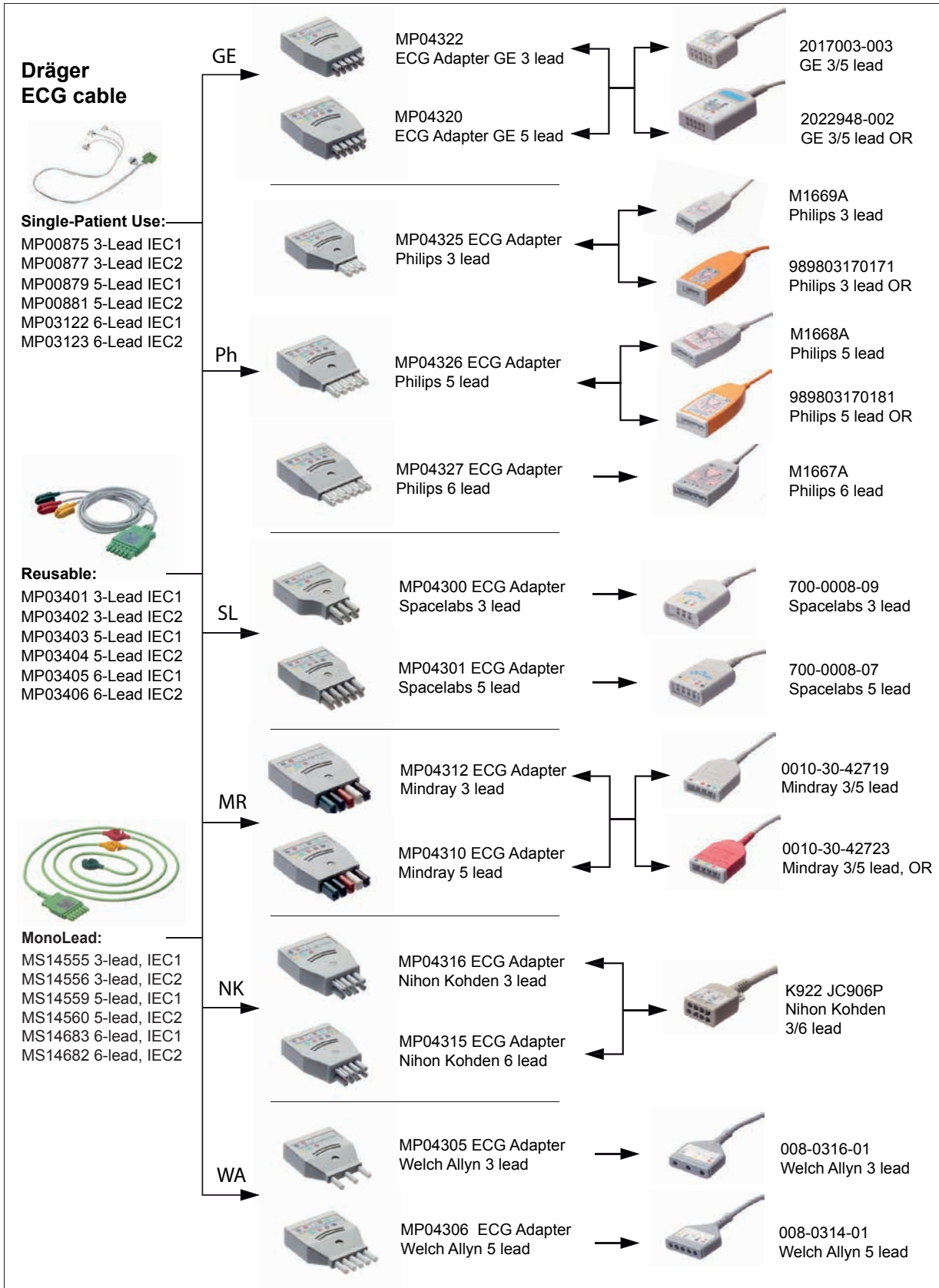

Not all articles are available worldwide.

MULTI VENDOR ADAPTER


6.2 ECG ADAPTER TELEMETRY*

A	<u>ECG adapter tel, GE-AP-CH, 3 Lead</u>	<u>MP04350</u>
	<u>ECG adapter tel, GE-AP-CH, 5 Lead</u>	<u>MP04351</u>
	<u>ECG adapter tel, GE-AP-CH, 6 Lead</u>	<u>MP04352</u>
B	<u>ECG adapter tel, GE-AP-FH, 5 Lead</u>	<u>MP04353</u>
C	<u>ECG adapter tel, PH MX40, 3 Lead</u>	<u>MP04355</u>
	<u>ECG adapter tel, PH MX40, 5 Lead</u>	<u>MP04356</u>
	<u>ECG adapter tel, PH MX40, 6 Lead</u>	<u>MP04357</u>
D	<u>ECG adapter tel, PH M, 5 Lead</u>	<u>MP04358</u>
	<u>ECG adapter tel, PH M, 6 Lead</u>	<u>MP04359</u>
E	<u>ECG adapter tel, SL AT, 5 Lead</u>	<u>MP04365</u>
F	<u>ECG adapter tel, SL UV, 5 Lead</u>	<u>MP04366</u>

6.3 TEMP ADAPTER

	<u>TEMP Adapter Jack</u>	<u>MP00943</u>
G	<u>TEMP Adapter 2-pin</u>	<u>MP00944</u>

*Tradeware

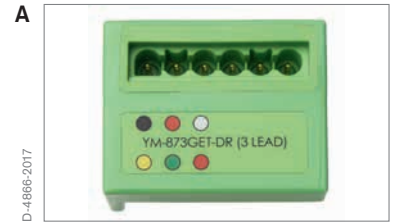

MULTI VENDOR ADAPTER


Dräger temperature sensors

Manufacturer	Model	Connector Type
GE	MP00943	GE 2016998-001, 2021700-001, 2021701-001
Philips	MP00943	Philips 989803163281
Spacelabs	MP00943	Spacelabs 700-0031-00
Welch Allyn	MP00943	Welch Allyn 008-0235-00
Nihon Kohden	MP00943	Nihon Kohden K961 JT 900P
Mindray	MP00944	Mindray 040-001235-00

Not all articles are available worldwide.

Dräger

Oxett

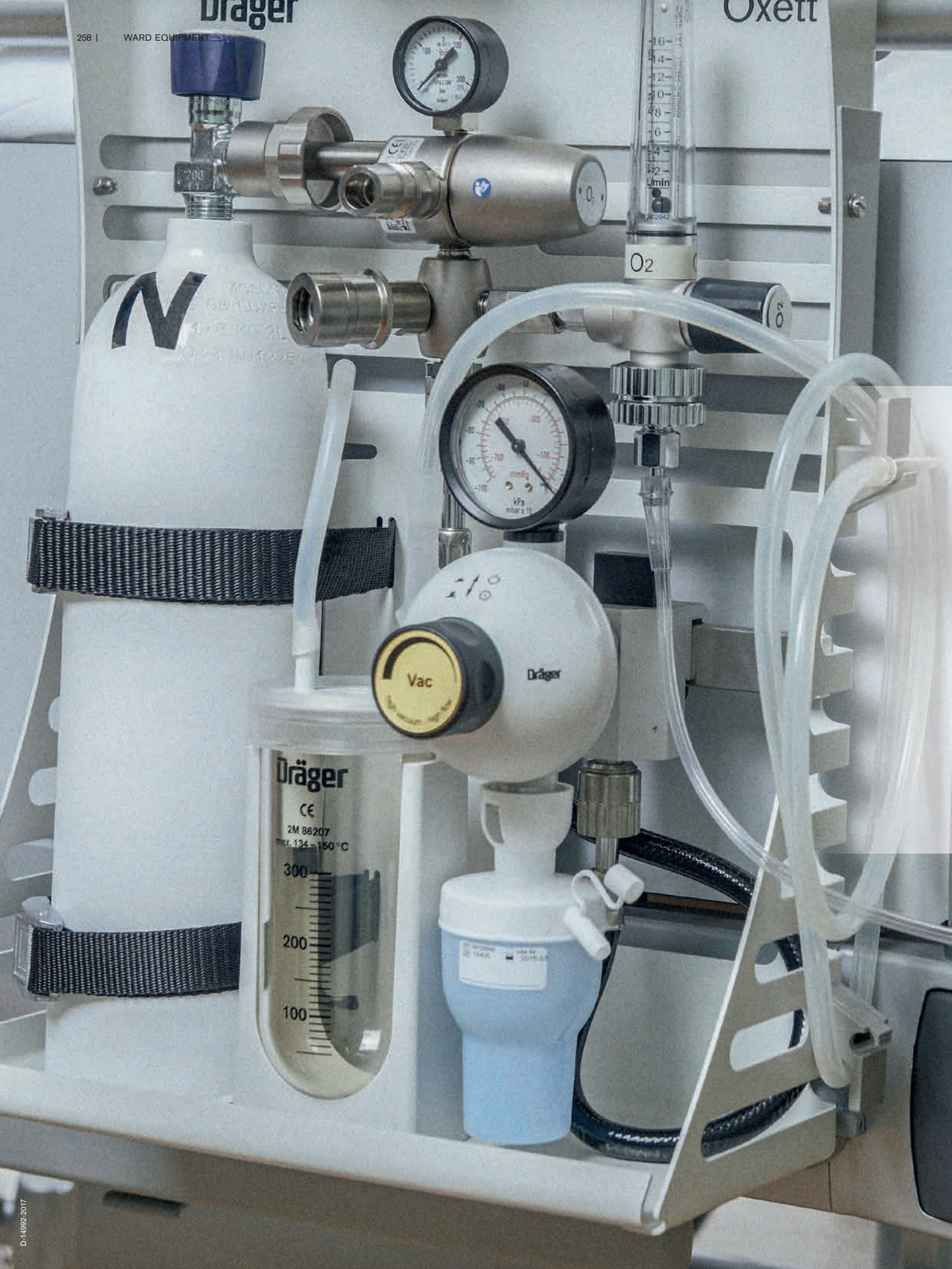

Dräger

CE

2M 86207
max. 134 - 150 °C

300

200

100

Vac

Dräger

O₂

O₂

Ward equipment

Adapting to the needs of your patients

Removing complexities and optimising your hospital's clinical workflows are key elements to ensure both patient and healthcare staff safety. Intelligent, user-friendly ward equipment does just that – promotes ergonomic workstations which reduce stress related to heavy lifting or uncomfortable access between equipment and the patient, thus improving staff satisfaction.

Our tailor-made workstations integrate optimised solutions so as to create a hygienic, healing environment that can help you deliver life-saving therapies to your patients.

CONTENT WARD EQUIPMENT

1. Oxygen therapy	261
1.1 Flowmeter	261
1.2 Oxymixer	271
1.3 Pressure reducer	273
1.4 Humidifier, adapter and application	281
2. Medical aspiration	283
2.1 Bronchial aspiration	283
2.2 OR aspiration	303
2.3 Continuous and thoracic drainage	307
3. Accessory for gas supply	319
3.1 CS hoses	319
3.1.1 German standard	319
3.1.2 French standard	321
3.1.3 British standard	323
3.1.4 Sweden standard	325
3.1.5 Swiss standard	326
3.1.6 Italian standard	327
3.1.7 Australian standard	328
3.1.8 Japanese standard	329
3.1.9 US american standard	330
3.1.10 Canadian standard	332
3.1.11 Chinese standard	333
3.2 Adapter	334

FLOWMETER

1. OXYGEN THERAPY

1.1 FLOWMETER

These O₂ Flowmeters are plugged directly into a wall terminal unit for oxygen and are used to control the flow of oxygen in inhalation and insufflation. Flow rates can be set with a high degree of accuracy. The large rotary knob makes setting easier.

Accuracy of flow rate at 4.5 bar ± 0.5 bar (1,013 hPa/20 °C)**
 ± 15% (corresponding to smallest setting value)
 ± 10% (all other setting values)

** Using the flowmeters outside the specified ambient conditions could affect flow accuracy.

Dimensions (H × W × D): 160 × 40 × 130 mm
 (6.3 × 1.6 × 5.1 in)
 Material: Housing: Aluminium
 Pressure dome: Polyamide
 Flow indicator: Glass
 Weight: Approx. 500 g (1.1 lb)

- | | |
|---|----------------|
| <p>A O₂ flowmeter, 4 l/min, CS, DIN, ISO and neutral color, M34×1.5 and 9/16" for accessories and sterile water systems</p> | MP04551 |
| <p>B O₂ flowmeter, 16 l/min, CS, DIN, ISO and neutral color, M34×1.5 and 9/16" for accessories and sterile water systems</p> | MP04558 |
| <p>C O₂ flowmeter, 16 l/min, CS, long CS plug (6 cm; 2.4 in), DIN, ISO and neutral color, M34×1.5 and 9/16" for accessories and sterile water systems</p> | MP04561 |



FLOWMETER

A <u>O₂ flowmeter, 16 l/min, CS, British Standard, ISO, M34×1.5, adapter 9/16"</u>	<u>MP04565</u>
B <u>O₂ flowmeter, 16 l/min, CS, French Standard, ISO, M34×1.5, adapter 9/16"</u>	<u>MP04566</u>
C <u>O₂ flowmeter, 16 l/min, CS, Swedish Standard, M34x1.5 and Adapter 9/16"</u> [*]	<u>MP13040</u>
D <u>O₂ flowmeter, 10 l/min, CS, CARBA, including adapter to 9/16" for sterile water</u>	<u>MP04554</u>
E <u>O₂ flowmeter, 16 l/min, CS, CARBA, including adapter to 9/16" for sterile water</u>	<u>MP04563</u>
F <u>O₂ flowmeter, 32 l/min, CS, DIN, ISO and neutral color, M34×1.5 and 9/16" for accessories and sterile water systems</u>	<u>MP04571</u>

*Available as of end 2022



FLOWMETER

TWIN FLOWMETER

The particular advantage of the twin flowmeter is, that two patients can be treated simultaneously and independently of one another, while occupying only one wall terminal.

Accuracy of flow rate at 4.5 bar \pm 0.5 bar (1,013 hPa/20 °C)**
 \pm 15% (corresponding to smallest setting value)
 \pm 10% (all other setting values)

** Using the flowmeters outside the specified ambient conditions could affect flow accuracy.

Dimensions (H \times W \times D): 160 \times 150 \times 190 mm
 (6.3 \times 5.9 \times 7.15 in)

Material: Housing: Aluminium
 Pressure dome: Polyamide
 Flow indicator: Glass

Weight: Approx. 900 g (1.98 lb)

A O ₂ twin flowmeter, 2 \times 16 l/min, CS, DIN, ISO, neutral color, M34 \times 1.5 and 9/16" for accessories and sterile water systems	MP04560
B O ₂ twin flowmeter, 2 \times 10 l/min, CS, CARBA, 2 \times M34 \times 1.5, including adapter to 9/16" for sterile water	MP04555
C O ₂ twin flowmeter, 2 \times 10 l/min, rail fitting with NIST gas inlet, 2 \times M34 \times 1.5, adapter 9/16"	MP04557
O ₂ twin flowmeter, 2 \times 16 l/min, rail fitting with NIST gas inlet, 2 \times M34 \times 1.5, adapter 9/16"	MP04569



FLOWMETER

These O₂ flowmeters for rail attachment can be used to control the flow of oxygen in inhalation and insufflation. Flow rates up to 4 and 16 l/min can be set with a high degree of accuracy. The large rotay knob makes setting easier. Using the rail clamp, the flowmeter can be quickly and safely attached to all rail profiles. The connection to the gas supply can be ensured by a CS hose with NIST coupling nut.

Accuracy of flow rate at 4.5 bar \pm 0.5 bar (1,013 hPa/20 °C)**
 \pm 15% (corresponding to smallest setting value)
 \pm 10% (all other setting values)

** Using the flowmeters outside the specified ambient conditions could affect flow accuracy.

Dimensions (H \times W \times D): 160 \times 40 \times 130 mm
 (6.3 \times 1.6 \times 5.1 in)

Material: Housing: Aluminium
 Pressure dome: Polyamide
 Flow indicator: Glass

Weight: Approx. 950 g (2.09 lbs)

A	O₂ flowmeter, 4 l/min, rail with NIST gas inlet, M34\times1.5 and 9/16" for accessories and sterile water systems	MP04552
B	O₂ flowmeter, 16 l/min, rail NIST, packed with hose 1.5 m, DIN, ISO, M34\times1.5 and 9/16" for accessories and sterile water systems	MP04573
C	O₂ flowmeter, 16 l/min, rail NIST, packed with hose 1.5 m, DIN, neutral color, M34\times1.5 and 9/16" for accessories and sterile water systems	MP04575
D	O₂ flowmeter, 16 l/min, rail, NIST gas inlet, M34\times1.5, adapter 9/16"	MP04567
E	O₂ flowmeter, 10 l/min, rail fitting with NIST gas inlet, M34\times1.5, adapter 9/16"	MP04556
F	O₂ flowmeter, 32 l/min, rail, NIST, M34\times1.5 and adapter 9/16"	MP04580



FLOWMETER

- A O₂ flowmeter, 4 l/min, rail, 90°, NIST gas inlet, M34×1.5 and 9/16" for accessories and sterile water systems MP04553
- B O₂ flowmeter, 16 l/min, rail, 90°, NIST gas inlet, M34×1.5 and 9/16" for accessories and sterile water systems MP04570


CS hoses

Several hoses to connect flowmeter with the central supply.
For CS hoses please refer to chapter "Accessory for gas supply".



FLOWMETER

<p>A O₂ flowmeter for anesthesia, 16 l/min, rail, with NIST ISO hose 1.2 m, connector 9/16" UNF to 6 mm</p>	<p>MP04576</p>
<p>B O₂ flowmeter for anesthesia, 16 l/min, rail, with NIST black hose 1 m, connector 9/16-18 UNF to 6 mm</p>	<p>MP04577</p>
<p>C O₂ flowmeter for anesthesia, 16 l/min, 90°, rail, with NIST ISO hose 1.2 m, connector 9/16-18 UNF to 6 mm</p>	<p>MP04578</p>
<p>D O₂ flowmeter for anesthesia, 16 l/min, 90°, rail, with NIST black hose 1 m, connector 9/16-18 UNF to 6 mm</p>	<p>MP04579</p>

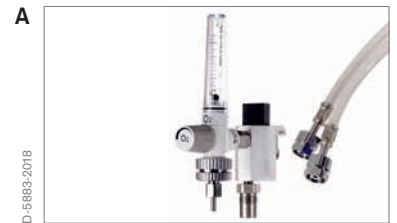

FLOWMETER

These flowmeters are plugged directly into a wall terminal unit for medical air and are used to control the flow of medical air for aerosoltherapy. Flow rates of up to 16 or 32 l/min can be set with a high degree of accuracy. The large rotary knob makes setting easier. With integrated sintered filter.

Accuracy of flow rate at 4.5 bar \pm 0.5 bar (1,013 hPa/20 °C)**
 \pm 15% (corresponding to smallest setting value)
 \pm 10% (all other setting values)

** Using the flowmeters outside the specified ambient conditions could affect flow accuracy.

Dimensions (H \times W \times D): 160 \times 40 \times 130 mm
 (6.3 \times 1.6 \times 5.1 in)

Material: Housing: Aluminium
 Pressure dome: Polyamide
 Flow indicator: Glass

Weight: Approx. 500 g (1.1 lb)

A	<u>AIR flowmeter, 16 l/min, CS, DIN, ISO and neutral color, M32x1.5</u>	<u>MP04559</u>
B	<u>AIR flowmeter, 16 l/min, CS, long cs plug (6 cm; 2.4 in), DIN, ISO and neutral color, M32x1.5</u>	<u>MP04562</u>
C	<u>AIR flowmeter, 16 l/min, CS, CARBA, ISO, M32x1.5</u>	<u>MP04564</u>
D	<u>AIR flowmeter, 16 l/min, CS, Swedish Standard, M32x1.5*</u>	<u>MP13041</u>
E	<u>AIR flowmeter, 32 l/min, CS, DIN, ISO and neutral color, M32x1.5</u>	<u>MP04572</u>

*Available as of end 2022



FLOWMETER

These Flowmeters for rail attachment can be used to control the flow of medical air for aerosol therapy. Flow rates of up to 16 l/min can be set with a high degree of accuracy. The large setting knob makes setting easier. Using the rail clamp, the flowmeter can be quickly and safely attached to standard rail profiles. The connection to the gas supply can be ensured by a CS hose with NIST coupling nut.

Accuracy of flow rate at 4.5 bar ± 0.5 bar (1,013 hPa/20 °C)**
 ± 15% (corresponding to smallest setting value)
 ± 10% (all other setting values)

** Using the flowmeters outside the specified ambient conditions could affect flow accuracy.

Dimensions (H × W × D): 160 × 40 × 130 mm
 (6.3 × 1.6 × 5.1 in)

Material: Housing: Aluminium
 Pressure dome: Polyamide
 Flow indicator: Glass

Weight: Approx. 950 g (2.09 lbs)

A AIR flowmeter, 16 l/min, rail, NIST, M32×1.5 **MP04568**

B AIR flowmeter, 32 l/min, rail, NIST, M32×1.5 **MP04581**

CS hoses

Several hoses to connect flowmeter with the central supply.
 For CS hoses please refer to chapter "Accessory for gas supply".



FLOWMETER

COMPACT FLOWMETER CLICKFLOW

Input pressure:	4.5 ± 0.5 bar*
Adjustable ranges:	0.1 – 1 l/min (0; 0.1; 0.2; 0.3; 0.4; 0.5; 0.6; 0.7; 0.8; 1 l/min) 0.5 – 5 l/min (0; 0.5; 1; 1.5; 2; 2.5; 3; 3.5; 4; 5 l/min) 1 – 16 l/min (0; 1; 2; 3; 4; 6; 8; 10; 12; 16 l/min) 1 – 32 l/min (0; 1; 5; 8; 12; 16; 20; 25; 28; 32 l/min)
Accuracy (at 4.5 ± 0.5 bar):	0.5 – 5 l/min, 1 – 16 l/min and 1 – 32 l/min: ± 20%; 0 – 1 l/min ± 30%
Gas inlet:	CS plug O ₂ or AIR; NIST (rail)
Gas outlet:	6 mm nozzle for application unit, thread 9/16" for sterile water or adapter for reusable humidifier
Weight:	approx. 420 g (CS); approx. 700 g (rail)
Housing material:	nickel-coated brass
Dimensions:	diameter about 38 mm, length about 92 mm (CS version); about 118 mm (rail version)

*1 bar = 1 kPa × 100



- | | |
|--|---|
| A Compact flowmeter Clickflow configurator, gives you the option of arranging the compact flowmeter according to your requirements, selectable options for configuration: | MP13000 |
| Gas type: | O ₂ ; AIR |
| Mounting: | CS; rail with NIST gas inlet |
| CS connection: | DIN; Carbamed; BS (England); NF (France); AGA (Sweden) |
| Version: | single compact flowmeter;
double compact flowmeter |
| Flow: | 0.1 – 1 l/min; 0.5 – 5 l/min;
1 – 16 l/min, 1 – 32 l/min |
| Gas outlet: | 6 mm nozzle, thread 9/16" for sterile water |
| B O₂ compact flowmeter Clickflow, DIN CS, 1 – 16 l/min, 9/16" | MP13011 |
| C O₂ compact flowmeter Clickflow, DIN CS, 1 – 16 l/min, 6 mm nozzle | MP13012 |

FLOWMETER

A	<u>O₂ compact flowmeter Clickflow, rail, NIST, 1 – 16 l/min, 9/16"</u>	<u>MP13013</u>
B	<u>O₂ twin compact flowmeter Clickflow, DIN, 2 × 1 – 16 l/min, 9/16"</u>	<u>MP13014</u>
C	<u>AIR compact flowmeter Clickflow, DIN, 1 – 16 l/min, 9/16"</u>	<u>MP13015</u>
D	<u>AIR compact flowmeter Clickflow, DIN, 1 – 16 l/min, 6 mm nozzle</u>	<u>MP13016</u>
E	<u>AIR compact flowmeter Clickflow, DIN, 0.1 – 1 l/min, 6 mm nozzle</u>	<u>MP13019</u>
F	<u>AIR compact flowmeter Clickflow, rail, 1 – 16 l/min, 9/16"</u>	<u>MP13017</u>

Accessories for Flowmeter in chapter 1.4
(Humidifier, adapter and application).

Several hoses to connect Clickflow with the central supply.
For CS hoses please refer to chapter "Accessory for gas supply".



OXYMIXER

1.2 OXYMIXER*

Oxymixer enables a continuous and individually adjusted mixture of medical air and medical oxygen to be administered to infants, children and adults. The O₂ concentration is easily adjusted to a highly accurate flow – even with flow rates of less than 3 l/min.

Oxymixer configuration

Free choice of Low- or High Flow Blender with: 1 or 2 flowmeter, NIST or DISS gas inlet, Monitor, Bleed swith, Nozzle, Sound damper. The rail fixation is included.

A Oxymixer High Flow	MP04201
Including: Rail clamp, Flowmeter 0 – 32 l/min right, Flowmeter 0 – 85 l/min left.	
B Oxymixer High Flow with monitor	MP04202
Including: Rail clamp, Flowmeter 0 – 32 l/min right, Flowmeter 0 – 85 l/min left, Monitor, Oxygen sensor, Oxygen sensor cable, Sensor receptacle.	
C Oxymixer Low Flow	MP04203
Including: Rail clamp, Flowmeter 0 – 3 l/min right, Flowmeter 0 – 15 l/min left, Bleed on/off switch.	
D Oxymixer Low Flow with monitor	MP04204
Including: Rail clamp, Flowmeter 0 – 3 l/min right, Flowmeter 0 – 15 l/min left, Bleed on/off switch, Monitor, Oxygen sensor, Oxygen sensor cable, Sensor receptacle.	
E Oxymixer Low Flow with 1 – 15 l flowmeter right, rail clamp	MP04229
F Multi-P trolley* equipped	MP13020
Including base with 5 castors, 38 mm Ø pole, 3 rail clamps, basket.	
G Multi-P trolley*	MP13021
Including base with 5 castors and 38 mm Ø pole for accessories.	

Several hoses to connect Oxymixer with the central supply.
For CS hoses please refer to chapter "Accessory for gas supply".

*Tradeware



OXYMIXER*

A	Oxymixer sound damper 9/16"	MP04218
B	Oxymixer sound damper 22 and 15 mm	MP04219
C	Adapter 9/16" to 15/22mm	MP05305
D	Adapter, with 6 mm nozzle, e.g., for disposable breathing masks, 9/16" connection	MP00416
E	Oxymixer oxygen sensor (spare part)	MP04221

*Tradeware



PRESSURE REDUCER

1.3 PRESSURE REDUCER

A Pressure reducer OxyLine, customized configuration of pressure reducers

MP00100

OxyLine's modular design means you can accurately choose the pressure reducer you need from the following components:

- O₂ or AIR pressure reducer with 200 bar technology
- Manometer, vertical or horizontal
- Flow regulator settings (1, 2, 3, 4, 6, 8, 12 and 15 l/min) and 9/16" sterile water connector
- Flowmeter for up to 16 l/min with connector for Dräger reusable humidifier, sterile water connector and 6 mm (0.2 in) nozzle for an oxygen mask
- Self-closing NIST connector
- Centralized gas supply coupling, e.g. to connect up a mains-operated ventilation device
- Ejector for short-time suction of eg. secret

Type:	Single-stage pre-set pressure reducer
Permissible inlet pressure range:	11 to 200 kPa × 100
Nominal outlet pressure (p2):	non-adjustable, output pressure as per rating plate
Weight:	1.2 kg to 2 kg (depending on connector type)
Material:	Brass with chemical nickel surface



PRESSURE REDUCER

OxyLine® Pressure reducer


MP00100

Order information international


Make your own configuration by setting a cross or a number.

☐ = One position must be chosen

☐☐ = One position can be chosen

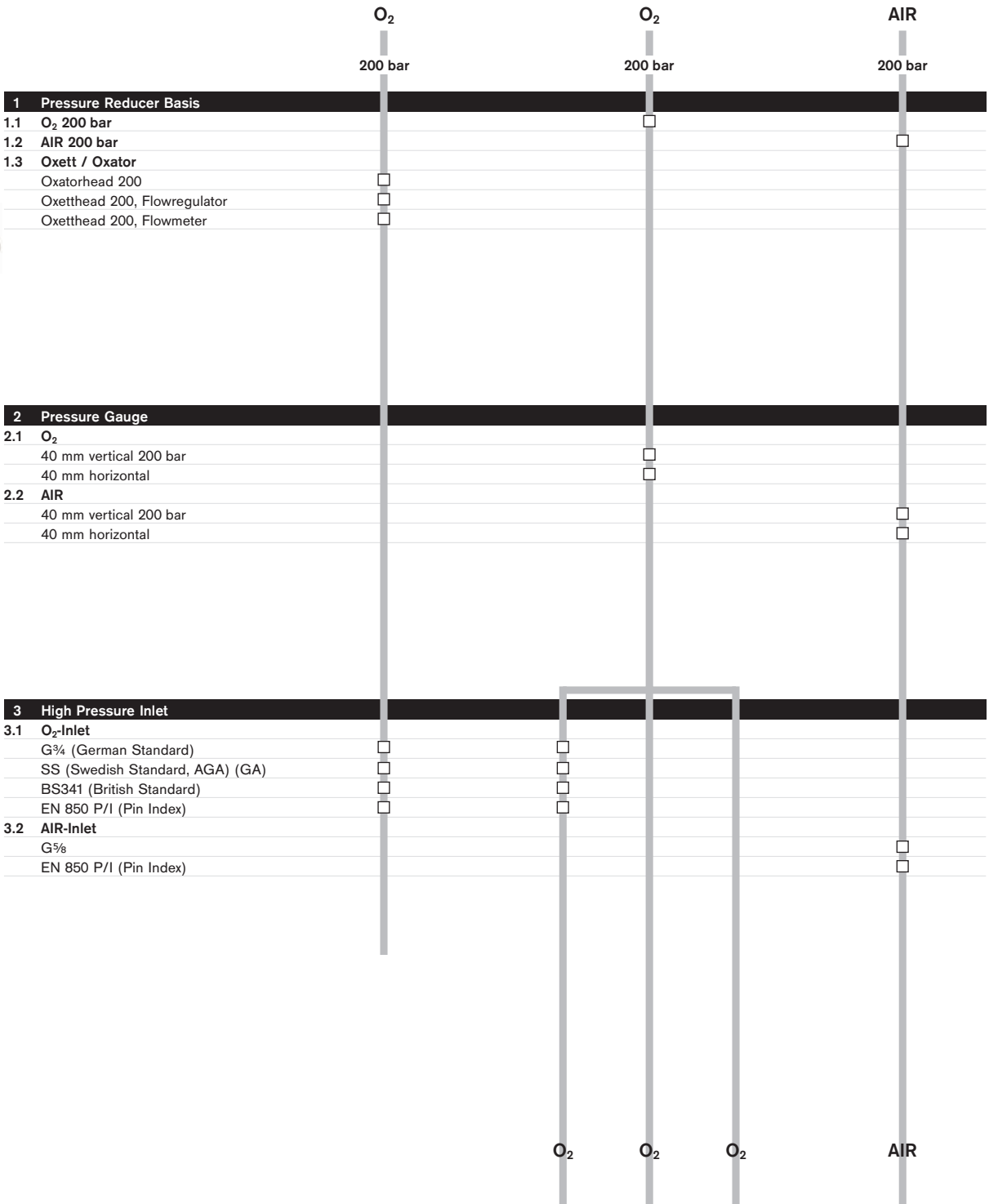

PRESSURE REDUCER

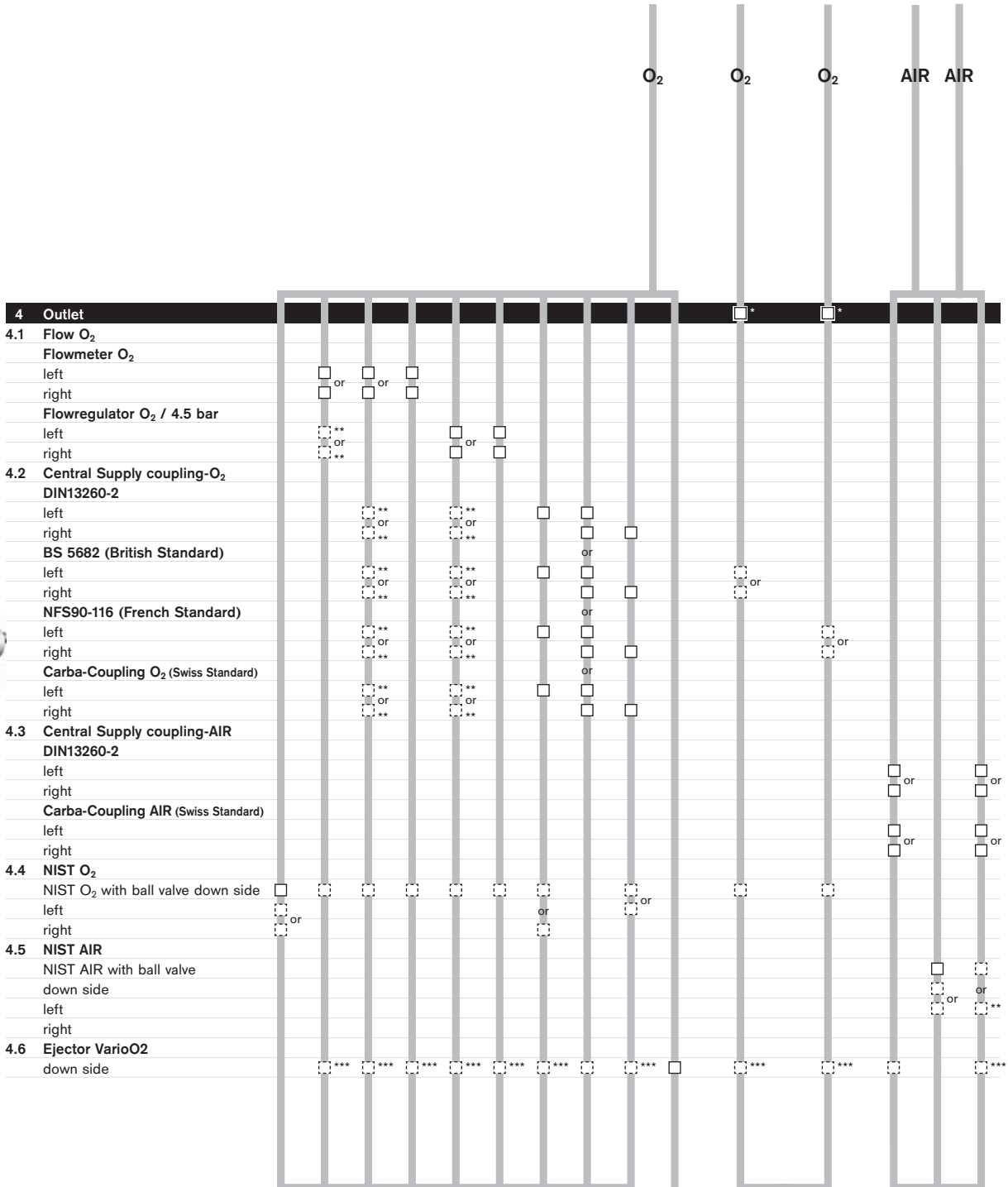

4 Outlet

4.1 Flow O₂

Flowmeter O₂
left
right

Flowregulator O₂ / 4.5 bar
left
right

4.2 Central Supply coupling-O₂

DIN13260-2
left
right

BS 5682 (British Standard)
left
right

NFS90-116 (French Standard)
left
right

Carba-Coupling O₂ (Swiss Standard)
left
right

4.3 Central Supply coupling-AIR
DIN13260-2
left
right

Carba-Coupling AIR (Swiss Standard)
left
right

4.4 NIST O₂

NIST O₂ with ball valve down side
left
right

4.5 NIST AIR

NIST AIR with ball valve
down side
left
right

4.6 Ejector VarioO₂

down side

* At least one outlet must be selected.
** Please check former selection (left/right).
*** Only possible if no NIST O₂ (4.4) or NIST AIR (4.5).

5 Accessories

Adapter NIST to DIN CS-Coupling O ₂	M35056	<input type="checkbox"/>
Adapter NIST to DIN CS-Coupling AIR	M35058	<input type="checkbox"/>
Adapter, with 6 mm nozzle, 9/16" connection	MP00416	<input type="checkbox"/>
Humidifier-bubbler, for oxygen humidification	2M85834	<input type="checkbox"/>
Humidifier-nebulizer, for oxygen humidification	2M85835	<input type="checkbox"/>

PRESSURE REDUCER

A Oxett head OxyLine, pressure reducer G $\frac{3}{4}$ " , 200 bar, NIST, flowmeter, DIN-CS, ejector VarioO2	MP01101
B Oxett head OxyLine, pressure reducer G $\frac{3}{4}$ " , 200 bar, NIST, flow regulator, DIN-CS, ejector VarioO2	MP01102
C Oxator head OxyLine, pressure reducer G $\frac{3}{4}$ " , 200 bar, NIST, 2 × DIN-CS, ejector VarioO2	MP01103
D Pressure reducer w/ flowmeter OxyLine, G $\frac{3}{4}$ " , 200 bar, 16 l/min, M34×1.5, 9/16", nozzle	MP01104
E Pressure reducer with NIST OxyLine, G $\frac{3}{4}$ " , for 200 bar, NIST	MP01105



PRESSURE REDUCER

Oxett

Oxett is a portable oxygen supply unit for mobile use. Oxett ensures consistent treatment, because the device can be mounted to bed, rail and stretcher. All devices for inhalation, insufflation and suction are available on a solid aluminum plate. Precise O₂ delivery is handled through the pressure compensated flow meter or the compact flowmeter. Manual or automatic short-term ventilation is even possible, because the Oxylog[®] can be connected via the O₂ standard coupling. The portable Oxett unit meets the requirements of daily hospital and emergency medical practice in a simple way. It removes all well-practiced temporary supply solutions.

A Oxett carrying frame, for Oxett head+O₂ cylinder up to 2.5 l, with secretion jar and rail for suction ejector	MP04110
Oxett head, configurable via OxyLine[®]	MP00100
B Oxett head OxyLine, pressure reducer G³/₄" , 200 bar, NIST, flowmeter, DIN-CS, ejector VarioO2	MP01101
B Oxett head OxyLine, pressure reducer G³/₄" , 200 bar, NIST, flow regulator, DIN-CS, ejector VarioO2	MP01102
C Oxett adapter plate, to mount Oxett on bed in higher position	2M86205



PRESSURE REDUCER

Small Oxator

Space-saving solution for transportation of large gas-cylinder.

The stable trolley can be used flexible.

Set up from components below:

A	Simple cylinder trolley, for one gas cylinders up to a volume of 11 L	MP00038
B	Oxator head, configurable via OxyLine® Oxator head OxyLine, pressure reducer G$\frac{3}{4}$" , 200 bar, NIST, 2 × DIN-CS, ejector VarioO2	MP00100 MP01103
C	Ejector VarioO2® B/OR Rail, for bronchial and OR suction, gas connection NIST O₂	MP00552
C	NIST-NIST Connecting hose 1.2 m	M36053



PRESSURE REDUCER

- A Compact rail,**
for 38 mm tube, standard rail profile 25 × 10 mm,
each side 10 cm, weight: 0.9 kg, max. load: 30 kg **2M85337**
-
- B Secretion jar set, with 700 ml secretion jars (2 pcs.)**
and 1.5 m suction hose,
maximum capacity up to the overflow safeguard 2 × 1000 ml,
scale in steps of 10 ml up to max. 700 ml **2M85056**
-
- C Rail clamp lever type,**
for attachment of devices to 25 × 10 mm standard rail,
for compact suction unit and secretion jar set **M26842**

For CS hoses please refer to chapter „Accessory for gas supply“.



PRESSURE REDUCER

A	Simple cylinder trolley, for one gas cylinders up to a volume of 11 L	<u>MP00038</u>
B	Rail holder, for 38 mm tube, standard rail profile 25 × 10 mm, 40 cm, length: 40 cm, weight: 0.9 kg, max. load: 5 kg	<u>2M86185</u>
C	Adapter NIST to DIN CS-Coupling O ₂	<u>M35056</u>
D	Adapter NIST to DIN CS-Coupling AIR	<u>M35058</u>



HUMIDIFIER, ADAPTER AND APPLICATION

1.4 HUMIDIFIER, ADAPTER AND APPLICATION

A Humidifier-bubbler, for oxygen humidification up to appr. 60% relative humidity and flows up to 5 l/min	2M85834
B Humidifier-nebulizer, for oxygen humidification up to 100% relative humidity and flows above 5 l/min	2M85835
C Adapter, to connect Dräger humidifier to mixing block, M15×1 to M34×1.5	2M85853
D Adapter 9/16" to M34×1.5 for Dräger reusable humidifier	MP03325



HUMIDIFIER, ADAPTER AND APPLICATION

A Adapter, for disposable sterilized water systems to O₂ flowmeters, M34×1.5 to 9/16"	<u>2M86680</u>
A Adapter, for disposable sterilized water systems to AIR flowmeter, M32×1.5 to 9/16"	<u>MP03366</u>
B Adapter, with 6 mm nozzle, e.g., for disposable breathing masks, 9/16" connection	<u>MP00416</u>

Oxygen and aerosol therapy masks, venturi mask kids, nasal cannulas and accessories in chapter 5.3, breathing accessories.




MEDICAL ASPIRATION

2.1 Bronchial aspiration

Dräger offers an extensive line of bronchial aspirators, while also providing individual solutions to meet your specific needs in the ICU and OR. Generally you have the choice between vacuum, medical air or oxygen as the drive gas. We differ by the driving units between

the ejectors and the vacuum regulators. In addition, we can supply the attachments and hose brackets you need to ensure an ergonomically efficient workplace. You could also choose the connection you require and the fitting jar system. We offer reusable solutions as well as disposables.

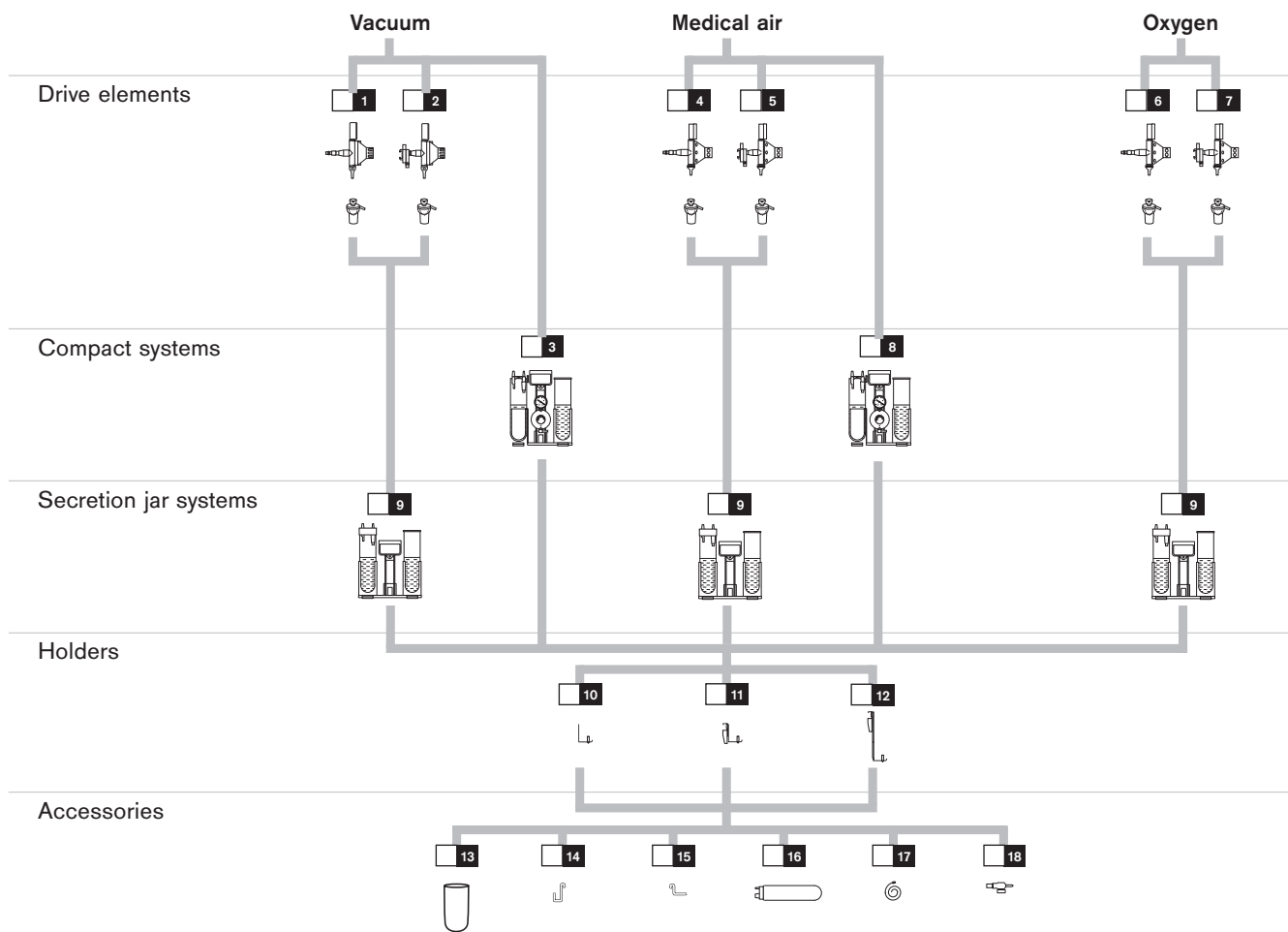

1	MP00501	VarioVac® B DIN	MP00532	VarioAir® B NF	10	M26859	Wall holder	
	MP00502	VarioVac® B BS	MP00533	VarioAir® B AGA	11	M26842	Rail clamp lever type	
	MP00503	VarioVac® B NF	MP00534	VarioAir® B UNI	12	2M85119	Ceiling holder	
	MP00504	VarioVac® B AGA	MP00548	VarioAir® B CARBA	13	2M86207	Secretion jar, 300 ml	
	MP00514	VarioVac® B CARBA	5	MP00535	VarioAir® B/OR Rail	2M85594	Secretion jar, 700 ml	
2	MP00505	VarioVac® B/OR Rail	MP00545	VarioAir® P Rail	14	2M85446	CS connection hose holder	
	MP00660	VarioVac® P Rail	MP03220	VarioAir® MRI B/OR Rail	15	2M85239	Suction hose holder	
	MP00516	Dual VarioVac® B/OR Rail	MP03221	VarioAir® MRI P Rail	16	MX23023	Secretion jar VacuSmart®, 700 ml	
	MP03222	VarioVac® MRI B/OR Rail	6	MP00550	VarioO2® B DIN	MX23025	Secretion jar VacuSmart® Gel, 700 ml	
	MP03223	VarioVac® MRI P Rail	7	MP00551	VarioO2® B CARBA	17	MX23026	Suction tube VacuSmart® Tube
3	MP00506	VarioVac® B compact	8	MP00552	VarioO2® B/OR Rail	18	MP00570	Shut-off valve for VarioVac®
	MP00661	VarioVac® P compact	8	MP00536	VarioAir® B compact			
4	MP00530	VarioAir® B DIN		MP00547	VarioAir® P compact			
	MP00531	VarioAir® B BS	9	2M85056	Secretion jar set			

Not all articles are available worldwide.

BRONCHIAL ASPIRATION

2. MEDICAL ASPIRATION

2.1 BRONCHIAL ASPIRATION

The vacuum regulator is connected directly to a central vacuum supply system or to rails. The vacuum is infinitely variable up to max. -90 kPa (-900 mbar) and is shown on a turnable, manually adjustable pressure gauge. Shut-off valve is included.

You can use the VarioSafe® disposable filter system to reliably protect the **VarioVac®** series and the central vacuum system against contamination. A VarioSafe® filter system is included with every VarioVac® vacuum regulator.

Technical data:

Maximum vacuum: -90 to 0 kPa (-900 to 0 mbar)
 Maximum suction flow: 55 l/min
 Weight: 600 to 800 g

800 g for rail version
 1,700 g for dual version
 (MP00516)

Dimensions (W/H/D) in mm: 80 × 210 × 160 (without VarioSafe)
 80 × 340 × 160 (with VarioSafe)
 MP00516 for dual version:
 200 × 210 × 210 (without VarioSafe)
 200 × 340 × 210 (with VarioSafe)

Maximum sound level: 43 dB(A)

Hose nozzle: Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention: 100% (with VarioSafe)*

Bacterial retention: 99.99997% (with VarioSafe)**

Viral retention: 99.99922% (with VarioSafe)**

NaCl retention: 99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA,
 Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc.,
 Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA,
 March 2021, Lübeck, Germany, according to EN ISO 23328-1

A	Vacuum regulator VarioVac® B DIN, for bronchial suction, gas connection VAC	MP00501
B	Vacuum regulator VarioVac® B BS, for bronchial suction, gas connection VAC	MP00502
C	Vacuum regulator VarioVac® B NF, for bronchial suction, gas connection VAC	MP00503
D	Vacuum regulator VarioVac® B AGA, for bronchial suction, gas connection VAC, Swedish standard	MP00504
E	Vacuum regulator VarioVac® B CARBA, for bronchial suction, gas connection VAC	MP00514

For CS hoses please refer to chapter „Accessory for gas supply“.



BRONCHIAL ASPIRATION

VarioVac® vacuum regulator for rail.

Technical data:

Maximum vacuum: -90 to 0 kPa (-900 to 0 mbar)
 Maximum suction flow: 55 l/min
 Weight: 800 g for rail version
 1,700 g for dual version (MP00516)

Dimensions (W/H/D) in mm: for rail version
 80 × 210 × 160 (without VarioSafe)
 80 × 340 × 160 (with VarioSafe)
 for dual version (MP00516):
 200 × 210 × 210 (without VarioSafe)
 200 × 340 × 210 (with VarioSafe)

Maximum sound level: 43 dB(A)
 Hose nozzle: Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention: 100% (with VarioSafe)*
 Bacterial retention: 99.99997% (with VarioSafe)**
 Viral retention: 99.99922% (with VarioSafe)**
 NaCl retention: 99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



A	Vacuum regulator VarioVac® B/OR Rail, for bronchial and OR suction, gas connection NIST VAC	MP00505
B	Vacuum regulator dual VarioVac® B/OR Rail, for bronchial and OR suction, gas connection NIST VAC	MP00516
C	Vacuum regulator VarioVac® B/OR Rail DISS, for bronchial and OR suction, gas connection DISS VAC	MP00631

For CS hoses please refer to chapter „Accessory for gas supply“.

BRONCHIAL ASPIRATION

VarioVac® vacuum regulator for rail.

The gas supply is ensured by VAC connection hose.

This rail version can be connected to all standard wall rails (10 mm width and 25 mm height).

You can use the VarioSafe® disposable filter system to reliably protect the VarioVac® series and the central vacuum system against contamination.

A VarioSafe® filter system is included with every VarioVac® vacuum regulator.

Technical data:	VarioVac P (MP00660, MP03223)
Maximum vacuum:	-55 to 0 kPa (-550 to 0 mbar)
Maximum suction flow:	20 l/min
Maximum sound level:	40 dB(A)

	VarioVac B/OP (MP03222)
Maximum vacuum:	-90 to 0 kPa (-900 to 0 mbar)
Maximum suction flow:	55 l/min
Maximum sound level:	43 dB(A)

Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 210 × 160 (without VarioSafe)
	80 × 340 × 160 (with VarioSafe)
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



A Vacuum regulator VarioVac® P Rail, for bronchial suction, gas connection NIST VAC, for connection to all standard wall rails MP00660

A Vacuum regulator VarioVac® P Rail DISS, for bronchial suction, gas connection DISS VAC, for connection to all standard wall rails MP00632

These bronchial aspirators are suited for use in MRI. The bronchial aspirators can be used in magnets with static magnetic fields of 1.5 tesla and 3 tesla. Minimum distance from the center of the magnetic field: Magnetic flux density of 40 millitesla (400 gauss).

B Vacuum regulator VarioVac® MRI P Rail, for bronchial suction, gas connection NIST VAC, for connection to all standard wall rails MP03223

C Vacuum regulator VarioVac® MRI B/OR Rail, for bronchial and OR suction, gas connection NIST VAC, for connection to all standard wall rails MP03222

For CS hoses please refer to chapter „Accessory for gas supply“.

Not all articles are available worldwide.

BRONCHIAL ASPIRATION

VACUUM SUCTION

For connecting via CS vacuum connecting hose.

This rail version can be connected to all standard wall rails (10 mm width and 25 mm height). The gas supply is ensured by VAC connection hose.

You can use the bacterial filter MK00514 to reliably protect the central vacuum system against contamination.

Two bacterial filters are included within the scope of delivery.

These bronchial aspirators are suited for use in MRI.

The bronchial aspirators can be used in magnets with static magnetic fields of 1.5 tesla and 3 tesla. Minimum distance from the center of the magnetic field: Magnetic flux density of 40 millitesla (400 gauss).

Technical data:

Maximum vacuum: -95 to 0 kPa (-950 to 0 mbar)

Maximum suction flow: 50 l/min

Dimensions (W/H/D) in mm: 130 × 190 × 120

Maximum sound level: 55 dB(A)

Hose nozzle: Ø 6 mm

Additional Information for Filter:

Bacterial retention: 99.99997%*

Viral retention: 99.99997%*

NaCl retention: 99.9982%**

* Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



A	Vacuum suction NIST, for bronchial suction, gas connection VAC	MK01422
A	Vacuum suction DISS, for bronchial suction, gas connection VAC	MK03140
A	Vacuum suction DISS Canada, for bronchial suction, gas connection VAC	MK03320
B	Vacuum aspirator DISS, version for Atlan® A300/A350 and Perseus® A500	MP00686
	Vacuum aspirator DISS Canada, version for Atlan® A300/A350 and Perseus® A500	MP00687
	Vacuum aspirator NIST, version for Atlan® A300/A350 and Perseus® A500	MP00689

For CS hoses please refer to chapter „Accessory for gas supply“.

BRONCHIAL ASPIRATION

The **VarioAir®** ejector is connected directly to a medical air supply system or to rails. The low-noise ejector generates a vacuum up to max. -90 kPa (-900 mbar) and shows the vacuum level on a turnable, manually adjustable pressure gauge.

All devices offer a pressure relief valve (ex 0.1 kPa) to avoid pressure increase at the patient side.

The vacuum adjusting knob is also the on/off switch to interrupt the suction without reset the values.

A new VarioSafe® disposable filter system reliably protects the VarioAir® series and the patient environment against contamination. One piece of the new VarioSafe® filter system is always included.

Technical data:

Maximum vacuum:	-90 to 0 kPa (-900 to 0 mbar)
Minimum suction flow:	40 l/min
Inlet pressure range:	4 to 5.5 bar
Maximum gas consumption:	55 l/min
Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 200 × 160 (without VarioSafe) 80 × 330 × 160 (with VarioSafe)
Maximum sound level:	57 dB(A)
Hose nozzle:	Ø 6 mm
Pressure relief valve:	Ex 0.1 kPa (1 mbar)

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1

A	Ejector VarioAir® B DIN, for bronchial suction, gas connection AIR	MP00530
B	Ejector VarioAir® B BS, for bronchial suction, gas connection AIR	MP00531
C	Ejector VarioAir® B NF, for bronchial suction, gas connection AIR	MP00532
D	Ejector VarioAir® B AGA, for bronchial suction, gas connection AIR, Swedish standard	MP00533
E	Ejector VarioAir® B UNI, for bronchial suction, gas connection AIR	MP00534
F	Ejector VarioAir® B CARBA, for bronchial suction, gas connection AIR	MP00548

For CS hoses please refer to chapter „Accessory for gas supply“.

Not all articles are available worldwide.



BRONCHIAL ASPIRATION

VarioAir® ejector for rail.

Technical data:

Maximum vacuum:	-90 to 0 kPa (-900 to 0 mbar)
Minimum suction flow:	40 l/min
Inlet pressure range:	4 to 5.5 bar
Maximum gas consumption:	55 l/min
Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 200 × 160 (without VarioSafe) 80 × 330 × 160 (with VarioSafe)
Maximum sound level:	57 dB(A)
Hose nozzle:	Ø 6 mm
Pressure relief valve:	Ex 0.1 kPa (1 mbar)


Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1

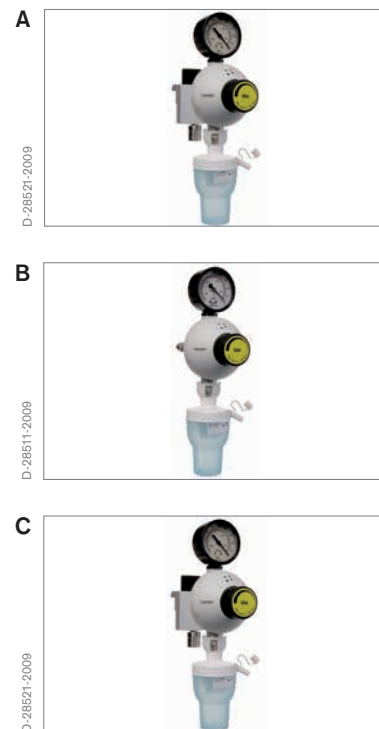

A	Ejector VarioAir® B/OR Rail, for bronchial and OR suction, gas connection NIST AIR	MP00535
A	Ejector VarioAir® B/OR Rail DISS, for bronchial and OR suction, gas connection DISS AIR	MP00638
B	Ejector VarioO2® B DIN, for bronchial suction, gas connection O₂	MP00550
B	Ejector VarioO2® B CARBA, for bronchial suction, gas connection O₂	MP00551
C	Ejector VarioO2® B/OR Rail, for bronchial and OR suction, gas connection NIST O₂	MP00552

For CS hoses please refer to chapter „Accessory for gas supply“.

BRONCHIAL ASPIRATION

VarioAir® ejector for rail.

The gas supply is ensured by AIR connection hose.

This rail version can be connected to all standard wall rails (10 mm width and 25 mm height).

You can use the VarioSafe® disposable filter system to reliably protect the VarioAir® series and the central vacuum system against contamination.

Technical data: VarioAir P (MP00545, MP03221)

Maximum vacuum: -48 to 0 kPa (-480 to 0 mbar)

Maximum suction flow: 20 l/min

Maximum gas consumption: 40 l/min

Maximum sound level: 47 dB(A)

VarioAir B/OP (MP03220)

Maximum vacuum: -90 to 0 kPa (-900 to 0 mbar)

Maximum suction flow: 40 l/min

Maximum gas consumption: 55 l/min

Maximum sound level: 57 dB(A)

Inlet pressure range: 4 to 5.5 bar

Weight: 800 g

Dimensions (W/H/D) in mm: 80 × 200 × 160 (without VarioSafe)

80 × 330 × 160 (with VarioSafe)

Hose nozzle: Ø 6 mm

Pressure relief valve: from 0.1 kPa on (1 mbar)

Additional Information for Filter:

Aerosole and Fluid retention: 100% (with VarioSafe)*

Bacterial retention: 99.99997% (with VarioSafe)**

Viral retention: 99.99922% (with VarioSafe)**

NaCl retention: 99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



**A Ejector VarioAir® P Rail,
for bronchial suction, gas connection NIST AIR,
for connection to all standard wall rails** MP00545

**A Ejector VarioAir® P Rail DISS,
for bronchial suction, gas connection DISS AIR,
for connection to all standard wall rails** MP00639

These bronchial aspirators are suited for use in MRI.

The bronchial aspirators can be used in magnets with static magnetic fields of 1.5 tesla and 3 tesla. Minimum distance from the center of the magnetic field: Magnetic flux density of 40 millitesla (400 gauss).

**Ejector VarioAir® MRI B/OR Rail,
for bronchial suction, gas connection NIST AIR,
for connection to all standard wall rails** MP03220

**Ejector VarioAir® MRI P Rail,
for bronchial suction, gas connection NIST AIR,
for connection to all standard wall rails** MP03221

For CS hoses please refer to chapter „Accessory for gas supply“.

Not all articles are available worldwide.

BRONCHIAL ASPIRATION

EJECTOR SUCTION

For connecting via CS vacuum connecting hose.
The gas supply is ensured by AIR connection hose.
You can use the bacterial filter MK00514 to reliably protect the central vacuum system against contamination.
Two bacterial filters are included within the scope of delivery.
These bronchial aspirators are suited for use in MRI.
The bronchial aspirators can be used in magnets with static magnetic fields of 1.5 tesla and 3 tesla. Minimum distance from the center of the magnetic field: Magnetic flux density of 40 millitesla (400 gauss).

Technical data:

Maximum vacuum:	-90 to 0 kPa (-900 to 0 mbar)
Maximum suction flow:	30 l/min
Pressure range:	2.7 to 6.9 bar
Maximal gas consumption:	60 l/min
Dimensions (W/H/D) in mm:	130 × 190 × 120
Maximum sound level:	55 dB(A)
Pressure relief valve:	ex 1 kPa (10 mbar)
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Bacterial retention:	99.99997%*
Viral retention:	99.99997%*
NaCl retention:	99.9982%**

* Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1



A	Ejector suction NIST AIR O₂ for bronchial suction, gas connection AIR/O₂, for connection to all standard wall rails	<u>MK01418</u>
A	Ejector suction NIST for bronchial suction, gas connection AIR, for connection to all standard wall rails	<u>MK01420</u>
A	Ejector suction NIST, for bronchial suction, gas connection AIR, version for Atlan® A300/A350 and Perseus® A500	<u>MP00688</u>

For CS hoses please refer to chapter „Accessory for gas supply“.

BRONCHIAL ASPIRATION

Drive unit (VarioVac® or VarioAir®) and jar set are joined together, so that you have a complete, convenient and easy solution. The ergonomically shaped handle is positioned exactly at the units centre of gravity, making it easy to transport the aspirator quickly. The gas supply is ensured by CS connection hose. You can use the VarioSafe® disposable filter system to reliably protect the device and the patient environment against contamination. This rail version can be connected to all standard wall rails (10 mm width and 25 mm height).

Technical data:

Weight: 1,600 g (empty)
 Dimensions (W/H/D) in mm: 300 × 440 × 160
 Hose nozzle: Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention: 100% (with VarioSafe)*
 Bacterial retention: 99.99997% (with VarioSafe)**
 Viral retention: 99.99922% (with VarioSafe)**
 NaCl retention: 99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



A Vacuum regulator VarioVac® B compact, for bronchial suction, gas connection NIST VAC, including secretion jar set and integrated vacuum regulator **MP00506**

Technical data:
 Maximum vacuum: -90 to 0 kPa (-900 to 0 mbar)
 Maximum suction flow: 55 l/min
 Maximum sound level: 43 dB(A)

A Vacuum regulator VarioVac® P compact, for bronchial suction, gas connection NIST VAC, including secretion jar set and integrated vacuum regulator **MP00661**

Technical data:
 Maximum vacuum: -55 to 0 kPa (-550 to 0 mbar)
 Maximum suction flow: 20 l/min
 Maximum sound level: 40 dB(A)

For CS hoses please refer to chapter „Accessory for gas supply“.

BRONCHIAL ASPIRATION

**A Ejector VarioAir® B compact,
for bronchial suction, gas connection NIST AIR,
including secretion jar set and integrated ejector** **MP00536**

Technical data:

Maximum vacuum:	-90 to 0 kPa (-900 to 0 mbar)
Minimum suction flow:	40 l/min
Inlet pressure range:	4 to 5.5 bar
Maximum gas consumption:	55 l/min
Maximum sound level:	57 dB(A)
Pressure relief valve:	Ex 0.1 kPa (1 mbar)



**A Ejector VarioAir® P compact,
for bronchial suction, gas connection NIST AIR,
including secretion jar set and integrated ejector** **MP00547**

Technical data:

Maximum vacuum:	-48 to 0 kPa (-480 to 0 mbar)
Maximum suction flow:	20 l/min
Inlet pressure range:	4 to 5.5 bar
Maximum gas consumption:	40 l/min
Maximum sound level:	47 dB(A)
Pressure relief valve:	Ex 0.1 kPa (1 mbar)

For CS hoses please refer to chapter „Accessory for gas supply“.

BRONCHIAL ASPIRATION

Extra lightweight portable secretion jar set for bronchial aspiration. The set is made of high quality and extremely robust plastic and can be autoclaved at up to 134 °C. The two 700 ml jars, which are used as secretion and rinsing jars, are identical, making hospital logistics considerably easier. The jars are made of high quality polysulfone, are reusable, extremely robust and implosion proof. The long service lives makes the jars both environmentally friendly and cost effective.

An automatic overflow safeguard is fitted in the secretion jar cap to prevent secretion from entering the drive unit and the piped medical gas supply. The ergonomically shaped handle is positioned exactly at the units centre of gravity, making it easy to transport the secretion jar set quickly.

- A** **Secretion jar set, with 700 ml secretion jars (2 pcs.) and 1.5 m suction hose, maximum capacity up to the overflow safeguard 2 × 1,000 ml, scale in steps of 10 ml up to max. 700 ml**

2M85056



BRONCHIAL ASPIRATION

A <u>Secretion jar set with clamp</u>	<u>MP00646</u>
B <u>Carrier rack CD</u>	<u>2M85380</u>
C <u>Carrying rack bronchial suction</u>	<u>2M85381</u>



BRONCHIAL ASPIRATION

Holders for attaching a compact bronchial aspirator or a secretion jar set to the wall/rail or to a rail on a ceiling supply unit. The devices are only latched and can be detached and moved easily if this is required.

- | | | |
|---|---|----------------|
| A | Wall holder, with fixing screws, for compact suction unit and secretion jar set | M26859 |
| B | Rail clamp lever type, for attachment of devices to 25 × 10 mm standard rail, for compact suction unit and secretion jar set | M26842 |
| C | Ceiling holder*, rail holder for compact suction systems and secretion jar sets | 2M85119 |
| <p> Holders for holding a central supply connecting hose or an aspiration hose. For attachment to an intensive care rail or compact rail.</p> | | |
| D | CS connection hose holder | 2M85446 |
| E | Suction hose holder, suitable for various diameters | 2M85239 |

*Tradeware



BRONCHIAL ASPIRATION

VacuSmart® disposable cartridge for bronchial suction

VacuSmart® is a disposable cartridge for 700 ml secretion jars. This product has been specifically designed and thoroughly tested to work with Dräger suction units. Therefore there is absolutely no need to modify the installed base. An integrated hydrophobic filter acts as an overflow safeguard by stopping suction of secretion once maximum capacity has been reached. The closed and leakproof system, optional with an integrated gelling agent, reduces the risk of cross contamination to a minimum. The optional gelling agent prevents leakages e.g. during the waste disposal in a waste press. The gelling agent is secured in a water soluble bag.

Technical data:

Maximum capacity:	700 ml
Maximum capacity up to the overflow safeguard:	1,000 ml
Overflow safeguard:	Integrated, hydrophobic filter
Recommended duration of use:	24 h
Gelling agent:	Sodium polyacrylate

A	Secretion jar VacuSmart®, PE (Polyethylene), disposable, 700 ml, 25 pcs.	MX23023
B	Secretion jar VacuSmart® Gel, PE (Polyethylene), disposable, 700 ml, 25 pcs.	MX23025
C	Suction tube VacuSmart® Tube, PE (Polyethylene), disposable, 2 m, w/ finger tip, 25 pcs.	MX23026
You can also enjoy some smart savings by choosing any of the following packing unit 25 pieces sets:		
D	Suction set, disposable, VacuSmart®, 700 ml and VacuSmart® Tube, 25 pcs. each	MP00061
D	Set VacuSmart® Gel and Tube, 25 pcs.	MP03211



BRONCHIAL ASPIRATION

A Bacterial filter 767, 5 pcs. 6723976

A Bacterial filter 767, disposable, for ejectors (model range until 2010) CH00102

The filter 767 reliably protects the ejector and the patient environment against contamination.

Technical data:

Material: PS (polystyrol)

Recommended duration of use: 7 days

Additional Information for Filter:

Bacterial retention: 99.992%****

Viral retention: 99.99%****

****According to EN 1822-1:1998 Filter for suspended particles (HEPA and ULPA)

B Bacterial filter disposable, 20 pcs. MK00514

You can use the bacterial filter to reliably protect the suction devices against contamination.

Technical data:

Material: Acryl

Recommended duration of use: 7 days

Additional Information for Filter:

Bacterial retention: 99.99997%*

Viral retention: 99.99997%*

NaCl retention: 99.9982%**

* Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1

C Bacterial and viral filter system VarioSafe, for VarioAir®, VarioVac®, VarioO2, disposable, 10 pcs. MP00555

The new VarioSafe® disposable filter system reliably protects the VarioAir®, VarioVac® and VarioO2 series and the patient environment against contamination.

Technical data:

Material: PS (polystyrol)

Recommended duration of use: 7 days

Maximum capacity: 90 ml

Additional Information for Filter:

Aerosole and Fluid retention: 100% (with VarioSafe)*

Bacterial retention: 99.99997% (with VarioSafe)**

Viral retention: 99.99922% (with VarioSafe)**

NaCl retention: 99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany


** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



BRONCHIAL ASPIRATION

A	Suction hose, silicone, reusable, 1.5 m, for medical suction	M25780
B	Connection hose, silicone, reusable, ID 5 mm, wall thickness 2 mm, by meter	1203606
	Silicone hose 6 × 2.5, 2 m	5001330
C	Cap, for 700 ml secretion jar, with overflow safeguard	2M85011
D	Cap, for secretion jar, for fistula indicator for continuous and thoracic drainage	2M85013
E	Cap, for secretion jar, w/ overflow safeguard and pressure relief valve (≥ 0.1 kPa)	2M85012

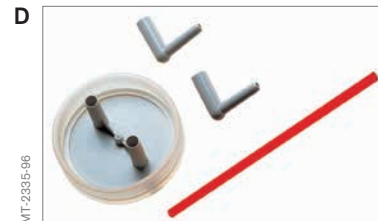

BRONCHIAL ASPIRATION

A	Secretion jar, PSU (Polysulfone), reusable, 300 ml, autoclavable up to 134°C	2M86207
B	Secretion jar, PSU (Polysulfone), reusable, 700 ml, autoclavable up to 134°C	2M85594
C	Elbow, angled hose connection for secretion jar caps	2M19063
D	Shut-off valve for VarioVac®, for temporary opening and closing of suction flow	MP00570
E	Secretion control glass	M07582



BRONCHIAL ASPIRATION

A Pressure gauge for VarioVac® + VarioAir®, double scale kPa/mmHg, -100 kPa (-1000 mbar)	MP00578
B Pressure gauge for VarioVac® + VarioAir®, double scale kPa/mmHg, -60 kPa (-600 mbar)	MP00580
C Pressure gauge -100 kPa MRI, double scale kPa/mmHg, -100 kPa (-1000 mbar)	MP03589
D Pressure gauge -60 kPa MRI, double scale kPa/mmHg, -60 kPa (-600 mbar)	MP03587



BRONCHIAL ASPIRATION

A Universal holder with groove 25 × 5 for wall rails 25 × 8 to 35 × 10 mm, for compact suction system and secretion jar set	2M86757
B Universal holder profile 30 × 5 for wall rails 25 × 8 to 35 × 10 mm	MP00634
C Universal holder for wall rails 25 × 8 to 35 × 10 mm with M5 drill hole 15 mm distance	MP00635
D Universal holder for wall rails 25 × 8 to 35 × 10 mm with M5 drill hole 22 mm distance	MP00636
E Universal holder flex 15 for wall rails 25 × 8 to 35 × 10 mm with bottle screw	MP00637




MEDICAL ASPIRATION

2.2 OR aspiration

During surgical, endoscopic and minimally invasive procedures you do not have to abandon the easy use of these products. A modular design facilitates the easy replacement of components and allows the system to be upgraded with any of Dräger's rail-mounted accessories. The negative pressure can be generated either by a vacuum regulator or an ejector.

An overflow safeguard in the secretion jar and an additional protection jar reduce the risk of contamination of the drive element and the pressure or vacuum generating system.

The OR aspirators can be connected either to the central vacuum system or, when used with an ejector drive unit, to the piped oxygen or medical air supply system.

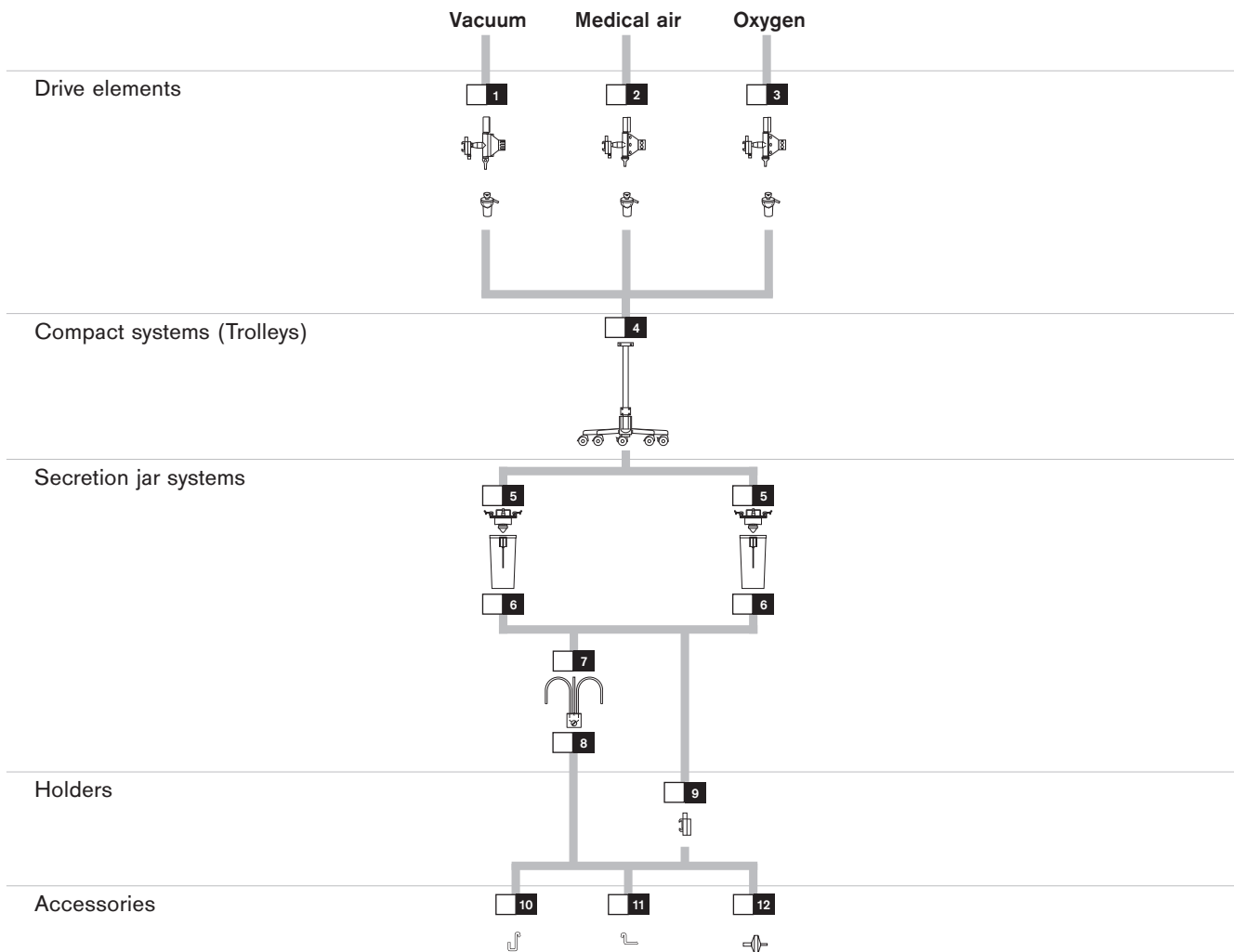

1	MP00505	VarioVac® B/OR Rail
	MP00516	Dual VarioVac® B/OR Rail
2	MP00535	VarioAir® B/OR Rail
3	MP00552	VarioO2® B/OR Rail
4	2M85441	Trolley, for OR suction unit
5	MP03231	Cap secretion jar Medela
6	MP03230	Secretion jar 3 L Medela

7	MP03233	Connection hose Medela
8	MP03232	Switchover valve Medela
9	MP00634	Universal holder profile 30 × 5
10	2M85446	CS connection hose holder
11	2M85239	Suction hose holder
12	MX37021	Bacterial filter, disposable, for Evac 40

OR ASPIRATION

2.2 OR ASPIRATION

A	Vacuum regulator VarioVac® B/OR Rail, for bronchial and OR suction, gas connection NIST VAC	MP00505
Technical data:		
Maximum vacuum:	-90 to 0 kPa (-900 to 0 mbar)	
Maximum suction flow:	55 l/min	
Weight:	600 to 800 g (800 g for rail versions)	
Dimensions (W/H/D) in mm:	80 × 200 × 160 (without VarioSafe) 80 × 330 × 160 (with VarioSafe)	
Maximum sound level:	43 dB(A)	
Hose nozzle:	Ø 6 mm	
<u>Additional Information for Filter:</u>		
Aerosole and Fluid retention:	100% (with VarioSafe)*	
Bacterial retention:	99.99997% (with VarioSafe)**	
Viral retention:	99.99922% (with VarioSafe)**	
NaCl retention:	99.8900% (with VarioSafe)***	
* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany		
** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA		
*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1		
B	Vacuum regulator dual VarioVac® B/OR Rail, for bronchial and OR suction, gas connection NIST VAC	MP00516
Weight: 1,700 g		
Dimensions (W/H/D) in mm: 200 × 210 × 210 (without VarioSafe) 200 × 340 × 210 (with VarioSafe)		
C	Ejector VarioAir® B/OR Rail, for bronchial and OR suction, gas connection NIST AIR	MP00535
C	Ejector VarioO2® B/OR Rail, for bronchial and OR suction, gas connection NIST O₂	MP00552
Technical data:		
Maximum vacuum:	-90 to 0 kPa (-900 to 0 mbar)	
Minimum suction flow:	40 l/min	
Inlet pressure range:	4 to 5.5 bar	
Maximum gas consumption:	55 l/min	
Weight:	800 g	
Dimensions (W/H/D) in mm:	80 × 210 × 160 (without VarioSafe) 80 × 340 × 160 (with VarioSafe)	
Maximum sound level:	57 dB(A)	
Hose nozzle:	Ø 6 mm	
Pressure relief valve:	Ex 0.1 kPa (1 mbar)	
<u>Additional Information for Filter:</u>		
Aerosole and Fluid retention:	100% (with VarioSafe)*	
Bacterial retention:	99.99997% (with VarioSafe)**	
Viral retention:	99.99922% (with VarioSafe)**	
NaCl retention:	99.8900% (with VarioSafe)***	
* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany		
** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA		
*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1		



For CS hoses please refer to chapter „Accessory for gas supply“.

Not all articles are available worldwide.

OR ASPIRATION

A	Trolley*, for OR suction unit, with protection jar, w/o 4 L secretion jars	2M85441
B	Secretion jar 3 L Medela*	MP03230
C	Cap secretion jar Medela*	MP03231
D	Switchover valve Medela*	MP03232
	Connection hose Medela*	MP03233

*Tradeware



OR ASPIRATION

A	Universal holder with groove 25 × 5 for wall rails 25 × 8 to 35 × 10 mm, for compact suction system and secretion jar set	<u>2M86757</u>
B	Universal holder profile 30 × 5 for wall rails 25 × 8 to 35 × 10 mm	<u>MP00634</u>
C	Bacterial filter*, disposable, for Evac 40, 10 pcs.	<u>MX37021</u>
D	CS connection hose holder	<u>2M85446</u>
E	Suction hose holder, suitable for various diameters	<u>2M85239</u>

*Tradeware




MEDICAL ASPIRATION

2.3 Continuous and thoracic drainages

Continuous drainages are used for drainage of secretion, blood and/or air from body cavities over long periods, following operations, particularly after open-thorax surgery.

All of Dräger's continuous drainages are made from high-quality materials and are designed with space in mind. Individual solutions allows you to choose between pneumatic or electrical driven drainages, while also assisting you in optimizing your specific needs and requirements.

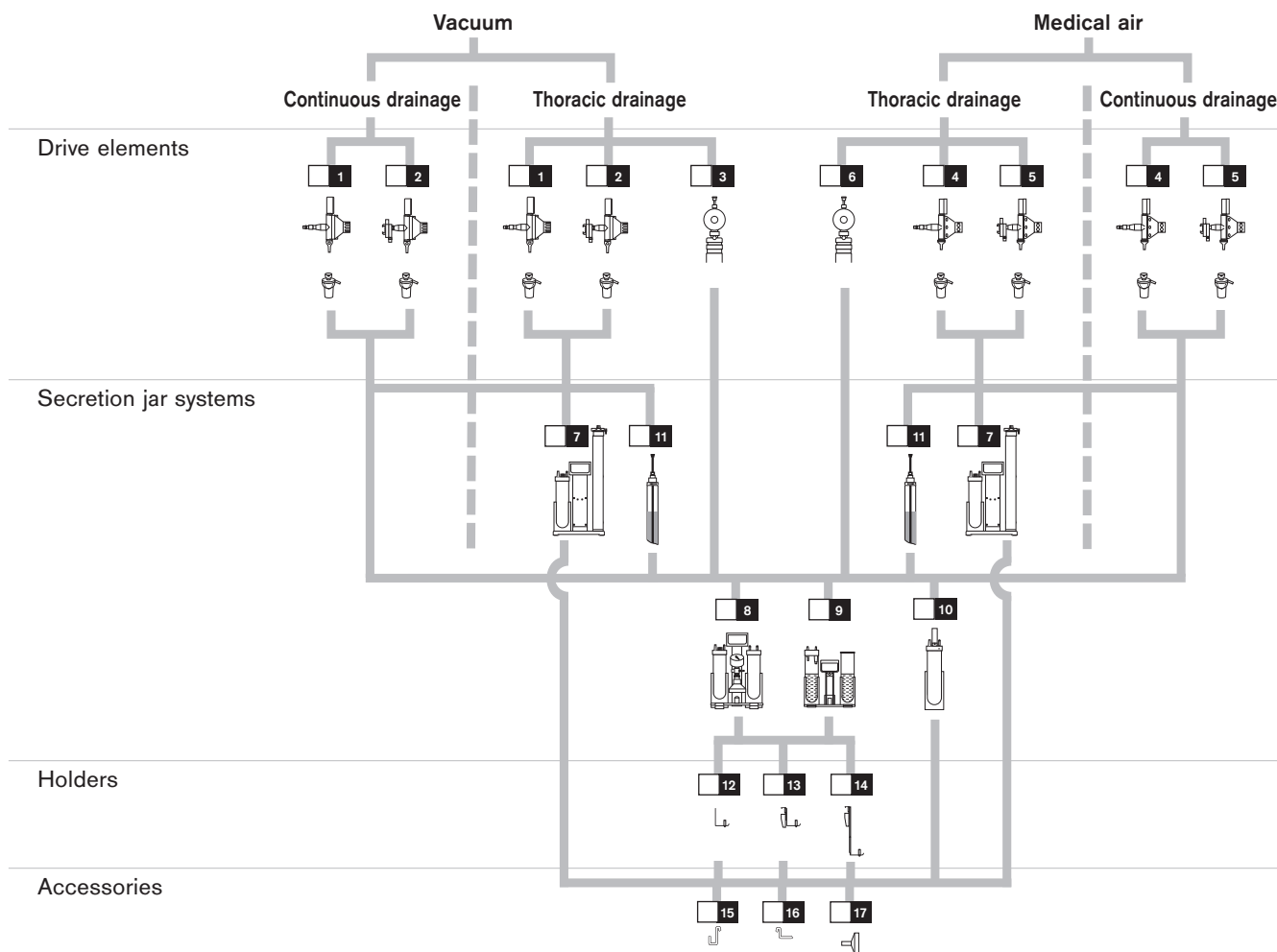

1	MP00507	VarioVac® D DIN
	MP00508	VarioVac® D BS
	MP00509	VarioVac® D NF
	MP00510	VarioVac® D AGA
	MP00515	VarioVac® D CARBA
2	MP00664	VarioVac® T DIN
	MP00511	VarioVac® D Rail
3	MP00663	VarioVac® T Rail
	MP00512	VarioVac® T Rail WV
4	MP00537	VarioAir® D DIN
	MP00538	VarioAir® D BS
	MP00539	VarioAir® D NF
	MP00540	VarioAir® D AGA
	MP00541	VarioAir® D UNI
	MP00549	VarioAir® D CARBA
	MP00553	VarioAir® T DIN

5	MP00542	VarioAir® D Rail
	MP00543	VarioAir® T Rail
6	MP00544	VarioAir® T Rail WV
7	2M85153	Continuous drainage 3
	2M85151	Continuous drainage 1
8	2M85056	Secretion jar set
9	2M85163	Secretion jar set, with rail holder and secretion jar 700 ml
10	MP00518	Water vacuum gauge
11	M26859	Wall holder
12	M26842	Rail clamp lever type
13	2M85119	Ceiling holder
14	2M85446	CS connection hose holder
15	2M85239	Suction hose holder
16	6723976	Bacterial filter 767

CONTINUOUS AND THORACIC DRAINAGES

2.3 CONTINUOUS AND THORACIC DRAINAGES

The vacuum regulator is directly connected to a central vacuum supply, or mounted on a rail. The smoothly adjustable vacuum, up to max. -160 mbar (-16 kPa), is displayed on a rotating manometer, legible from all sides.

All vacuum regulators include a shut-off valve.

You can use the VarioSafe® disposable filter system to reliably protect the **VarioVac®** series and the central vacuum system against contamination. A VarioSafe® filter system is included with every VarioVac® vacuum regulator.

Technical data:

Maximum vacuum: -16 to 0 kPa (-160 to 0 mbar)

Maximum suction flow: 18 l/min

Weight: 600 to 800 g

(800 g for rail versions)

Dimensions (W/H/D) in mm: 80 × 210 × 160 (without VarioSafe)

80 × 340 × 160 (with VarioSafe)

Maximum sound level: 50 dB(A)

Hose nozzle: Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention: 100% (with VarioSafe)*

Bacterial retention: 99.99997% (with VarioSafe)**

Viral retention: 99.99922% (with VarioSafe)**

NaCl retention: 99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA,
Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1

A	Vacuum regulator VarioVac® D DIN, for continuous and wound drainage, gas connection VAC	MP00507
B	Vacuum regulator VarioVac® D BS, for continuous and wound drainage, gas connection VAC	MP00508
C	Vacuum regulator VarioVac® D NF, for continuous and wound drainage, gas connection VAC	MP00509
D	Vacuum regulator VarioVac® D AGA, for continuous and wound drainage, gas connection VAC	MP00510
E	Vacuum regulator VarioVac® D CARBA, for continuous and wound drainage, gas connection VAC	MP00515
F	Vacuum regulator VarioVac® D Rail, for continuous and wound drainage, gas connection VAC	MP00511

For CS hoses please refer to chapter „Accessory for gas supply“.



CONTINUOUS AND THORACIC DRAINAGES

A Vacuum regulator VarioVac® T Rail, for thoracic drainage, gas connection VAC, for connection to all standard wall rails

MP00663

For connecting via CS vacuum connecting hose.
The gas supply is ensured by VAC connection hose.

Technical data:

Maximum vacuum:	-7 to 0 kPa (-70 to 0 mbar)
Minimum Flow:	15 l/min
Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 210 × 160 (without VarioSafe) 80 × 340 × 160 (with VarioSafe)
Maximum sound level:	50 dB(A)
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA,
Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1



B Vacuum regulator VarioVac® T DIN, for thoracic drainage, gas connection VAC, for connecting to the central vacuum supply

MP00664

Technical data:

Maximum vacuum:	-7 to 0 kPa (-70 to 0 mbar)
Minimum Flow:	15 l/min
Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 210 × 160 (without VarioSafe) 80 × 340 × 160 (with VarioSafe)
Maximum sound level:	50 dB(A)
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA,
Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1

For CS hoses please refer to chapter „Accessory for gas supply“.

CONTINUOUS AND THORACIC DRAINAGES

**A Vacuum regulator VarioVac® T Rail WV,
for thoracic drainage, gas connection VAC,
water vacuum gauge included,
for connection to all standard wall rails,
double-scaled immersion pipe (kPa and cmH₂O)**

MP00512

For connecting via CS vacuum connection hose.
The filter MK00514 protects the product and the
central vacuum supply against contamination.
Two filters are included in the scope of delivery.

Technical data:

Maximum vacuum:	-5 to 0 kPa (-50 to 0 mbar)
Minimum suction flow:	15 l/min
Weight:	1,100 g
Dimensions (W/H/D) in mm:	80 × 860 × 160
Maximum sound level:	40 dB(A)
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Bacterial retention:	99.99997%*
Viral retention:	99.99997%*
NaCl retention:	99.9982%**

* Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1



For CS hoses please refer to chapter „Accessory for gas supply“.

CONTINUOUS AND THORACIC DRAINAGES

The **VarioAir®** ejector is connected directly to a medical air supply system or to rails. The low-noise ejector generates a vacuum up to max. -16 kPa (-160 mbar) and shows the vacuum level on a turnable, manually adjustable pressure gauge. All devices offer a pressure relief valve (ex 0.1 kPa) to avoid pressure increase at the patient side.

A new **VarioSafe®** disposable filter system reliably protects the **VarioAir®** series and the patient environment against contamination. One piece of the new **VarioSafe®** filter system is always included.

Technical data:

Maximum vacuum:	-16 to 0 kPa (-160 to 0 mbar)
Maximum suction flow:	18 l/min
Inlet pressure range:	4 to 5.5 bar
Maximum gas consumption:	25 l/min
Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 200 × 160 (without VarioSafe) 80 × 330 × 160 (with VarioSafe)

Maximum sound level:	45 dB(A)
Hose nozzle:	Ø 6 mm
Pressure relief valve:	Ex 0.1 kPa (1 mbar)

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1

A	Ejector VarioAir® D DIN, for continuous and wound drainage, gas connection AIR	MP00537
B	Ejector VarioAir® D BS, for continuous and wound drainage, gas connection AIR	MP00538
C	Ejector VarioAir® D NF, for continuous and wound drainage, gas connection AIR	MP00539
D	Ejector VarioAir® D AGA, for continuous and wound drainage, gas connection AIR	MP00540
E	Ejector VarioAir® D UNI, for continuous and wound drainage, gas connection AIR	MP00541
F	Ejector VarioAir® D CARBA, for continuous and wound drainage, gas connection AIR	MP00549
G	Ejector VarioAir® D Rail, for continuous and wound drainage, gas connection AIR, for connection to all standard wall rails	MP00542
	For gas supply via CS vacuum connection hose. Connection thread: NIST.	

For CS hoses please refer to chapter „Accessory for gas supply“.

Not all articles are available worldwide.



CONTINUOUS AND THORACIC DRAINAGES

A Ejector VarioAir® T Rail, for thoracic drainage, gas connection NIST AIR, for connection to all standard wall rails

MP00543

The gas supply is ensured by AIR connection hose.

Technical data:

Maximum vacuum:	-7 to 0 kPa (-70 to 0 mbar)
Minimum suction flow:	15 l/min
Maximum gas consumption:	30 l/min
Maximum sound level:	40 dB(A)

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA,
Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1



B Ejector VarioAir® T DIN, for thoracic drainage, gas connection DIN, for connecting to the central AIR supply according to DIN standard

MP00553

Technical data:

Maximum vacuum:	-7 to 0 kPa (-70 to 0 mbar)
Minimum suction flow:	15 l/min
Maximum gas consumption:	30 l/min
Maximum sound level:	40 dB(A)

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA,
Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1

For CS hoses please refer to chapter „Accessory for gas supply“.

CONTINUOUS AND THORACIC DRAINAGES

A Ejector VarioAir® T Rail WV, for thoracic drainage, with water vacuum gauge, for connection to all standard wall rails

MP00544

For gas supply via CS vacuum connection hose.

Connection thread: NIST.

Double-scaled immersion pipe (kPa and cmH₂O).

The filter MK00514 protects the product and the central vacuum supply against contamination.

Two filters are included in the scope of delivery.

Technical data:

Maximum vacuum: -5 to 0 kPa (-50 to 0 mbar)

Minimum suction flow: 15 l/min

Maximum gas consumption: 30 l/min

Weight: 1,100 g

Dimensions (W/H/D) in mm: 80 × 860 × 160

Maximum sound level: 50 dB(A)

Hose nozzle: Ø 6 mm

Additional Information for Filter:

Bacterial retention: 99.99997%*

Viral retention: 99.99997%*

NaCl retention: 99.9982%**

* Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1



For CS hoses please refer to chapter „Accessory for gas supply“.

CONTINUOUS AND THORACIC DRAINAGES

The vacuum can either be generated by an ejector or via a vacuum regulator. When used with the ejector, a bacterial filter reliably protects against contamination of the environment. All these pneumatic continuous drainages can be used in potentially explosive areas.

All parts which may become contaminated can be autoclaved at up to 134 °C.

The continuous drainage systems 1 and 3 can be used for continuous- and wound drainage as well as for thoracic drainage in the combination with a regulator from the VarioVac® D and VarioAir® D series. An additional precision regulator make it easy to adjust lowest vacuum levels.

- A Continuous drainage 1,
for continuous and wound drainage,
precision regulator, fistula indicator
and bacterial filter 767 included, capacity 700 ml,
maximum capacity up to the overflow safeguard 1,000 ml,
scale in steps of 10 ml up to max. 700 ml** **2M85151**

Technical data:

Maximum vacuum:	-10 to 0 kPa (-100 to 0 mbar)
Maximum suction flow:	14 l/min
Weight:	3 kg (empty)
Dimensions (W/H/D) in mm:	285 × 340 × 160

- B Continuous drainage 3,
for continuous, wound, and thoracic drainage,
precision regulator, water vacuummeter, fistula indicator
and bacterial filter 767 included, capacity 700 ml,
maximum capacity up to the overflow safeguard 1,000 ml,
scale in steps of 10 ml up to max. 700 ml** **2M85153**

Technical data:

Maximum vacuum:	-5 to 0 kPa (-50 to 0 mbar)
Maximum suction flow:	14 l/min
Weight:	5.5 kg (empty)
Dimensions (W/H/D) in mm:	260 × 800 × 180



CONTINUOUS AND THORACIC DRAINAGES

- A Vacuum regulator VarioVac® T Rail WV, for thoracic drainage, gas connection VAC, water vacuum gauge included, for connection to all standard wall rails, double-scaled immersion pipe (kPa and cmH₂O)** MP00512

For connecting via CS vacuum connection hose.
The filter MK00514 protects the product and the central vacuum supply against contamination.
Two filters are included in the scope of delivery.

Technical data:

Maximum vacuum:	-5 to 0 kPa (-50 to 0 mbar)
Minimum suction flow:	15 l/min
Weight:	1,100 g
Dimensions (W/H/D) in mm:	80 × 860 × 160
Maximum sound level:	40 dB(A)
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Bacterial retention:	99.99997%*
Viral retention:	99.99997%*
NaCl retention:	99.9982%**

* Test has been performed by Nelson Laboratories Inc.,
Dec 2020, Salt Lake City, USA

** Test has been performed by Dräger Safety AG & Co. KGaA,
March 2021, Lübeck, Germany, according to EN ISO 23328-1
For CS hoses please refer to chapter „Accessory for gas supply“.

- B Water vacuum gauge, for thoracic drainage, combination w/ VarioVac® D and VarioAir® D, for connection to all standard wall rails, double-scaled immersion pipe (kPa and cmH₂O)** MP00518

Technical data:

Maximum vacuum:	-5 to 0 kPa (-50 to 0 mbar)
Maximum suction flow:	10 l/min
Pressure relief valve:	from 0.1 kPa on (1 mbar)
Hose nozzle:	2 × Ø 6 mm
Weight:	1 kg (empty)
Dimensions (W/H/D) in mm:	80 × 860 × 160

- C Secretion jar set, with rail holder and secretion jar 700 ml, for continuous drainage units** 2M85163
- D Cylinder support** 2M85180



CONTINUOUS AND THORACIC DRAINAGES

A Bacterial filter 767, 5 pcs. 6723976

A Bacterial filter 767, disposable, for ejectors (model range until 2010) CH00102

The filter 767 reliably protects the ejector and the patient environment against contamination.

Technical data:

Material: PS (polystyrol)

Recommended duration of use: 7 days

Additional Information for Filter:

Bacterial retention: 99.992%****

Viral retention: 99.99%****

****According to EN 1822-1:1998 Filter for suspended particles (HEPA and ULPA)

B Bacterial filter disposable, 20 pcs. MK00514

You can use the bacterial filter to reliably protect the suction devices against contamination.

Technical data:

Material: Acryl

Recommended duration of use: 7 days

Additional Information for Filter:

Bacterial retention: 99.99997%*

Viral retention: 99.99997%*

NaCl retention: 99.9982%**

* Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1

C Bacterial and viral filter system VarioSafe, for VarioAir®, VarioVac®, VarioO2, disposable, 10 pcs. MP00555

The new VarioSafe® disposable filter system reliably protects the VarioAir®, VarioVac® and VarioO2 series and the patient environment against contamination.

Technical data:

Material: PS (polystyrol)

Recommended duration of use: 7 days

Maximum capacity: 90 ml

Additional Information for Filter:

Aerosole and Fluid retention: 100% (with VarioSafe)*

Bacterial retention: 99.99997% (with VarioSafe)**

Viral retention: 99.99922% (with VarioSafe)**

NaCl retention: 99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany


** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1




CONTINUOUS AND THORACIC DRAINAGES

A	Immersion tube, for continuous drainages, not scaled	M16395
	Rubber plug, for immersion tube of water vacuum gauge	2M85139
B	Water container, for continuous drainage system 2 – 4, w/o screw thread	M20668
C	Pressure gauge for VarioVac® + VarioAir®, double scale kPa/mmHg, -16 kPa (-160 mbar)	MP00579
D	Pressure gauge for VarioVac® + VarioAir®, double scale kPa/mmHg, -8 kPa (-80 mbar)	MP03597

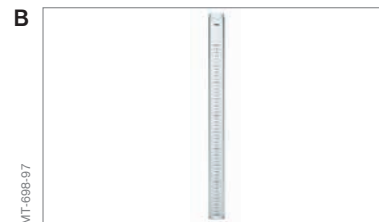

CONTINUOUS AND THORACIC DRAINAGES

A Immersion tube, for water vacuum gauge of for VarioVac® T + VarioAir® T, scaled	MP00569
B Locking screw, for immersion tube, with gaskets, for VarioVac® T Rail WV (MP00512), VarioAir® T Rail WV (MP00544) and water vacuum gauge (MP00518)	MP00596
C Water container with screw thread + gaskets, for VarioVac® T Rail WV (MP00512), VarioAir® T Rail WV (MP00544) and water vacuum gauge (MP00518)	MP00565
D Cap for water vacuum gauge, for VarioVac® T Rail WV (MP00512), VarioAir® T Rail WV (MP00544) and water vacuum gauge (MP00518), including gaskets, combination possible with filter MK00514 and plug nozzle MP00566	MP00595
E Plug nozzle water vacuum gauge, for VarioVac® T Rail WV (MP00512), VarioAir® T Rail WV (MP00544) and water vacuum gauge (MP00518)	MP00566
E Hose nozzle 8 mm small, connection nozzle double sided 8 mm	5000833

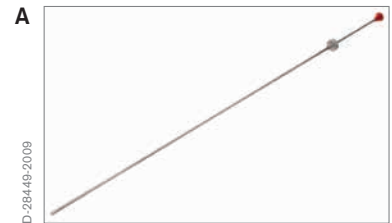

ACCESSORY FOR GAS SUPPLY

3. ACCESSORY FOR GAS SUPPLY

3.1 CS HOSES

CS hoses,
several hoses to connect flowmeter
with the central supply

8601697

3.1.1 GERMAN STANDARD

A	O ₂ CS hose 0.9 m (2.95 ft) NIST EN DIN	M36729
	O ₂ CS hose 1.5 m (4.9 ft) NIST EN DIN	M35489
	O ₂ CS hose 3 m (9.8 ft) NIST EN DIN	M35490
	O ₂ CS hose 5 m (16.4 ft) NIST EN DIN	M35491
B	AIR CS hose 0.9 m (2.95 ft) NIST EN DIN	M36730
	AIR CS hose 1.5 m (4.9 ft) NIST EN DIN	M35495
	AIR CS hose 3 m (9.8 ft) NIST EN DIN	M35496
	AIR CS hose 5 m (16.4 ft) NIST EN DIN	M35497
C	VAC CS hose 1.5 m (4.9 ft) NIST EN DIN	M35498
	VAC CS hose 3 m (9.8 ft) NIST EN DIN	M35499
	VAC CS hose 5 m (16.4 ft) NIST EN DIN	M36000
D	AIR/O ₂ CS hose 1.5 m (4.9 ft) NIST EN DIN	M36001
	AIR/O ₂ CS hose 3 m (9.8 ft) NIST EN DIN	M36002
	AIR/O ₂ CS hose 5 m (16.4 ft) NIST EN DIN	M36003
E	N ₂ O CS hose 0.9 m (2.95 ft) NIST EN DIN	M36731
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN DIN	M35492
	N ₂ O CS hose 3 m (9.8 ft) NIST EN DIN	M35493
	N ₂ O CS hose 5 m (16.4 ft) NIST EN DIN	M35494



ACCESSORY FOR GAS SUPPLY

A	O₂ CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36004
	O₂ CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36005
	O₂ CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36006
B	AIR CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36010
	AIR CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36011
	AIR CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36012
C	VAC CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36013
	VAC CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36014
	VAC CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36015
D	AIR/O₂ CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36016
	AIR/O₂ CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36017
	AIR/O₂ CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36018
E	N₂O CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36007
	N₂O CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36008
	N₂O CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36009



ACCESSORY FOR GAS SUPPLY

3.1.2 FRENCH STANDARD

A	O₂ CS hose 1.5 m (4.9 ft) NIST EN NF	M35416
	O₂ CS hose 3 m (9.8 ft) NIST EN NF	M35417
	O₂ CS hose 5 m (16.4 ft) NIST EN NF	M35330
B	AIR CS hose 1.5 m (4.9 ft) NIST EN NF	M35420
	AIR CS hose 3 m (9.8 ft) NIST EN NF	M35421
	AIR CS hose 5 m (16.4 ft) NIST EN NF	M35332
C	VAC CS hose 1.5 m (4.9 ft) NIST EN NF	M35422
	VAC CS hose 3 m (9.8 ft) NIST EN NF	M35423
	VAC CS hose 5 m (16.4 ft) NIST EN NF	M35424
D	N₂O CS hose 1.5 m (4.9 ft) NIST EN NF	M35418
	N₂O CS hose 3 m (9.8 ft) NIST EN NF	M35419
	N₂O CS hose 5 m (16.4 ft) NIST EN NF	M35331
E	O₂ CS hose 1.5 m (4.9 ft) NIST NEUTRAL NF	M35434
	O₂ CS hose 3 m (9.8 ft) NIST NEUTRAL NF	M35435
	O₂ CS hose 5 m (16.4 ft) NIST NEUTRAL NF	M35436
F	AIR CS hose 1.5 m (4.9 ft) NIST NEUTRAL NF	M35452
	AIR CS hose 3 m (9.8 ft) NIST NEUTRAL NF	M35453
	AIR CS hose 5 m (16.4 ft) NIST NEUTRAL NF	M35454
G	VAC CS hose 1.5 m (4.9 ft) NIST NEUTRAL NF	M35455
	VAC CS hose 3 m (9.8 ft) NIST NEUTRAL NF	M35456
	VAC CS hose 5 m (16.4 ft) NIST NEUTRAL NF	M35457



Not all articles are available worldwide.

ACCESSORY FOR GAS SUPPLY

A	<u>N₂O CS hose 1.5 m (4.9 ft) NIST NEUTRAL NF</u>	<u>M35437</u>
	<u>N₂O CS hose 3 m (9.8 ft) NIST NEUTRAL NF</u>	<u>M35438</u>
	<u>N₂O CS hose 5 m (16.4 ft) NIST NEUTRAL NF</u>	<u>M35439</u>
B	<u>O₂ CS hose 0.9 m (2.95 ft) NF EN NF</u>	<u>M36723</u>
	<u>O₂ CS hose 1.5 m (4.9 ft) NF EN NF</u>	<u>M35425</u>
	<u>O₂ CS hose 3 m (9.8 ft) NF EN NF</u>	<u>M35426</u>
	<u>O₂ CS hose 5 m (16.4 ft) NF EN NF</u>	<u>M35342</u>
C	<u>AIR CS hose 0.9 m (2.95 ft) NF EN NF</u>	<u>M36724</u>
	<u>AIR CS hose 1.5 m (4.9 ft) NF EN NF</u>	<u>M35429</u>
	<u>AIR CS hose 3 m (9.8 ft) NF EN NF</u>	<u>M35430</u>
	<u>AIR CS hose 5 m (16.4 ft) NF EN NF</u>	<u>M35343</u>
D	<u>VAC CS hose 1.5 m (4.9 ft) NF EN NF</u>	<u>M35431</u>
	<u>VAC CS hose 3 m (9.8 ft) NF EN NF</u>	<u>M35432</u>
	<u>VAC CS hose 5 m (16.4 ft) NF EN NF</u>	<u>M35433</u>
E	<u>N₂O CS hose 0.9 m (2.95 ft) NF EN NF</u>	<u>M36725</u>
	<u>N₂O CS hose 1.5 m (4.9 ft) NF EN NF</u>	<u>M35427</u>
	<u>N₂O CS hose 3 m (9.8 ft) NF EN NF</u>	<u>M35428</u>
	<u>N₂O CS hose 5 m (16.4 ft) NF EN NF</u>	<u>M35344</u>



ACCESSORY FOR GAS SUPPLY

3.1.3 BRITISH STANDARD

A	O ₂ CS hose 0.9 m (2.95 ft) NIST EN BS	M36726
	O ₂ CS hose 1.5 m (4.9 ft) NIST EN BS	M35389
	O ₂ CS hose 3 m (9.8 ft) NIST EN BS	M35390
	O ₂ CS hose 5 m (16.4 ft) NIST EN BS	M35333
B	AIR CS hose 0.9 m (2.95 ft) NIST EN BS	M36727
	AIR CS hose 1.5 m (4.9 ft) NIST EN BS	M35393
	AIR CS hose 3 m (9.8 ft) NIST EN BS	M35394
	AIR CS hose 5 m (16.4 ft) NIST EN BS	M35335
C	VAC CS hose 1.5 m (4.9 ft) NIST EN BS	M35395
	VAC CS hose 3 m (9.8 ft) NIST EN BS	M35396
	VAC CS hose 5 m (16.4 ft) NIST EN BS	M35397
D	N ₂ O CS hose 0.9 m (2.95 ft) NIST EN BS	M36728
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN BS	M35391
	N ₂ O CS hose 3 m (9.8 ft) NIST EN BS	M35392
	N ₂ O CS hose 5 m (16.4 ft) NIST EN BS	M35334



ACCESSORY FOR GAS SUPPLY

A	<u>O₂ CS hose 1.5 m (4.9 ft) NIST EN BS straight</u>	<u>M35398</u>
	<u>O₂ CS hose 3 m (9.8 ft) NIST EN BS straight</u>	<u>M35399</u>
	<u>O₂ CS hose 5 m (16.4 ft) NIST EN BS straight</u>	<u>M35339</u>
B	<u>AIR CS hose 1.5 m (4.9 ft) NIST BS straight</u>	<u>M35402</u>
	<u>AIR CS hose 3 m (9.8 ft) NIST BS straight</u>	<u>M35403</u>
	<u>AIR CS hose 5 m (16.4 ft) NIST BS straight</u>	<u>M35340</u>
C	<u>VAC CS hose 1.5 m (4.9 ft) NIST BS straight</u>	<u>M35404</u>
	<u>VAC CS hose 3 m (9.8 ft) NIST BS straight</u>	<u>M35405</u>
	<u>VAC CS hose 5 m (16.4 ft) NIST BS straight</u>	<u>M35406</u>
D	<u>N₂O CS hose 1.5 m (4.9 ft) NIST EN BS straight</u>	<u>M35400</u>
	<u>N₂O CS hose 3 m (9.8 ft) NIST EN BS straight</u>	<u>M35401</u>
	<u>N₂O CS hose 5 m (16.4 ft) NIST EN BS straight</u>	<u>M35341</u>



ACCESSORY FOR GAS SUPPLY

3.1.4 SWEDEN STANDARD

A	O₂ CS hose 1.5 m (4.9 ft) NIST EN SS	M35375
	O₂ CS hose 3 m (9.8 ft) NIST EN SS	M35376
	O₂ CS hose 5 m (16.4 ft) NIST EN SS	M35377
B	AIR CS hose 1.5 m (4.9 ft) NIST EN SS	M35383
	AIR CS hose 3 m (9.8 ft) NIST EN SS	M35384
	AIR CS hose 5 m (16.4 ft) NIST EN SS	M35385
C	VAC CS hose 1.5 m (4.9 ft) NIST EN SS	M35386
	VAC CS hose 3 m (9.8 ft) NIST EN SS	M35387
	VAC CS hose 5 m (16.4 ft) NIST EN SS	M35388
D	N₂O CS hose 1.5 m (4.9 ft) NIST EN SS	M35378
	N₂O CS hose 3 m (9.8 ft) NIST EN SS	M35379
	N₂O CS hose 5 m (16.4 ft) NIST EN SS	M35382



ACCESSORY FOR GAS SUPPLY

3.1.5 SWISS STANDARD

A	<u>O₂ CS hose 1.5 m (4.9 ft) NIST EN Carba</u>	<u>M35469</u>
	<u>O₂ CS hose 3 m (9.8 ft) NIST EN Carba</u>	<u>M35472</u>
	<u>O₂ CS hose 5 m (16.4 ft) NIST EN Carba</u>	<u>M35345</u>
B	<u>AIR CS hose 1.5 m (4.9 ft) NIST EN Carba</u>	<u>M35475</u>
	<u>AIR CS hose 3 m (9.8 ft) NIST EN Carba</u>	<u>M35476</u>
	<u>AIR CS hose 5 m (16.4 ft) NIST EN Carba</u>	<u>M35346</u>
C	<u>VAC CS hose 1.5 m (4.9 ft) NIST EN Carba</u>	<u>M35477</u>
	<u>VAC CS hose 3 m (9.8 ft) NIST EN Carba</u>	<u>M35478</u>
	<u>VAC CS hose 5 m (16.4 ft) NIST EN Carba</u>	<u>M35479</u>
D	<u>N₂O CS hose 1.5 m (4.9 ft) NIST EN Carba</u>	<u>M35473</u>
	<u>N₂O CS hose 3 m (9.8 ft) NIST EN Carba</u>	<u>M35474</u>
	<u>N₂O CS hose 5 m (16.4 ft) NIST EN Carba</u>	<u>M35347</u>



ACCESSORY FOR GAS SUPPLY

3.1.6 ITALIAN STANDARD


A	O₂ CS hose 1.5 m (4.9 ft) NIST EN UNI	M35407
	O₂ CS hose 3 m (9.8 ft) NIST EN UNI	M35408
	O₂ CS hose 5 m (16.4 ft) NIST EN UNI	M35357
B	AIR CS hose 1.5 m (4.9 ft) NIST EN UNI	M35411
	AIR CS hose 3 m (9.8 ft) NIST EN UNI	M35412
	AIR CS hose 5 m (16.4 ft) NIST EN UNI	M35358
C	VAC CS hose 1.5 m (4.9 ft) NIST EN UNI	M35413
	VAC CS hose 3 m (9.8 ft) NIST EN UNI	M35414
	VAC CS hose 5 m (16.4 ft) NIST EN UNI	M35415
D	N₂O CS hose 1.5 m (4.9 ft) NIST EN UNI	M35409
	N₂O CS hose 3 m (9.8 ft) NIST EN UNI	M35410
	N₂O CS hose 5 m (16.4 ft) NIST EN UNI	M35359



ACCESSORY FOR GAS SUPPLY

3.1.7 AUSTRALIAN STANDARD


A	O₂ CS hose 1.5 m (4.9 ft) NIST EN AS	M36153
	O₂ CS hose 3 m (9.8 ft) NIST EN AS	M36152
	O₂ CS hose 5 m (16.4 ft) NIST EN AS	M36151
B	AIR CS hose 1.5 m (4.9 ft) NIST EN AS	M36150
	AIR CS hose 3 m (9.8 ft) NIST EN AS	M36149
	AIR CS hose 5 m (16.4 ft) NIST EN AS	M36148
C	N₂O CS hose 1.5 m (4.9 ft) NIST EN AS	M36156
	N₂O CS hose 3 m (9.8 ft) NIST EN AS	M36155
	N₂O CS hose 5 m (16.4 ft) NIST EN AS	M36154
D	VAC CS hose 1.5 m (4.9 ft) NIST EN AS	M36826
	VAC CS hose 3 m (9.8 ft) NIST EN AS	M36827
	VAC CS hose 5 m (16.4 ft) NIST EN AS	M36828

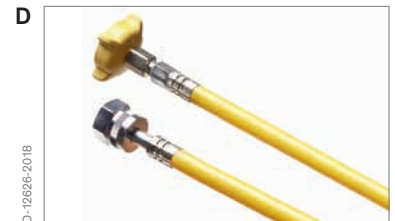

ACCESSORY FOR GAS SUPPLY

3.1.8 JAPANESE STANDARD

A	<u>O₂ CS hose 1.5 m (4.9 ft) DISS US KAWA</u>	<u>M35806</u>
	<u>O₂ CS hose 3 m (9.8 ft) DISS US KAWA</u>	<u>M35805</u>
	<u>O₂ CS hose 5 m (16.4 ft) DISS US KAWA</u>	<u>M35804</u>
	<u>O₂ CS hose 8 m (26,25 ft) DISS US KAWA</u>	<u>M35866</u>
B	<u>AIR CS hose 1.5 m (4.9 ft) DISS US KAWA</u>	<u>M35803</u>
	<u>AIR CS hose 3 m (9.8 ft) DISS US KAWA</u>	<u>M35802</u>
	<u>AIR CS hose 5 m (16.4 ft) DISS US KAWA</u>	<u>M35801</u>
	<u>AIR CS hose 8 m (26,25 ft) DISS US KAWA</u>	<u>M35865</u>
C	<u>N₂O CS hose 1.5 m (4.9 ft) DISS US KAWA</u>	<u>M35809</u>
	<u>N₂O CS hose 3 m (9.8 ft) DISS US KAWA</u>	<u>M35808</u>
	<u>N₂O CS hose 5 m (16.4 ft) DISS US KAWA</u>	<u>M35807</u>
	<u>N₂O CS hose 8 m (26,25 ft) DISS US KAWA</u>	<u>M35867</u>
D	<u>VAC CS hose 5 m (16.4 ft) DISS JIS KAWA</u>	<u>M35820</u>

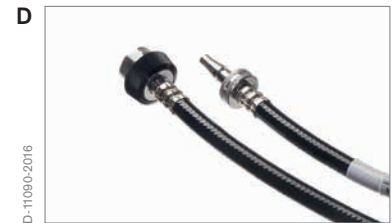


ACCESSORY FOR GAS SUPPLY

3.1.9 US AMERICAN STANDARD

A	O₂ CS hose DISS N/DISS N 10L*	4185140
	O₂ CS hose DISS N/DISS N 15L*	4185141
	O₂ CS hose DISS N/DISS N 20L*	4199714
B	O₂ CS hose DISS Hand 10L*	4199592
	O₂ CS hose DISS Hand 15L*	4199591
C	AIR CS hose DISS N/DISS N 10L*	4199605
	AIR CS hose DISS N/DISS N 15L*	4199604
	AIR CS hose DISS N/DISS N 20L*	4199712
D	AIR CS hose DISS Hand 10L*	4199594
	AIR CS hose DISS Hand 15L*	4199593
E	VAC CS hose DISS N/DISS N 10L*	4199602
	VAC CS hose DISS N/DISS N 15L*	4199601
	VAC CS hose DISS N/DISS N 20L*	4199600
F	VAC CS hose DISS Hand 10L*	4199590
	VAC CS hose DISS Hand 15L*	4199589

*Tradeware

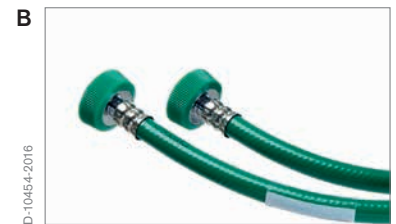

ACCESSORY FOR GAS SUPPLY

A	<u>N₂O CS hose DISS N/DISS N 10L*</u>	<u>4185142</u>
	<u>N₂O CS hose DISS N/DISS N 15L*</u>	<u>4185143</u>
	<u>N₂O CS hose DISS N/DISS N 20L*</u>	<u>4199713</u>
B	<u>N₂O CS hose DISS Hand 10L*</u>	<u>4199596</u>
	<u>N₂O CS hose DISS Hand 15L*</u>	<u>4199595</u>
C	<u>Hose ASM-WAGD DISS N/ DISS N 10L*</u>	<u>MP20022</u>
	<u>Hose ASM-WAGD DISS N/ DISS N 15L*</u>	<u>MP20023</u>
	<u>Hose ASM-WAGD DISS N/ DISS N 20L*</u>	<u>MP20024</u>
D	<u>Hose ASM-WAGD HAND DISS N 10L*</u>	<u>MP20029</u>
	<u>Hose ASM-WAGD HAND DISS N 15L*</u>	<u>MP20030</u>

*Tradeware



ACCESSORY FOR GAS SUPPLY

3.1.10 CANADIAN STANDARD

Now with check valve as standard.

A	<u>O₂ CS hose 1.5 m (4.9 ft) DISSB EN HAND</u>	<u>M36290</u>
	<u>O₂ CS hose 3 m (9.8 ft) DISSB EN HAND</u>	<u>M36291</u>
	<u>O₂ CS hose 5 m (16.4 ft) DISSB EN HAND</u>	<u>M36292</u>
B	<u>AIR CS hose 1.5 m (4.9 ft) DISSB EN HAND</u>	<u>M36296</u>
	<u>AIR CS hose 3 m (9.8 ft) DISSB EN HAND</u>	<u>M36297</u>
	<u>AIR CS hose 5 m (16.4 ft) DISSB EN HAND</u>	<u>M36298</u>
C	<u>VAC CS hose 1.5 m (4.9 ft) DISSB EN HAND</u>	<u>M36299</u>
	<u>VAC CS hose 3 m (9.8 ft) DISSB EN HAND</u>	<u>M36300</u>
	<u>VAC CS hose 5 m (16.4 ft) DISSB EN HAND</u>	<u>M36301</u>
D	<u>N₂O CS hose 1.5 m (4.9 ft) DISSB EN HAND</u>	<u>M36293</u>
	<u>N₂O CS hose 3 m (9.8 ft) DISSB EN HAND</u>	<u>M36294</u>
	<u>N₂O CS hose 5 m (16.4 ft) DISSB EN HAND</u>	<u>M36295</u>
E	<u>Hose ASM-WAGD HAND DISS N 10L</u>	<u>MP20029</u>
	<u>Hose ASM-WAGD HAND DISS N 15L</u>	<u>MP20030</u>



ACCESSORY FOR GAS SUPPLY

3.1.11 CHINESE STANDARD

	<u>O₂ CS hose 3 m (9.8 ft) NIST GB BS</u>	<u>M36732</u>
	<u>O₂ CS hose 5 m (16.4 ft) NIST GB BS</u>	<u>M36733</u>
	<u>AIR CS hose 3 m (9.8 ft) NIST GB BS</u>	<u>M36736</u>
	<u>AIR CS hose 5 m (16.4 ft) NIST GB BS</u>	<u>M36737</u>
	<u>N₂O CS hose 3 m (9.8 ft) NIST GB BS</u>	<u>M36734</u>
	<u>N₂O CS hose 5 m (16.4 ft) NIST GB BS</u>	<u>M36735</u>
	<u>O₂ CS hose 0.9 m (2.95 ft) NIST GB DIN</u>	<u>M36784</u>
	<u>O₂ CS hose 1.5 m (4.9 ft) NIST GB DIN</u>	<u>M36785</u>
	<u>O₂ CS hose 3 m (9.8 ft) NIST GB DIN</u>	<u>M36738</u>
A	<u>O₂ CS hose 5 m (16.4 ft) NIST GB DIN</u>	<u>M36739</u>
	<u>N₂O CS hose 0.9 m (2.95 ft) NIST GB DIN</u>	<u>M36788</u>
	<u>N₂O CS hose 1.5 m (4.9 ft) NIST GB DIN</u>	<u>M36789</u>
	<u>N₂O CS hose 3 m (9.8 ft) NIST GB DIN</u>	<u>M36740</u>
B	<u>N₂O CS hose 5 m (16.4 ft) NIST GB DIN</u>	<u>M36741</u>
	<u>AIR CS hose 0.9 m (2.95 ft) NIST GB DIN</u>	<u>M36786</u>
	<u>AIR CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN</u>	<u>M36010</u>
	<u>AIR CS hose 3 m (9.8 ft) NIST NEUTRAL DIN</u>	<u>M36011</u>
	<u>AIR CS hose 5 m (16.4 ft) NIST NEUTRAL DIN</u>	<u>M36012</u>



ACCESSORY FOR GAS SUPPLY

3.2 ADAPTER

A	Adapter DISS M/QC-N₂O-Chemetron*	4119065
	Adapter DISS M/QC-AIR-Chemetron*	4119066
	Adapter DISS M/QC-VAC-Chemetron*	4119067
	Adapter DISS M/QC-O₂-Chemetron*	4119068
	Adapter DISS M/QC-EVAC-Chemetron*	4119069
B	Adapter DISS M/QC-N₂O-Ohmeda*	4119070
	Adapter DISS M/QC-AIR-Ohmeda*	4119071
	Adapter DISS M/QC-VAC-Ohmeda*	4119072
	Adapter DISS M/QC-O₂-Ohmeda*	4119073
	Adapter DISS M/QC-EVAC-Ohmeda*	4119074
C	Adapter DISS M/QC-N₂O-Pur Ben*	4119075
	Adapter DISS M/QC-AIR-Pur Ben*	4119076
	Adapter DISS M/QC-VAC-Pur Ben*	4119077
	Adapter DISS M/QC-O₂-Pur Ben*	4119078
	Adapter DISS M/QC-EVAC-Pur Ben*	4119079
	Adapter DISS M/QC-N₂O-Schrader*	4119080
	Adapter DISS M/QC-AIR-Schrader*	4119081
	Adapter DISS M/QC-VAC-Schrader*	4119082
	Adapter DISS M/QC-O₂-Schrader*	4119083
D	Double-Plug coupling AIR	G01166
	Double-Plug coupling O₂	G01168
	Double-Plug coupling VAC	G01169

*Tradeware



ACCESSORY FOR GAS SUPPLY

A	Adapter O₂ NIST/DIN	M32493
	Adapter N₂O NIST/DIN	M32494
	Adapter AIR NIST/DIN	M32495
	Adapter VAC NIST/DIN	M32496
	Adapter AIR/O₂ NIST/DIN	M32497
B	Adapter O₂, DISS-NIST	M34875
	Adapter N₂O, DISS-NIST	M34876
	Adapter AIR, DISS-NIST	M34877
	Adapter VAC, DISS-NIST	M34878
C	Adapter NIST to DIN CS-Coupling O₂	M35056
	Adapter N₂O NIST/DIN	M35057
	Adapter NIST to DIN CS-Coupling AIR	M35058
	Adapter AIR-O₂/O₂ NIST/DIN	M36042
D	T-piece AIR NIST*	M36056
E	Adapter O₂ NF/NIST*	M35336
	Adapter N₂O NF/NIST*	M35337
	Adapter AIR NF/NIST*	M35338

*Tradeware





Rail accessories and examination lights

Workstations designed with healing in mind

Hospital lighting has a vital impact on the health and well-being of your patients. Our smart solutions equip you with ideally suited lightning, giving you excellent illumination and working conditions.

Sophisticated without being complicated, our Quickstop hinged arms are easy to handle, allowing you more flexibility for improved patient care.

We can help you conceptualise your hospital installations with intelligent workstations integrated with complementing accessory products for safer and more optimised patient care.

**CONTENT RAIL ACCESSORIES
AND EXAMINATION LIGHT**


1. Rail accessories	339
2. Examination light	340
3. Hinged arms	342
4. Device related accessories	343

RAIL ACCESSORIES AND EXAMINATION LIGHT

1. RAIL ACCESSORIES

A	Rail clamp, for attachment to poles of different diameters, rotatable by 90° to left and right, length of standard rail (profile 25 × 10 mm): 20 cm	2M85274
B	Rail clamp lever type, for attachment of devices to 25 × 10 mm standard rail, for compact suction unit and secretion jar set	M26842
C	CS connection hose holder	2M85446
D	Suction hose holder, suitable for various diameters	2M85239
E	Ceiling holder, rail holder for compact suction systems and secretion jar sets	2M85119
Wire baskets*		
The wire baskets are used for storage of disposable catheters or disposable gloves. The catheters can be simply stored in bundles and are then easy to take out one by one. The baskets are chromeplated and are thus easy to clean and can be autoclaved at up to 134 °C. The amount of space for a basket on the rails is around 150 mm (5.9 in).		
F	Basket 150, 15 × 20 × 10 cm, e.g., for disposable gloves	M26146
F	Basket 300, 15 × 28 × 10 cm, e.g., for 300 mm catheters	M26145
F	Basket 600, 15 × 48 × 10 cm, e.g., for 600 mm catheters	M25121
G	Litter bag holder, stainless steel, with base plate and litter bags (100 pcs.)	M24695
	Litter bags, disposable, for litter bag holder, 100 pcs.	M26240

*Tradeware

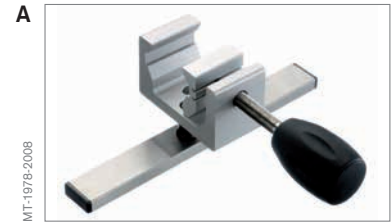

RAIL ACCESSORIES AND EXAMINATION LIGHT

2. EXAMINATION LIGHT

The Dräger **VarioLux®** offers cool, highly efficient, and variable intensity examination light and was specially designed for use in neonatology and pediatric settings. Made for single-handed operation, it features adjustable color temperature settings as well as a dimmer function.



Operating mode:	Continuous operation
Weight:	1.5 kg (3.3 lbs)
Electrical values:	100 – 240 V, 50/60 Hz
Max. illumination intensity E_c at 0.5 m:	30,000 lx
Dimming levels:	30,000 lx, 24,000 lx, 18,000 lx, 12,000 lx, 7,000 lx
Light field diameter d_{10} at 0.5 m:	20 cm (7.9 in)
Light field diameter d_{50} at 0.5 m:	9 cm (3.5 in)
Color temperature adjustable:	Cool white (blue) 4,700 K, Neutral white (white) 4,100 K, Warm white (red) 3,500 K
Min. service life of bulb (LEDs):	> 50,000 h
Color reproduction index:	$R_a \geq 96$ (at 3,500 K) $R_9 = 89$ (at 3,500 K)
Max. radiant power at 0.5 m:	105 W/m ²

A	Medical examination light VarioLux® with LED technology, EU	MP00601
	Medical examination light VarioLux® with LED technology, US	MP00602
	Medical examination light VarioLux® with LED technology, GB	MP00603
	Medical examination light VarioLux® with LED technology, AUS	MP00606
	Medical examination light VarioLux® with LED technology, JPN	MP00608
	Power cable EU, 1.8 m, 2 pin	MP00591
	Power cable US/JPN, 1.8 m, 2 pin	MP00592
	Power cable GB, 1.8 m, 2 pin	MP00593
	Power cable AUS, 1.8 m, 2 pin	MP00626



RAIL ACCESSORIES AND EXAMINATION LIGHT

A	Rail claw for lamp	MP00615
B	Turn knob, compl.	MP00587
C	Control panel	MP00589
D	Rail clamp, for attachment to poles of different diameters, rotatable by 90° to left and right, length of standard rail (profile 25 × 10 mm): 20 cm	2M85274
E	Compact rail, for 38 mm tube, standard rail profile 25 × 10 mm, each side 10 cm, weight: 0.9 kg, max. load: 30 kg	2M85337

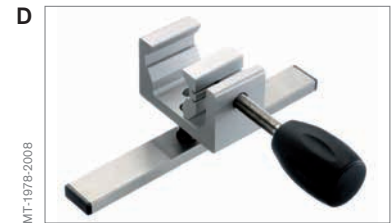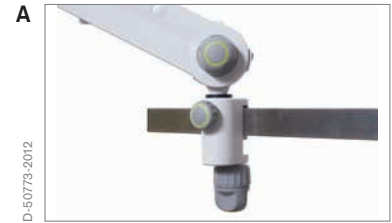

RAIL ACCESSORIES AND EXAMINATION LIGHT

3. HINGED ARMS

Keeping everything close at hand. To leave your hands free to care for the patient, all the joints can be fixed into place and released with just one quick hand movement. All cables, lines and hose systems can be arranged in such a way as to more or less completely prevent accidental disconnection.

- | | | |
|----------|---|----------------|
| A | Hinged arm Infinity® ACS, for positioning and fixation of breathing hoses, this arm can be fixed into place without a hand movement through a patent registered friction torque-breaking system, length 1,080 mm, max. load 0.5 kg, includes hose clamps: MP00693 (1×), 8417396 (2×) | MP00690 |
| B | Hinged arm 240°, 3 hinges, fixable, ball holder system, length 1,155 mm, max. load 0.3 kg, includes hose clamps: MP00493 (1×), 8409841 (1×) | 8409609 |
| C | Hinged arm Quickstop 2, with universal ball holder system for optional accessories length 858 mm, max. load 0.4 kg, includes hose clamps: MP00493 (1×), 8404761 (2×), 8409841 (1×) | 2M85706 |
| D | Hinged arm Quickstop short, with mounted universal tube holder plate, for various tube dimensions/coaxial hoses, specifically designed for OR tables, length 323 mm, max. load 1 kg | 2M86463 |



DEVICE RELATED ACCESSORIES

4. DEVICE RELATED ACCESSORIES

A	Disposable Handle Polaris®, Dual, for Polaris® 100/200, 600 and Multimedia, adapter required, 60 pcs.	MP05700
A	Disposable Handle Polaris®, Single, for Polaris® 100/200, 600 and Multimedia, adapter required, 120 pcs.	MP05701

CONVERSION KITS

	Polaris® 100 with sterilizable handle (G92099), 1 pc.	G35332
	Polaris® 200 with sterilizable handle (G92099), 1 pc.	G35614
	Polaris® 100/200 with sterilizable handle E (G36993), 1 pc.	G35325
	Polaris® 600 with sterilizable handles (G92099 or G36993), 1 pc.	G35324
	Polaris® Multimedia (display adapter) with sterilizable handles (G92099 or G36993), 1 pc.	G35333





Neonatal care and thermoregulation

Improving care for your smallest of patients

The majority of all neonatal deaths (75%) occurs during the first week of life¹⁾. This sobering reality drives us to develop gentle and efficient solutions for improving acute care in neonatal intensive care units. Together with you, our goal is to lower morbidity and mortality rates while reducing stress for newborns and clinical staff.

We understand that “small” doesn’t mean less important. And that’s why we give our accessories and consumables the highest level of attention – so that when you need to give your tiniest patients the best possible care, you can rely on our smallest products for:

- Thermoregulation,
- neonatal respiratory support,
- vital sign monitoring,
- jaundice treatment,
- and so much more.

To help you do just that – save lives.

¹⁾<https://www.who.int/news-room/fact-sheets/detail/newborns-reducing-mortality>

CONTENT NEONATAL CARE AND THERMOREGULATION

1. Warming therapy	347
Babyleo® TN500	347
Babyroo® TN300	349
Caleo®	351
Isolette® 8000/8000 plus	353
Air-Shields® Isolette® C2000/8000	355
Babytherm® 8000 WB/8004/8010/4200	357
Air-Shields® Resuscitaire®	360
Globe-Trotter GT5400	363
Air-Shields® TI 500 Globe-Trotter	366
Air-Shields® Isolette® TI 500	369
Photo Therapy	371
Developmental Care	372
Miscellaneous	373
2. Ventilation accessories	374
3. Sensors	400
4. Monitoring accessories	410
5. Ward equipment	418
6. Rail accessories and examination light	433

BABYLEO® TN500

1. WARMING THERAPY

BABYLEO® TN500

A	Skin temperature probe ThermoTrace, disposable, yellow, central, 5 pcs.	MX11000
B	Skin temperature probe ThermoTrace, disposable, white, peripheral, 5 pcs.	MX11001
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
D	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
E	Mattress SoftBed, for Babyleo® TN500, reusable	2M60390
	Heating mattress, 1 pc.	2M60170
F	Incubator cover for Babyleo® TN500, reusable	MP03730
G	Cuddling toy Babyleo® TN500*, lion, orange	MP03737

*Tradeware



Not all articles are available worldwide.

BABYLEO® TN500

A	O₂ sensor Oxycell, for Babyleo® TN500	6873065
B	Ventilation hose holder, reusable, 40 cm	8411075
C	Air intake filter for Caleo® and Babyleo® TN500, reusable, 20 pcs.	MX17015
D	Transition Set, disposable, 20 pcs.	MX17018
E	Gastroschisis arm, complete	2M60180
F	Holder for cylinder packet for GCX profile	4118453
G	Rail, 100 mm for GCX profile	4118555
	Clamp for GCX rail, 5 pcs.	2M60401



BABYROO® TN300**BABYROO® TN300**

A	Skin temperature probe ThermoTrace, disposable, yellow, central, 5 pcs.	MX11000
B	Skin temperature probe ThermoTrace, disposable, white, peripheral, 5 pcs.	MX11001
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
D	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
E	Mattress SoftBed, for BabyRoo TN300, reusable	MP16297
F	Gel Mattress for BabyTherm 8000/8010 and Babyroo TN300, reusable	2M20827
G	Cuddling toy BabyRoo, kangaroo	MP16090



BABYROO® TN300

A	VentStar Autobreath N for Resuscitaire® and TN300, N, disposable, 25 pcs.	MU10841
B	VentStar® Resus Neo, disposable, N, 1 m, 25 pcs., not made with natural rubber latex	MP00310
C	VentStar® Resus Neo with pressure measurement line, disposable, N, 1 m, 25 pcs., not made with natural rubber latex	MP00311
	F&P configuration for MR850 and F&P950	8414144
D	Suction bottle with filter and tube, 800 ml (27.05 fl. oz.), 20 pcs.	MU10918
	Suction canister holder for GCX	3702751



CALEO®

CALEO®

A	Skin temperature probe ThermoTrace, disposable, yellow, central, 5 pcs.	MX11000
B	Skin temperature probe ThermoTrace, disposable, white, peripheral, 5 pcs.	MX11001
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
D	Mattress SoftBed, for Caleo®, reusable	MX17012
E	Incubator cover for Caleo®, reusable	MP01423
F	Tubing grommet for Caleo®, reusable, angle	2M50385
G	Tubing grommet for Caleo®, reusable, flap	2M50412



Not all articles are available worldwide.

CALEO®

A	Ventilation hose holder, reusable, 40 cm	8411075
B	O₂ sensor OxyTrace INCU, for Caleo®	MX01050
C	Air intake filter for Caleo® and Babyleo® TN500, reusable, 20 pcs.	MX17015
D	Transition Set, disposable, 20 pcs.	MX17018
E	Water tank and water transfer hose, for Caleo®, reusable	2M50040
F	Infusion holder, reusable, 3 hooks, for poles with Ø 38 mm	2M21514
G	Swivel table, reusable, for poles with Ø 38 mm	2M21186




ISOLETTE® 8000/8000 PLUS**ISOLETTE® 8000/8000 PLUS**

A	Skin temperature probe Premi-Probe® for Isolette® 8000 plus, disposable, yellow, central, 10 pcs.	MU26041
B	Skin temperature probe Premi-Probe® for Isolette® 8000 plus, disposable, white, peripheral, 10 pcs.	MU26042
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
D	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
E	Mattress SoftBed, for Isolette® C2000, 8000 and 8000 plus, reusable	MP01401
F	Incubator cover for Isolette® C2000, 8000 and 8000 plus, reusable	MP01424



ISOLETTE® 8000/8000 PLUS

A	Tubing grommet for Isolette® C2000, 8000 and 8000 plus, reusable	MU12609
B	O₂ sensor for Isolette® C2000, 8000 and 8000 plus, reusable, 2 pcs.	MU24903
C	Air intake filter for Isolette® C2000, 8000 and 8000 plus, reusable, 4 pcs.	MU12504
D	Condensate management hoses for Isolette® 8000 and 8000 plus, disposable, 20 pcs.	MU21120
E	Suction bottle with filter and tube for Isolette® 8000, disposable, 800 ml, 20 pcs.	MU10918
F	Iris port sleeve, disposable, 100 pcs.	MU03876

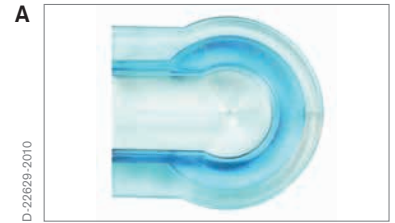

AIR-SHIELDS® ISOLETTE® C2000/8000**AIR-SHIELDS® ISOLETTE® C2000/8000**

A	Skin temperature probe Premi-Probe® 5, reusable	MU12533
B	Skin temperature probe Premi-Probe® 5, disposable, 10 pcs.	MU12551
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
D	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
E	Mattress SoftBed, for Isolette® C2000, 8000 and 8000 plus, reusable	MP01401
F	Incubator cover for Isolette® C2000, 8000 and 8000 plus, reusable	MP01424
G	Cuddling toy Isi*, deer, green	MP03738


*Tradeware



Not all articles are available worldwide.

AIR-SHIELDS® ISOLETTE® C2000/8000

A	Tubing grommet for Isolette® C2000, 8000 and 8000 plus, reusable	MU12609
B	O₂ sensor for Isolette® C2000, 8000 and 8000 plus, reusable, 2 pcs.	MU24903
C	Air intake filter for Isolette® C2000, 8000 and 8000 plus, reusable, 4 pcs.	MU12504
	Humidifier chamber for Isolette® C2000, reusable	MU17191
D	Condensate management hoses for Isolette® 8000 and 8000 plus, disposable, 20 pcs.	MU21120
E	Suction bottle with filter and tube for Isolette® 8000, disposable, 800 ml, 20 pcs.	MU10918
F	Iris port sleeve, disposable, 100 pcs.	MU03876
G	Hand port cuff, disposable, oval, 100 pcs.	MU06570

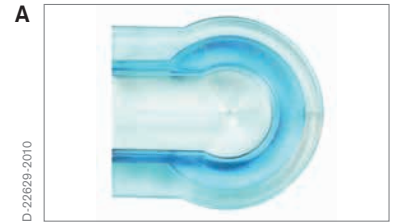

BABYTHERM® 8000 WB/8004/8010/4200**BABYTHERM® 8000 WB/8004/8010/4200**

A	Skin temperature probe ThermoTrace, disposable, yellow, central, 5 pcs.	MX11000
B	Skin temperature probe ThermoTrace, disposable, white, peripheral, 5 pcs.	MX11001
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
D	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
E	Gel Mattress for BabyTherm 8000/8010 and Babyroo TN300, reusable	2M20827
F	Mattress SoftBed, for Babytherm®, reusable	2M21012
F	Mattress SoftBed, for BabyRoo TN300, reusable	MP16297



BABYTHERM® 8000 WB/8004/8010/4200

A	Holder for hinged hood for Babytherm® 8000 WB, 8004 and 8010, reusable, neonatal, mounts to hand rail	2M21342
B	Deposit plate for Babytherm® 8000	2M21158
C	Handle, complete for Babytherm® 8000	2M21161
D	Breathing hose holder for Babytherm®, neonatal, for 10 mm breathing hose	2M21191
	Tubing grommet, reusable	2M20434



BABYTHERM® 8000 WB/8004/8010/4200

A <u>Infusion holder, reusable, 3 hooks, for poles with Ø 38 mm</u>	<u>2M21514</u>
B <u>Infusion bottle holder for Babytherm®, reusable, neonatal</u>	<u>2M16520</u>
C <u>Swivel table, reusable, for poles with Ø 38 mm</u>	<u>2M21186</u>



AIR-SHIELDS® RESUSCITAIRE®

AIR-SHIELDS® RESUSCITAIRE®

A	Skin temperature probe Premi-Probe® 3, reusable	MU11225
B	Skin temperature probe Premi-Probe® 3, disposable, 10 pcs.	MU11228
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
D	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
E	Mattress SoftBed, for Resuscitaire®, reusable, with scale	MP01404
E	Mattress SoftBed, for Resuscitaire®, reusable, without scale	MP01402



AIR-SHIELDS® RESUSCITAIRE®

A	VentStar Autobreath N for Resuscitaire® and TN300, N, disposable, 25 pcs.	MU10841
B	Breathing hose for breathing bag, disposable, 15 mm, 25 pcs.	MU10900
C	VentStar® Resus Neo, disposable, N, 1 m, 25 pcs., not made with natural rubber latex	MP00310
D	VentStar® Resus Neo with pressure measurement line, disposable, N, 1 m, 25 pcs., not made with natural rubber latex	MP00311
E	T-piece adapter for MX 300/MX 300i and Resuscitaire® breathing circuit, disposable, 20 pcs.	MP01880
F	Silicone mask, reusable, size 00, premature baby	MP02900
F	Silicone mask, reusable, size 0, small neonate	MP01910
F	Silicone mask, reusable, size 1, large neonate	MP01911



AIR-SHIELDS® RESUSCITAIRE®

A	T-connector, disposable, Ø 15 mm, 25 pcs.	MU10902
B	Suction bottle with holder for Resuscitaire, disposable, 40 ml	MU10916
C	Suction bottle, disposable, 40 ml, 10 pcs.	MU10917
D	Filter for suction bottle 40 ml for Resuscitaire®, disposable, 10 pcs.	MU10915
E	Suction bottle with filter and tube, disposable, 800 ml, 20 pcs.	MU10918



GLOBE-TROTTER GT5400**GLOBE-TROTTER GT5400**

A	Skin temperature probe Premi-Probe® 1, disposable, 10 pcs.	MU06933
B	Skin temperature probe Premi-Probe® 1, reusable	MU06951
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
D	Mattress SoftBed, for TI 500 and TI 500 Globe-Trotter, reusable	MP01419
E	Restraining strap, disposable, 10 pcs.	MU06011
	Tubing grommet for TI 500 and TI 500 Globe-Trotter, reusable	MU04366
F	Iris port sleeve, disposable, 100 pcs.	MU03876
G	Hand port cuff, disposable, oval, 100 pcs.	MU06570




Not all articles are available worldwide.

GLOBE-TROTTER GT5400

A <u>Air intake filter for TI 500 and TI 500 Globe-Trotter, disposable, 6 pcs.</u>	<u>MU04472</u>
B <u>Humidity pad, disposable, antibacterial, 25 pcs.</u>	<u>MU06008</u>
C <u>O₂ sensor* for Crossvent 2i+</u>	<u>MU24402</u>
D <u>Flow sensor* for Crossvent 2i+, disposable, infant, 10 pcs.</u>	<u>MU24401</u>
E <u>Breathing circuit* for active humidification, disposable, neonate, 20 pcs., not made with natural rubber latex</u>	<u>MU24409</u>

*Tradeware

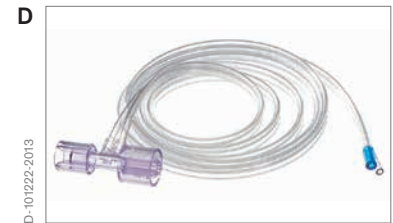

GLOBE-TROTTER GT5400

- | | |
|---|----------------|
| A SpO ₂ adhesive sensor* Masimo® LNCS Neo, disposable, neonate, 20 pcs. | MP00792 |
| B SpO ₂ adhesive sensor* Masimo® LNCS NeoPt, disposable, premature baby, 20 pcs. | MP00794 |

*Tradeware



AIR-SHIELDS® TI 500 GLOBE-TROTTER

AIR-SHIELDS® TI 500 GLOBE-TROTTER

A	Skin temperature probe Premi-Probe® 1, disposable, 10 pcs.	MU06933
B	Skin temperature probe Premi-Probe® 1, reusable	MU06951
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
D	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
E	Mattress SoftBed, for TI 500 and TI 500 Globe-Trotter, reusable	MP01419
	Tubing grommet for TI 500 and TI 500 Globe-Trotter, reusable	MU04366
	Breathing hose for TI 500 Globe-Trotter, disposable, 20 pcs.	MU04771



AIR-SHIELDS® TI 500 GLOBE-TROTTER

A Air intake filter for TI 500 and TI 500 Globe-Trotter, disposable, 6 pcs.	MU04472
B Iris port sleeve, disposable, 100 pcs.	MU03876
C Hand port cuff, disposable, oval, 100 pcs.	MU06570



AIR-SHIELDS® TI 500 GLOBE-TROTTER

- A Humidity pad, disposable, antibacterial, 25 pcs. MU06008
- B Restraining strap, disposable, 10 pcs. MU06011



AIR-SHIELDS® ISOLETTE® TI 500**AIR-SHIELDS® ISOLETTE® TI 500**

A	Skin temperature probe Premi-Probe® 1, disposable, 10 pcs.	MU06933
B	Skin temperature probe Premi-Probe® 1, reusable	MU06951
C	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944
	Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943
D	Mattress SoftBed, for TI 500 and TI 500 Globe-Trotter, reusable	MP01419
	Tubing grommet for TI 500 and TI 500 Globe-Trotter, reusable	MU04366
E	Air intake filter for TI 500 and TI 500 Globe-Trotter, disposable, 6 pcs.	MU04472
F	Iris port sleeve, disposable, 100 pcs.	MU03876
G	Hand port cuff, disposable, oval, 100 pcs.	MU06570



AIR-SHIELDS® ISOLETTE® TI 500

- A Humidity pad, disposable, antibacterial, 25 pcs. MU06008
- B Restraining strap, disposable, 10 pcs. MU06011



PHOTO THERAPY

PHOTO THERAPY

A	Light shield curtain for BiliLux, reusable	MP03760
B	Phototherapy mask Eyemax 2*, disposable, micro, head circumference 20-25 cm, 20 pcs.	MP03770
C	Phototherapy mask Eyemax 2*, disposable, premie, head circumference 26-32 cm, 20 pcs.	MP03771
D	Phototherapy mask Eyemax 2*, disposable, regular, head circumference 33-38 cm, 20 pcs.	MP03772

*Tradeware



DEVELOPMENTAL CARE

DEVELOPMENTAL CARE

A	Cuddling toy Cally*, reusable, light blue or pink	2M30462
B	Cuddling toy Babyleo® TN500*, lion, orange	MP03737
C	Cuddling toy Isi*, deer, green	MP03738
D	Cuddling toy BabyRoo, kangaroo	MP16090
E	Incubator cover for Babyleo® TN500, reusable	MP03730
F	Incubator cover for Caleo®, reusable	MP01423
G	Incubator cover for Isolette® C2000, 8000 and 8000 plus, reusable	MP01424

*Tradeware



MISCELLANEOUS

MISCELLANEOUS

<p>A Warming mattress TransWarmer®*, disposable, self warming for approx. 2 h, 6 pcs.</p>	<p>MP01414</p>
<p>B Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.</p>	<p>MU06944</p>
<p>C Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.</p>	<p>MU06943</p>

*Tradeware



VENTILATION ACCESSORIES

2. VENTILATION ACCESSORIES

OXYMIXER*

Oxymixer enables a continuous and individually adjusted mixture of medical air and medical oxygen to be administered to infants, children and adults. The O₂ concentration is easily adjusted to a highly accurate flow – even with flow rates of less than 3 l/min.

Oxymixer configuration

MP04200

Free choice of Low- or High Flow Blender with:
1 or 2 flowmeter, NIST or DISS gas inlet,
Monitor, Bleed swith, Nozzle, Sound damper.
The rail fixation is included.

A Oxymixer Low Flow

MP04203

Including: Rail clamp, Flowmeter 0 – 3 l/min right,
Flowmeter 0 – 15 l/min left, Bleed on/off switch.

B Oxymixer Low Flow with monitor

MP04204

Including: Rail clamp, Flowmeter 0 – 3 l/min right,
Flowmeter 0 – 15 l/min left, Bleed on/off switch,
Monitor, Oxygen sensor, Oxygen sensor cable,
Sensor receptacle.

C Oxymixer Low Flow with 1 – 15 l flowmeter right, rail clamp

MP04229

Several hoses to connect Oxymixer with the central supply.
For CS hoses please refer to chapter "Accessory for gas supply".

*Tradeware



VENTILATION ACCESSORIES

OXYMIXER*

A	Oxymixer sound damper 9/16"	MP04218
B	Oxymixer sound damper 22 and 15 mm	MP04219
C	Adapter, with 6 mm nozzle, e.g., for disposable breathing masks, 9/16" connection	MP00416
D	Oxymixer oxygen sensor (spare part)	MP04221

*Tradeware



VENTILATION ACCESSORIES

SEATTLE PAP PLUS – BUBBLE CPAP SYSTEM

Helping infants with respiratory distress breathe easier. The Seattle-Positive Airway Pressure (PAP) system is an innovation, which uses the proven advantages of Bubble CPAP therapy, such as oscillatory effects similar to high frequency ventilation and combines it with a unique design.

Seattle-PAP plus is safe, effective and easy to use.

It was designed to:

- Promote infant safety with a uniquely designed pressure manifold with pressure relief valve.
- Safety tube lock to avoid unintended changes of the PEEP level.
- Provide comfort and optimum CPAP support to the infant due to disposable BabyFlow® plus interfaces complementing the system.
- Include funnel with the filling tube which minimises spillages.
- The swivel connector to the expiratory limb prevents unnecessary disconnections and pressure loss.
- The water trap in the expiratory limb prevents unwanted pressure peaks due to condensation in the expiratory limb.

It requires only oxygen and compressed air, both to be administered by an O₂ blender like the Oxymixer from Dräger and a humidification source for the circuit.



A	Seattle PAP plus	MP02260
	Trolley	5790148
B	Multi-P trolley* equipped	MP13020
	Including base with 5 castors, 38 mm Ø pole, 3 rail clamps, basket.	
C	Multi-P trolley*	MP13021
	Including base with 5 castors and 38 mm Ø pole for accessories.	

*Tradeware

VENTILATION ACCESSORIES

A Compact rail, for 38 mm tube, standard rail profile 25 × 10 mm, each side 10 cm, weight: 0.9 kg, max. load: 30 kg	2M85337
B Universal holder profile 30 × 5 for wall rails 25 × 8 to 35 × 10 mm	MP00634
C Adapter, with 6 mm nozzle, e.g., for disposable breathing masks, 9/16" connection	MP00416
D Rail clamp for F&P humidifier and Seattle PAP bubble box*	8411074
E Rail clamp, for attachment to poles of different diameters, rotatable by 90° to left and right, length of standard rail (profile 25 × 10 mm): 20 cm	2M85274



*Tradeware



VENTILATION ACCESSORIES

ACCESSORIES FOR SEATTLE PAP PLUS


A	BabyFlow® plus, disposable, size S, 20 pcs.	MP03700
	BabyFlow® plus, disposable, size M, 20 pcs.	MP03701
	BabyFlow® plus, disposable, size L, 20 pcs.	MP03702
B	Demo package BabyFlow® plus, disposable, size S	MP03703
	Demo package BabyFlow® plus, disposable, size M	MP03704
C	Soft connector N, disposable, 20 pcs.	MP03826
D	① Prong BabyFlow®, standard, disposable, size XS, 10 pcs.	8418415
	② Prong BabyFlow®, standard, disposable, size S, 10 pcs.	8418605
	③ Prong BabyFlow®, standard, disposable, size M, 10 pcs.	8418416
	④ Prong BabyFlow®, standard, disposable, size L, 10 pcs.	8418531
	⑤ Prong BabyFlow®, standard, disposable, size XL, 10 pcs.	8418417
E	Prong BabyFlow®, standard, disposable, size XXL, 10 pcs.	8418617
F	① Prong BabyFlow®, wide, disposable, size 3, 10 pcs.	8418603
	② Prong BabyFlow®, wide, disposable, size 4, 10 pcs.	8418607
	③ Prong BabyFlow®, wide, disposable, size 5, 10 pcs.	8418609
	④ Prong BabyFlow®, wide, disposable, size 6, 10 pcs.	8418611

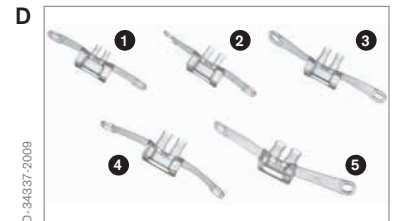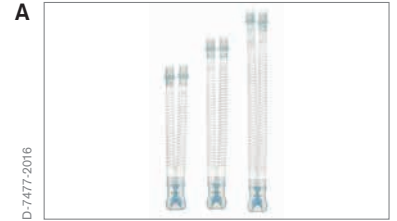

VENTILATION ACCESSORIES

A	① Neo-mask, silicone, disposable, size S, 10 pcs.	8418491
	② Neo-mask, silicone, disposable, size M, 10 pcs.	8418490
	③ Neo-mask, silicone, disposable, size L, 10 pcs.	8418619
B	Neo-mask, TPE, disposable, size 4, 10 pcs.	MP01494
C	Neo-mask, TPE, disposable, size 5, 10 pcs.	MP01495
D	Neo-mask, TPE, disposable, size 6, 10 pcs.	MP01496
E	Headgear, disposable, size XXS, 17 – 20 cm (6.7 – 7.9 in) head circumference, 5 pcs.	MP03880
	Headgear, disposable, size XS, 20 – 23 cm (7.9 – 9.1 in) head circumference, 5 pcs.	MP03881
F	Headgear, disposable, size S, 23 – 27 cm (9.1 – 10.6 in) head circumference, 5 pcs.	MP03882
	Headgear, disposable, size M, 27 – 31 cm (10.6 – 12.2 in) head circumference, 5 pcs.	MP03883
	Headgear, disposable, size L, 31 – 37 cm (12.2 – 14.6 in) head circumference, 5 pcs.	MP03884
	Headgear, disposable, size XL, 37 – 42 cm (14.6 – 16.5 in) head circumference, 5 pcs.	MP03885
	Cap BabyFlow®, disposable, size XS, white, 18 – 20.5 cm (7.1 – 8.1 in) head circumference, 5 pcs.	8418533
	Cap BabyFlow®, disposable, size S, yellow, 20.5 – 23 cm (8.1 – 9.1 in) head circumference, 5 pcs.	8418534
	Cap BabyFlow®, disposable, size M, red, 23 – 25.5 cm (9.1 – 10 in) head circumference, 5 pcs.	8418535
	Cap BabyFlow®, disposable, size L, green, 25.5 – 28 cm (10 – 11 in) head circumference, 5 pcs.	8418536
	Cap BabyFlow®, disposable, size XL, orange, 28 – 30.5 cm (11 – 12 in) head circumference, 5 pcs.	8418537
	Cap BabyFlow®, disposable, size XXL, blue, 30.5 – 33 cm (12 – 13 in) head circumference, 5 pcs.	8418538
G	Cap BabyFlow®, disposable, size XXL+, white/black, 33 – 36 cm (13 – 14.2 in) head circumference, 5 pcs.	8418539

Not all articles are available worldwide.

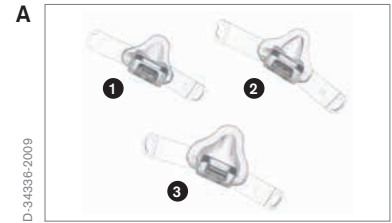


VENTILATION ACCESSORIES

BABYFLOW® PLUS

BabyFlow® plus – The superior non-invasive respiratory support system, designed to maximise performance and patient comfort.

A	BabyFlow® plus, disposable, size S, 20 pcs.	MP03700
	BabyFlow® plus, disposable, size M, 20 pcs.	MP03701
	BabyFlow® plus, disposable, size L, 20 pcs.	MP03702
B	Demo package BabyFlow® plus, disposable, size S	MP03703
	Demo package BabyFlow® plus, disposable, size M	MP03704
	Demo package BabyFlow® plus, disposable, size L	MP03705
C	Soft connector N, disposable, 20 pcs.	MP03826

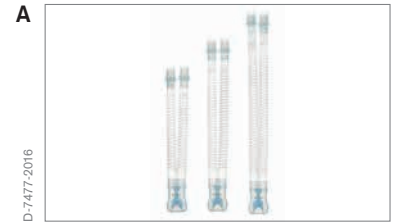

VENTILATION ACCESSORIES



STANDARD PRONGS

Soft anatomically formed nasale silicone prongs, disposable.

A	① Prong BabyFlow®, standard, disposable, size XS, 10 pcs.	8418415
	② Prong BabyFlow®, standard, disposable, size S, 10 pcs.	8418605
	③ Prong BabyFlow®, standard, disposable, size M, 10 pcs.	8418416
	④ Prong BabyFlow®, standard, disposable, size L, 10 pcs.	8418531
	⑤ Prong BabyFlow®, standard, disposable, size XL, 10 pcs.	8418417
B	Prong BabyFlow®, standard, disposable, size XXL, 10 pcs.	8418617

WIDE PRONGS

C	① Prong BabyFlow®, wide, disposable, size 3, 10 pcs.	8418603
	② Prong BabyFlow®, wide, disposable, size 4, 10 pcs.	8418607
	③ Prong BabyFlow®, wide, disposable, size 5, 10 pcs.	8418609
	④ Prong BabyFlow®, wide, disposable, size 6, 10 pcs.	8418611

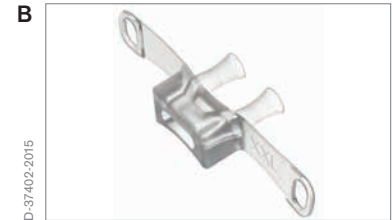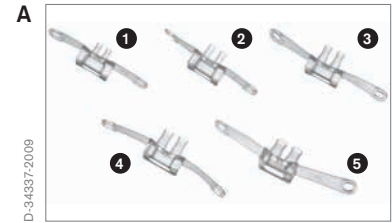


VENTILATION ACCESSORIES

MASKS

Nasal masks of silicone and TPE complete the portfolio and increase the comfort for the patient!

A	① Neo-mask, silicone, disposable, size S, 10 pcs.	8418491
	② Neo-mask, silicone, disposable, size M, 10 pcs.	8418490
	③ Neo-mask, silicone, disposable, size L, 10 pcs.	8418619
B	Neo-mask, TPE, disposable, size 4, 10 pcs.	MP01494
C	Neo-mask, TPE, disposable, size 5, 10 pcs.	MP01495
D	Neo-mask, TPE, disposable, size 6, 10 pcs.	MP01496

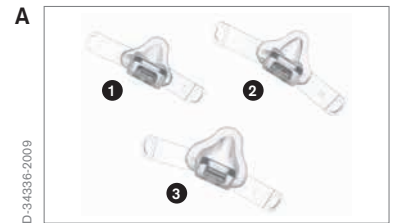

VENTILATION ACCESSORIES

HEADGEARS

The new headgear enable the customer to choose the most comfortable fit and combination for each of their tiny patients.

A	Headgear, disposable, size XXS, 17 – 20 cm (6.7 – 7.9 in) head circumference, 5 pcs.	MP03880
	Headgear, disposable, size XS, 20 – 23 cm (7.9 – 9.1 in) head circumference, 5 pcs.	MP03881
B	Headgear, disposable, size S, 23 – 27 cm (9.1 – 10.6 in) head circumference, 5 pcs.	MP03882
	Headgear, disposable, size M, 27 – 31 cm (10.6 – 12.2 in) head circumference, 5 pcs.	MP03883
C	Headgear, disposable, size L, 31 – 37 cm (12.2 – 14.6 in) head circumference, 5 pcs.	MP03884
	Headgear, disposable, size XL, 37 – 42 cm (14.6 – 16.5 in) head circumference, 5 pcs.	MP03885



VENTILATION ACCESSORIES

CAPS

The disposable and breathable microfiber caps (95% cotton, 5% spandex) are available in seven sizes (color-coded).

A	① Cap BabyFlow [®] , disposable, size XS, white, 18 – 20.5 cm (7.1 – 8.1 in) head circumference, 5 pcs.	8418533
	Cap BabyFlow [®] , disposable, size S, yellow, 20.5 – 23 cm (8.1 – 9.1 in) head circumference, 5 pcs.	8418534
	Cap BabyFlow [®] , disposable, size M, red, 23 – 25.5 cm (9.1 – 10 in) head circumference, 5 pcs.	8418535
	② Cap BabyFlow [®] , disposable, size L, green, 25.5 – 28 cm (10 – 11 in) head circumference, 5 pcs.	8418536
	③ Cap BabyFlow [®] , disposable, size XL, orange, 28 – 30.5 cm (11 – 12 in) head circumference, 5 pcs.	8418537
	Cap BabyFlow [®] , disposable, size XXL, blue, 30.5 – 33 cm (12 – 13 in) head circumference, 5 pcs.	8418538
B	Cap BabyFlow [®] , disposable, size XXL+, white/black, 33 – 36 cm (13 – 14.2 in) head circumference, 5 pcs.	8418539



VENTILATION ACCESSORIES

The **MiniMe2*** nasal mask is designed for children ages 2 – 12. It features an ergonomic silicone-gel cushion and a shapeable shell mask to fit individual needs and thus ensuring an effective seal and a comfortable fit.

A	Pediatric nasal mask MiniMe2 non-vented, disposable, size S	MP01555
	Pediatric nasal mask MiniMe2 non-vented, disposable, size L	MP01556
B	Pediatric nasal mask MiniMe2 vented, disposable, size S	MP01553
	Pediatric nasal mask MiniMe2 vented, disposable, size L	MP01554

*Tradeware



VENTILATION ACCESSORIES

A	VentStar® Resus Neo, disposable, N, 1 m, 25 pcs., not made with natural rubber latex	MP00310
B	VentStar® Resus Neo with pressure measurement line, disposable, N, 1 m, 25 pcs., not made with natural rubber latex	MP00311
C	Breathing circuit VentStar®, disposable, N, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP00353
D	Breathing circuit VentStar®, disposable, N, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	MP00363
	Breathing circuit VentStar®, disposable, N, MRI, 3 m, 10 pcs., not made with natural rubber latex	MP04802
E	Breathing bag set VentStar®, disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP00383
	Breathing bag set VentStar® (Cone 22OD), disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP16021
F	Anesthesia set VentStar®, disposable, N, 1.8 m/1.1 m, 25 pcs., not made with natural rubber latex	MP00333
	Anesthesia set VentStar®, disposable, N, 1.8 m/1.1 m, Bag cone 22OD, 25 pcs., not made with natural rubber latex	MP16032
G	VentStar® Breathing Bag 0.5 L, 25 pcs.	MP02710



VENTILATION ACCESSORIES

A Breathing circuit Infinity® ID, disposable, N, basic, 1.8 m, 25 pcs., not made with natural rubber latex	<u>MP01327</u>
Breathing circuit Infinity® ID, disposable, N, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	<u>MP01328</u>
B Breathing bag set Infinity® ID, disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex	<u>MP01329</u>



VENTILATION ACCESSORIES

A Breathing circuit VentStar®, disposable, N, heated, 1.2 m, 10 pcs., not made with natural rubber latex	MP00308
B Breathing circuit VentStar®, disposable, N, heated, w/o chamber, 1.2 m, 10 pcs., not made with natural rubber latex	MP00309
C Breathing circuit Infinity® ID, disposable, N, heated, 1.2 m, 10 pcs., not made with natural rubber latex	MP00343
D Breathing circuit VentStar® Helix, disposable, N, dual heated, 1.7 m, 10 pcs., not made with natural rubber latex	MP02650
E Breathing circuit VentStar® Helix, disposable, N, heated, 1.7 m, 10 pcs., not made with natural rubber latex	MP02608
VentStar® Aquapor, dual heated, N plus, 10 pcs.	MP07150
VentStar® Aquapor, heated, N plus, 10 pcs.	MP07151
F T-adapter neonatal for NO therapy, 10 pcs.	8422215



VENTILATION ACCESSORIES

A	Humidifier chamber VentStar®, disposable, 10 pcs.	MP00234
B	Humidifier chamber VentStar®, disposable, N, 10 pcs.	MP02590
C	Double connector 11ID/15OD, disposable, 10 pcs.	MP02655
ACCESSORIES FOR F&P 950*		
D	F&P 950 sensor cartridge	MP16002
E	F&P 950 heater wire adaptor	MP16003
F	F&P 950 A80 breathing circuit (A)	MP16004
G	F&P 950 neonatal disposable breathing circuit	MP16005

*Tradeware



Not all articles are available worldwide.

VENTILATION ACCESSORIES

FISHER & PAYKEL HUMIDIFIER AND RELATED ACCESSORIES

A	Humidifier chamber* for F&P MR 290, disposable, 10 pcs.	8418282
B	Humidifier chamber* for F&P MR 340, reusable, N/P	8411047
C	Rail clamp* for F&P humidifier and Seattle PAP bubble box	8411074
D	Holder* for F&P humidifier U-shaped, for trolley, with rail clamp	MP01002
E	Water bag holder* 900 MR 290	MP01008
	Holder, lateral (F&P 900 MR 309), including 8411074	MP01009
	Universal holder with DIN rail	G93140
	Humidifier holder for side rail	8416325
	Humidifier holder chassis, swiveling	G93111
F	Hose heater* for F&P MR 700/730 and MR 850, for reusable breathing circuit, 0.7 m	8411154
G	Hose heater* for F&P MR 700/730 and MR 850, for reusable breathing circuit, 1.1 m	8414993

*Tradeware



VENTILATION ACCESSORIES

	Temperature probe/flow sensor* for F&P MR 850 for HFV reusable breathing circuits, 0.7 m	8418507
A	Temperature probe/flow sensor* for F&P MR 850 for all breathing circuits, 1.1 m	8414989
A	Temperature probe/flow sensor* for F&P MR 850, for all breathing circuits, 1.5 m	8414966
B	Hose heater adapter* for F&P MR 850, for reusable breathing circuit*	8414968
C	Hose heater adapter for F&P MR 850, for dual-heated disposable breathing circuit	MP01006
D	Hose heater adapter for F&P MR 850, for single-heated disposable breathing circuit	8414992
E	Draw wire* for hose heater 1.5 m, for reusable breathing circuit	8411050

*Tradeware



VENTILATION ACCESSORIES

REUSABLE BREATHING CIRCUITS/SETS

Reusable Breathing circuit (N), 900MR782* MP17128

*Tradeware

REUSABLE BREATHING SYSTEMS: BLUESET VENTILATION


A	BlueSet Ventilation (N)	8412344
	Consists of:	
	1 Y-piece, reusable, N, 90°, with suction opening	8403075
	2 Corrugated hose, reusable, N/P, 14 cm, connector 10 mm/10 mm, Ø 10 mm	8403333
	3 Temperature sensor mount, reusable, N	8403056
	4 Water trap with outlet, reusable N/P	8404760
	5 Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/10 mm, Ø 10 mm	2165848
	6 Water trap, reusable, N/P	8409627
	7 Red cap, 7 mm	8401645
	8 Breathing hose, silicone, reusable, N/P, 40 cm, connector 10 mm/22 mm, Ø 10 mm	2165856
	9 Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/22 mm, Ø 10 mm	2165821
	Dummy Plug	8410049



VENTILATION ACCESSORIES

REUSABLE ANESTHESIA SYSTEMS:
BLUESET ANESTHESIA

A	BlueSet Anesthesia (N)	M27542
Consists of:		
1	Breathing hose, silicone, reusable, N/P, 1.5 m, connector 10 mm/22 mm, Ø 10 mm	2165961
2	Breathing hose, silicone, reusable, N/P, 1.1 m, connector 22 mm/22 mm, Ø 10 mm	2165678
3	Breathing bag, silicone, reusable, 0.5 L, not made with natural rubber latex	2166054
4	Adapter, reusable, 22 mm/22 mm	M25647
5	Y-piece, reusable, N, 90°, with suction opening	8403075
6	Elbow, reusable, N/P, metal, Luer Lock	M34568
7	Set of 5 red caps, 7 mm	8402953

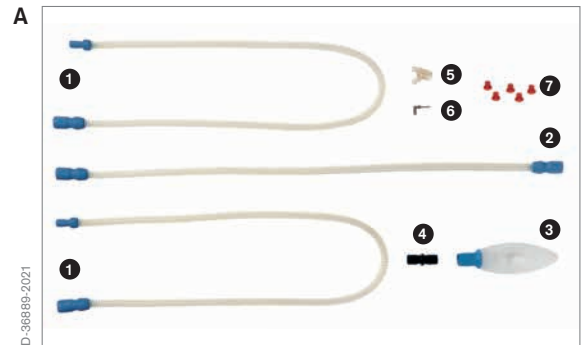

VENTILATION ACCESSORIES


Silicone breathing hose (pediatric) with standard sleeves 10 mm (0.4 in) on both sides and 10 mm (0.4 in) inner diameter

A Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/10 mm, Ø 10 mm

2165848

A

D-36216-2011

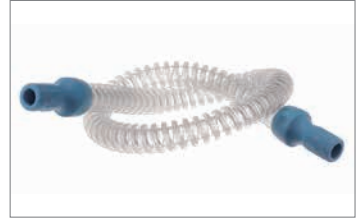

VENTILATION ACCESSORIES

Silicone breathing hose (pediatric) with standard sleeves with one 10 mm (0.4 in) and one 22 mm (0.9 in) sleeve and 10 mm (0.4 in) inner diameter

A	Breathing hose, silicone, reusable, N/P, 40 cm, connector 10 mm/22 mm, Ø 10 mm	2165856
B	Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/22 mm, Ø 10 mm	2165821
C	Breathing hose, silicone, reusable, N/P, 1.1 m, connector 10 mm/22 mm, Ø 10 mm	2165651
D	Breathing hose, silicone, reusable, N/P, 1.5 m, connector 10 mm/22 mm, Ø 10 mm	2165961



VENTILATION ACCESSORIES

A <u>Y-piece straight, reusable, N/P</u>	<u>M27077</u>
B <u>Y-piece, reusable, N, 90°, with suction opening</u>	<u>8403075</u>
C <u>Elbow, reusable, N/P, metal, Luer Lock</u>	<u>M34568</u>



VENTILATION ACCESSORIES

Transparent breathing bags of silicone elastomer

A	Breathing bag, silicone, reusable, 0.5 L, not made with natural rubber latex	2166054
B	Breathing bag, silicone, reusable, 1.5 L, not made with natural rubber latex	2166127
C	Demo lung to simulate neonatal lungs, with suitcase	8418077
D	Neonatal Test Lung, reusable, 50 ml, not made with natural rubber latex	MP20050



VENTILATION ACCESSORIES

The TwinStar® Plus filter/HME from Dräger is a combination of an efficient HME and a high-performance breathing system filter, thereby significantly contributing to infection prophylaxis in ventilation treatment.

- Combination of efficient HME and high-performance breathing system filter
- High retention rates for bacteria and viruses
- High return of moisture
- Sampling connector with tethered cap
- Transparent housing for visual control
- Standardized connectors for safe connection to other components

A Filter + HME TwinStar® 9 Plus*, disposable, ~9 ml deadspace, 100 pcs.

MP05820

HumidStar® Plus moisture exchangers from Dräger partially take over the function of the upper airways and efficiently humidify the inspired air.

- High return moisture
- Sampling connector with tethered cap
- Transparent housing for visual control
- Standardized connectors for safe connection to other components

B HME HumidStar® 2 Plus, disposable, 2 ml deadspace, 100 pcs.

MP05845

For local availability please get in touch with your local Dräger contact.



VENTILATION ACCESSORIES

CO₂ CUVETTES

CO₂ cuvettes are designed for the monitoring of CO₂ concentration during patient ventilation.

A CO₂ cuvette disposable, pediatric, 10 pcs. **MP01063**

B CO₂ cuvette reusable, pediatric **6870280**

CO₂ CABLES

C CO₂ mainstream sensor, for Evita® 2 dura, Evita® 4 Edition and Evita® XL **6871500**

D MCable Mainstream CO₂, for Evita V-family and Babylog VN-family, Infinity® M540, Vista120, Vista120 S, and Delta VF9.1 and higher **6871950**



VENTILATION ACCESSORIES

SENSORS

A	Neonatal flow sensor ISO 15, PSU, sterilizable	8411130
B	Insert for neonatal flow sensor, PSU, sterilizable, 5 pcs.	8410179
C	Flow sensor cable, for neonatal flow sensor, e.g., Babylog	8409626



VENTILATION ACCESSORIES

A	<u>O₂ sensor capsule, for mainstream measurement</u>	<u>6850645</u>
B	<u>O₂ sensor Oxycell, for Babyleo® TN500</u>	<u>6873065</u>
C	<u>O₂ sensor housing Oxydig, for O₂ sensor capsule, with spiral cable</u>	<u>6850250</u>



VENTILATION ACCESSORIES

ANESTHESIA MASKS SEDISTAR®

A	Anesthesia mask SediStar®, disposable, w/o hook ring, size 1, infant S, 30 pcs.	MP01701
B	Anesthesia mask SediStar®, disposable, w/o hook ring, size 2, infant L, 30 pcs.	MP01702

ANESTHESIA MASKS LITESTAR®

C	Anesthesia mask LiteStar®, disposable, w/o hook ring, size 1, infant S, 30 pcs.	MP11001
D	Anesthesia mask LiteStar®, disposable, w/o hook ring, size 2, infant L, 30 pcs.	MP11002



VENTILATION ACCESSORIES

ANESTHESIA MASKS COMFORTSTAR®

A	Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 0, premature baby, 20 pcs.	MP01510
B	Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 1, neonatal, 20 pcs.	MP01511
C	Anesthesia mask ComfortStar®, disposable, hook ring, standard, size 2, baby, 20 pcs.	MP01512
D	Anesthesia mask ComfortStar®, disposable, hook ring, strawberry, size 1, neonatal, 20 pcs.	MP01521
E	Anesthesia mask ComfortStar®, disposable, hook ring, strawberry, size 2, baby, 20 pcs.	MP01522
F	Anesthesia mask ComfortStar®, disposable, hook ring, bubble gum, size 1, neonatal, 20 pcs.	MP01531
G	Anesthesia mask ComfortStar®, disposable, hook ring, bubble gum, size 2, baby, 20 pcs.	MP01532

Pumpball as well as different headgears for fixation are available (cf. Mask accessories).



VENTILATION ACCESSORIES

SILICONE MASKS

A	Silicone mask, reusable, size 00, premature baby	MP02900
B	Silicone mask, reusable, size 0, small neonate	MP01910
C	Silicone mask, reusable, size 1, large neonate	MP01911
D	Hook ring for silicone masks, reusable, 15 mm OD, for sizes 00 – 1, 10 pcs.	MP01916


Additional mask accessories, e. g. reusable hook rings as well as reusable headgears are available for fixation (cf. Mask accessories).



VENTILATION ACCESSORIES

ACCESSORIES FOR ANESTHESIA MASKS


A	Headgear, disposable, soft fabric, small	MP01559
B	Headgear, disposable, soft rubber, small	MP01561
C	Two-sided pump ball for air management	MP01590

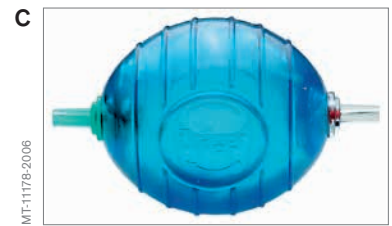

VENTILATION ACCESSORIES

NASAL CANNULAS

A	O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size S, 50 pcs.	MP02920
	O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size M, 50 pcs.	MP02921
B	O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size S, 50 pcs.	MP02924
	O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size M, 50 pcs.	MP02925

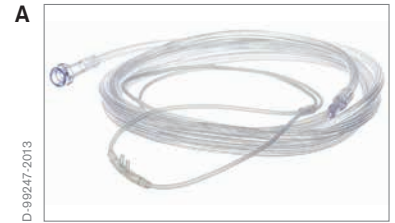

VENTILATION ACCESSORIES

SET2GO

- A** With Dräger's Set2go solutions we offer you totally customised sets consisting of an individual selection of single-patient-use accessories. Dräger Set2go makes using and ordering consumables and accessories much simpler for you. You just need one article number to get all the items in a single package and accomplish your specific task.

Please contact your local Dräger sales organisation to create your individual accessory solution.



VENTILATION ACCESSORIES

AEROGEN PRO*

A	Medication nebulizer Aerogen Pro	MP01010
B	Nebulizer unit with filler cap for Aerogen Pro	MP01030
C	T-piece set with plug for Aerogen Pro, N, 15 mm (0.6 in), w/ adapter 10 mm (0.4 in), 5 pcs.	MP01032
D	T-connector for Aerogen Pro, 12 mm, with plug, 5 pcs.	MP01037

*Tradeware



VENTILATION ACCESSORIES

AEROGEN SOLO*

A	Aerogen Solo Starter Kit, Int.	MP01071
	Aerogen Solo Starter Kit, FR, NL, IT, GE	MP01072
	Aerogen Solo Starter Kit, ES, IT, PT, GR, TK	MP01073
	Aerogen Solo Starter Kit, SW, NO, FI, DK	MP01074
	Aerogen Solo Starter Kit, UK	MP01075
	Aerogen Solo Starter Kit, US	MP01076
B	Nebulizer head, 5 pcs.	MP01078
C	Nebulizer head, 10 pcs.	MP01079
	T-Adapter, pediatric, 10 pcs.	MP01086
	T-Adapter**, neonate, 10 pcs.	MP02689
	Kit T-Adapter, Pediatric with Neo connectors, 10 pcs.	MP01087
	Aerogen Solo USB Controller Starter Kit NE	MP01089
	Aerogen USB Controller Replacement w/o AC/DC adapter	MP01095

*Tradeware

**Made by Drägerwerk AG & Co. KGaA



MONITORING ACCESSORIES

4. MONITORING ACCESSORIES

Infinity® Monitoring: Gamma, Gamma XL, Gamma X XL, Delta, Delta XL, Kappa, Kappa XLT.

ELECTROCARDIOGRAPHY (ECG)

A	Multiparameter cable NeoMed Pod, for ECG, resp., SpO₂, temp., and FiO₂, 2.5 m	5590539
	ECG electrodes*, single-patient use, neonate, 300 pcs.	5195024
B	Electrode adapter, to connect neonatal electrodes to MultiMed® pod, 10 pcs.	5194779
C	ECG intermediate cable, 3-lead, for NeoMed Pod, 1.5 m	5592162

PULSE OXIMETRY (SPO₂)

DRÄGER/NELLCOR® CABLES

D	SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 1 m	3368433
	SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 2 m	3375834
E	SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 1.2 m	MS18683
F	SpO₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 3 m	MS17330

*Tradeware



MONITORING ACCESSORIES

DRÄGER DISPOSABLE SENSORS

A SpO ₂ adhesive sensor Dräger, cloth, disposable, baby, 24 pcs.	MS16447
B SpO ₂ adhesive sensor Dräger, cloth, disposable, neonate, 24 pcs.	MS16446

NELLCOR® SENSORS*


C SpO ₂ adhesive sensor Nellcor® Oximax® Max-I, disposable, child, 24 pcs.	MX50067
D SpO ₂ adhesive sensor Nellcor® Oximax® Max-N, disposable, neonate/adult, 24 pcs.	MX50068

NELLCOR® ACCESSORIES*

E Adhesive wrap, for Nellcor® Oxiband™ A/N, 100 pcs.	8201012
---	----------------

*Tradeware

INFINITY
MONITORING
MICRO₂+

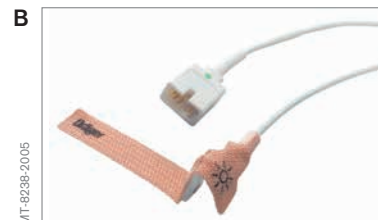

MONITORING ACCESSORIES

INFINITY®
MONITORING
MICRO₂+

MASIMO® CABLES*

**SpO₂ intermediate cables
for connection to Masimo® SET Pod.**

- A** SpO₂ intermediate cable Masimo® LNCS, 3 m MS17522
- A** SpO₂ intermediate cable Masimo® LNCS, 1.2 m MS24303

MASIMO® RD SET ADAPTER CABLES*

- B** Adapter cable RD SET to LNOP,
for connection of RD SET sensors to
Masimo® SET Pod via the
LNOP intermediate cables MS18680
and MS17041, as well as to MultiMed® 5/6/12,
NeoMed Pods and MultiMed® Plus/Plus OR
via the according intermediate cables
MS13926/7492601 or MS20163* MS33731
- C** Adapter cable RD SET to LNCS,
for connection of RD SET sensors to
Masimo® SET Pod via the
LNCS intermediate cables MS24303
or MS17522 MS33732
- C** Adapter cable RD SET to MultiMed®,
for connection of RD SET sensors to
MultiMed® 5/6/12, NeoMed Pods as well as
MultiMed® Plus/Plus OR via the according
intermediate cables 3368433/3375834
or MS18683/MS17330 MS33733

*Tradeware

**Products are not available any more.

Adapter is to be used to switch from sensors
to the installed base.



MONITORING ACCESSORIES

INFINITY
MONITORING

MASIMO® LNCS SENSORS*

For Masimo® accessories, specific part numbers for USA and Japan apply. Please contact your local Dräger Sales representative for further information.

- A** SpO₂ adhesive sensor Masimo® LNCS Neo, disposable, neonate, 20 pcs. **MP00792**
- B** SpO₂ adhesive sensor Masimo® LNCS NeoPt, disposable, premature baby, 20 pcs. **MP00794**
- SpO₂ Masimo® sensor LNCS-YI Multisite, reusable, neonate/pediatric/adult **MP00789**

MASIMO® LNCS TRAUMA DISPOSABLE SENSORS*


Masimo® LNCS Trauma sensors are for use with Masimo® SET technology only.

- C** Masimo® Trauma sensor LNCS, neonate, 20 pcs. **MP03192**

MASIMO® RD SET DISPOSABLE SENSORS*

- D** Masimo® RD SET sensor, neonatal, disposable **MS33729**
- E** Masimo® RD SET sensor, neonatal preterm, disposable **MS33730**

*Tradeware

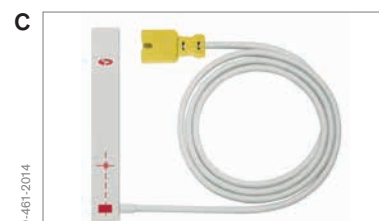


MONITORING ACCESSORIES

INFINITY®
MONITORING
MICRO₂+

MASIMO® ACCESSORIES*

A	SpO₂ adhesive wrap Masimo® Posey-Wrap, for NeoPt SpO₂ sensor, 12 pcs.	7496933
B	SpO₂ adhesive wrap Masimo®, for NeoPt SpO₂ sensor, 100 pcs.	7496925
C	SpO₂ adhesive wrap Masimo® Multisite-Wrap, for LNOP YI SpO₂ sensor, 100 pcs.	7496917
D	SpO₂ adhesive wrap Masimo® Petite-Wrap, for LNOP YI SpO₂ sensor, 100 pcs.	7496909
E	SpO₂ foam wrap Masimo® Foam-Wrap, neonate, LNOP YI SpO₂ sensor, 12 pcs.	7496891

*Tradeware

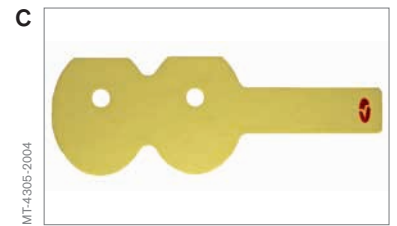

MONITORING ACCESSORIES

TEMPERATURE

SINGLE-PATIENT USE TEMPERATURE PROBES

A	Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 1.6 m	MP00997
B	Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 3 m	MP00998

REUSABLE TEMPERATURE PROBES

C	Skin temperature probe, reusable, with 7-pin connector, 1.5 m	4329822
C	Skin temperature probe, reusable, with 7-pin connector, 3 m	5204669

INFINITY®
MONITORING



Not all articles are available worldwide.

MONITORING ACCESSORIES

NON-INVASIVE BLOOD PRESSURE CUFFS (NIBP)

A	NIBP hose, neonate, 2.4 m	2870298
B	NIBP cuff, disposable, neonate, size 1, 3.1 – 5.7 cm, 10 pcs.	2870181
C	NIBP cuff, disposable, neonate, size 2, 4.3 – 8 cm, 10 pcs.	2870199
D	NIBP cuff, disposable, neonate, size 3, 5.8 – 10.9 cm, 10 pcs.	2870207
E	NIBP cuff, disposable, neonate, size 4, 7.1 – 13.1 cm, 10 pcs.	2870215
F	NIBP cuff, disposable, neonate, size 5, 8.3 – 15 cm, 10 pcs.	2870173

INFINITY
MONITORING

MONITORING ACCESSORIES

ETCO₂ AND RESPIRATORY MECHANICS

For Delta, Delta XL, Kappa, Kappa XLT.

A etCO₂ Microstr. Intubated Infant/Neonatal – Extended Duration, 25 pcs. MP15003

B etCO₂ Microstr. Nasal Infant/Neonatal – Extended Duration, 25 pcs. MP15007

For Gamma X XL, Delta, Delta XL, Kappa, Kappa XLT.

C etCO₂ airway adapter, reusable, neonate 4721788


D CO₂ airway adapter*, neonatal (infant/pediatric), disposable, 10 pcs. 2606228

RECORDER PAPER

E Recording paper for Infinity® R 50, 50 mm, approx. 30 m, 10 rolls 4711201

*Tradeware

INFINITY
MONITORING

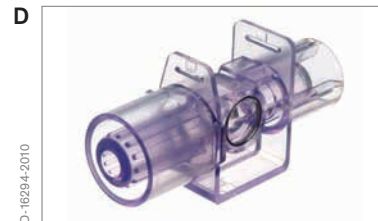

WARD EQUIPMENT

5. WARD EQUIPMENT

Please find more information in chapter 5 „Ward equipment“.

FLOWMETER

These O₂ Flowmeters are plugged directly into a wall terminal unit for oxygen and are used to control the flow of oxygen in inhalation and insufflation. Flow rates can be set with a high degree of accuracy. The large rotary knob makes setting easier.

Accuracy of flow rate at 4.5 bar ± 0.5 bar (1,013 hPa/20 °C)**
 ± 15% (corresponding to smallest setting value)
 ± 10% (all other setting values)

** Using the flowmeters outside the specified ambient conditions could affect flow accuracy.

Dimensions (H × W × D):	160 × 40 × 130 mm (6.3 × 1.6 × 5.1 in)
Material:	Housing: Aluminium Pressure dome: Polyamide Flow indicator: Glass
Weight:	Approx. 500 g (1.1 lb)

A O₂ flowmeter, 4 l/min, CS, DIN, ISO and neutral color, M34×1.5 and 9/16“ for accessories and sterile water systems

MP04551



WARD EQUIPMENT

TWIN FLOWMETER

The particular advantage of the twin flowmeter is, that two patients can be treated simultaneously and independently of one another, while occupying only one wall terminal.

Accuracy of flow rate at 4.5 bar \pm 0.5 bar (1,013 hPa/20 °C)**
 \pm 15% (corresponding to smallest setting value)
 \pm 10% (all other setting values)

** Using the flowmeters outside the specified ambient conditions could affect flow accuracy.

Dimensions (H \times W \times D): 160 \times 150 \times 190 mm
 (6.3 \times 5.9 \times 7.15 in)

Material: Housing: Aluminium
 Pressure dome: Polyamide
 Flow indicator: Glass

Weight: Approx. 900 g (1.98 lb)



A O₂ twin flowmeter, CS, 2 \times 16 l/min, CS, DIN, ISO, neutral color, M34\times1.5 and 9/16" for accessories and sterile water systems	MP04560
O₂ twin flowmeter, 2 \times 16 l/min, rail fitting with NIST gas inlet, 2 \times M34\times1.5, adapter 9/16"	MP04569

WARD EQUIPMENT

These O₂ flowmeters for rail attachment can be used to control the flow of oxygen in inhalation and insufflation. Flow rates up to 4 and 16 l/min can be set with a high degree of accuracy. The large rotay knob makes setting easier. Using the rail clamp, the flowmeter can be quickly and safely attached to all rail profiles. The connection to the gas supply can be ensured by a CS hose with NIST coupling nut.

Accuracy of flow rate at 4.5 bar \pm 0.5 bar (1,013 hPa/20 °C)**
 \pm 15% (corresponding to smallest setting value)
 \pm 10% (all other setting values)

** Using the flowmeters outside the specified ambient conditions could affect flow accuracy.

Dimensions (H \times W \times D): 160 \times 40 \times 130 mm
 (6.3 \times 1.6 \times 5.1 in)

Material: Housing: Aluminium
 Pressure dome: Polyamide
 Flow indicator: Glass

Weight: Approx. 950 g (2.09 lbs)

- | | | |
|----------|---|----------------|
| A | O₂ flowmeter, 4 l/min, rail with NIST gas inlet, M34\times1.5 and 9/16" for accessories and sterile water systems | MP04552 |
| B | O₂ flowmeter, 4 l/min, rail, 90°, NIST gas inlet, M34\times1.5 and 9/16" for accessories and sterile water systems | MP04553 |
| C | O₂ flowmeter, 16 l/min, rail, NIST gas inlet, M34\times1.5, adapter 9/16" | MP04567 |

CS hoses

Several hoses to connect flowmeter with the central supply.
 For CS hoses please refer to chapter "Accessory for gas supply".




WARD EQUIPMENT

COMPACT FLOWMETER CLICKFLOW

Input pressure:	4.5 ± 0.5 bar*
Adjustable ranges:	0.1 – 1 l/min (0; 0.1; 0.2; 0.3; 0.4; 0.5; 0.6; 0.7; 0.8; 1 l/min) 0.5 – 5 l/min (0; 0.5; 1; 1.5; 2; 2.5; 3; 3.5; 4; 5 l/min) 1 – 16 l/min (0; 1; 2; 3; 4; 6; 8; 10; 12; 16 l/min) 1 – 32 l/min (0; 1; 5; 8; 12; 16; 20; 25; 28; 32 l/min)
Accuracy (at 4.5 ± 0.5 bar):	0.5 – 5 l/min, 1 – 16 l/min and 1 – 32 l/min: ± 20%; 0 – 1 l/min ± 30%
Gas inlet:	CS plug O ₂ or AIR; NIST (rail)
Gas outlet:	6 mm nozzle for application unit, thread 9/16" for sterile water or adapter for reusable humidifier
Weight:	approx. 420 g (CS); approx. 700 g (rail)
Housing material:	nickel-coated brass
Dimensions:	diameter about 38 mm, length about 92 mm (CS version); about 118 mm (rail version)

*1 bar = 1 kPa × 100

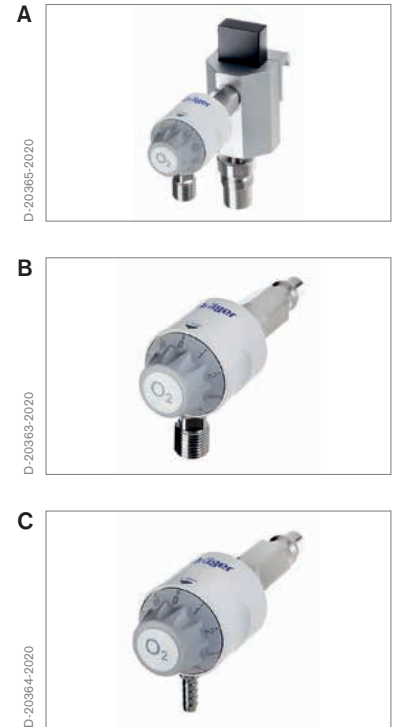

- A Compact flowmeter Clickflow configurator, gives you the option of arranging the compact flowmeter according to your requirements, selectable options for configuration: MP13000**
- | | |
|----------------|---|
| Gas type: | O ₂ ; AIR |
| Mounting: | CS; rail with NIST gas inlet |
| CS connection: | DIN; Carbamed; BS (England); NF (France); AGA (Sweden) |
| Version: | single compact flowmeter;
double compact flowmeter |
| Flow: | 0.1 – 1 l/min; 0.5 – 5 l/min;
1 – 16 l/min, 1 – 32 l/min |
| Gas outlet: | 6 mm nozzle, thread 9/16" for sterile water |
- B O₂ compact flowmeter Clickflow, DIN CS, 1 – 16 l/min, 9/16" MP13011**
- C O₂ compact flowmeter Clickflow, DIN CS, 1 – 16 l/min, 6 mm nozzle MP13012**

Accessories for flowmeter: See section Ward Equipment, chapter 1.4 (humidifier, adapter and application).

WARD EQUIPMENT

A	Humidifier-bubbler, for oxygen humidification up to appr. 60% relative humidity and flows up to 5 l/min	2M85834
B	Adapter, to connect e.g a tube or nasal cannula to O₂ flowmeter, M34×1.5 to 6 mm nozzle	M06258
C	Adapter, for disposable sterilized water systems to O₂ flowmeters, M34×1.5 to 9/16"	2M86680
C	Adapter, for disposable sterilized water systems to AIR flowmeter, M32×1.5 to 9/16"	MP03366
D	Adapter, with 6 mm nozzle, e.g., for disposable breathing masks, 9/16" connection	MP00416
E	Adapter 9/16" to M34×1.5 for Dräger reusable humidifier	MP03325



WARD EQUIPMENT

MEDICAL ASPIRATION

The vacuum regulator is directly connected to a central vacuum supply, or mounted on a rail. The smoothly adjustable vacuum, up to max. -160 mbar (-16 kPa), is displayed on a rotating manometer, legible from all sides.

All vacuum regulators include a shut-off valve.

You can use the VarioSafe® disposable filter system to reliably protect the **VarioVac®** series and the central vacuum system against contamination. A VarioSafe® filter system is included with every VarioVac® vacuum regulator.

Technical data:

Maximum vacuum:	-16 to 0 kPa (-160 to 0 mbar)
Maximum suction flow:	18 l/min
Weight:	600 to 800 g (800 g for rail versions)
Dimensions (W/H/D) in mm:	80 × 210 × 160 (without VarioSafe) 80 × 340 × 160 (with VarioSafe)
Maximum sound level:	50 dB (A)
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



A Vacuum regulator VarioVac® D Rail, for continuous and wound drainage, gas connection VAC MP00511

For CS hoses please refer to chapter „Accessory for gas supply“.

WARD EQUIPMENT

VarioVac® vacuum regulator for rail.

The gas supply is ensured by VAC connection hose.

This rail version can be connected to all standard wall rails (10 mm width and 25 mm height).

You can use the VarioSafe® disposable filter system to reliably protect the VarioVac® series and the central vacuum system against contamination.

A VarioSafe® filter system is included with every VarioVac® vacuum regulator.

Technical data:

VarioVac P (MP00660, MP03223)

Maximum vacuum:	-55 to 0 kPa (-550 to 0 mbar)
Maximum suction flow:	20 l/min
Maximum sound level:	40 dB(A)
Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 210 × 160 (without VarioSafe) 80 × 340 × 160 (with VarioSafe)
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



- A Vacuum regulator VarioVac® P Rail,
for bronchial suction, gas connection NIST VAC,
for connection to all standard wall rails**

MP00660

The bronchial aspirator can be used in magnets with static magnetic fields of 1.5 tesla and 3 tesla.

Minimum distance from the center of the magnetic field:

Magnetic flux density of 40 millitesla (400 gauss).

- B Vacuum regulator VarioVac® MRI P Rail,
for bronchial suction, gas connection NIST VAC,
for connection to all standard wall rails**

MP03223

For CS hoses please refer to chapter „Accessory for gas supply“.

WARD EQUIPMENT

Drive unit VarioVac® and jar set are joined together, so that you have a complete, convenient and easy solution. The ergonomically shaped handle is positioned exactly at the units centre of gravity, making it easy to transport the aspirator quickly. The gas supply is ensured by CS connection hose. You can use the VarioSafe® disposable filter system to reliably protect the device and the patient environment against contamination. This rail version can be connected to all standard wall rails (10 mm width and 25 mm height).

Technical data:

Weight: 1,600 g (empty)
 Dimensions (W/H/D) in mm: 300 × 440 × 160
 Hose nozzle: Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention: 100% (with VarioSafe)*
 Bacterial retention: 99.99997% (with VarioSafe)**
 Viral retention: 99.99922% (with VarioSafe)**
 NaCl retention: 99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



A Vacuum regulator VarioVac® P compact, for bronchial suction, gas connection NIST VAC, including secretion jar set and integrated vacuum regulator **MP00661**

Technical data:

Maximum vacuum: -55 to 0 kPa (-550 to 0 mbar)
 Maximum suction flow: 20 l/min
 Maximum sound level: 40 dB(A)

For CS hoses please refer to chapter „Accessory for gas supply“.

WARD EQUIPMENT

The **VarioAir®** ejector is connected directly to a medical air supply system or to rails. The low-noise ejector generates a vacuum up to max. -16 kPa (-160 mbar) and shows the vacuum level on a turnable, manually adjustable pressure gauge. All devices offer a pressure relief valve (ex 0.1 kPa) to avoid pressure increase at the patient side.

A new VarioSafe® disposable filter system reliably protects the VarioAir® series and the patient environment against contamination. One piece of the new VarioSafe® filter system is always included.

Technical data:

Maximum vacuum:	-16 to 0 kPa (-160 to 0 mbar)
Maximum suction flow:	18 l/min
Inlet pressure range:	4 to 5.5 bar
Maximum gas consumption:	25 l/min
Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 200 × 160 (without VarioSafe) 80 × 330 × 160 (with VarioSafe)
Maximum sound level:	45 dB (A)
Hose nozzle:	Ø 6 mm
Pressure relief valve:	Ex 0.1 kPa (1 mbar)

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



A Ejector VarioAir® D Rail, for continuous and wound drainage, gas connection AIR, for connection to all standard wall rails

MP00542

For gas supply via CS vacuum connection hose.

Connection thread: NIST.

For CS hoses please refer to chapter „Accessory for gas supply“.

WARD EQUIPMENT

VarioAir® ejector for rail.

The gas supply is ensured by AIR connection hose.

This rail version can be connected to all standard wall rails (10 mm width and 25 mm height).

You can use the VarioSafe® disposable filter system to reliably protect the VarioAir® series and the central vacuum system against contamination.

Technical data:

VarioAir P (MP00545, MP03221)

Maximum vacuum:	-48 to 0 kPa (-480 to 0 mbar)
Maximum suction flow:	20 l/min
Maximum gas consumption:	40 l/min
Maximum sound level:	47 dB(A)
Inlet pressure range:	4 to 5.5 bar
Weight:	800 g
Dimensions (W/H/D) in mm:	80 × 200 × 160 (without VarioSafe) 80 × 330 × 160 (with VarioSafe)
Hose nozzle:	Ø 6 mm
Pressure relief valve:	from 0.1 kPa on (1 mbar)

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



- A Ejector VarioAir® P Rail,
for bronchial suction, gas connection NIST AIR,
for connection to all standard wall rails** MP00545

The bronchial aspirator can be used in magnets with static magnetic fields of 1.5 tesla and 3 tesla.

Minimum distance from the center of the magnetic field:
Magnetic flux density of 40 millitesla (400 gauss).

- A Ejector VarioAir® MRI P Rail,
for bronchial suction, gas connection NIST AIR,
for connection to all standard wall rails** MP03221

For CS hoses please refer to chapter „Accessory for gas supply“.

WARD EQUIPMENT

Drive unit VarioAir® and jar set are joined together, so that you have a complete, convenient and easy solution. The ergonomically shaped handle is positioned exactly at the units centre of gravity, making it easy to transport the aspirator quickly.

The gas supply is ensured by CS connection hose. You can use the VarioSafe® disposable filter system to reliably protect the device and the patient environment against contamination.

This rail version can be connected to all standard wall rails (10 mm width and 25 mm height).

Technical data:

Weight:	1,600 g (empty)
Dimensions (W/H/D) in mm:	300 × 440 × 160
Hose nozzle:	Ø 6 mm

Additional Information for Filter:

Aerosole and Fluid retention:	100% (with VarioSafe)*
Bacterial retention:	99.99997% (with VarioSafe)**
Viral retention:	99.99922% (with VarioSafe)**
NaCl retention:	99.8900% (with VarioSafe)***

* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany

** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA

*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1



A Ejector VarioAir® P compact, for bronchial suction, gas connection NIST AIR, including secretion jar set and integrated ejector

MP00547

Technical data:

Maximum vacuum:	-48 to 0 kPa (-480 to 0 mbar)
Maximum suction flow:	20 l/min
Inlet pressure range:	4 to 5.5 bar
Maximum gas consumption:	40 l/min
Maximum sound level:	47 dB(A)
Pressure relief valve:	Ex 0.1 kPa (1 mbar)

For CS hoses please refer to chapter „Accessory for gas supply“.

WARD EQUIPMENT

Extra lightweight portable secretion jar set for bronchial aspiration. The set is made of high quality and extremely robust plastic and can be autoclaved at up to 134 °C. The two 700 ml jars, which are used as secretion and rinsing jars, are identical, making hospital logistics considerably easier. The jars are made of high quality polysulfone, are reusable, extremely robust and implosion proof. The long service lives makes the jars both environmentally friendly and cost effective.

An automatic overflow safeguard is fitted in the secretion jar cap to prevent secretion from entering the drive unit and the piped medical gas supply. The ergonomically shaped handle is positioned exactly at the units centre of gravity, making it easy to transport the secretion jar set quickly.



- A Secretion jar set, with 700 ml secretion jars (2 pcs.) and 1.5 m suction hose, maximum capacity up to the overflow safeguard 2 × 1,000 ml, scale in steps of 10 ml up to max. 700 ml**

2M85056

WARD EQUIPMENT

VacuSmart disposable cartridge for bronchial suction

VacuSmart is a disposable cartridge for 700 ml secretion jars. This product has been specifically designed and thoroughly tested to work with Dräger suction units. Therefore there is absolutely no need to modify the installed base. An integrated hydrophobic filter acts as an overflow safeguard by stopping suction of secretion once maximum capacity has been reached. The closed and leakproof system, optional with an integrated gelling agent, reduces the risk of cross contamination to a minimum. The optional gelling agent prevent leakages e.g. during the waste disposal in a waste press. The gelling agent is secured in a water soluble bag.



Technical data:

Maximum capacity:	700 ml
Maximum capacity up to the overflow safeguard:	1,000 ml
Overflow safeguard:	Integrated, hydrophobic filter
Recommended duration of use:	24 h
Gelling agent:	Sodium polyacrylate

A	Secretion jar VacuSmart®, PE (Polyethylene), disposable, 700 ml, 25 pcs.	MX23023
B	Secretion jar VacuSmart® Gel, PE (Polyethylene), disposable, 700 ml, 25 pcs.	MX23025
C	Suction tube VacuSmart® Tube, PE (Polyethylene), disposable, 2 m, w/ finger tip, 25 pcs.	MX23026

You can also enjoy some smart savings by choosing any of the following packing unit 25 pieces sets:

D	Suction set, disposable, VacuSmart®, 700 ml and VacuSmart® Tube, 25 pcs. each	MP00061
D	Set VacuSmart® Gel and Tube, 25 pcs.	MP03211

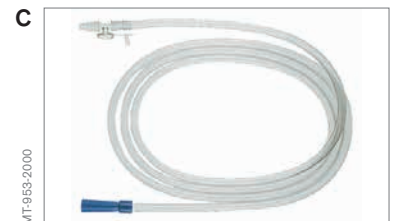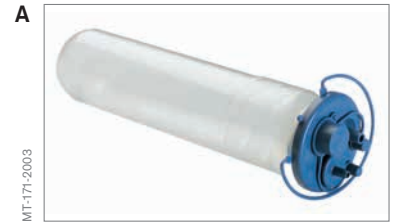

WARD EQUIPMENT

A	Bacterial filter 767, 5 pcs.	6723976
B	Bacterial and viral filter system VarioSafe, for VarioAir®, VarioVac®, VarioO2, disposable, 10 pcs.	MP00555
	The new VarioSafe® disposable filter system reliably protects the VarioAir®, VarioVac® and VarioO2 series and the patient environment against contamination.	
	Technical data:	
	Material:	PS (polystyrol)
	Recommended duration of use:	7 days
	Maximum capacity:	90 ml
	<u>Additional Information for Filter:</u>	
	Aerosole and Fluid retention:	100% (with VarioSafe)*
	Bacterial retention:	99.99997% (with VarioSafe)**
	Viral retention:	99.99922% (with VarioSafe)**
	NaCl retention:	99.8900% (with VarioSafe)***
	* Test has been performed by Drägerwerk AG & Co. KGaA, Dec 2020, Lübeck, Germany	
	** Test has been performed by Nelson Laboratories Inc., Dec 2020, Salt Lake City, USA	
	*** Test has been performed by Dräger Safety AG & Co. KGaA, March 2021, Lübeck, Germany, according to EN ISO 23328-1	
C	Suction hose, silicone, reusable, 1.5 m, for medical suction	M25780
D	Connection hose, silicone, reusable, ID 5 mm, wall thickness 2 mm, by meter	1203606
E	Secretion jar, PSU (Polysulfone), reusable, 300 ml, autoclavable up to 134°C	2M86207
F	Secretion jar, PSU (Polysulfone), reusable, 700 ml, autoclavable up to 134°C	2M85594
G	Elbow, angled hose connection for secretion jar caps	2M19063



WARD EQUIPMENT

Holders for attaching a compact bronchial aspirator or a secretion jar set to the wall/ rail or to a rail on a ceiling supply unit. The devices are only latched and can be detached and moved easily if this is required.

A Wall holder, with fixing screws, for compact suction unit and secretion jar set **M26859**

B Rail clamp lever type, for attachment of devices to 25 × 10 mm standard rail, for compact suction unit and secretion jar set **M26842**

C Ceiling holder, rail holder for compact suction systems and secretion jar sets **2M85119**

Holders for holding a central supply connecting hose or an aspiration hose. For attachment to an intensive care rail or compact rail.

D CS connection hose holder **2M85446**

E Suction hose holder, suitable for various diameters **2M85239**



RAIL ACCESSORIES AND EXAMINATION LIGHT

6. RAIL ACCESSORIES AND EXAMINATION LIGHT

RAIL ACCESSORIES

Wire baskets*

The wire baskets are used for storage of disposable catheters or disposable gloves. The catheters can be simply stored in bundles and are then easy to take out one by one. The baskets are chromeplated and are thus easy to clean and can be autoclaved at up to 134 °C. The amount of space for a basket on the rails is around 150 mm (5.9 in).

A	Basket 150, 15 × 20 × 10 cm, e.g., for disposable gloves	M26146
A	Basket 300, 15 × 28 × 10 cm, e.g., for 300 mm catheters	M26145
A	Basket 600, 15 × 48 × 10 cm, e.g., for 600 mm catheters	M25121
B	Litter bag holder, stainless steel, with base plate and litter bags (100 pcs.)	M24695
	Litter bags, disposable, for litter bag holder, 100 pcs.	M26240

*Tradeware



RAIL ACCESSORIES AND EXAMINATION LIGHT

EXAMINATION LIGHT

The Dräger **VarioLux®** offers cool, highly efficient, and variable intensity examination light and was specially designed for use in neonatology and pediatric settings. Made for single-handed operation, it features adjustable color temperature settings as well as a dimmer function.


Operating mode:	Continuous operation
Weight:	1.5 kg (3.3 lbs)
Electrical values:	100 – 240 V, 50/60 Hz
Max. illumination intensity E_c at 0.5 m:	30,000 lx
Dimming levels:	30,000 lx, 24,000 lx, 18,000 lx, 12,000 lx, 7,000 lx
Light field diameter d_{10} at 0.5 m:	20 cm (7.9 in)
Light field diameter d_{50} at 0.5 m:	9 cm (3.5 in)
Color temperature adjustable:	Cool white (blue) 4,700 K, Neutral white (white) 4,100 K, Warm white (red) 3,500 K
Min. service life of bulb (LEDs):	> 50,000 h
Color reproduction index:	$R_a \geq 96$ (at 3,500 K) $R_9 = 89$ (at 3,500 K)
Max. radiant power at 0.5 m:	105 W/m ²

A	Medical examination light VarioLux® with LED technology, EU	MP00601
	Medical examination light VarioLux® with LED technology, US	MP00602
	Medical examination light VarioLux® with LED technology, GB	MP00603
	Medical examination light VarioLux® with LED technology, AUS	MP00606
	Medical examination light VarioLux® with LED technology, JPN	MP00608
	Power cable EU, 1.8 m, 2 pin	MP00591
	Power cable US/JPN, 1.8 m, 2 pin	MP00592
	Power cable GB, 1.8 m, 2 pin	MP00593
	Power cable AUS, 1.8 m, 2 pin	MP00626



RAIL ACCESSORIES AND EXAMINATION LIGHT

A	Rail claw for lamp	MP00615
B	Turn knob, compl.	MP00587
C	Control panel	MP00589
D	Rail clamp, for attachment to poles of different diameters, rotatable by 90° to left and right, length of standard rail (profile 25 × 10 mm): 20 cm	2M85274
E	Compact rail, for 38 mm tube, standard rail profile 25 × 10 mm, each side 10 cm, weight: 0.9 kg, max. load: 30 kg	2M85337

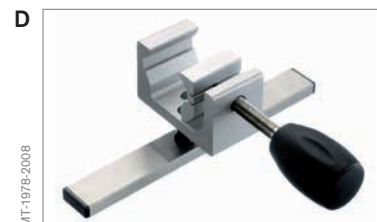

RAIL ACCESSORIES AND EXAMINATION LIGHT

HINGED ARMS

Keeping everything close at hand. To leave your hands free to care for the patient, all the joints can be fixed into place and released with just one quick hand movement. All cables, lines and hose systems can be arranged in such a way as to more or less completely prevent accidental disconnection.

- | | | |
|----------|--|----------------|
| A | Hinged arm Quickstop 2, with universal ball holder system for optional accessories length 858 mm, max. load 0.4 kg, includes hose clamps: MP00493 (1×), 8404761 (2×), 8409841 (1×) | 2M85706 |
| B | Hinged arm Quickstop short, with mounted universal tube holder plate, for various tube dimensions/coaxial hoses, specifically designed for OR tables, length 323 mm, max. load 1 kg | 2M86463 |





NOTICE

Size classification to adults (A), pediatrics (P) and neonats (N)

Not made with natural rubber latex: product is produced without latex

INDEX BY ORDER NUMBER

1	1203606	299, 431	2601227	152, 202	2607559	232	
			2601228	151, 201	2607560	237	
2	2111764	95	2601229	151, 201	2607561	237	
	2112760	45, 46	2601230	151, 201	2612000	221	
	2120461	95	2601231	151, 201	2612001	221	
	2120488	95	2601232	150, 200	2612002	221	
	2120496	95	2601233	150, 200	2612003	221	
	2165619	40, 47	2601234	150, 200	2612020	220	
	2165627	40, 43, 47	2601235	150, 200	2612021	220	
	2165651	49, 395	2601236	153, 203	2612024	220	
	2165678	44, 48, 393	2605124	234	2612025	220	
	2165821	39, 49, 392, 395	2606151	227, 248	2612026	220	
	2165848	39, 48, 392, 394	2606152	227, 248	2612027	220	
	2165856	39, 49, 392, 395	2606153	227, 248	2612028	220	
	2165961	44, 49, 393, 395	2606154	227, 248	2612029	220	
	2166038	42, 43, 47	2606155	227, 248	2870173	145, 197, 416	
	2166046	41, 42, 45, 46, 47	2606156	227, 248	2870181	145, 197, 416	
	2166054	44, 52, 393, 397	2606208	225, 238	2870199	145, 197, 416	
	2166062	42, 43, 52	2606209	225, 238	2870207	145, 197, 416	
	2166070	52	2606210	225, 238	2870215	145, 197, 416	
	2166127	52, 397	2606211	225, 238	2870298	142, 194, 416	
	2170221	135	2606225	232			
	2600146	251	2606226	236	3	3368383	156
	2600147	251	2606227	167, 236		3368391	111
	2600930	249	2606228	167, 236, 417		3368425	112, 147, 198
	2600932	249	2606229	237		3368433	121, 124, 134, 410
	2600933	249	2606231	237		3368458	159, 204
	2600934	249	2606483	225, 238		3375834	121, 124, 134, 410
	2600935	249	2606484	225, 238		3375933	158
	2600936	249	2606485	225, 238		3375941	154
	2600937	249	2606486	225, 238		3375958	154, 161
	2600938	249	2606487	225, 238		3702751	350
	2600939	251	2606488	232		3702831	185
	2600940	251	2606489	232		3702834	185
	2600942	251	2606490	232		3702837	185
	2600943	251	2606491	232		3702841	185
	2600945	251	2606492	224		3718728	129, 183, 216
	2600946	250	2606493	224		3722467	185
	2600948	250	2606494	224		3722471	185
	2600949	250	2606495	224			
	2601190	233	2606496	224	4	4118453	348
	2601191	233	2606497	224		4118555	348
	2601192	233	2606498	224		4119065	334
	2601195	227, 248	2607032	235		4119066	334
	2601197	229	2607034	228		4119067	334
	2601198	229	2607035	228		4119068	334
	2601199	229	2607036	228		4119069	334
	2601200	229	2607037	228		4119070	334
	2601201	233	2607038	228		4119071	334
	2601220	153, 203	2607550	234		4119072	334
	2601221	153, 203	2607551	234		4119073	334
	2601222	153, 203	2607552	234		4119074	334
	2601223	153, 203	2607553	234		4119075	334
	2601224	152, 202	2607554	234		4119076	334
	2601225	152, 202	2607555	234		4119077	334
	2601226	152, 202	2607558	232		4119078	334

INDEX BY ORDER NUMBER

4119079	334	5590539	111, 410	8290274	40
4119080	334	5591925	148	8290286	25, 172
4119081	334	5591933	148	8306370	106
4119082	334	5592147	156	8401645	39, 51, 392
4119083	334	5592154	139, 191	8402868	25, 172
4185140	330	5592162	119, 410	8402953	44, 393
4185141	330	5598128	112	8403056	39, 51, 392
4185142	331	5703041	32	8403075	39, 44, 50, 392, 393, 396
4185143	331	5704962	34	8403076	51
4199589	330	5704964	33	8403201	53
4199590	330	5731281	156	8403333	39, 392
4199591	330	5741975	159, 204	8403566	40
4199592	330	5790124	95	8403735	105
4199593	330	5790148	376	8404760	39, 392
4199594	330	5947226	120	8404985	40, 43, 51
4199595	331	5947804	165	8405435	40, 50
4199596	331	5950196	111	8407475	45
4199600	330	5954859	165	8408950	85
4199601	330			8409609	342
4199602	330	6 6723976	298, 316, 431	8409626	106, 400
4199604	330	6850250	107, 401	8409627	39, 51, 392
4199605	330	6850435	169, 205	8409897	51
4199712	330	6850645	107, 401	8410049	39, 40, 392
4199713	331	6870279	62, 168, 205, 235	8410179	105, 400
4199714	330	6870280	62, 168, 205, 235, 399	8410547	95
4321563	157	6870281	169, 205	8410580	85
4329822	140, 192, 231, 415	6871500	62, 399	8410929	46
4329848	140, 192, 231	6871950	62, 168, 205, 235, 399	8411030	92
4329889	140, 192, 231	6871980	105	8411044	51
4528741	158	6872020	25, 172	8411045	29
4711201	173, 417	6872130	25, 172	8411047	28, 390
4721416	155	6872666	107	8411050	31, 391
4721424	155	6873065	107, 348, 401	8411074	29, 377, 390
4721614	155			8411075	348, 352
4721788	167, 417	7 7014616	140, 192, 229, 231	8411130	105, 400
4721796	167	7257988	112, 198	8411154	29, 390
		7262764	135, 190, 218, 226, 239	8412001	46
5 5000833	318	7265155	204	8412015	85
5001330	299	7265163	204	8412034	46, 105
5191221	111	7270460	155	8412068	46
5192112	155	7485621	112, 147, 148	8412217	28
5194779	119, 410	7486140	120	8412235	46, 51
5195024	119, 178, 224, 410	7489417	157	8412344	39, 392
5195180	157	7489433	157	8412710	169, 205
5195198	112, 198	7496875	132, 187, 217	8412860	41
5196980	154	7496891	132, 187, 217, 414	8412888	50
5196998	154	7496909	132, 187, 217, 414	8412913	46
5198333	139	7496917	132, 187, 217, 414	8412935	92
5202747	112	7496925	414	8412985	92
5204644	140, 192, 231	7496933	414	8413125	85
5204651	140, 192, 231	7498947	157	8413146	40
5204669	140, 192, 231, 415	7864155	157	8413293	46
5206581	158	7864908	157	8413660	86
5206607	156	7872174	166	8414144	28, 350
5588095	156			8414776	85
5588822	148	8 8201012	135, 411	8414966	30, 391

INDEX BY ORDER NUMBER

8414968	31, 391	2M50412	351	M07582	300
8414969	29	2M60170	347	M16395	317
8414989	30, 391	2M60180	348	M20668	317
8414992	31, 391	2M60390	347	M24695	339, 433
8414993	29, 390	2M60401	348	M25121	339, 433
8415270	85	2M85011	299	M25647	40, 42, 43, 44, 51, 393
8415296	85	2M85012	299	M25649	51
8416325	29, 390	2M85013	299	M25650	41, 50
8416352	169, 205	2M85056	279, 294, 429	M25780	299, 431
8416600	106	2M85119	296, 339, 432	M26145	339, 433
8416750	86	2M85139	317	M26146	339, 433
8417050	86	2M85151	314	M26240	339, 433
8418077	397	2M85153	314	M26842	279, 296, 339, 432
8418282	28, 390	2M85163	315	M26859	296, 432
8418289	51	2M85180	315	M27077	50, 396
8418415	378, 381	2M85239	296, 306, 339, 432	M27542	44, 393
8418416	378, 381	2M85274	339, 341, 377, 435	M30146	42
8418417	378, 381	2M85337	279, 341, 377, 435	M32493	335
8418490	379, 382	2M85380	295	M32494	335
8418491	379, 382	2M85381	295	M32495	335
8418507	30, 391	2M85441	305	M32496	335
8418531	378, 381	2M85446	296, 306, 339, 432	M32497	335
8418533	379, 384	2M85594	300, 431	M32692	25, 172
8418534	379, 384	2M85706	342, 436	M33278	42, 43, 50
8418535	379, 384	2M85834	281, 422	M34568	44, 51, 393, 396
8418536	379, 384	2M85835	281	M34875	335
8418537	379, 384	2M85853	281	M34876	335
8418538	379, 384	2M86185	280	M34877	335
8418539	379, 384	2M86205	277	M34878	335
8418603	378, 381	2M86207	300, 431	M35056	280, 335
8418605	378, 381	2M86291	95	M35057	335
8418607	378, 381	2M86292	95	M35058	280, 335
8418609	378, 381	2M86293	95	M35330	321
8418611	378, 381	2M86463	342, 436	M35331	321
8418617	378, 381	2M86680	282, 422	M35332	321
8418619	379, 382	2M86757	302, 306	M35333	323
8419160	159, 204	2M86800	45	M35334	323
8420077	159, 204	2M86832	45	M35335	323
8422215	23, 388	2M86837	32	M35336	335
8539983	159, 204	2M86913	95	M35337	335
8601697	319			M35338	335
		C CH00102	298, 316	M35339	324
2M 2M16520	359			M35340	324
2M19063	300, 431	G G01166	334	M35341	324
2M20434	358	G01168	334	M35342	322
2M20827	349, 357	G01169	334	M35343	322
2M21012	357	G35324	343	M35344	322
2M21158	358	G35325	343	M35345	326
2M21161	358	G35332	343	M35346	326
2M21186	352, 359	G35333	343	M35347	326
2M21191	358	G35614	343	M35357	327
2M21342	358	G93111	29, 390	M35358	327
2M21514	352, 359	G93140	29, 390	M35359	327
2M30462	372			M35375	325
2M50040	352	M M06258	422	M35376	325
2M50385	351	M07152	95	M35377	325

INDEX BY ORDER NUMBER

M35378	325	M35436	321	M36013	320
M35379	325	M35437	322	M36014	320
M35382	325	M35438	322	M36015	320
M35383	325	M35439	322	M36016	320
M35384	325	M35452	321	M36017	320
M35385	325	M35453	321	M36018	320
M35386	325	M35454	321	M36042	335
M35387	325	M35455	321	M36053	278
M35388	325	M35456	321	M36056	335
M35389	323	M35457	321	M36148	328
M35390	323	M35469	326	M36149	328
M35391	323	M35472	326	M36150	328
M35392	323	M35473	326	M36151	328
M35393	323	M35474	326	M36152	328
M35394	323	M35475	326	M36153	328
M35395	323	M35476	326	M36154	328
M35396	323	M35477	326	M36155	328
M35397	323	M35478	326	M36156	328
M35398	324	M35479	326	M36290	332
M35399	324	M35489	319	M36291	332
M35400	324	M35490	319	M36292	332
M35401	324	M35491	319	M36293	332
M35402	324	M35492	319	M36294	332
M35403	324	M35493	319	M36295	332
M35404	324	M35494	319	M36296	332
M35405	324	M35495	319	M36297	332
M35406	324	M35496	319	M36298	332
M35407	327	M35497	319	M36299	332
M35408	327	M35498	319	M36300	332
M35409	327	M35499	319	M36301	332
M35410	327	M35801	329	M36723	322
M35411	327	M35802	329	M36724	322
M35412	327	M35803	329	M36725	322
M35413	327	M35804	329	M36726	323
M35414	327	M35805	329	M36727	323
M35415	327	M35806	329	M36728	323
M35416	321	M35807	329	M36729	319
M35417	321	M35808	329	M36730	319
M35418	321	M35809	329	M36731	319
M35419	321	M35820	329	M36732	333
M35420	321	M35865	329	M36733	333
M35421	321	M35866	329	M36734	333
M35422	321	M35867	329	M36735	333
M35423	321	M36000	319	M36736	333
M35424	321	M36001	319	M36737	333
M35425	322	M36002	319	M36738	333
M35426	322	M36003	319	M36739	333
M35427	322	M36004	320	M36740	333
M35428	322	M36005	320	M36741	333
M35429	322	M36006	320	M36784	333
M35430	322	M36007	320	M36785	333
M35431	322	M36008	320	M36786	333
M35432	322	M36009	320	M36788	333
M35433	322	M36010	320, 333	M36789	333
M35434	321	M36011	320, 333	M36826	328
M35435	321	M36012	320, 333	M36827	328

INDEX BY ORDER NUMBER

	M36828	328		MP00372	19		MP00579	317
				MP00373	12		MP00580	301
MK	MK00514	298, 316		MP00374	20		MP00587	435, 341
	MK01418	291		MP00375	20		MP00589	435, 341
	MK01420	291		MP00378	36		MP00591	434, 340
	MK01422	287		MP00379	12		MP00592	434, 340
	MK01900	105		MP00380	14		MP00593	434, 340
	MK02588	101		MP00381	14		MP00595	318
	MK03140	287		MP00383	15, 386		MP00596	318
	MK03320	287		MP00384	15		MP00601	340, 434
	MK04248	43		MP00416	84, 272, 282, 375, 377, 422		MP00602	340, 434
				MP00501	284		MP00603	340, 434
MP	MP00038	278, 280		MP00502	284		MP00606	340, 434
	MP00061	297, 430		MP00503	284		MP00608	340, 434
	MP00100	273, 277, 278		MP00504	284		MP00615	341, 435
	MP00220	86		MP00505	285, 304		MP00626	340, 434
	MP00222	14		MP00506	292		MP00631	285
	MP00224	86		MP00507	308		MP00632	286
	MP00234	22, 389		MP00508	308		MP00634	302, 306, 377
	MP00300	9		MP00509	308		MP00635	302
	MP00301	16		MP00510	308		MP00636	302
	MP00303	16		MP00511	308, 423		MP00637	302
	MP00305	13		MP00512	310, 315		MP00638	289
	MP00308	22, 388		MP00514	284		MP00639	290
	MP00309	22, 388		MP00515	308		MP00646	295
	MP00310	350, 361, 386		MP00516	285, 304		MP00660	286, 424
	MP00311	350, 361, 386		MP00518	315		MP00661	292, 425
	MP00312	34		MP00530	288		MP00663	309
	MP00313	34		MP00531	288		MP00664	309
	MP00330	16		MP00532	288		MP00686	287
	MP00331	17		MP00533	288		MP00687	287
	MP00333	17, 386		MP00534	288		MP00688	291
	MP00334	11		MP00535	289, 304		MP00689	287
	MP00335	32		MP00536	293		MP00690	342
	MP00337	11		MP00537	311		MP00721	112, 141, 193
	MP00338	10		MP00538	311		MP00748	190, 226, 239
	MP00341	38		MP00539	311		MP00788	125, 180, 215, 241
	MP00343	388		MP00540	311		MP00789	125, 180, 215, 241, 413
	MP00349	9		MP00541	311		MP00790	126, 181, 215, 242
	MP00350	9		MP00542	311, 426		MP00791	126, 181, 215, 242
	MP00351	10		MP00543	312		MP00792	126, 181, 242, 365, 413
	MP00352	10		MP00544	313		MP00793	126, 181, 215, 242
	MP00353	10, 386		MP00545	290, 427		MP00794	126, 181, 242, 365, 413
	MP00355	13		MP00547	293, 428		MP00795	125, 180, 215, 241
	MP00357	36		MP00548	288		MP00796	125, 180, 215, 241
	MP00358	21		MP00549	311		MP00797	179, 240
	MP00359	21		MP00550	289		MP00799	125, 180, 215, 241
	MP00361	11		MP00551	289		MP00800	179, 240
	MP00362	11		MP00552	278, 289, 304		MP00875	116, 177, 212, 222
	MP00363	11, 386		MP00553	312		MP00877	116, 177, 212, 222
	MP00365	18		MP00555	298, 316, 431		MP00879	116, 177, 212, 222
	MP00366	18		MP00565	318		MP00881	116, 177, 212, 222
	MP00367	36		MP00566	318		MP00911	143, 195
	MP00368	36		MP00569	318		MP00912	143, 195
	MP00370	36		MP00570	300		MP00913	143, 195
	MP00371	18		MP00578	301		MP00915	143, 195

INDEX BY ORDER NUMBER

MP00916	143, 195	MP01086	94, 409	MP01553	71, 385
MP00917	143, 195	MP01087	94, 409	MP01554	71, 385
MP00918	143, 195	MP01089	94, 409	MP01555	71, 385
MP00919	143, 195	MP01095	94, 409	MP01556	71, 385
MP00921	143, 195	MP01101	276, 277	MP01559	78, 405
MP00924	144, 196	MP01102	276, 277	MP01560	78
MP00925	144, 196	MP01103	276, 278	MP01561	78, 405
MP00926	144, 196	MP01104	276	MP01562	78
MP00928	144, 196	MP01105	276	MP01570	65
MP00929	144, 196	MP01326	35	MP01571	65
MP00930	144, 196	MP01327	37, 387	MP01572	65
MP00931	144, 196	MP01328	37, 387	MP01573	65
MP00932	144, 196	MP01329	38, 387	MP01574	65
MP00934	144, 196	MP01330	38	MP01575	65
MP00943	256	MP01340	37	MP01590	72, 78, 405
MP00944	256	MP01341	37	MP01620	66
MP00953	142, 194	MP01343	38	MP01621	66
MP00989	141, 193	MP01347	35	MP01622	66
MP00990	219	MP01348	35	MP01623	66
MP00991	139, 191, 230	MP01350	35	MP01624	66
MP00992	139, 191, 230	MP01370	33	MP01625	66
MP00993	139, 191, 230	MP01371	33	MP01701	73, 402
MP00994	139, 191, 230	MP01372	46	MP01702	73, 402
MP00995	139, 191, 230	MP01390	46	MP01703	73
MP00996	139, 191, 230	MP01391	46	MP01704	73
MP00997	139, 191, 230, 415	MP01401	353, 355	MP01705	73
MP00998	139, 191, 230, 415	MP01402	360	MP01706	73
MP00999	141, 193	MP01404	360	MP01707	73
MP01000	27	MP01414	373	MP01708	73
MP01002	29, 390	MP01419	363, 366, 369	MP01781	87
MP01006	31, 391	MP01423	351, 372	MP01840	63
MP01007	31, 50	MP01424	353, 355, 372	MP01845	63
MP01008	29, 390	MP01442	88, 107	MP01850	63
MP01009	29, 390	MP01445	88	MP01855	63
MP01010	93, 408	MP01494	379, 382	MP01860	63
MP01030	93, 408	MP01495	379, 382	MP01880	361
MP01031	93	MP01496	379, 382	MP01890	63
MP01032	93, 408	MP01510	75, 403	MP01895	63
MP01035	93	MP01511	75, 403	MP01910	77, 361, 404
MP01036	93	MP01512	75, 403	MP01911	77, 361, 404
MP01037	93, 408	MP01513	75	MP01912	77
MP01060/		MP01514	75	MP01913	77
MP20108	86	MP01515	75	MP01914	77
MP01061/		MP01516	75	MP01915	77
MP20109	86	MP01521	76, 403	MP01916	77, 404
MP01062	62, 168, 205, 235	MP01522	76, 403	MP01917	77
MP01063	62, 168, 205, 235, 399	MP01523	76	MP01920	79
MP01071	94, 409	MP01524	76	MP01921	79
MP01072	94, 409	MP01531	76, 403	MP01922	79
MP01073	94, 409	MP01532	76, 403	MP01924	79
MP01074	94, 409	MP01533	76	MP01925	79
MP01075	94, 409	MP01534	76	MP01926	79
MP01076	94, 409	MP01543	76	MP01927	81
MP01078	94, 409	MP01544	76	MP01930	79
MP01079	94, 409	MP01545	76	MP01931	79
MP01085	94	MP01546	76	MP01932	79

INDEX BY ORDER NUMBER

MP01933	79	MP03196	179	MP03705	380
MP01934	79	MP03197	133, 188	MP03730	347, 372
MP01935	79	MP03198	133, 188	MP03737	347, 372
MP01940	81	MP03203	188	MP03738	355, 372
MP01942	81	MP03211	297, 430	MP03760	371
MP01943	81	MP03220	290	MP03770	371
MP01944	81	MP03221	290, 427	MP03771	371
MP01945	81	MP03222	286	MP03772	371
MP02201	141, 193	MP03223	286, 424	MP03826	378, 380
MP02260	376	MP03230	305	MP03880	379, 383
MP02400	53	MP03231	305	MP03881	379, 383
MP02405	53	MP03232	305	MP03882	379, 383
MP02590	23, 389	MP03233	305	MP03883	379, 383
MP02600	85	MP03325	281, 422	MP03884	379, 383
MP02605	24	MP03366	282, 422	MP03885	379, 383
MP02606	24	MP03373	12	MP04110	277
MP02607	24	MP03374	12	MP04200	83, 271, 374
MP02608	23, 388	MP03375	36	MP04201	83, 271
MP02650	23, 388	MP03376	12	MP04202	83, 271
MP02655	23, 389	MP03377	36	MP04203	271, 374
MP02689	23, 94, 409	MP03378	12	MP04204	271, 374
MP02710	15, 386	MP03379	12	MP04218	84, 272, 375
MP02711	15	MP03381	21	MP04219	84, 272, 375
MP02712	15	MP03383	21	MP04221	84, 272, 375
MP02713	15	MP03384	21	MP04229	271, 374
MP02720	19	MP03388	36	MP04300	254
MP02721	19	MP03401	117, 176, 211, 222	MP04301	254
MP02900	77, 361, 404	MP03402	117, 176, 211, 222	MP04305	254
MP02910	67	MP03403	117, 176, 211, 222	MP04306	254
MP02911	67	MP03404	117, 176, 211, 222	MP04310	254
MP02912	67	MP03405	117, 176, 211	MP04312	254
MP02913	67	MP03406	117, 176, 211	MP04315	254
MP02914	67	MP03407	176	MP04316	254
MP02915	67	MP03408	176	MP04320	254
MP02916	67	MP03411	118	MP04322	254
MP02917	67	MP03412	118	MP04325	254
MP02918	67	MP03413	118	MP04326	254
MP02920	80, 406	MP03414	118	MP04327	254
MP02921	80, 406	MP03415	118	MP04350	256
MP02922	80	MP03416	118	MP04351	256
MP02923	80	MP03417	118	MP04352	256
MP02924	80, 406	MP03418	118	MP04353	256
MP02925	80, 406	MP03421	211	MP04355	256
MP02926	80	MP03422	211	MP04356	256
MP02927	80	MP03423	211	MP04357	256
MP02991	179	MP03424	211	MP04358	256
MP02992	179	MP03425	211	MP04359	256
MP02993	125, 180, 215, 241	MP03426	211	MP04365	256
MP03120	177	MP03587	301	MP04366	256
MP03121	177	MP03589	301	MP04551	261, 418
MP03122	116, 177, 212	MP03597	317	MP04552	264, 420
MP03123	116, 177, 212	MP03700	378, 380	MP04553	265, 420
MP03191	127, 182, 243	MP03701	378, 380	MP04554	262
MP03192	127, 182, 243, 413	MP03702	378, 380	MP04555	263
MP03193	127, 182, 243	MP03703	378, 380	MP04556	264
MP03194	179, 240	MP03704	378, 380	MP04557	263

INDEX BY ORDER NUMBER

MP04558	261	MP05072	70	MP11002	74, 402
MP04559	267	MP05073	70	MP11003	74
MP04560	263, 419	MP05074	70	MP11004	74
MP04561	261	MP05080	72	MP11005	74
MP04562	267	MP05114	178	MP11006	74
MP04563	262	MP05115	120, 178	MP11007	74
MP04564	267	MP05120	114, 174, 210, 223	MP11008	74
MP04565	262	MP05121	114, 174, 210, 223	MP12060	64
MP04566	262	MP05122	114, 174, 210, 223	MP12061	64
MP04567	264, 420	MP05123	114, 174, 210, 223	MP12062	64
MP04568	268	MP05124	114, 174, 210	MP12063	64
MP04569	263, 419	MP05125	114, 174, 210	MP13000	269, 421
MP04570	265	MP05126	175	MP13011	269, 421
MP04571	262	MP05127	175	MP13012	269, 421
MP04572	267	MP05128	175	MP13013	270
MP04573	264	MP05129	175	MP13014	270
MP04575	264	MP05130	115, 174, 223	MP13015	270
MP04576	266	MP05131	115, 174, 223	MP13016	270
MP04577	266	MP05132	115, 174, 223	MP13017	270
MP04578	266	MP05133	115, 174, 223	MP13019	270
MP04579	266	MP05134	115, 174	MP13020	83, 271, 376
MP04580	264	MP05135	115, 174	MP13021	83, 271, 376
MP04581	268	MP05305	84, 272	MP13040	262
MP04711	68	MP05507	27, 82	MP13041	267
MP04712	68	MP05511	82	MP15001	170, 206, 252
MP04713	68	MP05512	82	MP15002	170, 206, 252
MP04714	68	MP05513	82	MP15003	170, 206, 252, 417
MP04721	68	MP05601	82	MP15004	170, 206, 252
MP04722	68	MP05700	343	MP15005	170, 206, 252
MP04723	68	MP05701	343	MP15006	170, 206, 252
MP04724	68	MP05730	58	MP15007	170, 206, 252, 417
MP04731	68	MP05735	58	MP15008	170, 206, 252
MP04732	68	MP05750	58	MP15009	170, 206, 252
MP04733	68	MP05755	56	MP15010	171, 207, 253
MP04734	68	MP05770	56	MP15011	171, 207, 253
MP04741	69	MP05785	57	MP15012	171, 207, 253
MP04742	69	MP05790	57	MP15013	171, 207, 253
MP04743	69	MP05795	57	MP16002	28, 389
MP04744	69	MP05800	55	MP16003	28, 389
MP04751	69	MP05801	55	MP16004	28, 389
MP04752	69	MP05805	55	MP16005	28, 389
MP04753	69	MP05810	55	MP16017	14
MP04754	69	MP05815	55	MP16018	14
MP04800	13	MP05820	55, 398	MP16019	14
MP04802	13, 386	MP05845	58, 398	MP16020	15
MP04803	35	MP07101	27	MP16021	15, 386
MP04804	35	MP07102	27	MP16028	16
MP05051	70	MP07106	27	MP16029	16
MP05052	70	MP07107	27	MP16030	16
MP05053	70	MP07108	27	MP16031	17
MP05054	70	MP07109	27	MP16032	17, 386
MP05061	70	MP07130	27	MP16033	18
MP05062	70	MP07131	27	MP16034	18
MP05063	70	MP07150	27, 388	MP16035	18
MP05064	70	MP07151	27, 388	MP16036	19
MP05071	70	MP11001	74, 402	MP16037	20

INDEX BY ORDER NUMBER

MP16038	20	MS16230	113	MS33730	130, 184, 245, 413
MP16039	19	MS16232	113	MS33731	128, 216, 412
MP16040	19	MS16233	113	MS33732	128, 216, 412
MP16041	21	MS16444	122, 214	MS33733	128, 412
MP16042	21	MS16445	122, 214	MS33736	183, 244
MP16044	21	MS16446	123, 411	MS33737	183, 244
MP16045	36	MS16447	123, 411	MS33738	183, 185, 186
MP16049	21	MS16448	122, 214	MS33739	183, 185, 186
MP16050	21	MS16449	122, 214	MS33740	129, 183, 216, 246
MP16051	36	MS16492	120	MS33741	129, 183, 216, 246
MP16052	36	MS16916	161	MS34209	208
MP16090	349, 372	MS16918	161	MS34436	124, 216
MP16297	349, 357	MS16919	161	MS34437	124, 216
MP17008	11	MS16920	161	MS34677	137, 190
MP17009	9	MS17330	121, 124, 134, 213, 218, 410	MS34678	137, 190
MP17010	9	MS17522	124, 215, 412	MS40021	209
MP17011	10, 11	MS18683	121, 124, 134, 213, 218, 410	MS40067	129, 183, 246
MP17012	10	MS20093	111	MS40070	129, 183, 246
MP17090	27, 82	MS20094	111	MS40073	129, 183, 246
MP17128	28, 392	MS20725	148, 198	MS40076	186
MP17129	28	MS20783	198	MS40079	186
MP20020	53	MS20979	190, 226, 239	MS40082	186
MP20022	331	MS22147	149, 199	MS40085	186
MP20023	331	MS22148	149, 199	MS40088	131, 184, 247
MP20024	331	MS22259	168, 205	MS40091	131, 184, 247
MP20029	331, 332	MS22532	149, 161, 199	MS40094	131, 184, 247
MP20030	331, 332	MS22533	149, 199	MS40243	140, 192, 229, 231
MP20050	53, 397	MS22534	149, 199	MS40317	208
MP20110	88	MS22535	149, 199	MS40318	208
MP20370	89	MS22832	178	MS40319	208
MP20371	89	MS22905	219	MS40320	208
MP20372	90	MS23133	161	MS40321	208
MP20373	90	MS23598	154	MS40322	208
MP20374	90	MS23897	198	MS40323	208
MP20375	90	MS23898	198		
MP20376	90	MS24303	124, 215, 412	MU MU03876	354, 356, 363, 367, 369
MP20377	89	MS25951	178	MU04366	363, 366, 369
MP20378	89	MS26122	146	MU04472	364, 367, 369
MP20379	90	MS26123	146	MU04771	366
		MS26125	146	MU06008	364, 368, 370
MS MS13218	166	MS29468	165	MU06011	363, 368, 370
MS13235	122, 214	MS32150	209	MU06570	356, 363, 367, 369
MS14679	120, 178	MS32151	209	MU06933	363, 366, 369
MS14680	120, 178	MS32153	209	MU06943	347, 349, 351, 353,
MS14796	164	MS32154	208		355, 357, 360, 363,
MS15053	166	MS32156	209		366, 369, 373
MS15085	166	MS32282	161	MU06944	347, 349, 351, 353,
MS15086	166	MS32289	208		355, 357, 360, 363,
MS15087	166	MS32290	208		366, 369, 373
MS15893	146	MS32291	208	MU06951	363, 366, 369
MS15894	147	MS32318	208	MU10841	350, 361
MS15895	147	MS33239	193	MU10900	361
MS15896	147	MS33726	130, 184, 216, 245	MU10902	362
MS16160	113	MS33727	130, 184, 216, 245	MU10915	362
MS16161	113	MS33728	130, 184, 245	MU10916	362
MS16229	113	MS33729	130, 184, 245, 413	MU10917	362

INDEX BY ORDER NUMBER

MU10918	350, 354, 356, 362
MU11225	360
MU11228	360
MU12504	354, 356
MU12533	355
MU12551	355
MU12609	354, 356
MU17191	356
MU21120	354, 356
MU24401	364
MU24402	364
MU24409	364
MU24903	354, 356
MU26041	353
MU26042	353

MX MX00001	99
MX00004	100
MX01004	135, 226, 239
MX01005	135, 226, 239
MX01049	107
MX01050	107, 352
MX10209	99
MX11000	347, 349, 351, 357
MX11001	347, 349, 351, 357
MX17012	351
MX17015	348, 352
MX17018	348, 352
MX23023	297, 430
MX23025	297, 430
MX23026	297, 430
MX37021	306
MX50004	100
MX50050	99
MX50065	136, 190, 218
MX50066	136, 190, 218
MX50067	136, 190, 411
MX50068	136, 190, 218, 411
MX50069	137
MX50070	137, 226, 239
MX50071	136, 190, 218
MX50081	101
MX50090	101
MX50100	100
MX50115	101
MX50120	100

INDEX BY PRODUCT NAME

1	1.5 m spiral hose and M12×1 screw connector	2M86291	95
	1.5 m spiral hose and probe for central supply system	2M86292	95
3	3.0 m spiral hose and probe for central supply system	2M86293	95
A	Absorber adapter CLIC, for Dräger anesthesia workstations	MX50090	101
	Adapter 9/16" to 15/22mm	MP05305	84, 272
	Adapter 9/16" to M34×1.5 for Dräger reusable humidifier	MP03325	281, 422
	Adapter AIR NF/NIST	M35338	335
	Adapter AIR NIST/DIN	M32495	335
	Adapter AIR-O ₂ /O ₂ NIST/DIN	M36042	335
	Adapter AIR, DISS-NIST	M34877	335
	Adapter AIR/O ₂ NIST/DIN	M32497	335
	Adapter cable for Infinity® PiCCO SmartPod®, for 10-pin connector	MS16920	161
	Adapter cable RD SET to LNCS, for connection of RD SET sensors to Infinity® M300 and Masimo® SET Pod via the LNCS intermediate cables MS24303 or MS17522	MS33732	128, 216, 412
	Adapter cable RD SET to LNOP, for connection of RD SET sensors to Infinity® M300 and Masimo® SET Pod via the LNOP intermediate cables MS18680 and MS17041, as well as to MultiMed® 5/6/12, NeoMed Pods and MultiMed® Plus/Plus OR via the according intermediate cables MS13926/7492601 or MS20163	MS33731	128, 216, 412
	Adapter cable RD SET to MultiMed®, for connection of RD SET sensors to MultiMed® 5/6/12, NeoMed Pods as well as MultiMed® Plus/Plus OR via the according intermediate cables 3368433/3375834 or MS18683/MS17330	MS33733	128, 412
	Adapter cable, 7-pin connector for temperature probe with ¼" connector	5198333	139
	Adapter DISS M/QC-AIR-Chemetron	4119066	334
	Adapter DISS M/QC-AIR-Ohmeda	4119071	334
	Adapter DISS M/QC-AIR-Pur Ben	4119076	334
	Adapter DISS M/QC-AIR-Schrader	4119081	334
	Adapter DISS M/QC-EVAC-Chemetron	4119069	334
	Adapter DISS M/QC-EVAC-Ohmeda	4119074	334
	Adapter DISS M/QC-EVAC-Pur Ben	4119079	334
	Adapter DISS M/QC-N ₂ O-Chemetron	4119065	334
	Adapter DISS M/QC-N ₂ O-Ohmeda	4119070	334
	Adapter DISS M/QC-N ₂ O-Pur Ben	4119075	334
	Adapter DISS M/QC-N ₂ O-Schrader	4119080	334
	Adapter DISS M/QC-O ₂ -Chemetron	4119068	334
	Adapter DISS M/QC-O ₂ -Ohmeda	4119073	334
	Adapter DISS M/QC-O ₂ -Pur Ben	4119078	334
	Adapter DISS M/QC-O ₂ -Schrader	4119083	334
	Adapter DISS M/QC-VAC-Chemetron	4119067	334
	Adapter DISS M/QC-VAC-Ohmeda	4119072	334
	Adapter DISS M/QC-VAC-Pur Ben	4119077	334
	Adapter DISS M/QC-VAC-Schrader	4119082	334
	Adapter for Ambu Mark III	2120488	95
	Adapter for Combibag	2M86913	95
	Adapter for Laerdal	2120496	95
	Adapter for Laerdal	5790124	95
	Adapter for Oxidem for M12×1 (for existing hose systems)	8410547	95
	Adapter for Resutator 2000/MR-100	2120461	95
	Adapter N ₂ O NF/NIST	M35337	335
	Adapter N ₂ O NIST/DIN	M32494	335
	Adapter N ₂ O NIST/DIN	M35057	335
Adapter N ₂ O, DISS-NIST	M34876	335	
Adapter NIST to DIN CS-Coupling AIR	M35058	280, 335	

INDEX BY PRODUCT NAME

Adapter NIST to DIN CS-Coupling O ₂	M35056	280, 335
Adapter O ₂ NF/NIST	M35336	335
Adapter O ₂ NIST/DIN	M32493	335
Adapter O ₂ , DISS-NIST	M34875	335
Adapter straight, reusable, A, 18 mm/11 mm	8403076	51
Adapter VAC NIST/DIN	M32496	335
Adapter VAC, DISS-NIST	M34878	335
Adapter, for disposable sterilized water systems to AIR flowmeter, M32×1.5 to 9/16"	MP03366	282, 422
Adapter, for disposable sterilized water systems to O ₂ flowmeters, M34×1.5 to 9/16"	2M86680	282, 422
Adapter, reusable, 22 mm/11 mm, N/P	8409897	51
Adapter, reusable, 22 mm/15 mm, 90°	8290274	40
Adapter, reusable, 22 mm/22 mm	M25647	40, 42, 43, 44, 51, 393
Adapter, reusable, 22 mm/22 mm, 90°	M25649	51
Adapter, to connect Dräger humidifier to mixing block, M15×1 to M34×1.5	2M85853	281
Adapter, to connect e.g a tube or nasal cannula to O ₂ flowmeter, M34×1.5 to 6 mm nozzle	M06258	422
Adapter, with 6 mm nozzle, e.g., for disposable breathing masks, 9/16" connection	MP00416	84, 272, 282, 375, 377, 422
Adhesive Pad LNCS TF-I, 20 pcs.	MP03198	133, 188
Adhesive wrap for Nellcor® Oxiband™ A/N, 100 pcs.	8201012	135, 411
Adhesive wrap for Nellcor® Oxiband™ P/I, 100 pcs.	2170221	135
Aerogen Solo Starter Kit, ES, IT, PT, GR, TK	MP01073	94, 409
Aerogen Solo Starter Kit, FR, NL, IT, GE	MP01072	94, 409
Aerogen Solo Starter Kit, Int.	MP01071	94, 409
Aerogen Solo Starter Kit, SW, NO, FI, DK	MP01074	94, 409
Aerogen Solo Starter Kit, UK	MP01075	94, 409
Aerogen Solo Starter Kit, US	MP01076	94, 409
Aerogen Solo USB Controller Starter Kit NE	MP01089	94, 409
Aerogen T-Adapter, 11.5 mm, neonates, 10 pcs.	MP02689	23, 94, 409
Aerogen USB Controller Replacement w/o AC/DC adapter	MP01095	94, 409
AIR compact flowmeter Clickflow, DIN, 0.1 – 1 l/min, 6 mm nozzle	MP13019	270
AIR compact flowmeter Clickflow, DIN, 1 – 16 l/min, 6 mm nozzle	MP13016	270
AIR compact flowmeter Clickflow, DIN, 1 – 16 l/min, 9/16"	MP13015	270
AIR compact flowmeter Clickflow, rail, 1 – 16 l/min, 9/16"	MP13017	270
AIR CS hose 0.9 m (2.95 ft) NF EN NF	M36724	322
AIR CS hose 0.9 m (2.95 ft) NIST EN BS	M36727	323
AIR CS hose 0.9 m (2.95 ft) NIST EN DIN	M36730	319
AIR CS hose 0.9 m (2.95 ft) NIST GB DIN	M36786	333
AIR CS hose 1.5 m (4.9 ft) DISS US KAWA	M35803	329
AIR CS hose 1.5 m (4.9 ft) DISSB EN HAND	M36296	332
AIR CS hose 1.5 m (4.9 ft) NF EN NF	M35429	322
AIR CS hose 1.5 m (4.9 ft) NIST BS straight	M35402	324
AIR CS hose 1.5 m (4.9 ft) NIST EN AS	M36150	328
AIR CS hose 1.5 m (4.9 ft) NIST EN BS	M35393	323
AIR CS hose 1.5 m (4.9 ft) NIST EN Carba	M35475	326
AIR CS hose 1.5 m (4.9 ft) NIST EN DIN	M35495	319
AIR CS hose 1.5 m (4.9 ft) NIST EN NF	M35420	321
AIR CS hose 1.5 m (4.9 ft) NIST EN SS	M35383	325
AIR CS hose 1.5 m (4.9 ft) NIST EN UNI	M35411	327
AIR CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36010	320, 333
AIR CS hose 1.5 m (4.9 ft) NIST NEUTRAL NF	M35452	321
AIR CS hose 3 m (9.8 ft) DISS US KAWA	M35802	329

INDEX BY PRODUCT NAME

AIR CS hose 3 m (9.8 ft) DISSB EN HAND	M36297	332
AIR CS hose 3 m (9.8 ft) NF EN NF	M35430	322
AIR CS hose 3 m (9.8 ft) NIST BS straight	M35403	324
AIR CS hose 3 m (9.8 ft) NIST EN AS	M36149	328
AIR CS hose 3 m (9.8 ft) NIST EN BS	M35394	323
AIR CS hose 3 m (9.8 ft) NIST EN Carba	M35476	326
AIR CS hose 3 m (9.8 ft) NIST EN DIN	M35496	319
AIR CS hose 3 m (9.8 ft) NIST EN NF	M35421	321
AIR CS hose 3 m (9.8 ft) NIST EN SS	M35384	325
AIR CS hose 3 m (9.8 ft) NIST EN UNI	M35412	327
AIR CS hose 3 m (9.8 ft) NIST GB BS	M36736	333
AIR CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36011	320, 333
AIR CS hose 3 m (9.8 ft) NIST NEUTRAL NF	M35453	321
AIR CS hose 5 m (16.4 ft) DISS US KAWA	M35801	329
AIR CS hose 5 m (16.4 ft) DISSB EN HAND	M36298	332
AIR CS hose 5 m (16.4 ft) NF EN NF	M35343	322
AIR CS hose 5 m (16.4 ft) NIST BS straight	M35340	324
AIR CS hose 5 m (16.4 ft) NIST EN AS	M36148	328
AIR CS hose 5 m (16.4 ft) NIST EN BS	M35335	323
AIR CS hose 5 m (16.4 ft) NIST EN Carba	M35346	326
AIR CS hose 5 m (16.4 ft) NIST EN DIN	M35497	319
AIR CS hose 5 m (16.4 ft) NIST EN NF	M35332	321
AIR CS hose 5 m (16.4 ft) NIST EN SS	M35385	325
AIR CS hose 5 m (16.4 ft) NIST EN UNI	M35358	327
AIR CS hose 5 m (16.4 ft) NIST GB BS	M36737	333
AIR CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36012	320, 333
AIR CS hose 5 m (16.4 ft) NIST NEUTRAL NF	M35454	321
AIR CS hose 8 m (26,25 ft) DISS US KAWA	M35865	329
AIR CS hose DISS Hand 10L	4199594	330
AIR CS hose DISS Hand 15L	4199593	330
AIR CS hose DISS N/DISS N 10L	4199605	330
AIR CS hose DISS N/DISS N 15L	4199604	330
AIR CS hose DISS N/DISS N 20L	4199712	330
AIR flowmeter, 16 l/min, CS, CARBA, ISO, M32×1.5	MP04564	267
AIR flowmeter, 16 l/min, CS, DIN, ISO and neutral color, M32×1.5	MP04559	267
AIR flowmeter, 16 l/min, CS, long cs plug (6 cm; 2.4 in), DIN, ISO and neutral color, M32×1.5	MP04562	267
AIR flowmeter, 16 l/min, CS, Swedish Standard, M32x1.5	MP13041	267
AIR flowmeter, 16 l/min, rail, NIST, M32×1.5	MP04568	268
AIR flowmeter, 32 l/min, CS, DIN, ISO and neutral color, M32×1.5	MP04572	267
AIR flowmeter, 32 l/min, rail, NIST, M32×1.5	MP04581	268
Air intake filter for Caleo® and Babyleo® TN500, reusable, 20 pcs.	MX17015	348, 352
Air intake filter for Isolette® C2000, 8000 and 8000 plus, reusable, 4 pcs.	MU12504	354, 356
Air intake filter for TI 500 and TI 500 Globe-Trotter, disposable, 6 pcs.	MU04472	364, 367, 369
AIR/O ₂ CS hose 1.5 m (4.9 ft) NIST EN DIN	M36001	319
AIR/O ₂ CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36016	320
AIR/O ₂ CS hose 3 m (9.8 ft) NIST EN DIN	M36002	319
AIR/O ₂ CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36017	320
AIR/O ₂ CS hose 5 m (16.4 ft) NIST EN DIN	M36003	319
AIR/O ₂ CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36018	320
Airway connector ErgoStar® AC 90, disposable, double-swivel elbow w/ double cap, 100 pcs.	MP01890	63
Airway connector ErgoStar® AC 95, disposable, w/ Luer Lock sampling port, 100 pcs.	MP01895	63

INDEX BY PRODUCT NAME

Anesthesia mask ComfortStar [®] , disposable, hook ring, bubble gum, size 1, neonatal, 20 pcs.	MP01531	76, 403
Anesthesia mask ComfortStar [®] , disposable, hook ring, bubble gum, size 2, baby, 20 pcs.	MP01532	76, 403
Anesthesia mask ComfortStar [®] , disposable, hook ring, bubble gum, size 3, infant, 20 pcs.	MP01533	76
Anesthesia mask ComfortStar [®] , disposable, hook ring, bubble gum, size 4, adult S, 20 pcs.	MP01534	76
Anesthesia mask ComfortStar [®] , disposable, hook ring, mint, size 3, infant, 20 pcs.	MP01543	76
Anesthesia mask ComfortStar [®] , disposable, hook ring, mint, size 4, adult S, 20 pcs.	MP01544	76
Anesthesia mask ComfortStar [®] , disposable, hook ring, mint, size 5, adult M, 20 pcs.	MP01545	76
Anesthesia mask ComfortStar [®] , disposable, hook ring, mint, size 6, adult L, 20 pcs.	MP01546	76
Anesthesia mask ComfortStar [®] , disposable, hook ring, standard, size 0, premature baby, 20 pcs.	MP01510	75, 403
Anesthesia mask ComfortStar [®] , disposable, hook ring, standard, size 1, neonatal, 20 pcs.	MP01511	75, 403
Anesthesia mask ComfortStar [®] , disposable, hook ring, standard, size 2, baby, 20 pcs.	MP01512	75, 403
Anesthesia mask ComfortStar [®] , disposable, hook ring, standard, size 3, infant, 20 pcs.	MP01513	75
Anesthesia mask ComfortStar [®] , disposable, hook ring, standard, size 4, adult S, 20 pcs.	MP01514	75
Anesthesia mask ComfortStar [®] , disposable, hook ring, standard, size 5, adult M, 20 pcs.	MP01515	75
Anesthesia mask ComfortStar [®] , disposable, hook ring, standard, size 6, adult L, 20 pcs.	MP01516	75
Anesthesia mask ComfortStar [®] , disposable, hook ring, strawberry, size 1, neonatal, 20 pcs.	MP01521	76, 403
Anesthesia mask ComfortStar [®] , disposable, hook ring, strawberry, size 2, baby, 20 pcs.	MP01522	76, 403
Anesthesia mask ComfortStar [®] , disposable, hook ring, strawberry, size 3, infant, 20 pcs.	MP01523	76
Anesthesia mask ComfortStar [®] , disposable, hook ring, strawberry, size 4, adult S, 20 pcs.	MP01524	76
Anesthesia mask LiteStar [®] , disposable, w/o hook ring, size 1, infant S, 30 pcs.	MP11001	74, 402
Anesthesia mask LiteStar [®] , disposable, w/o hook ring, size 2, infant L, 30 pcs.	MP11002	74, 402
Anesthesia mask LiteStar [®] , disposable, w/o hook ring, size 3, child, 30 pcs.	MP11003	74
Anesthesia mask LiteStar [®] , disposable, w/o hook ring, size 4, adult S, 30 pcs.	MP11004	74
Anesthesia mask LiteStar [®] , disposable, w/o hook ring, size 5, adult M, 30 pcs.	MP11005	74
Anesthesia mask LiteStar [®] , disposable, w/o hook ring, size 6, adult L, 30 pcs.	MP11006	74
Anesthesia mask LiteStar [®] , disposable, w/o hook ring, size 7, adult XL, 30 pcs.	MP11007	74
Anesthesia mask SediStar [®] , disposable, w/o hook ring, size 1, infant S, 30 pcs.	MP01701	73, 402
Anesthesia mask SediStar [®] , disposable, w/o hook ring, size 2, infant L, 30 pcs.	MP01702	73, 402
Anesthesia mask SediStar [®] , disposable, w/o hook ring, size 3, child, 30 pcs.	MP01703	73
Anesthesia mask SediStar [®] , disposable, w/o hook ring, size 4, adult S, 30 pcs.	MP01704	73
Anesthesia mask SediStar [®] , disposable, w/o hook ring, size 5, adult M, 30 pcs.	MP01705	73
Anesthesia mask SediStar [®] , disposable, w/o hook ring, size 6, adult L, 30 pcs.	MP01706	73
Anesthesia mask SediStar [®] , disposable, w/o hook ring, size 7, adult XL, 30 pcs.	MP01707	73
Anesthesia set Infinity ID [®] , disposable, coaxial, GSL, 2.3 m/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16051	36
Anesthesia set Infinity ID [®] , disposable, coaxial, water trap, GSL, 2.1 m/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16052	36

INDEX BY PRODUCT NAME

Anesthesia set Infinity® ID, disposable, coaxial, GSL, 2.3 m/1.2 m, 10 pcs., not made with natural rubber latex	MP00368	36
Anesthesia set Infinity® ID, disposable, coaxial, water trap, 2 L, 2.1/1.2 m, 10 pcs., not made with natural rubber latex	MP03388	36
Anesthesia set Infinity® ID, disposable, coaxial, water trap, 2 L, 2.1/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16045	36
Anesthesia set Infinity® ID, disposable, coaxial, water trap, GSL, 2.1 m/1.2 m, 10 pcs., not made with natural rubber latex	MP00370	36
Anesthesia set VentStar, disposable, coaxial, 2 L, 1.9 m/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16041	21
Anesthesia set VentStar, disposable, coaxial, 2 L, 2.3 m/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16042	21
Anesthesia set VentStar®, coaxial, water trap, GSL, 2 L, 2.1 m /1.2 m, 10 pcs., not made with natural rubber latex	MP00359	21
Anesthesia set VentStar®, coaxial, water trap, GSL, 2 L, 2.1 m /1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16050	21
Anesthesia set VentStar®, disposable, basic, 2 L, 1.8 m/1.5 m, 25 pcs., not made with natural rubber latex	MP00365	18
Anesthesia set VentStar®, disposable, basic, 2 L, 1.8 m/1.5 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex	MP16033	18
Anesthesia set VentStar®, disposable, basic, 2 L, 2.5 m/1.5 m, 10 pcs., not made with natural rubber latex	MP00366	18
Anesthesia set VentStar®, disposable, basic, 2 L, 2.5 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16034	18
Anesthesia set VentStar®, disposable, coaxial, 2 L, 1.9 m/1.2 m, 10 pcs., not made with natural rubber latex	MP03384	21
Anesthesia set VentStar®, disposable, coaxial, 2 L, 2.3 m/1.2 m, 10 pcs., not made with natural rubber latex	MP03383	21
Anesthesia set VentStar®, disposable, coaxial, GSL, 2 L, 2.3 m /1.2 m, 10 pcs., not made with natural rubber latex	MP00358	21
Anesthesia set VentStar®, disposable, coaxial, GSL, 2 L, 2.3 m /1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16049	21
Anesthesia set VentStar®, disposable, coaxial, water trap, 2 L, 2.1 m/1.2 m, 10 pcs., not made with natural rubber latex	MP03381	21
Anesthesia set VentStar®, disposable, coaxial, water trap, 2 L, 2.1/1.2 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16044	21
Anesthesia set VentStar®, disposable, N, 1.8 m/1.1 m, 25 pcs., not made with natural rubber latex	MP00333	17, 386
Anesthesia set VentStar®, disposable, N, 1.8 m/1.1 m, Bag cone 22OD, 25 pcs., not made with natural rubber latex	MP16032	17, 386
Anesthesia set VentStar®, disposable, P, w/o Luer Lock, 1.5 m/1.1 m, 25 pcs., not made with natural rubber latex	MP00331	17
Anesthesia set VentStar®, disposable, P, w/o Luer Lock, 1.5 m/1.1 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex	MP16031	17
Anesthesia set VentStar®, disposable, P, water trap, 1 L, 1.8 m/1.1 m, 10 pcs., not made with natural rubber latex	MP00374	20
Anesthesia set VentStar®, disposable, P, water trap, 1 L, 1.8 m/1.1 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16037	20
Anesthesia set VentStar®, disposable, P, water trap, w/o Luer Lock, 1 L, 1.8 m/1.1 m, 10 pcs., not made with natural rubber latex	MP00375	20
Anesthesia set VentStar®, disposable, P, water trap, w/o Luer Lock, 1 L, 1.8 m/1.1 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16038	20
Anesthesia set VentStar®, disposable, w/o Luer Lock, 1.8 m/0.8 m, 25 pcs., not made with natural rubber latex	MP00330	16
Anesthesia set VentStar®, disposable, w/o Luer Lock, 1.8 m/0.8 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex	MP16030	16
Anesthesia set VentStar®, disposable, water trap, 2 L, 1.8 m/1.5 m, 10 pcs., not made with natural rubber latex	MP00372	19

INDEX BY PRODUCT NAME

Anesthesia set VentStar®, disposable, water trap, 2 L, 1.8 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16036	19
Anesthesia set VentStar®, disposable, water trap, 2 L, 2.8 m/1.5 m, 10 pcs., not made with natural rubber latex	MP02721	19
Anesthesia set VentStar®, disposable, water trap, 2 L, 2.8 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16040	19
Anesthesia set VentStar®, disposable, water trap, w/o Luer Lock, 2 L, 1.8 m/1.5 m, 10 pcs., not made with natural rubber latex	MP00371	18
Anesthesia set VentStar®, disposable, water trap, w/o Luer Lock, 2 L, 1.8 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16035	18
Anesthesia set VentStar®, disposable, water trap, w/o Luer Lock, 2 L, 2.8 m/1.5 m, 10 pcs., not made with natural rubber latex	MP02720	19
Anesthesia set VentStar®, disposable, water trap, w/o Luer Lock, 2 L, 2.8 m/1.5 m, Bag Cone 22OD, 10 pcs., not made with natural rubber latex	MP16039	19
Anesthesia set, disposable, basic, 1.5 m/0.8 m, 25 pcs., not made with natural rubber latex	MP00301	16
Anesthesia set, disposable, basic, 1.5 m/0.8 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex	MP16028	16
Anesthesia set, disposable, extendable, 1.8 m/0.8 m, 25 pcs., not made with natural rubber latex	MP00303	16
Anesthesia set, disposable, extendable, 1.8 m/0.8 m, Bag Cone 22OD, 25 pcs., not made with natural rubber latex	MP16029	16
APC NIBP Cuff, adult, 23 cm – 33 cm, reusable, Vista 120 SC	2600933	249
APC NIBP Cuff, child, 12 cm – 19 cm, reusable, Vista 120 SC	2600930	249
APC NIBP Cuff, large adult, 31 cm – 40 cm, reusable, Vista 120	2600934	249
APC NIBP Cuff, small adult, 17 cm – 25 cm, reusable, Vista 120 SC	2600932	249
Aquapor H300 – Configuration	MP01000	27
Aquapor H300 Basic unit 115V	MP07102	27
Aquapor H300 Basic unit 230V	MP07101	27
B BabyFlow® plus, disposable, size L, 20 pcs.	MP03702	378, 380
BabyFlow® plus, disposable, size M, 20 pcs.	MP03701	378, 380
BabyFlow® plus, disposable, size S, 20 pcs.	MP03700	378, 380
Bacteria filter IBF absorber, only use with Zeus®/Zeus® IE, 5 pcs.	MK02588	101
Bacterial and viral filter system VarioSafe, for VarioAir®, VarioVac®, VarioO2, disposable, 10 pcs.	MP00555	298, 316, 431
Bacterial filter 767, 5 pcs.	6723976	298, 316, 431
Bacterial filter 767, disposable, for ejectors (model range until 2010)	CH00102	298, 316
Bacterial filter disposable, 20 pcs.	MK00514	298, 316
Bacterial filter for Dräger patient monitors, disposable	8402868	25, 172
Bacterial filter, disposable, for Evac 40, 10 pcs.	MX37021	306
Basket 150, 15 × 20 × 10 cm, e.g., for disposable gloves	M26146	339, 433
Basket 300, 15 × 28 × 10 cm, e.g., for 300 mm catheters	M26145	339, 433
Basket 600, 15 × 48 × 10 cm, e.g., for 600 mm catheters	M25121	339, 433
Bed clip for cables	MS15087	166
BISx adapter cable Vista 120	2607560	237
BISx pod Vista 120	2607561	237
Block for EEG referential lead conversion	5954859	165
BlueSet Anesthesia (N)	M27542	44, 393
BlueSet Anesthesia (P/A)	M30146	42
BlueSet Anesthesia water trap (P/A)	MK04248	43
BlueSet Ventilation (N)	8412344	39, 392
BlueSet Ventilation (P/A)	8413146	40
BlueSet Ventilation HME (P/A)	8412860	41

INDEX BY PRODUCT NAME

Breathing bag set Infinity® ID, disposable, 2 L, 1.2 m, 25 pcs., not made with natural rubber latex	MP00341	38
Breathing bag set Infinity® ID, disposable, 2 L, 1.5 m, 25 pcs., not made with natural rubber latex	MP01330	38
Breathing bag set Infinity® ID, disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP01329	38, 387
Breathing bag set Infinity® ID, disposable, P, 1 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP01343	38
Breathing bag set VentStar® (Cone 22OD), disposable, 2 L, 0.8 m, 25 pcs., not made with natural rubber latex	MP16017	14
Breathing bag set VentStar® (Cone 22OD), disposable, 2 L, 1.2 m, 25 pcs., not made with natural rubber latex	MP16018	14
Breathing bag set VentStar® (Cone 22OD), disposable, 2 L, 1.5 m, 25 pcs., not made with natural rubber latex	MP16019	14
Breathing bag set VentStar® (Cone 22OD), disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP16021	15, 386
Breathing bag set VentStar® (Cone 22OD), disposable, P, 1 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP16020	15
Breathing bag set VentStar®, disposable, 2 L, 0.8 m, 25 pcs., not made with natural rubber latex	MP00222	14
Breathing bag set VentStar®, disposable, 2 L, 1.2 m, 25 pcs., not made with natural rubber latex	MP00380	14
Breathing bag set VentStar®, disposable, 2 L, 1.5 m, 25 pcs., not made with natural rubber latex	MP00381	14
Breathing bag set VentStar®, disposable, N, 0.5 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP00383	15, 386
Breathing bag set VentStar®, disposable, P, 1 L, 1.1 m, 25 pcs., not made with natural rubber latex	MP00384	15
Breathing bag, silicone, reusable, 0.5 L, not made with natural rubber latex	2166054	44, 52, 393, 397
Breathing bag, silicone, reusable, 1.5 L, not made with natural rubber latex	2166127	52, 397
Breathing bag, silicone, reusable, 2.3 L, not made with natural rubber latex	2166062	42, 43, 52
Breathing bag, silicone, reusable, 3 L, not made with natural rubber latex	2166070	52
Breathing circuit for active humidification, disposable, neonate, 20 pcs., not made with natural rubber latex	MU24409	364
Breathing circuit Infinity® ID, disposable, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP01350	35
Breathing circuit Infinity® ID, disposable, basic, 2.5 m, 25 pcs., not made with natural rubber latex	MP01348	35
Breathing circuit Infinity® ID, disposable, coaxial, 1.9 m, 20 pcs., not made with natural rubber latex	MP03375	36
Breathing circuit Infinity® ID, disposable, coaxial, 2.3 m, 20 pcs., not made with natural rubber latex	MP03377	36
Breathing circuit Infinity® ID, disposable, coaxial, GSL, 2.3 m, 20 pcs., not made with natural rubber latex	MP00357	36
Breathing circuit Infinity® ID, disposable, extendable, 2.2 m, 25 pcs., not made with natural rubber latex	MP01347	35
Breathing circuit Infinity® ID, disposable, N, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP01327	37, 387
Breathing circuit Infinity® ID, disposable, N, heated, 1.2 m, 10 pcs., not made with natural rubber latex	MP00343	388
Breathing circuit Infinity® ID, disposable, N, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	MP01328	37, 387
Breathing circuit Infinity® ID, disposable, P, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP01340	37
Breathing circuit Infinity® ID, disposable, P, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	MP01341	37
Breathing circuit Infinity® ID, disposable, smoothbore, 1.6 m, 25 pcs., not made with natural rubber latex	MP04803	35

INDEX BY PRODUCT NAME

Breathing circuit Infinity® ID, disposable, smoothbore, 3 m, 10 pcs., not made with natural rubber latex	MP04804	35
Breathing circuit Infinity® ID, disposable, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	MP01326	35
Breathing circuit Infinity® ID, disposable, water trap, 2.1 m, 10 pcs., not made with natural rubber latex	MP00378	36
Breathing circuit Infinity® ID, disposable, water trap, GSL, 2.1 m, 10 pcs., not made with natural rubber latex	MP00367	36
Breathing circuit Vent Set, disposable, basic, 1.5 m, 25 pcs., not made with natural rubber latex	MP00300	9
Breathing circuit Vent Set, disposable, extendable, 1.8 m, 25 pcs., not made with natural rubber latex	MP00305	13
Breathing circuit VentStar® Carina® ExpV, disposable, 1.5 m, 10 pcs., not made with natural rubber latex	MP00313	34
Breathing circuit VentStar® Carina® ExpV, disposable, 2 m, 10 pcs., not made with natural rubber latex	5704962	34
Breathing circuit VentStar® Carina® LeakV, disposable, 1.5 m, 25 pcs., not made with natural rubber latex	MP00312	34
Breathing circuit VentStar® Helix, disposable, N, dual heated, 1.7 m, 10 pcs., not made with natural rubber latex	MP02650	23, 388
Breathing circuit VentStar® Helix, disposable, N, heated, 1.7 m, 10 pcs., not made with natural rubber latex	MP02608	23, 388
Breathing circuit VentStar® ICU (P) 180, disposable, P, basic, y-piece detachable, 1.8 m, 25 pcs., not made with natural rubber latex	MP17011	10, 11
Breathing circuit VentStar® ICU (P) 250, disposable, P, basic, y-piece detachable, 2.5 m, 25 pcs., not made with natural rubber latex	MP17012	10
Breathing circuit VentStar® ICU 180, disposable, basic, y-piece detachable 1.8 m, 25 pcs., not made with natural rubber latex	MP17009	9
Breathing circuit VentStar® ICU 250, disposable, basic, y-piece detachable, 2.5 m, 10 pcs., not made with natural rubber latex	MP17010	9
Breathing circuit VentStar® Oxylog 1000, disposable, 1.8 m, 25 pcs., not made with natural rubber latex	2M86837	32
Breathing circuit VentStar® Oxylog 3000, 3000+, disposable, P, 1.9 m, 5 pcs., Oxylog 3000: ≥ SW 1.20, VE300 (Software ≥ 2.0, Pediatric Option), not made with natural rubber latex	5704964	33
Breathing circuit VentStar® Oxylog 3000F, 2000p, VE300, disposable, 1.5 m (4.92 ft), 5 pcs., for Oxylog VE300 only with Option Plus	5703041	32
Breathing circuit VentStar® Oxylog 3000F, 2000p, VE300, disposable, 3 m (9.8 ft), 5 pcs., for Oxylog VE300 only with Option Plus	MP00335	32
Breathing circuit VentStar® Oxylog VE300, disposable, 3 m, 5 pcs., not made with natural rubber latex	MP01371	33
Breathing circuit VentStar® Oxylog VE300, not for VE300 with Option plus, disposable, 1.5 m, 5 pcs., not made with natural rubber latex	MP01370	33
Breathing circuit VentStar® WT ICU 180, disposable, water trap, y-piece detachable, w/o Luer Lock, 1.8 m, 25 pcs., not made with natural rubber latex	MP17008	11
Breathing circuit VentStar®, disposable, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP00349	9
Breathing circuit VentStar®, disposable, basic, 2.5 m, 10 pcs., not made with natural rubber latex	MP00350	9
Breathing circuit VentStar®, disposable, coaxial, 1.7 m, 20 pcs., not made with natural rubber latex	MP03373	12
Breathing circuit VentStar®, disposable, coaxial, 1.9 m, 20 pcs., not made with natural rubber latex	MP03374	12
Breathing circuit VentStar®, disposable, coaxial, 2.3 m, 20 pcs., not made with natural rubber latex	MP03376	12
Breathing circuit VentStar®, disposable, coaxial, GSL, 2.3 m, 20 pcs., not made with natural rubber latex	MP03378	12

INDEX BY PRODUCT NAME

Breathing circuit VentStar®, disposable, coaxial, water trap, 2.1 m, 10 pcs., not made with natural rubber latex	MP00373	12
Breathing circuit VentStar®, disposable, coaxial, water trap, GSL, 2.1 m, 10 pcs., not made with natural rubber latex	MP03379	12
Breathing circuit VentStar®, disposable, extendable, 2.2 m, 25 pcs., not made with natural rubber latex	MP00355	13
Breathing circuit VentStar®, disposable, N, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP00353	10, 386
Breathing circuit VentStar®, disposable, N, heated, 1.2 m, 10 pcs., not made with natural rubber latex	MP00308	22, 388
Breathing circuit VentStar®, disposable, N, heated, w/o chamber, 1.2 m, 10 pcs., not made with natural rubber latex	MP00309	22, 388
Breathing circuit VentStar®, disposable, N, MRI, 3 m, 10 pcs., not made with natural rubber latex	MP04802	13, 386
Breathing circuit VentStar®, disposable, N, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	MP00363	11, 386
Breathing circuit VentStar®, disposable, P, basic, 1.8 m, 25 pcs., not made with natural rubber latex	MP00351	10
Breathing circuit VentStar®, disposable, P, basic, 2.5 m, 25 pcs., not made with natural rubber latex	MP00352	10
Breathing circuit VentStar®, disposable, P, basic, w/o Luer Lock, 1.5 m, 25 pcs., not made with natural rubber latex	MP00338	10
Breathing circuit VentStar®, disposable, P, coaxial, 1.5 m, 20 pcs., not made with natural rubber latex	MP00379	12
Breathing circuit VentStar®, disposable, P, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	MP00361	11
Breathing circuit VentStar®, disposable, P, water trap, w/o Luer Lock, 1.8 m, 25 pcs., not made with natural rubber latex	MP00362	11
Breathing circuit VentStar®, disposable, smoothbore, for MRI, 3 m, 10 pcs., not made with natural rubber latex	MP04800	13
Breathing circuit VentStar®, disposable, water trap, 1.8 m, 25 pcs., not made with natural rubber latex	MP00337	11
Breathing circuit VentStar®, disposable, water trap, w/o Luer Lock, 1.8 m, 25 pcs., not made with natural rubber latex	MP00334	11
Breathing hose for breathing bag, disposable, 15 mm, 25 pcs.	MU10900	361
Breathing hose for Oxylog® 1000, reusable, silicone, 3 m	2112760	45, 46
Breathing hose for TI 500 Globe-Trotter, disposable, 20 pcs.	MU04771	366
Breathing hose holder for Babytherm®, neonatal, for 10 mm breathing hose	2M21191	358
Breathing hose set for Oxylog® 2000/3000 family, reusable, w/ flow measurement line, 1.5 m	8412068	46
Breathing hose set for Oxylog® 2000/3000 family, reusable, w/ flow measurement line, 3 m	8412913	46
Breathing hose, silicone, reusable, A, 1.1 m, connector 22 mm/22 mm, Ø 22 mm	2166038	42, 43, 47
Breathing hose, silicone, reusable, A, 1.5 m, connector 22 mm/22 mm, Ø 22 mm	2166046	41, 42, 45, 46, 47
Breathing hose, silicone, reusable, A, 35 cm, connector 22 mm/22 mm, Ø 22 mm	2165619	40, 47
Breathing hose, silicone, reusable, A, 60 cm, connector 22 mm/22 mm, Ø 22 mm	2165627	40, 43, 47
Breathing hose, silicone, reusable, N/P, 1.1 m, connector 10 mm/22 mm, Ø 10 mm	2165651	49, 395
Breathing hose, silicone, reusable, N/P, 1.1 m, connector 22 mm/22 mm, Ø 10 mm	2165678	44, 48, 393
Breathing hose, silicone, reusable, N/P, 1.5 m, connector 10 mm/22 mm, Ø 10 mm	2165961	44, 49, 393, 395
Breathing hose, silicone, reusable, N/P, 40 cm, connector 10 mm/22 mm, Ø 10 mm	2165856	39, 49, 392, 395
Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/10 mm, Ø 10 mm	2165848	39, 48, 392, 394
Breathing hose, silicone, reusable, N/P, 60 cm, connector 10 mm/22 mm, Ø 10 mm	2165821	39, 49, 392, 395
C Calibration set, with test connection and test gas cylinder	8412710	169, 205
Canister opener, for canisters MX00001 and MX50050, 5 pcs.	MX10209	99

INDEX BY PRODUCT NAME

Cap BabyFlow [®] , disposable, size L, green, 25.5 – 28 cm (10 – 11 in) head circumference, 5 pcs.	8418536	379, 384
Cap BabyFlow [®] , disposable, size M, red, 23 – 25.5 cm (9.1 – 10 in) head circumference, 5 pcs.	8418535	379, 384
Cap BabyFlow [®] , disposable, size S, yellow, 20.5 – 23 cm (8.1 – 9.1 in) head circumference, 5 pcs.	8418534	379, 384
Cap BabyFlow [®] , disposable, size XL, orange, 28 – 30.5 cm (11 – 12 in) head circumference, 5 pcs.	8418537	379, 384
Cap BabyFlow [®] , disposable, size XS, white, 18 – 20.5 cm (7.1 – 8.1 in) head circumference, 5 pcs.	8418533	379, 384
Cap BabyFlow [®] , disposable, size XXL, blue, 30.5 – 33 cm (12 – 13 in) head circumference, 5 pcs.	8418538	379, 384
Cap BabyFlow [®] , disposable, size XXL+, white/black, 33 – 36 cm (13 – 14.2 in) head circumference, 5 pcs.	8418539	379, 384
Cap for water vacuum gauge, for VarioVac [®] T Rail WV (MP00512), VarioAir [®] T Rail WV (MP00544) and water vacuum gauge (MP00518), including gaskets, combination possible with filter MK00514 and plug nozzle MP00566	MP00595	318
Cap secretion jar Medela	MP03231	305
Cap, for 700 ml secretion jar, with overflow safeguard	2M85011	299
Cap, for secretion jar, for fistula indicator for continuous and thoracic drainage	2M85013	299
Cap, for secretion jar, w/ overflow safeguard and pressure relief valve (≥ 0.1 kPa)	2M85012	299
Cardiac Output cable, 3 m	2601201	233
Cardiac output kit, Argon/Merit Medical	7265163	204
Cardiac output kit, Baxter	7265155	204
Carrier rack CD	2M85380	295
Carrying rack bronchial suction	2M85381	295
Catheter mount ErgoStar [®] CM 40, disposable, 10 cm, double swivel elbow, 50 pcs.	MP01840	63
Catheter mount ErgoStar [®] CM 45, disposable, 15 cm, double swivel elbow, 50 pcs.	MP01845	63
Catheter mount ErgoStar [®] CM 50, disposable, 10 – 16 cm, extendable, 50 pcs.	MP01850	63
Catheter mount ErgoStar [®] CM 55, disposable, 15 cm, straight, 50 pcs.	MP01855	63
Catheter mount ErgoStar [®] CM 60, disposable, 10 cm, small, double swivel elbow, 50 pcs.	MP01860	63
Ceiling holder, rail holder for compact suction systems and secretion jar sets	2M85119	296, 339, 432
Clamp for GCX rail, 5 pcs.	2M60401	348
Cleaning plate for CLIC adapter	MX50081	101
CLIC Absorber 800+, disposable absorber, 1.2 l, 6 pcs.	MX00004	100
CLIC Absorber Free, disposable absorber, 1.2 l, 6 pcs.	MX50100	100
CNAP [™] cuff controller and cradle	MS26125	146
CNAP [™] cuff controller forearm strap	MS26122	146
CNAP [™] intermediate cable for Infinity [®] CNAP [™] SmartPod [®] , for CNAP [™] controller	MS15893	146
CNAP [™] sensor cuff, size L	MS15896	147
CNAP [™] sensor cuff, size M	MS15895	147
CNAP [™] sensor cuff, size S	MS15894	147
CO cable for Infinity [®] PiCCO Pod	MS16916	161
CO catheter cable	8419160	159, 204
CO catheter cable for Infinity [®] PiCCO Pod	MS16918	161
CO connection cable, 1 m	3368458	159, 204
CO Control Syringe	2601192	233
CO In-line Injection temperature probe Becton Dickinson	2601190	233
CO In-line Injection temperature probe housing	2601191	233
CO thermistor cable Argon/Merit Medical	8420077	159, 204
CO thermistor cable Baxter/Edwards	8539983	159, 204
CO thermistor cable for Infinity [®] PiCCO Pod	MS16919	161
CO thermistor T-piece Argon/Merit Medical	5741975	159, 204

INDEX BY PRODUCT NAME

CO ₂ airway adapter, adult, disposable, 10 pcs.	2606227	167, 236
CO ₂ airway adapter, neonatal (infant/pediatric), disposable, 10 pcs.	2606228	167, 236, 417
CO ₂ cuvette disposable, adult, 10 pcs.	MP01062	62, 168, 205, 235
CO ₂ cuvette disposable, pediatric, 10 pcs.	MP01063	62, 168, 205, 235, 399
CO ₂ cuvette reusable, adult	6870279	62, 168, 205, 235
CO ₂ cuvette reusable, pediatric	6870280	62, 168, 205, 235, 399
CO ₂ mainstream sensor, for Evita® 2 dura, Evita® 4 Edition and Evita® XL	6871500	62, 399
CO ₂ Microstr. O ₂ Nasal Adult - Extended Duration, 25 pcs.	MP15010	171, 207, 253
CO ₂ Microstr. O ₂ Nasal Pediatric - Extended Duration, 25 pcs.	MP15011	171, 207, 253
CO ₂ Microstr. O ₂ Nasal/Oral Adult - Short-term use, 25 pcs.	MP15012	171, 207, 253
CO ₂ Microstr. O ₂ Nasal/Oral Pediatric - Short-term use, 25 pcs.	MP15013	171, 207, 253
CO ₂ nasal cannula + O ₂ , adult, disposable, for Vista 120, 50 pcs.	2607554	234
CO ₂ nasal cannula + O ₂ , child, disposable, for Vista 120, 50 pcs.	2607555	234
CO ₂ nasal cannula, adult, disposable, for Vista 120, 100 pcs.	2607551	234
CO ₂ nasal cannula, infant, disposable, for Vista 120, 100 pcs.	2607553	234
CO ₂ nasal cannula, pediatric, disposable, for Vista 120, 100 pcs.	2607552	234
CO ₂ sampling line, disposable, for Vista 120, 100 pcs.	2607550	234
CO ₂ sidestream water trap, for Vista 120, 20 pcs.	2605124	234
CombiStar A Filter HME, combination of airway connector MP01845 and filter TwinStar® 55 Plus (MP05805), 25 pcs.	MP12060	64
CombiStar F-HME HEPA flex, combination of flexible tube extension MP01850 and filter TwinStar® HEPA Plus (MP05801), 25 pcs.	MP12062	64
CombiStar Filter HME flex, combination of flexible tube extension MP01850 and filter TwinStar® 55 Plus (MP05805), 25 pcs.	MP12061	64
CombiStar mechanical Filter flex, combination of flexible tube extension MP01850 and filter SafeStar® 55 Plus (MP05790), 25 pcs.	MP12063	64
Compact flowmeter Clickflow configurator, gives you the option of arranging the compact flowmeter according to your requirements, selectable options for configuration	MP13000	269, 421
Compact rail, for 38 mm tube, standard rail profile 25 × 10 mm, each side 10 cm, weight: 0.9 kg, max. load: 30 kg	2M85337	279, 341, 377, 435
Condensate management hoses for Isolette® 8000 and 8000 plus, disposable, 20 pcs.	MU21120	354, 356
Connection cable, connects HemoMed Pod to monitor, 3 m	5591925	148
Connection cable, connects HemoMed Pod to monitor, 5 m	5591933	148
Connection hose Medela	MP03233	305
Connection hose, silicone, reusable, ID 5 mm, wall thickness 2 mm, by meter	1203606	299, 431
Connector pressure Oxylog® VE300	MP01372	46
Continuous drainage 1, for continuous and wound drainage, precision regulator, fistula indicator and bacterial filter 767 included, capacity 700 ml, maximum capacity up to the overflow safeguard 1,000 ml, scale in steps of 10 ml up to max. 700 ml	2M85151	314
Continuous drainage 3, for continuous, wound, and thoracic drainage, precision regulator, water vacuummeter, fistula indicator and bacterial filter 767 included, capacity 700 ml, maximum capacity up to the overflow safeguard 1,000 ml, scale in steps of 10 ml up to max. 700 ml	2M85153	314
Control panel	MP00589	435, 341
Controller cable for Aerogen Pro	MP01036	93
Corrugated hose, reusable, N/P, 14 cm, connector 10 mm/10 mm, Ø 10 mm	8403333	39, 392
Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, large, 100 pcs.	MU06943	347, 349, 351, 353, 355, 357, 360, 363, 366, 369, 373
Cover for skin temperature probe CARE-FOR-ME, hydrogel, disposable, small, 100 pcs.	MU06944	347, 349, 351, 353, 355, 357, 360, 363, 366, 369, 373
Covidien Quick Temp Pod Oral	2600146	251
Covidien Quick Temp Pod Rectal	2600147	251
CS connection hose holder	2M85446	296, 306, 339, 432
CS hoses, several hoses to connect flowmeter with the central supply	8601697	319

INDEX BY PRODUCT NAME

Cuddling toy Babyleo® TN500, lion, orange	MP03737	347, 372
Cuddling toy BabyRoo, kangaroo	MP16090	349, 372
Cuddling toy Cally, reusable, light blue or pink	2M30462	372
Cuddling toy Isi, deer, green	MP03738	355, 372
Cuff controller cradle allows CNAP™ cuff controller to “slide” with movement of patient’s arm	MS26123	146
Cylinder support	2M85180	315
D Demo lung SmartLung, max.1 L, reusable, not made with natural rubber latex	MP02405	53
Demo lung to simulate neonatal lungs, with suitcase	8418077	397
Demo package BabyFlow® plus, disposable, size L	MP03705	380
Demo package BabyFlow® plus, disposable, size M	MP03704	378, 380
Demo package BabyFlow® plus, disposable, size S	MP03703	378, 380
Deposit plate for Babytherm® 8000	2M21158	358
Disposable hand sensor, adult, as of Software 1.8.5	MS40320	208
Disposable Handle Polaris®, Dual, for Polaris® 100/200, 600 and Multimedia, adapter required, 60 pcs.	MP05700	343
Disposable Handle Polaris®, Single, for Polaris® 100/200, 600 and Multimedia, adapter required, 120 pcs.	MP05701	343
Double connector 11ID/15OD, disposable, 10 pcs.	MP02655	23, 389
Double-Plug coupling AIR	G01166	334
Double-Plug coupling O ₂	G01168	334
Double-Plug coupling VAC	G01169	334
Draw wire for hose heater 1.5 m, for reusable breathing circuit	8411050	31, 391
Dummy Plug	8410049	39, 40, 392
Dust filter Perseus®, 5 pcs.	MX50115	101
E Ear clip for Nellcor® DURA-Y™	MX01005	135, 226, 239
ECG "4-lead" spacer, 10 pcs.	MS22832	178
ECG adapter for single-pin ECG cable to MultiMed® 5	MS14679	120, 178
ECG adapter for single-pin ECG cable to MultiMed® 6	MS14680	120, 178
ECG adapter tel, GE-AP-CH, 3 Lead	MP04350	256
ECG adapter tel, GE-AP-CH, 5 Lead	MP04351	256
ECG adapter tel, GE-AP-CH, 6 Lead	MP04352	256
ECG adapter tel, GE-AP-FH, 5 Lead	MP04353	256
ECG adapter tel, PH M, 5 Lead	MP04358	256
ECG adapter tel, PH M, 6 Lead	MP04359	256
ECG adapter tel, PH MX40, 3 Lead	MP04355	256
ECG adapter tel, PH MX40, 5 Lead	MP04356	256
ECG adapter tel, PH MX40, 6 Lead	MP04357	256
ECG adapter tel, SL AT, 5 Lead	MP04365	256
ECG adapter tel, SL UV, 5 Lead	MP04366	256
ECG Adapter, 3 Lead GE	MP04322	254
ECG Adapter, 3 Lead Mindray	MP04312	254
ECG Adapter, 3 Lead NK	MP04316	254
ECG Adapter, 3 Lead Philips	MP04325	254
ECG Adapter, 3 Lead Spacelabs	MP04300	254
ECG Adapter, 3 Lead Welch Allyn	MP04305	254
ECG Adapter, 5 Lead GE	MP04320	254
ECG Adapter, 5 Lead Mindray	MP04310	254
ECG Adapter, 5 Lead Philips	MP04326	254
ECG Adapter, 5 Lead Spacelabs	MP04301	254
ECG Adapter, 5 Lead Welch Allyn	MP04306	254

INDEX BY PRODUCT NAME

ECG Adapter, 6 Lead NK	MP04315	254
ECG Adapter, 6 Lead Philips	MP04327	254
ECG cable for C lead MonoLead® 5, IEC1 (European color code), 1.5 m	MS16232	113
ECG cable for C lead, IEC1 (European color code), 1 m	MP03417	118
ECG cable for V lead MonoLead® 5, IEC2 (AHA/US color code), 1.5 m	MS16230	113
ECG cable for V lead, IEC2 (AHA/US color code), 1 m	MP03418	118
ECG cable, 3-lead Infinity® TruST™, IEC1 (European color code), 70 cm	MP03421	211
ECG cable, 3-lead Infinity® TruST™, IEC2 (AHA/US color code), 70 cm	MP03422	211
ECG cable, 3-lead, dual-pin connector MonoLead® 3, IEC1 (European color code), 2 m	MS16160	113
ECG cable, 3-lead, dual-pin connector MonoLead® 3, IEC2 (AHA/US color code), 2 m	MS16233	113
ECG cable, 3-lead, dual-pin connector, IEC1 (European color code), 1 m	MP03411	118
ECG cable, 3-lead, dual-pin connector, IEC2 (AHA/US color code), 1 m	MP03412	118
ECG cable, 3-lead, single-pin connector, IEC1 (European color code), 1 m	MP03401	117, 176, 211, 222
ECG cable, 3-lead, single-pin connector, single-patient use, IEC1 (European color code), 1 m, 20 pcs.	MP00875	116, 177, 212, 222
ECG cable, 3-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1 m, 20 pcs.	MP00877	116, 177, 212, 222
ECG cable, 3-lead, single-pin connector, IEC2 (AHA/US color code), 1 m	MP03402	117, 176, 211, 222
ECG cable, 4-lead, single-pin connector, IEC1 (European color code), 1 m	MP03407	176
ECG cable, 4-lead, single-pin connector, IEC2 (AHA/US color code), 1 m	MP03408	176
ECG cable, 4-lead, single-pin connector, single-patient use, IEC1 (European color code), 1 m, 20 pcs.	MP03120	177
ECG cable, 4-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1 m, 20 pcs.	MP03121	177
ECG cable, 5-lead Infinity® TruST™, IEC1 (European color code), 70 cm	MP03423	211
ECG cable, 5-lead Infinity® TruST™, IEC2 (AHA/US color code), 70 cm	MP03424	211
ECG cable, 5-lead, dual-pin connector MonoLead® 5, IEC1 (European color code), 2.5 m	MS16161	113
ECG cable, 5-lead, dual-pin connector MonoLead® 5, IEC2 (AHA/US color code), 2.5 m	MS16229	113
ECG cable, 5-lead, dual-pin connector, IEC1 (European color code), 1/1.5 m	MP03413	118
ECG cable, 5-lead, dual-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03414	118
ECG cable, 5-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03403	117, 176, 211, 222
ECG cable, 5-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03404	117, 176, 211, 222
ECG cable, 5-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs.	MP00879	116, 177, 212, 222
ECG cable, 5-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs.	MP00881	116, 177, 212, 222
ECG cable, 6-lead Infinity® TruST™, IEC1 (European color code), 70 cm	MP03425	211
ECG cable, 6-lead Infinity® TruST™, IEC2 (AHA/US color code), 70 cm	MP03426	211
ECG cable, 6-lead, dual-pin connector, IEC1 (European color code), 1/1.5 m	MP03415	118
ECG cable, 6-lead, dual-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03416	118
ECG cable, 6-lead, single-pin connector, IEC1 (European color code), 1/1.5 m	MP03405	117, 176, 211
ECG cable, 6-lead, single-pin connector, IEC2 (AHA/US color code), 1/1.5 m	MP03406	117, 176, 211
ECG cable, 6-lead, single-pin connector, single-patient use, IEC1 (European color code), 1/1.5 m, 20 pcs.	MP03122	116, 177, 212
ECG cable, 6-lead, single-pin connector, single-patient use, IEC2 (AHA/US color code), 1/1.5 m, 20 pcs.	MP03123	116, 177, 212
ECG defib trunk cable, 3/5-lead single-pin, IEC/AHA, 1.5 m	2612020	220
ECG defib trunk cable, 3/5-lead single-pin, IEC/AHA, 2.5 m	2612021	220
ECG defib trunk cable, Neonate, single-pin, IEC/AHA, 1.5 m	2612026	220
ECG defib trunk cable, Neonate, single-pin, IEC/AHA, 2.5 m	2612027	220
ECG electrodes, single-patient use, neonate, 300 pcs.	5195024	119, 178, 224, 410

INDEX BY PRODUCT NAME

ECG ESU block, for OR use of 3- or 5-lead ECG cables with MultiMed® 5	5947226	120
ECG ESU block, for OR use of 3-, 5-, or 6-lead ECG cables with MultiMed® 6	7486140	120
ECG ESU trunk cable, 3/5-lead single-pin, IEC/AHA, 1.5 m	2612024	220
ECG ESU trunk cable, 3/5-lead single-pin, IEC/AHA, 2.5 m	2612025	220
ECG ESU trunk cable, Neonate, single-pin, IEC/AHA, 1.5 m	2612028	220
ECG ESU trunk cable, Neonate, single-pin, IEC/AHA, 2.5 m	2612029	220
ECG extension cable, 4-lead, 2 m	MP05114	178
ECG extension cable, 6-lead, 2 m	MP05115	120, 178
ECG extension cable, dual-pin connector, 2 m	MS16492	120
ECG intermediate cable, 3-lead, for NeoMed Pod, 1.5 m	5592162	119, 410
ECG leadwire, 3-lead single-pin, AHA, 1 m	2612001	221
ECG leadwire, 3-lead single-pin, IEC, 1 m	2612000	221
ECG leadwire, 5-lead single-pin, AHA, 1.5 m	2612003	221
ECG leadwire, 5-lead single-pin, IEC, 1.5 m	2612002	221
ECG limb wires, 3-lead, Grabber, AHA colour code	2606497	224
ECG limb wires, 3-lead, Grabber, IEC colour code	2606493	224
ECG limb wires, 5-lead, Grabber, AHA colour code	2606498	224
ECG limb wires, 5-lead, Grabber, IEC colour code	2606494	224
ECG MonoLead®, 3-lead, AHA/US, 2 m	MP05121	114, 174, 210, 223
ECG MonoLead®, 3-lead, AHA/US, 4.1 m	MP05131	115, 174, 223
ECG MonoLead®, 3-lead, EUR, 2 m	MP05120	114, 174, 210, 223
ECG MonoLead®, 3-lead, EUR, 4.1 m	MP05130	115, 174, 223
ECG MonoLead®, 4-lead, AHA/US, 1.7 m	MP05127	175
ECG MonoLead®, 4-lead, AHA/US, 3.7 m	MP05129	175
ECG MonoLead®, 4-lead, EUR, 1.7 m	MP05126	175
ECG MonoLead®, 4-lead, EUR, 3.7 m	MP05128	175
ECG MonoLead®, 5-lead, AHA/US, 2.5 m	MP05123	114, 174, 210, 223
ECG MonoLead®, 5-lead, AHA/US, 5.3 m	MP05133	115, 174, 223
ECG MonoLead®, 5-lead, EUR, 2.5 m	MP05122	114, 174, 210, 223
ECG MonoLead®, 5-lead, EUR, 5.3 m	MP05132	115, 174, 223
ECG MonoLead®, 6-lead, AHA/US, 3 m	MP05125	114, 174, 210
ECG MonoLead®, 6-lead, AHA/US, 5.7 m	MP05135	115, 174
ECG MonoLead®, 6-lead, EUR, 3 m	MP05124	114, 174, 210
ECG MonoLead®, 6-lead, EUR, 5.7 m	MP05134	115, 174
ECG neonatal 3-lead adapter cable, 2 m	MS25951	178
ECG trunk cable, 3-lead, for neonates, IEC & AHA	2606492	224
ECG trunk cable, 3-lead, IEC/AHA colour code	2606495	224
ECG trunk cable, 5-lead, IEC/AHA colour code	2606496	224
EEG electrode cable, color-coded, 9 pcs.	5947804	165
EEG electrodes, disposable, 25 pcs.	MS29468	165
Ejector suction NIST AIR O ₂ for bronchial suction, gas connection AIR/O ₂ , for connection to all standard wall rails	MK01418	291
Ejector suction NIST for bronchial suction, gas connection AIR, for connection to all standard wall rails	MK01420	291
Ejector suction NIST, for bronchial suction, gas connection AIR, version for Atlan® A300/A350 and Perseus® A500	MP00688	291
Ejector VarioAir® B AGA, for bronchial suction, gas connection AIR, Swedish standard	MP00533	288
Ejector VarioAir® B BS, for bronchial suction, gas connection AIR	MP00531	288
Ejector VarioAir® B CARBA, for bronchial suction, gas connection AIR	MP00548	288
Ejector VarioAir® B compact, for bronchial suction, gas connection NIST AIR, including secretion jar set and integrated ejector	MP00536	293
Ejector VarioAir® B DIN, for bronchial suction, gas connection AIR	MP00530	288
Ejector VarioAir® B NF, for bronchial suction, gas connection AIR	MP00532	288

INDEX BY PRODUCT NAME

Ejector VarioAir® B UNI, for bronchial suction, gas connection AIR	MP00534	288
Ejector VarioAir® B/OR Rail DISS, for bronchial and OR suction, gas connection DISS AIR	MP00638	289
Ejector VarioAir® B/OR Rail, for bronchial and OR suction, gas connection NIST AIR	MP00535	289, 304
Ejector VarioAir® D AGA, for continuous and wound drainage, gas connection AIR	MP00540	311
Ejector VarioAir® D BS, for continuous and wound drainage, gas connection AIR	MP00538	311
Ejector VarioAir® D CARBA, for continuous and wound drainage, gas connection AIR	MP00549	311
Ejector VarioAir® D DIN, for continuous and wound drainage, gas connection AIR	MP00537	311
Ejector VarioAir® D NF, for continuous and wound drainage, gas connection AIR	MP00539	311
Ejector VarioAir® D Rail, for continuous and wound drainage, gas connection AIR, for connection to all standard wall rails	MP00542	311, 426
Ejector VarioAir® D UNI, for continuous and wound drainage, gas connection AIR	MP00541	311
Ejector VarioAir® MRI B/OR Rail, for bronchial suction, gas connection NIST AIR, for connection to all standard wall rails	MP03220	290
Ejector VarioAir® MRI P Rail, for bronchial suction, gas connection NIST AIR, for connection to all standard wall rails	MP03221	290, 427
Ejector VarioAir® P compact, for bronchial suction, gas connection NIST AIR, including secretion jar set and integrated ejector	MP00547	293, 428
Ejector VarioAir® P Rail DISS, for bronchial suction, gas connection DISS AIR, for connection to all standard wall rails	MP00639	290
Ejector VarioAir® P Rail, for bronchial suction, gas connection NIST AIR, for connection to all standard wall rails	MP00545	290, 427
Ejector VarioAir® T DIN, for thoracic drainage, gas connection DIN, for connecting to the central AIR supply according to DIN standard	MP00553	312
Ejector VarioAir® T Rail WV, for thoracic drainage, with water vacuum gauge, for connection to all standard wall rails	MP00544	313
Ejector VarioAir® T Rail, for thoracic drainage, gas connection NIST AIR, for connection to all standard wall rails	MP00543	312
Ejector VarioO2® B CARBA, for bronchial suction, gas connection O ₂	MP00551	289
Ejector VarioO2® B DIN, for bronchial suction, gas connection O ₂	MP00550	289
Ejector VarioO2® B/OR Rail, for bronchial and OR suction, gas connection NIST O ₂	MP00552	278, 289, 304
Elbow for F&P, reusable	8418289	51
Elbow, angled hose connection for secretion jar caps	2M19063	300, 431
Elbow, reusable	8412235	46, 51
Elbow, reusable, N/P, metal, Luer Lock	M34568	44, 51, 393, 396
Electrode adapter, to connect neonatal electrodes to MultiMed® pod, 10 pcs.	5194779	119, 410
Electrostatic filter CareStar® 20 Plus, disposable, ~20 ml deadspace, 100 pcs.	MP05770	56
Electrostatic filter CareStar® 35 Plus, disposable, ~35 ml deadspace, 100 pcs.	MP05755	56
etCO ₂ adapter cable M11.1	2607032	235
etCO ₂ airway adapter, reusable, adult	4721796	167
etCO ₂ airway adapter, reusable, neonate	4721788	167, 417
etCO ₂ Microstr. Intubated Adult/Pediatric - Extended Duration, 25 pcs.	MP15001	170, 206, 252
etCO ₂ Microstr. Intubated Adult/Pediatric - Short-term use, 25 pcs.	MP15002	170, 206, 252
etCO ₂ Microstr. Intubated Infant/Neonatal - Extended Duration, 25 pcs.	MP15003	170, 206, 252, 417
etCO ₂ Microstr. Nasal Adult - Extended Duration, 25 pcs.	MP15004	170, 206, 252
etCO ₂ Microstr. Nasal Adult - Short-term use, 25 pcs.	MP15005	170, 206, 252
etCO ₂ Microstr. Nasal Infant/Neonatal - Extended Duration, 25 pcs.	MP15007	170, 206, 252, 417

INDEX BY PRODUCT NAME

etCO ₂ Microstr. Nasal Pediatric - Short-term use, 25 pcs.	MP15006	170, 206, 252
etCO ₂ Microstr. Nasal/Oral Adult - Short-term use, 25 pcs.	MP15008	170, 206, 252
etCO ₂ Microstr. Nasal/Oral Pediatric - Short-term use, 25 pcs.	MP15009	170, 206, 252
etCO ₂ sensor Capnostat 5, reusable	2606226	236
Expiratory filter, disposable, 20 pcs.	MP01781	87
Expiratory valve for Babylog® 8000	8408950	85
Expiratory valve for Babylog® VN 500, reusable (RFID)	8415270	85
Expiratory valve for Carina®, disposable, 25 pcs.	MP00220	86
Expiratory valve for Evita® 2 dura, 4 edition, XL, disposable, 10 pcs.	8414776	85
Expiratory valve for Evita® 2 dura, 4 edition, XL, disposable, 10 pcs.	MP02600	85
Expiratory valve for Evita® Infinity® V500, reusable	8416750	86
Expiratory valve for Evita®, reusable	8410580	85
Expiratory valve for Savina®, reusable	8413660	86
Expiratory valve, disposable, 10 pcs.	MP01061/ MP20109	86
Expiratory valve, for neonatal/pediatric patients (Non-RFID), reusable	8415296	85
Expiratory valve, reusable	8417050	86
Extension cable	MS32153	209
Eyebrow sensor, reusable	MS32291	208
Eyebrow sensor, reusable, as of Software 1.8	MS40321	208
F F&P 950 A80 breathing circuit (A)	MP16004	28, 389
F&P 950 heater wire adaptor	MP16003	28, 389
F&P 950 neonatal disposable breathing circuit	MP16005	28, 389
F&P 950 sensor cartridge	MP16002	28, 389
F&P configuration for MR810, MR850 and F&P 950	8414144	28, 350
Filac 3000 Probe Covers, 20 pcs. Vista 120 SC	2600945	251
Filac 3000 Probe, oral, reusable, 1.2 m, Vista 120 SC	2600939	251
Filac 3000 Probe, oral, reusable, 2.7 m, Vista 120 SC	2600940	251
Filac 3000 Probe, rectal, reusable, 1.2 m, Vista 120 SC	2600942	251
Filac 3000 Probe, rectal, reusable, 2.7 m, Vista 120 SC	2600943	251
Filler cap for Aerogen Pro, silicone, 5 pcs.	MP01035	93
Filter + HME TwinStar® 25 Plus, disposable, ~25 ml deadspace, 100 pcs.	MP05815	55
Filter + HME TwinStar® 55 Plus, disposable, ~55 ml deadspace, 100 pcs.	MP05805	55
Filter + HME TwinStar® 60A Plus, disposable, ~60 ml deadspace, angled, 100 pcs.	MP05810	55
Filter + HME TwinStar® 9 Plus, disposable, ~9 ml deadspace, 100 pcs.	MP05820	55, 398
Filter + HME TwinStar® 90 Plus, disposable, ~90 ml deadspace, 100 pcs.	MP05800	55
Filter for suction bottle 40 ml for Resuscitaire®, disposable, 10 pcs.	MU10915	362
Flow measuring line for Oxylog® 2000/3000 family, 1.5 m	8410929	46
Flow measuring line for Oxylog® 2000/3000 family, 3 m	8413293	46
Flow sensor cable for V-series with Neo option and VN-series	8416600	106
Flow sensor cable, for neonatal flow sensor, e.g., Babylog®	8409626	106, 400
Flow sensor for Crossvent 2i+, disposable, infant, 10 pcs.	MU24401	364
Flow sensor Infinity® ID, PC/ABS, 5 pcs.	6871980	105
Flow sensor SpiroLife®, PPSU, sterilizable	MK01900	105
Flow sensor Spirolog®, PC/ABS, 5 pcs.	8403735	105
Fluxmed TM, Pressure Probe with esophageal balloon, disposable	MP20110	88
Foot sensor, reusable	MS32290	208
Foot sensor, reusable, as of Software 1.8	MS40322	208
G Gastroschisis arm, complete	2M60180	348
Gel Mattress for BabyTherm 8000/8010 and Babyroo TN300, reusable	2M20827	349, 357
General-purpose temperature probe, reusable, adult, 7-pin connector, 1.5 m	4329889	140, 192, 231

INDEX BY PRODUCT NAME

	General-purpose temperature probe, reusable, adult, 7-pin connector, 3 m	5204644	140, 192, 231
	General-purpose temperature probe, reusable, child, 7-pin connector, 1.5 m	4329848	140, 192, 231
	General-purpose temperature probe, reusable, child, 7-pin connector, 3 m	5204651	140, 192, 231
	General-purpose temperature probe, single-patient use, adult, 7-pin connector, 1.6 m	MP00991	139, 191, 230
	General-purpose temperature probe, single-patient use, adult, 7-pin connector, 3 m	MP00992	139, 191, 230
	General-purpose temperature probe, single-patient use, child, 7-pin connector, 1.6 m	MP00993	139, 191, 230
	General-purpose temperature probe, single-patient use, child, 7-pin connector, 3 m	MP00994	139, 191, 230
H	H300 Heater wire adapter, dual	MP07107	27
	H300 Heater wire adapter, insp	MP07106	27
	H300 Temp. sensor cable	MP07109	27
	H300 Temp. sensor cable, 160 cm	MP07108	27
	Hand port cuff, disposable, oval, 100 pcs.	MU06570	356, 363, 367, 369
	Hand sensor, reusable, adult	MS32154	208
	Hand sensor, reusable, adult, as of Software 1.8	MS40317	208
	Hand sensor, reusable, pediatric	MS32289	208
	Hand sensor, reusable, pediatric, as of Software 1.8	MS40318	208
	Hand sensor, reusable, small pediatric, as of Software 1.8	MS40319	208
	Handle, complete for Babytherm® 8000	2M21161	358
	Headband LNCS TF-I, 5 pcs.	MP03197	133, 188
	Headgear, disposable, size L, 31 – 37 cm (12.2 – 14.6 in) head circumference, 5 pcs.	MP03884	379, 383
	Headgear, disposable, size M, 27 – 31 cm (10.6 – 12.2 in) head circumference, 5 pcs.	MP03883	379, 383
	Headgear, disposable, size S, 23 – 27 cm (9.1 – 10.6 in) head circumference, 5 pcs.	MP03882	379, 383
	Headgear, disposable, size XL, 37 – 42 cm (14.6 – 16.5 in) head circumference, 5 pcs.	MP03885	379, 383
	Headgear, disposable, size XS, 20 – 23 cm (7.9 – 9.1 in) head circumference, 5 pcs.	MP03881	379, 383
	Headgear, disposable, size XXS, 17 – 20 cm (6.7 – 7.9 in) head circumference, 5 pcs.	MP03880	379, 383
	Headgear, disposable, soft fabric, large	MP01560	78
	Headgear, disposable, soft fabric, small	MP01559	78, 405
	Headgear, disposable, soft rubber, large	MP01562	78
	Headgear, disposable, soft rubber, small	MP01561	78, 405
	Heating mattress, 1 pc.	2M60170	347
	HI-Flow Star Kit Aquapor A, heated inspiratory breathing circuit for the H300 humidifier, heated inspiratory hose, humidifier chamber, breathing hose when used with a ventilator, O ₂ hose when used with an O ₂ blender, bed clip, disposable, 10 pcs.	MP17090	27, 82
	HI-Flow Star Nasal Cannula, disposable, adult L, Flow Range < 60 l/min, 25 pcs.	MP05513	82
	HI-Flow Star Nasal Cannula, disposable, adult M, Flow Range < 60 l/min, 25 pcs.	MP05512	82
	HI-Flow Star Nasal Cannula, disposable, adult S, Flow Range < 60 l/min, 25 pcs.	MP05511	82
	HI-Flow Star System, heated inspiratory breathing circuit, consists of: heated inspiratory hose, humidifier chamber, breathing hose when used with a ventilator, O ₂ hose when used with an O ₂ blender, bed clip, disposable, 25 pcs.	MP05601	82
	HI-Flow Star Valve Kit w/ Pop-Up Valve and Connector 22 mm, disposable, 15 pcs.	MP05507	27, 82
	Hinged arm 240°, 3 hinges, fixable, ball holder system, length 1,155 mm, max. load 0.3 kg, includes hose clamps: MP00493 (1×), 8409841 (1×)	8409609	342
	Hinged arm Infinity® ACS, for positioning and fixation of breathing hoses, this arm can be fixed into place without a hand movement through a patent registered friction torque-breaking system, length 1,080 mm, max. load 0.5 kg, includes hose clamps: MP00693 (1×), 8417396 (2×)	MP00690	342

INDEX BY PRODUCT NAME

Hinged arm Quickstop 2, with universal ball holder system for optional accessories length 858 mm, max. load 0.4 kg, includes hose clamps: MP00493 (1×), 8404761 (2×), 8409841 (1×)	2M85706	342, 436
Hinged arm Quickstop short, with mounted universal tube holder plate, for various tube dimensions/coaxial hoses, specifically designed for OR tables, length 323 mm, max. load 1 kg	2M86463	342, 436
HME HumidStar® 2 Plus, disposable, 2 ml deadspace, 100 pcs.	MP05845	58, 398
HME HumidStar® 25 Plus, disposable, ~25 ml deadspace, 100 pcs.	MP05735	58
HME HumidStar® 55 Plus, disposable, ~55 ml deadspace, 100 pcs.	MP05730	58
HME HumidStar® Trach Plus, disposable, 6 ml deadspace, 100 pcs.	MP05750	58
Holder for cylinder packet for GCX profile	4118453	348
Holder for F&P humidifier U-shaped, for trolley, with rail clamp	MP01002	29, 390
Holder for hinged hood for Babytherm® 8000 WB, 8004 and 8010, reusable, neonatal, mounts to hand rail	2M21342	358
Holder for Infinity® M300 bedside charger, swiveling 360°	MP00990	219
Holder for MultiMed® 12 Pod	5202747	112
Holder for MultiMed® 5/6, MultiMed® Plus/OR, NeoMed Pods and Tcore®, swiveling 360°	MP00721	112, 141, 193
Holder for MultiMed® 5/6, MultiMed® Plus/OR, and NeoMed Pods	5598128	112
Holder, lateral (F&P 900 MR 309), including 8411074	MP01009	29, 390
Hook ring for LiteStar® anesthesia masks, disposable, 22 mm ID, for sizes 3 – 7, 30 pcs.	MP11008	74
Hook ring for SediStar® anesthesia masks, disposable, 22 mm ID, for sizes 3 – 7, 30 pcs.	MP01708	73
Hook ring for silicone masks, reusable, 15 mm OD, for sizes 00 – 1, 10 pcs.	MP01916	77, 404
Hook ring for silicone masks, reusable, 22 mm ID, for sizes 2 – 5, 10 pcs.	MP01917	77
Hose ASM-WAGD DISS N/ DISS N 10L	MP20022	331
Hose ASM-WAGD DISS N/ DISS N 15L	MP20023	331
Hose ASM-WAGD DISS N/ DISS N 20L	MP20024	331
Hose ASM-WAGD HAND DISS N 10L	MP20029	331, 332
Hose ASM-WAGD HAND DISS N 15L	MP20030	331, 332
Hose clamp, 10 pcs.	8403566	40
Hose heater adapter for F&P MR 850, for dual-heated disposable breathing circuit	MP01006	31, 391
Hose heater adapter for F&P MR 850, for reusable breathing circuit	8414968	31, 391
Hose heater adapter for F&P MR 850, for single-heated disposable breathing circuit	8414992	31, 391
Hose heater for F&P MR 700/730 and MR 850, for reusable breathing circuit, 0.7 m	8411154	29, 390
Hose heater for F&P MR 700/730 and MR 850, for reusable breathing circuit, 1.1 m	8414993	29, 390
Hose heater for F&P MR 700/730 and MR 850, for reusable breathing circuit, 1.1 m (Low -Flow)	8411045	29
Hose heater for F&P MR 700/730 and MR 850, for reusable breathing circuit, 1.5 m	8414969	29
Hose nozzle 8 mm small, connection nozzle double sided 8 mm	5000833	318
Hose Oxylog® VE300, 150, reusable	MP01390	46
Hose Oxylog® VE300, 300, reusable	MP01391	46
Humidifier chamber for F&P MR 290, disposable, 10 pcs.	8418282	28, 390
Humidifier chamber for F&P MR 340, reusable, N/P	8411047	28, 390
Humidifier chamber for F&P MR 370, reusable, A	8412217	28
Humidifier chamber for Isolette® C2000, reusable	MU17191	356
Humidifier chamber VentStar®, disposable, 10 pcs.	MP00234	22, 389
Humidifier chamber VentStar®, disposable, N, 10 pcs.	MP02590	23, 389

INDEX BY PRODUCT NAME

Humidifier chamber, 10 pcs.	MP02605	24
Humidifier holder chassis, swiveling	G93111	29, 390
Humidifier holder for side rail	8416325	29, 390
Humidifier-bubbler, for oxygen humidification up to appr. 60% relative humidity and flows up to 5 l/min	2M85834	281, 422
Humidifier-nebulizer, for oxygen humidification up to 100% relative humidity and flows above 5 l/min	2M85835	281
Humidity pad, disposable, antibacterial, 25 pcs.	MU06008	364, 368, 370
I IBP adapter, 10-pin to 7-pin	3368383	156
IBP cable 10pin	2607559	232
IBP cable 7pin	2607558	232
IBP cable Abbott, Medex	2606490	232
IBP cable Becton Dickinson/Argon	2606488	232
IBP cable Edwards	2606489	232
IBP cable Utah	2606491	232
IBP connection cable Argon/Merit Medical, 3.7 m (12.1 ft)	3375933	158
IBP holder for pressure transducer Memscap (Capto/SensoNor) 844, 4 pcs.	7864155	157
IBP pressure transducer cable for Memscap (Capto/SensoNor), 3.7 m	4321563	157
IBP pressure transducer for Memscap (Capto/SensoNor) 844, disposable, sterile, 50 pcs.	7864908	157
IBP pressure transducer for Memscap (Capto/SensoNor) 844, reusable	7489417	157
IBP pressure transducer kit, BD/Argon, disposable, 5 kits	2606225	232
IBP pressure transducer plate, for Hemo2 Pod, for Memscap (Capto/SensoNor), 5 pcs.	4721614	155
IBP pressure transducer plate, for Hemo2/4 Pods, for Abbott Transpac IV, 5 pcs.	7270460	155
IBP pressure transducer plate, for Hemo2/4 Pods, for Abbott TranspacII/III/Braun, 5 pcs.	5192112	155
IBP pressure transducer plate, for Hemo4 Pod, for Memscap (Capto/SensoNor), 5 pcs.	4721416	155
IBP pressure transducer plate, for Hemo4 Pod, universal, 5 pcs.	4721424	155
IBP pressure transducer set Argon/Merit Medical, disposable, sterile, 5 pcs.	4528741	158
IBP pressure transducer set for Memscap (Capto/SensoNor) 844	7498947	157
IBP pressure transducer set for Memscap (Capto/SensoNor) 844, disposable, 10 pcs.	7489433	157
IBP transducer adapter Abbott/Medex/ICU Medical	MS22535	149, 199
IBP transducer adapter Argon/Merit Medical	MS22148	149, 199
IBP transducer adapter Baxter/Edwards	MS22147	149, 199
IBP transducer adapter cable, Abbott/Medex, Smiths/TranStar, 0.3 m	2601230	151, 201
IBP transducer adapter cable, Abbott/Medex, Smiths/TranStar, 2 m	2601231	151, 201
IBP transducer adapter cable, Argon/BD/Ohmeda, 0.3 m	2601222	153, 203
IBP transducer adapter cable, Argon/BD/Ohmeda, 2 m	2601223	153, 203
IBP transducer adapter cable, B.Braun, 0.3 m	2601232	150, 200
IBP transducer adapter cable, B.Braun, 2 m	2601233	150, 200
IBP transducer adapter cable, Baxter/Edwards, 0.3 m	2601220	153, 203
IBP transducer adapter cable, Baxter/Edwards, 2 m	2601221	153, 203
IBP transducer adapter cable, Codan/DPT, 0.3 m	2601234	150, 200
IBP transducer adapter cable, Codan/DPT, 2 m	2601235	150, 200
IBP transducer adapter cable, Codan/Xtrans, 0.3 m	2601236	153, 203
IBP transducer adapter cable, Dräger 10-pin	MS22532	149, 161, 199
IBP transducer adapter cable, Dräger 10-pin, 0.3 m	2601224	152, 202
IBP transducer adapter cable, Dräger 10-pin, 2 m	2601225	152, 202
IBP transducer adapter cable, Dräger 7-pin	MS22533	149, 199
IBP transducer adapter cable, Dräger 7-pin, 0.3 m	2601226	152, 202

INDEX BY PRODUCT NAME

IBP transducer adapter cable, Dräger 7-pin, 2 m	2601227	152, 202
IBP transducer adapter cable, Utah Medical, 0.3 m	2601228	151, 201
IBP transducer adapter cable, Utah Medical, 2 m	2601229	151, 201
IBP transducer adapter Utah Medical	MS22534	149, 199
IBP Y-adapter, 10-pin	5731281	156
IBP Y-adapter, 10-pin to 7-pin	5588095	156
IBP Y-adapter, 7-pin	5592147	156
IBP Y-cable, Baxter/Edwards, 3.7 m	5206607	156
IBP Y-cable, for Memscap (Capto/SensoNor), 3.7 m	5195180	157
IBP Y-cableArgon/Merit Medical 3.7 m (12.1 ft)	5206581	158
Immersion tube, for continuous drainages, not scaled	M16395	317
Immersion tube, for water vacuum gauge of for VarioVac® T + VarioAir® T, scaled	MP00569	318
Incubator cover for Babyleo® TN500, reusable	MP03730	347, 372
Incubator cover for Caleo®, reusable	MP01423	351, 372
Incubator cover for Isolette® C2000, 8000 and 8000 plus, reusable	MP01424	353, 355, 372
Infinity® ID CLIC Absorber 800+, disposable absorber, 1.2 l, 6 pcs.	MX50004	100
Infinity® ID CLIC Absorber Free, disposable absorber, 1.2 l, 6 pcs.	MX50120	100
Infinity® ID, expiratory valve, disposable, 10 pcs.	MP01060/ MP20108	86
Infinity® MPod-QuadHemo, cardiac output, 2 temperatures, up to 4 invasive pressures	MS20725	148, 198
Infinity® PiCCO Gen II SmartPod®	MS23133	161
Infinity® PiCCO Gen II SmartPod® kit	MS32282	161
Infusion bottle holder for Babytherm®, reusable, neonatal	2M16520	359
Infusion holder, reusable, 3 hooks, for poles with Ø 38 mm	2M21514	352, 359
Insert for neonatal flow sensor, PSU, sterilizable, 5 pcs.	8410179	105, 400
Intermediate cable for disposable sensor 3 m, reusable	MS40323	208
Iris port sleeve, disposable, 100 pcs.	MU03876	354, 356, 363, 367, 369
K Kit T-Adapter, Pediatric with Neo connectors, 10 pcs.	MP01087	94, 409
L Leakage valve for Carina® SyncVent, disposable, 5 pcs.	MP00224	86
Light shield curtain for BiliLux, reusable	MP03760	371
Litter bag holder, stainless steel, with base plate and litter bags (100 pcs.)	M24695	339, 433
Litter bags, disposable, for litter bag holder, 100 pcs.	M26240	339, 433
Locking screw, for immersion tube, with gaskets, for VarioVac® T Rail WV (MP00512), VarioAir® T Rail WV (MP00544) and water vacuum gauge (MP00518)	MP00596	318
M M540 Hemo connection cable, 1.5 m	MS23897	198
M540 Hemo connection cable, 3 m	MS23898	198
Masimo rainbow® sensor 8 Lambda (SpO ₂ , SpCO, SpMet), disposable, adult, 10 pcs.	3702841	185
Masimo rainbow® sensor 8 Lambda (SpO ₂ , SpCO, SpMet), disposable, adult/neonate, 10 pcs.	3702831	185
Masimo rainbow® sensor 8 Lambda (SpO ₂ , SpCO, SpMet), disposable, child, 10 pcs.	3702837	185
Masimo rainbow® sensor 8 Lambda (SpO ₂ , SpCO, SpMet), disposable, infant, 10 pcs.	3702834	185
Masimo rainbow® sensor DCI (SpO ₂ , SpCO, SpMet), resusable, adult	3722471	185
Masimo rainbow® sensor DCIP (SpO ₂ , SpCO, SpMet), resusable, pediatric	3722467	185
Masimo® LNCS Trauma sensor, adult, 20 pcs.	MP03191	127, 182, 243
Masimo® LNCS Trauma sensor, infant/pedi, 20 pcs.	MP03193	127, 182, 243
Masimo® LNCS Trauma sensor, neonate, 20 pcs.	MP03192	127, 182, 243, 413
Masimo® rainbow sensor cover, 10 pcs.	MP03203	188

INDEX BY PRODUCT NAME

Masimo® RD rainbow SET sensor, disposable, adult	MS40076	186
Masimo® RD rainbow SET sensor, disposable, infant	MS40082	186
Masimo® RD rainbow SET sensor, disposable, neonate	MS40085	186
Masimo® RD rainbow SET sensor, disposable, pediatric	MS40079	186
Masimo® RD SET DC-IP sensor, pediatric, reusable	MS33741	129, 183, 216, 246
Masimo® RD SET DCI sensor, adult, reusable	MS33740	129, 183, 216, 246
Masimo® RD SET sensor DB-I, adult, reusable	3718728	129, 183, 216
Masimo® RD SET sensor TC-I Ear, pediatric/adult, reusable	MS40067	129, 183, 246
Masimo® RD SET sensor TF-I, pediatric/adult, reusable	MS40073	129, 183, 246
Masimo® RD SET sensor YI Multisite, neonate/pediatric/adult, reusable	MS40070	129, 183, 246
Masimo® RD SET sensor, adult, disposable	MS33726	130, 184, 216, 245
Masimo® RD SET sensor, infant, disposable	MS33728	130, 184, 245
Masimo® RD SET sensor, neonatal preterm, disposable	MS33730	130, 184, 245, 413
Masimo® RD SET sensor, neonatal, disposable	MS33729	130, 184, 245, 413
Masimo® RD SET sensor, pediatric, disposable	MS33727	130, 184, 216, 245
Masimo® Trauma sensor RD SET, adult, disposable	MS40088	131, 184, 247
Masimo® Trauma sensor RD SET, infant/pedi, disposable	MS40091	131, 184, 247
Masimo® Trauma sensor RD SET, neonate, disposable	MS40094	131, 184, 247
Mattress SoftBed, for TI 500 and TI 500 Globe-Trotter, reusable	MP01419	363, 366, 369
Mattress SoftBed, for Babyleo® TN500, reusable	2M60390	347
Mattress SoftBed, for BabyRoo TN300, reusable	MP16297	349, 357
Mattress SoftBed, for Babytherm®, reusable	2M21012	357
Mattress SoftBed, for Caleo®, reusable	MX17012	351
Mattress SoftBed, for Isolette® C2000, 8000 and 8000 plus, reusable	MP01401	353, 355
Mattress SoftBed, for Resuscitaire®, reusable, with scale	MP01404	360
Mattress SoftBed, for Resuscitaire®, reusable, without scale	MP01402	360
MCable Mainstream CO ₂ , for Evita V-family and Babylog VN-family, Infinity® M540, Vista120, Vista120 S, and Delta VF9.1 and higher	6871950	62, 168, 205, 235, 399
Mechanical filter + HME TwinStar® HEPA Plus, disposable, ~90 ml deadspace, 100 pcs.	MP05801	55
Mechanical filter SafeStar® 55 Plus, disposable, ~55 ml deadspace, 100 pcs.	MP05790	57
Mechanical filter SafeStar® 60A Plus, disposable, ~60 ml deadspace, angled, 100 pcs.	MP05795	57
Mechanical filter SafeStar® 90 Plus, disposable, ~90 ml deadspace, 100 pcs.	MP05785	57
Medical examination light VarioLux® with LED technology, AUS	MP00606	340, 434
Medical examination light VarioLux® with LED technology, EU	MP00601	340, 434
Medical examination light VarioLux® with LED technology, GB	MP00603	340, 434
Medical examination light VarioLux® with LED technology, JPN	MP00608	340, 434
Medical examination light VarioLux® with LED technology, US	MP00602	340, 434
Medication nebulizer Aerogen Pro	MP01010	93, 408
Membrane for Evita® 2 dura, Evita® 4, Evita® XL, for expiratory valve	8412015	85
Mount clamp 10 – 40 mm (0.4 – 1.57 in)	MS32150	209
Mount clamp 20 – 60 mm (0.79 – 2.36 in)	MS32151	209
MPod-Dual Hemo,2 invasive blood pressures for M540	MS20783	198
MRI flow sensor extension cable, for MRI patient monitor	8306370	106
Multi-P trolley	MP13021	83, 271, 376
Multi-P trolley equipped	MP13020	83, 271, 376
Multiparameter cable MultiMed® 5 Pod, for ECG, resp., SpO ₂ , and temp., 1.5 m	5950196	111
Multiparameter cable MultiMed® 5 Pod, for ECG, resp., SpO ₂ , and temp., 2.5 m	3368391	111
Multiparameter cable MultiMed® 6 Pod, for ECG, resp., SpO ₂ , and temp., 2.5 m	5191221	111
Multiparameter cable MultiMed® plus OR, for ECG, SpO ₂ , and temp., 2.5 m	MS20094	111
Multiparameter cable MultiMed® plus, for ECG, resp., SpO ₂ , and temp., 2.5 m	MS20093	111

INDEX BY PRODUCT NAME

	Multiparameter cable NeoMed Pod, for ECG, resp., SpO ₂ , temp., and FiO ₂ , 2.5 m	5590539	111, 410
N	N ₂ O CS hose 0.9 m (2.95 ft) NF EN NF	M36725	322
	N ₂ O CS hose 0.9 m (2.95 ft) NIST EN BS	M36728	323
	N ₂ O CS hose 0.9 m (2.95 ft) NIST EN DIN	M36731	319
	N ₂ O CS hose 0.9 m (2.95 ft) NIST GB DIN	M36788	333
	N ₂ O CS hose 1.5 m (4.9 ft) DISS US KAWA	M35809	329
	N ₂ O CS hose 1.5 m (4.9 ft) DISSB EN HAND	M36293	332
	N ₂ O CS hose 1.5 m (4.9 ft) NF EN NF	M35427	322
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN AS	M36156	328
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN BS	M35391	323
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN BS straight	M35400	324
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN Carba	M35473	326
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN DIN	M35492	319
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN NF	M35418	321
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN SS	M35378	325
	N ₂ O CS hose 1.5 m (4.9 ft) NIST EN UNI	M35409	327
	N ₂ O CS hose 1.5 m (4.9 ft) NIST GB DIN	M36789	333
	N ₂ O CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36007	320
	N ₂ O CS hose 1.5 m (4.9 ft) NIST NEUTRAL NF	M35437	322
	N ₂ O CS hose 3 m (9.8 ft) DISS US KAWA	M35808	329
	N ₂ O CS hose 3 m (9.8 ft) DISSB EN HAND	M36294	332
	N ₂ O CS hose 3 m (9.8 ft) NF EN NF	M35428	322
	N ₂ O CS hose 3 m (9.8 ft) NIST EN AS	M36155	328
	N ₂ O CS hose 3 m (9.8 ft) NIST EN BS	M35392	323
	N ₂ O CS hose 3 m (9.8 ft) NIST EN BS straight	M35401	324
	N ₂ O CS hose 3 m (9.8 ft) NIST EN Carba	M35474	326
	N ₂ O CS hose 3 m (9.8 ft) NIST EN DIN	M35493	319
	N ₂ O CS hose 3 m (9.8 ft) NIST EN NF	M35419	321
	N ₂ O CS hose 3 m (9.8 ft) NIST EN SS	M35379	325
	N ₂ O CS hose 3 m (9.8 ft) NIST EN UNI	M35410	327
	N ₂ O CS hose 3 m (9.8 ft) NIST GB BS	M36734	333
	N ₂ O CS hose 3 m (9.8 ft) NIST GB DIN	M36740	333
	N ₂ O CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36008	320
	N ₂ O CS hose 3 m (9.8 ft) NIST NEUTRAL NF	M35438	322
	N ₂ O CS hose 5 m (16.4 ft) DISS US KAWA	M35807	329
	N ₂ O CS hose 5 m (16.4 ft) DISSB EN HAND	M36295	332
	N ₂ O CS hose 5 m (16.4 ft) NF EN NF	M35344	322
	N ₂ O CS hose 5 m (16.4 ft) NIST EN AS	M36154	328
	N ₂ O CS hose 5 m (16.4 ft) NIST EN BS	M35334	323
	N ₂ O CS hose 5 m (16.4 ft) NIST EN BS straight	M35341	324
	N ₂ O CS hose 5 m (16.4 ft) NIST EN Carba	M35347	326
	N ₂ O CS hose 5 m (16.4 ft) NIST EN DIN	M35494	319
	N ₂ O CS hose 5 m (16.4 ft) NIST EN NF	M35331	321
	N ₂ O CS hose 5 m (16.4 ft) NIST EN SS	M35382	325
	N ₂ O CS hose 5 m (16.4 ft) NIST EN UNI	M35359	327
	N ₂ O CS hose 5 m (16.4 ft) NIST GB BS	M36735	333
	N ₂ O CS hose 5 m (16.4 ft) NIST GB DIN	M36741	333
	N ₂ O CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36009	320
	N ₂ O CS hose 5 m (16.4 ft) NIST NEUTRAL NF	M35439	322
	N ₂ O CS hose 8 m (26,25 ft) DISS US KAWA	M35867	329
	N ₂ O CS hose DISS Hand 10L	4199596	331
	N ₂ O CS hose DISS Hand 15L	4199595	331

INDEX BY PRODUCT NAME

N ₂ O CS hose DISS N/DISS N 10L	4185142	331
N ₂ O CS hose DISS N/DISS N 15L	4185143	331
N ₂ O CS hose DISS N/DISS N 20L	4199713	331
Nebulizer head, 10 pcs.	MP01079	94, 409
Nebulizer head, 5 pcs.	MP01078	94, 409
Nebulizer hose	8412985	92
Nebulizer unit with filler cap for Aerogen Pro	MP01030	93, 408
Nellcor® intermediate cable, 1.2 m	MP00748	190, 226, 239
Nellcor® intermediate cable, 3 m	MS20979	190, 226, 239
Neo-mask, silicone, disposable, size L, 10 pcs.	8418619	379, 382
Neo-mask, silicone, disposable, size M, 10 pcs.	8418490	379, 382
Neo-mask, silicone, disposable, size S, 10 pcs.	8418491	379, 382
Neo-mask, TPE, disposable, size 4, 10 pcs.	MP01494	379, 382
Neo-mask, TPE, disposable, size 5, 10 pcs.	MP01495	379, 382
Neo-mask, TPE, disposable, size 6, 10 pcs.	MP01496	379, 382
Neonatal flow sensor ISO 15, PSU, sterilizable	8411130	105, 400
Neonatal Test Lung, reusable, 50 ml, not made with natural rubber latex	MP20050	53, 397
NIBP connection tube, adult, 3.7 m, for Infinity® monitors	MP00953	142, 194
NIBP cuff L, reusable, large adult, 31 – 40 cm/40 cm	MP00918	143, 195
NIBP cuff L, single-patient use, large adult, 31 – 40 cm/40 cm, 10 pcs.	MP00931	144, 196
NIBP cuff L+, reusable, large adult, long, 31 – 40 cm/55 cm	MP00919	143, 195
NIBP cuff L+, single-patient use, large adult long, 31 – 40 cm/55 cm, 10 pcs.	MP00932	144, 196
NIBP cuff M, reusable, adult, 23 – 33 cm/33 cm	MP00915	143, 195
NIBP cuff M, single-patient use, adult, 23 – 33 cm/33 cm, 10 pcs.	MP00928	144, 196
NIBP cuff M+, reusable, adult, long, 23 – 33 cm/43 cm	MP00916	143, 195
NIBP cuff M+, single-patient use, adult long, 23 – 33 cm/43 cm, 10 pcs.	MP00929	144, 196
NIBP cuff M++, reusable, adult, extra long, 23 – 33 cm/53 cm	MP00917	143, 195
NIBP cuff M++, single-patient use, adult extra long, 23 – 33 cm/53 cm, 10 pcs.	MP00930	144, 196
NIBP cuff S, reusable, small adult, 17 – 25 cm/29 cm	MP00913	143, 195
NIBP cuff S, single-patient use, small adult, 17 – 25 cm/29 cm, 10 pcs.	MP00926	144, 196
NIBP cuff XL, reusable, thigh, 38 – 50 cm/50 cm	MP00921	143, 195
NIBP cuff XL, single-patient use, thigh, 38 – 50 cm/50 cm, 10 pcs.	MP00934	144, 196
NIBP cuff XS, reusable, child, 12 – 19 cm/19 cm	MP00912	143, 195
NIBP cuff XS, single-patient use, child, 12 – 19 cm/19 cm, 10 pcs.	MP00925	144, 196
NIBP cuff XXS, reusable, infant, 8 – 13 cm/13 cm	MP00911	143, 195
NIBP cuff XXS, single-patient use, infant, 8 – 13 cm/13 cm, 10 pcs.	MP00924	144, 196
NIBP cuff, adult, 27 cm – 35 cm, reusable, Vista 120	2606155	227, 248
NIBP cuff, child, 16 cm - 21 cm, reusable, Vista 120	2606153	227, 248
NIBP cuff, disposable, neonate, size 1, 3.1 – 5.7 cm, 10 pcs.	2870181	145, 197, 416
NIBP cuff, disposable, neonate, size 2, 4.3 – 8 cm, 10 pcs.	2870199	145, 197, 416
NIBP cuff, disposable, neonate, size 3, 5.8 – 10.9 cm, 10 pcs.	2870207	145, 197, 416
NIBP cuff, disposable, neonate, size 4, 7.1 – 13.1 cm, 10 pcs.	2870215	145, 197, 416
NIBP cuff, disposable, neonate, size 5, 8.3 – 15 cm, 10 pcs.	2870173	145, 197, 416
NIBP cuff, infant, 10 cm – 15 cm, reusable, Vista 120	2606151	227, 248
NIBP cuff, larger adult, 34 cm – 43 cm, reusable, Vista 120	2606156	227, 248
NIBP cuff, neonatal #1, 3 – 6 cm, disposable, 20 pcs., Vista 120	2607034	228
NIBP cuff, neonatal #2, 4 – 8 cm, disposable, 20 pcs., Vista 120	2607035	228
NIBP cuff, neonatal #3, 6 – 11 cm, disposable, 20 pcs., Vista 120	2607036	228
NIBP cuff, neonatal #4, 7 – 13 cm, disposable, 20 pcs., Vista 120	2607037	228
NIBP cuff, neonatal #5, 8 – 15 cm, disposable, 20 pcs., Vista 120	2607038	228
NIBP cuff, small adult, 20.5 cm – 28 cm, reusable, Vista 120	2606154	227, 248
NIBP cuff, small child, 13 cm – 17 cm, reusable, Vista 120	2606152	227, 248
NIBP hose, neonate, 2.4 m	2870298	142, 194, 416

INDEX BY PRODUCT NAME

NIBP tube, 3 m	2601195	227, 248
NIST-NIST Connecting hose 1.2 m	M36053	278
NIV nasal mask ClassicStar® non-vented, disposable, size L	MP01625	66
NIV nasal mask ClassicStar® non-vented, disposable, size M	MP01624	66
NIV nasal mask ClassicStar® non-vented, disposable, size S	MP01623	66
NIV nasal mask ClassicStar® plus non-vented, disposable, size L, 5 pcs.	MP04743	69
NIV nasal mask ClassicStar® plus non-vented, disposable, size M, 5 pcs.	MP04742	69
NIV nasal mask ClassicStar® plus non-vented, disposable, size S, 5 pcs.	MP04741	69
NIV nasal mask ClassicStar® plus non-vented, disposable, size XL, 5 pcs.	MP04744	69
NIV nasal mask ClassicStar® plus vented, disposable, size L, 5 pcs.	MP04753	69
NIV nasal mask ClassicStar® plus vented, disposable, size M, 5 pcs.	MP04752	69
NIV nasal mask ClassicStar® plus vented, disposable, size S, 5 pcs.	MP04751	69
NIV nasal mask ClassicStar® plus vented, disposable, size XL, 5 pcs.	MP04754	69
NIV nasal mask ClassicStar® vented, disposable, size L	MP01622	66
NIV nasal mask ClassicStar® vented, disposable, size M	MP01621	66
NIV nasal mask ClassicStar® vented, disposable, size S	MP01620	66
NIV oronasal mask ClassicStar® AAV, disposable, size L	MP01572	65
NIV oronasal mask ClassicStar® AAV, disposable, size M	MP01571	65
NIV oronasal mask ClassicStar® AAV, disposable, size S	MP01570	65
NIV oronasal mask ClassicStar® plus AAV without exhalation ports, disposable, size L, 5 pcs.	MP04733	68
NIV oronasal mask ClassicStar® plus AAV without exhalation ports, disposable, size M, 5 pcs.	MP04732	68
NIV oronasal mask ClassicStar® plus AAV without exhalation ports, disposable, size S, 5 pcs.	MP04731	68
NIV oronasal mask ClassicStar® plus AAV without exhalation ports, disposable, size XL, 5 pcs.	MP04734	68
NIV oronasal mask ClassicStar® plus AAV, disposable, size L, 5 pcs.	MP04723	68
NIV oronasal mask ClassicStar® plus AAV, disposable, size M, 5 pcs.	MP04722	68
NIV oronasal mask ClassicStar® plus AAV, disposable, size S, 5 pcs.	MP04721	68
NIV oronasal mask ClassicStar® plus AAV, disposable, size XL, 5 pcs.	MP04724	68
NIV oronasal mask ClassicStar® plus SE, disposable, size L, 5 pcs.	MP04713	68
NIV oronasal mask ClassicStar® plus SE, disposable, size M, 5 pcs.	MP04712	68
NIV oronasal mask ClassicStar® plus SE, disposable, size S, 5 pcs.	MP04711	68
NIV oronasal mask ClassicStar® plus SE, disposable, size XL, 5 pcs.	MP04714	68
NIV oronasal mask ClassicStar® SE, disposable, size L	MP01575	65
NIV oronasal mask ClassicStar® SE, disposable, size M	MP01574	65
NIV oronasal mask ClassicStar® SE, disposable, size S	MP01573	65
NIV oronasal mask NovaStar® plus, AAV, reusable, size L, 3 pcs.	MP05063	70
NIV oronasal mask NovaStar® plus, AAV, reusable, size M, 3 pcs.	MP05062	70
NIV oronasal mask NovaStar® plus, AAV, reusable, size S, 3 pcs.	MP05061	70
NIV oronasal mask NovaStar® plus, AAV, reusable, size XL, 3 pcs.	MP05064	70
NIV oronasal mask NovaStar® plus, AAV, without exhalation ports, reusable, size L, 3 pcs.	MP05073	70
NIV oronasal mask NovaStar® plus, AAV, without exhalation ports, reusable, size M, 3 pcs.	MP05072	70
NIV oronasal mask NovaStar® plus, AAV, without exhalation ports, reusable, size S, 3 pcs.	MP05071	70
NIV oronasal mask NovaStar® plus, AAV, without exhalation ports, reusable, size XL, 3 pcs.	MP05074	70
NIV oronasal mask NovaStar® plus, SE, reusable, size L, 3 pcs.	MP05053	70
NIV oronasal mask NovaStar® plus, SE, reusable, size M, 3 pcs.	MP05052	70
NIV oronasal mask NovaStar® plus, SE, reusable, size S, 3 pcs.	MP05051	70
NIV oronasal mask NovaStar® plus, SE, reusable, size XL, 3 pcs.	MP05054	70

INDEX BY PRODUCT NAME

NIV total-face mask ClassicStar® AAV without exhalation ports, disposable, size L	MP02918	67
NIV total-face mask ClassicStar® AAV without exhalation ports, disposable, size M	MP02917	67
NIV total-face mask ClassicStar® AAV without exhalation ports, disposable, size S	MP02916	67
NIV total-face mask ClassicStar® AAV, disposable, size L	MP02915	67
NIV total-face mask ClassicStar® AAV, disposable, size M	MP02914	67
NIV total-face mask ClassicStar® AAV, disposable, size S	MP02913	67
NIV total-face mask ClassicStar® SE, disposable, size L	MP02912	67
NIV total-face mask ClassicStar® SE, disposable, size M	MP02911	67
NIV total-face mask ClassicStar® SE, disposable, size S	MP02910	67
NMT accelerometer Trident®	MS15085	166
NMT electrode cable Trident®, 2-lead, 0.3 m	MS13218	166
NMT hand adapter Trident®	MS15086	166
NMT intermediate cable for Infinity® Trident (NMT) SmartPod®, 1.9 m	7872174	166
NMT thermistor Trident®	MS15053	166
NovaStar® plus Headgear inclusive magnetic clips, disposable, 3 pcs.	MP05080	72
O O ₂ compact flowmeter Clickflow, DIN CS, 1 – 16 l/min, 6 mm nozzle	MP13012	269, 421
O ₂ compact flowmeter Clickflow, DIN CS, 1 – 16 l/min, 9/16"	MP13011	269, 421
O ₂ compact flowmeter Clickflow, rail, NIST, 1 – 16 l/min, 9/16"	MP13013	270
O ₂ CS hose 0.9 m (2.95 ft) NF EN NF	M36723	322
O ₂ CS hose 0.9 m (2.95 ft) NIST EN BS	M36726	323
O ₂ CS hose 0.9 m (2.95 ft) NIST EN DIN	M36729	319
O ₂ CS hose 0.9 m (2.95 ft) NIST GB DIN	M36784	333
O ₂ CS hose 1.5 m (4.9 ft) DISS US KAWA	M35806	329
O ₂ CS hose 1.5 m (4.9 ft) DISSB EN HAND	M36290	332
O ₂ CS hose 1.5 m (4.9 ft) NF EN NF	M35425	322
O ₂ CS hose 1.5 m (4.9 ft) NIST EN AS	M36153	328
O ₂ CS hose 1.5 m (4.9 ft) NIST EN BS	M35389	323
O ₂ CS hose 1.5 m (4.9 ft) NIST EN BS straight	M35398	324
O ₂ CS hose 1.5 m (4.9 ft) NIST EN Carba	M35469	326
O ₂ CS hose 1.5 m (4.9 ft) NIST EN DIN	M35489	319
O ₂ CS hose 1.5 m (4.9 ft) NIST EN NF	M35416	321
O ₂ CS hose 1.5 m (4.9 ft) NIST EN SS	M35375	325
O ₂ CS hose 1.5 m (4.9 ft) NIST EN UNI	M35407	327
O ₂ CS hose 1.5 m (4.9 ft) NIST GB DIN	M36785	333
O ₂ CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36004	320
O ₂ CS hose 1.5 m (4.9 ft) NIST NEUTRAL NF	M35434	321
O ₂ CS hose 3 m (9.8 ft) DISS US KAWA	M35805	329
O ₂ CS hose 3 m (9.8 ft) DISSB EN HAND	M36291	332
O ₂ CS hose 3 m (9.8 ft) NF EN NF	M35426	322
O ₂ CS hose 3 m (9.8 ft) NIST EN AS	M36152	328
O ₂ CS hose 3 m (9.8 ft) NIST EN BS	M35390	323
O ₂ CS hose 3 m (9.8 ft) NIST EN BS straight	M35399	324
O ₂ CS hose 3 m (9.8 ft) NIST EN Carba	M35472	326
O ₂ CS hose 3 m (9.8 ft) NIST EN DIN	M35490	319
O ₂ CS hose 3 m (9.8 ft) NIST EN NF	M35417	321
O ₂ CS hose 3 m (9.8 ft) NIST EN SS	M35376	325
O ₂ CS hose 3 m (9.8 ft) NIST EN UNI	M35408	327
O ₂ CS hose 3 m (9.8 ft) NIST GB BS	M36732	333
O ₂ CS hose 3 m (9.8 ft) NIST GB DIN	M36738	333

INDEX BY PRODUCT NAME

O ₂ CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36005	320
O ₂ CS hose 3 m (9.8 ft) NIST NEUTRAL NF	M35435	321
O ₂ CS hose 5 m (16.4 ft) DISS US KAWA	M35804	329
O ₂ CS hose 5 m (16.4 ft) DISSB EN HAND	M36292	332
O ₂ CS hose 5 m (16.4 ft) NF EN NF	M35342	322
O ₂ CS hose 5 m (16.4 ft) NIST EN AS	M36151	328
O ₂ CS hose 5 m (16.4 ft) NIST EN BS	M35333	323
O ₂ CS hose 5 m (16.4 ft) NIST EN BS straight	M35339	324
O ₂ CS hose 5 m (16.4 ft) NIST EN Carba	M35345	326
O ₂ CS hose 5 m (16.4 ft) NIST EN DIN	M35491	319
O ₂ CS hose 5 m (16.4 ft) NIST EN NF	M35330	321
O ₂ CS hose 5 m (16.4 ft) NIST EN SS	M35377	325
O ₂ CS hose 5 m (16.4 ft) NIST EN UNI	M35357	327
O ₂ CS hose 5 m (16.4 ft) NIST GB BS	M36733	333
O ₂ CS hose 5 m (16.4 ft) NIST GB DIN	M36739	333
O ₂ CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36006	320
O ₂ CS hose 5 m (16.4 ft) NIST NEUTRAL NF	M35436	321
O ₂ CS hose 8 m (26,25 ft) DISS US KAWA	M35866	329
O ₂ CS hose DISS Hand 10L	4199592	330
O ₂ CS hose DISS Hand 15L	4199591	330
O ₂ CS hose DISS N/DISS N 10L	4185140	330
O ₂ CS hose DISS N/DISS N 15L	4185141	330
O ₂ CS hose DISS N/DISS N 20L	4199714	330
O ₂ flowmeter for anesthesia, 16 l/min, 90°, rail, with NIST black hose 1 m, connector 9/16-18 UNF to 6 mm	MP04579	266
O ₂ flowmeter for anesthesia, 16 l/min, 90°, rail, with NIST ISO hose 1.2 m, connector 9/16-18 UNF to 6 mm	MP04578	266
O ₂ flowmeter for anesthesia, 16 l/min, rail, with NIST black hose 1 m, connector 9/16-18 UNF to 6 mm	MP04577	266
O ₂ flowmeter for anesthesia, 16 l/min, rail, with NIST ISO hose 1.2 m, connector 9/16" UNF to 6 mm	MP04576	266
O ₂ flowmeter, 10 l/min, CS, CARBA, including adapter to 9/16" for sterile water	MP04554	262
O ₂ flowmeter, 10 l/min, rail fitting with NIST gas inlet, M34×1.5, adapter 9/16"	MP04556	264
O ₂ flowmeter, 16 l/min, CS, British Standard, ISO, M34×1.5, adapter 9/16"	MP04565	262
O ₂ flowmeter, 16 l/min, CS, CARBA, including adapter to 9/16" for sterile water	MP04563	262
O ₂ flowmeter, 16 l/min, CS, DIN, ISO and neutral color, M34×1.5 and 9/16" for accessories and sterile water systems	MP04558	261
O ₂ flowmeter, 16 l/min, CS, French Standard, ISO, M34×1.5, adapter 9/16"	MP04566	262
O ₂ flowmeter, 16 l/min, CS, long CS plug (6 cm; 2.4 in), DIN, ISO and neutral color, M34×1.5 and 9/16" for accessories and sterile water systems	MP04561	261
O ₂ flowmeter, 16 l/min, CS, Swedish Standard, M34x1.5 and Adapter 9/16"	MP13040	262
O ₂ flowmeter, 16 l/min, rail NIST, packed with hose 1.5 m, DIN, ISO, M34×1.5 and 9/16" for accessories and sterile water systems	MP04573	264
O ₂ flowmeter, 16 l/min, rail NIST, packed with hose 1.5 m, DIN, neutral color, M34×1.5 and 9/16" for accessories and sterile water systems	MP04575	264
O ₂ flowmeter, 16 l/min, rail, 90°, NIST gas inlet, M34×1.5 and 9/16" for accessories and sterile water systems	MP04570	265
O ₂ flowmeter, 16 l/min, rail, NIST gas inlet, M34×1.5, adapter 9/16"	MP04567	264, 420
O ₂ flowmeter, 32 l/min, CS, DIN, ISO and neutral color, M34×1.5 and 9/16" for accessories and sterile water systems	MP04571	262
O ₂ flowmeter, 32 l/min, rail, NIST, M34×1.5 and adapter 9/16"	MP04580	264
O ₂ flowmeter, 4 l/min, CS, DIN, ISO and neutral color, M34×1.5 and 9/16" for accessories and sterile water systems	MP04551	261, 418

INDEX BY PRODUCT NAME

O ₂ flowmeter, 4 l/min, rail with NIST gas inlet, M34×1.5 and 9/16" for accessories and sterile water systems	MP04552	264, 420
O ₂ flowmeter, 4 l/min, rail, 90°, NIST gas inlet, M34×1.5 and 9/16" for accessories and sterile water systems	MP04553	265, 420
O ₂ sensor "Oxycell A" for Atlan	6872666	107
O ₂ sensor capsule, for mainstream measurement	6850645	107, 401
O ₂ sensor for Crossvent 2i+	MU24402	364
O ₂ sensor for Isolette® C2000, 8000 and 8000 plus, reusable, 2 pcs.	MU24903	354, 356
O ₂ sensor for MX 300 and MX 300-i, reusable	MP01442	88, 107
O ₂ sensor housing Oxydig, for O ₂ sensor capsule, with spiral cable	6850250	107, 401
O ₂ sensor Oxycell, for Babyleo® TN500	6873065	107, 348, 401
O ₂ sensor OxyTrace INCU, for Caleo®	MX01050	107, 352
O ₂ sensor OxyTrace VE, for Savina® ventilator	MX01049	107
O ₂ twin compact flowmeter Clickflow, DIN, 2 × 1 – 16 l/min, 9/16"	MP13014	270
O ₂ twin flowmeter, 2 × 10 l/min, CS, CARBA, 2 × M34×1.5, including adapter to 9/16" for sterile water	MP04555	263
O ₂ twin flowmeter, 2 × 10 l/min, rail fitting with NIST gas inlet, 2 × M34×1.5, adapter 9/16"	MP04557	263
O ₂ twin flowmeter, 2 × 16 l/min, CS, DIN, ISO, neutral color, M34×1.5 and 9/16" for accessories and sterile water systems	MP04560	263, 419
O ₂ twin flowmeter, 2 × 16 l/min, rail fitting with NIST gas inlet, 2 × M34×1.5, adapter 9/16"	MP04569	263, 419
O2Star™ Connector 22 mm ID / 22 mm ID, disposable, 5 pcs.	MP01942	81
O2Star™ Connector kit, 6 mm OD/22 mm OD + 6 mm OD/22 mm ID with connection tube, disposable, 5 pcs.	MP01940	81
O2Star™ Face tent, disposable, size M/L, 50 pcs.	MP01925	79
O2Star™ Face tent, disposable, size S, 50 pcs.	MP01933	79
O2Star™ Medium concentration oxygen mask with connection tube 2.1 m, disposable, size M/L, 50 pcs.	MP01921	79
O2Star™ Medium concentration oxygen mask with connection tube 2.1 m, disposable, size S, 50 pcs.	MP01931	79
O2Star™ Non-rebreather oxygen mask with connection tube 2.1 m, disposable, size M/L, 50 pcs.	MP01920	79
O2Star™ Non-rebreather oxygen mask with connection tube 2.1 m, disposable, size S, 50 pcs.	MP01930	79
O2Star™ Oxygen and aerosol mask with 22 mm OD, disposable, size M/L, 50 pcs.	MP01926	79
O2Star™ Oxygen and aerosol mask with 22 mm OD, disposable, size S, 50 pcs.	MP01934	79
O2Star™ Oxygen connection tube, 2.1 m, disposable, 25 pcs.	MP01927	81
O2Star™ Oxygen connection tube, 2.1 m, with 9/16" connector, disposable, 25 pcs.	MP01944	81
O2Star™ Oxygen connection tube, 4 m, disposable, 25 pcs.	MP01943	81
O2Star™ Oxygen connection tube, 4 m, with 9/16" connector, disposable, 25 pcs.	MP01945	81
O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size L, 50 pcs.	MP02926	80
O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size M, 50 pcs.	MP02925	80, 406
O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size S, 50 pcs.	MP02924	80, 406
O2Star™ Oxygen nasal cannula curved with connection tube 2.1 m, disposable, size XL, 50 pcs.	MP02927	80
O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size L, 50 pcs.	MP02922	80
O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size M, 50 pcs.	MP02921	80, 406
O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size S, 50 pcs.	MP02920	80, 406

INDEX BY PRODUCT NAME

O2Star™ Oxygen nasal cannula straight with connection tube 2.1 m, disposable, size XL, 50 pcs.	MP02923	80
O2Star™ Tracheostomy mask, disposable, size M/L, 50 pcs.	MP01924	79
O2Star™ Tracheostomy mask, disposable, size S, 50 pcs.	MP01932	79
O2Star™ Venturi mask kit, adjustable venturi valve, disposable, size M/L, 25 pcs.	MP01922	79
O2Star™ Venturi mask kit, adjustable venturi valve, disposable, size S, 25 pcs.	MP01935	79
OPC NIBP Cuff, adult, 23 cm – 33 cm, reusable, 5 pcs., Vista 120 SC	2600937	249
OPC NIBP Cuff, child, 12 cm – 19 cm, reusable, 5 pcs., Vista 120 SC	2600935	249
OPC NIBP Cuff, large adult, 31 cm – 40 cm, reusable, 5 pcs., Vista 120 SC	2600938	249
OPC NIBP Cuff, small adult, 17 cm – 25 cm, reusable, 5 pcs., Vista 120 SC	2600936	249
Optic serial RS 232 cable (1 meter) to connect the device to a Dräger monitorAC	MS32156	209
Optic serial RS 232 cable (2.5 meter) to connect the device to a Dräger monitorAC	MS40021	209
Oxator head OxyLine, pressure reducer G ³ / ₄ " , 200 bar, NIST, 2 × DIN-CS, ejector VarioO2	MP01103	276, 278
Oxett adapter plate, to mount Oxett on bed in higher position	2M86205	277
Oxett carrying frame, for Oxett head+O ₂ cylinder up to 2.5 l, with secretion jar and rail for suction ejector	MP04110	277
Oxett head OxyLine, pressure reducer G ³ / ₄ " , 200 bar, NIST, flow regulator, DIN-CS, ejector VarioO2	MP01102	276, 277
Oxett head OxyLine, pressure reducer G ³ / ₄ " , 200 bar, NIST, flowmeter, DIN-CS, ejector VarioO2	MP01101	276, 277
Oxidem 3000 demand valve incl. adapter (2120496, 2120488)	2111764	95
Oxymixer configuration	MP04200	83, 271, 374
Oxymixer High Flow	MP04201	83, 271
Oxymixer High Flow with monitor	MP04202	83, 271
Oxymixer Low Flow	MP04203	271, 374
Oxymixer Low Flow with 1 – 15 l flowmeter right, rail clamp	MP04229	271, 374
Oxymixer Low Flow with monitor	MP04204	271, 374
Oxymixer oxygen sensor (spare part)	MP04221	84, 272, 375
Oxymixer sound damper 22 and 15 mm	MP04219	84, 272, 375
Oxymixer sound damper 9/16"	MP04218	84, 272, 375
P Pack2Go® anesthesia pack, A, with anesthesia set + watertrap, 10 pcs.	MP20379	90
Pack2Go® Carina® NIV ventilation pack with oronasal mask ClassicStar® SE, size L, 10 pcs.	MP20374	90
Pack2Go® Carina® NIV ventilation pack with oronasal mask ClassicStar® SE, size M, 10 pcs.	MP20373	90
Pack2Go® Carina® NIV ventilation pack with oronasal mask ClassicStar® SE, size S, 10 pcs.	MP20372	90
Pack2Go® Evita® invasive ventilation pack with breathing circuit VentStar®, basic, + cuvette, 10 pcs.	MP20370	89
Pack2Go® Evita® invasive ventilation pack with breathing circuit VentStar®, basic, 1.8 m, 10 pcs.	MP20371	89
Pack2Go® Evita® invasive ventilation pack with breathing circuit VentStar®, coaxial, + cuvette, 10 pcs.	MP20377	89
Pack2Go® Evita® invasive ventilation pack with breathing circuit VentStar®, coaxial, 1.8 m, 10 pcs.	MP20378	89
Pack2Go® invasive ventilation pack with breathing circuit VentStar®, basic, 1.8 m, 10 pcs.	MP20376	90
Pack2Go® invasive ventilation pack with breathing circuit VentStar®, coaxial, 1.8 m, 10 pcs.	MP20375	90
Parking holder for CO ₂ sensor	8416352	169, 205
Patient parameter unit Infinity® BISx SmartPod®	MS14796	164
Patient parameter unit Infinity® HemoMed Pod	5588822	148

INDEX BY PRODUCT NAME

Pediatric nasal mask MiniMe2 non-vented, disposable, size L	MP01556	71, 385
Pediatric nasal mask MiniMe2 non-vented, disposable, size S	MP01555	71, 385
Pediatric nasal mask MiniMe2 vented, disposable, size L	MP01554	71, 385
Pediatric nasal mask MiniMe2 vented, disposable, size S	MP01553	71, 385
PEEP valve for Oxylog® 1000, disposable, adjustable from 5 – 20 mbar	2M86832	45
PEEP valve for Oxylog® 1000, reusable, adjustable from 0 – 10 mbar	8407475	45
Phototherapy mask Eyemax 2, disposable, micro, head circumference 20-25 cm, 20 pcs.	MP03770	371
Phototherapy mask Eyemax 2, disposable, preemie, head circumference 26-32 cm, 20 pcs.	MP03771	371
Phototherapy mask Eyemax 2, disposable, regular, head circumference 33-38 cm, 20 pcs.	MP03772	371
Plug nozzle water vacuum gauge, for VarioVac® T Rail WV (MP00512), VarioAir® T Rail WV (MP00544) and water vacuum gauge (MP00518)	MP00566	318
Pneum. Medication Nebulizer	8412935	92
Pneum. Medication Nebulizer N	8411030	92
Pod adapter SensoNor, for connecting pressure transducer cables, 7-pin	MS23598	154
Pod adapter, 10-pin, for pressure transducer	3375958	154, 161
Pod adapter, for Edwards pressure transducer	5196980	154
Pod adapter, for pressure transducer Argon/Merit Medical (former Becton Dickinson/Ohmeda)	3375941	154
Pod adapter, for pressure transducer Abbott/Medex/ICU Medical	5196998	154
Pod connection cable, 3 m	3368425	112, 147, 198
Pod connection cable, 30 cm	7257988	112, 198
Pod connection cable, 5 m	5195198	112, 198
Polaris® 100 with sterilizable handle (G92099), 1 pc.	G35332	343
Polaris® 100/200 with sterilizable handle E (G36993), 1 pc.	G35325	343
Polaris® 200 with sterilizable handle (G92099), 1 pc.	G35614	343
Polaris® 600 with sterilizable handles (G92099 or G36993), 1 pc.	G35324	343
Polaris® Multimedia (display adapter) with sterilizable handles (G92099 or G36993), 1 pc.	G35333	343
Power cable AUS, 1.8 m, 2 pin	MP00626	340, 434
Power cable EU, 1.8 m, 2 pin	MP00591	434, 340
Power cable GB, 1.8 m, 2 pin	MP00593	434, 340
Power cable US/JPN, 1.8 m, 2 pin	MP00592	434, 340
Pressure difference flow sensor, PPSU, sterilizable	8412034	46, 105
Pressure gauge -100 kPa MRI, double scale kPa/mmHg, -100 kPa (-1000 mbar)	MP03589	301
Pressure gauge -60 kPa MRI, double scale kPa/mmHg, -60 kPa (-600 mbar)	MP03587	301
Pressure gauge for VarioVac® + VarioAir®, double scale kPa/mmHg, -100 kPa (-1000 mbar)	MP00578	301
Pressure gauge for VarioVac® + VarioAir®, double scale kPa/mmHg, -16 kPa (-160 mbar)	MP00579	317
Pressure gauge for VarioVac® + VarioAir®, double scale kPa/mmHg, -60 kPa (-600 mbar)	MP00580	301
Pressure gauge for VarioVac® + VarioAir®, double scale kPa/mmHg, -8 kPa (-80 mbar)	MP03597	317
Pressure reducer OxyLine, customized configuration of pressure reducers	MP00100	273, 277, 278
Pressure reducer w/ flowmeter OxyLine, G¾", 200 bar, 16 l/min, M34×1.5, 9/16", nozzle	MP01104	276
Pressure reducer with NIST OxyLine, G¾", for 200 bar, NIST	MP01105	276
Prong BabyFlow®, standard, disposable, size L, 10 pcs.	8418531	378, 381
Prong BabyFlow®, standard, disposable, size M, 10 pcs.	8418416	378, 381
Prong BabyFlow®, standard, disposable, size S, 10 pcs.	8418605	378, 381
Prong BabyFlow®, standard, disposable, size XL, 10 pcs.	8418417	378, 381
Prong BabyFlow®, standard, disposable, size XS, 10 pcs.	8418415	378, 381

INDEX BY PRODUCT NAME

Prong BabyFlow [®] , standard, disposable, size XXL, 10 pcs.	8418617	378, 381
Prong BabyFlow [®] , wide, disposable, size 3, 10 pcs.	8418603	378, 381
Prong BabyFlow [®] , wide, disposable, size 4, 10 pcs.	8418607	378, 381
Prong BabyFlow [®] , wide, disposable, size 5, 10 pcs.	8418609	378, 381
Prong BabyFlow [®] , wide, disposable, size 6, 10 pcs.	8418611	378, 381
Protective cover for general-purpose temperature probe, contains latex, 10 pcs.	7014616	140, 192, 229, 231
Protective cover for general-purpose temperature probe, latexfree, 100 pcs.	MS40243	140, 192, 229, 231
R Rail clamp for F&P humidifier and Seattle PAP bubble box	8411074	29, 377, 390
Rail clamp lever type, for attachment of devices to 25 × 10 mm standard rail, for compact suction unit and secretion jar set	M26842	279, 296, 339, 432
Rail clamp, for attachment to poles of different diameters, rotatable by 90° to left and right, length of standard rail (profile 25 × 10 mm): 20 cm	2M85274	339, 341, 377, 435
Rail claw for lamp	MP00615	341, 435
Rail holder, for 38 mm tube, standard rail profile 25 × 10 mm, 40 cm, length: 40 cm, weight: 0.9 kg, max. load: 5 kg	2M86185	280
Rail, 100 mm for GCX profile	4118555	348
RD SET intermediate cable for Masimo [®] rainbow SET Mcable, 1.5 m	MS33738	183, 185, 186
RD SET intermediate cable for Masimo [®] rainbow SET Mcable, 3.6 m	MS33739	183, 185, 186
RD SET intermediate cable for Masimo [®] SET Mcable, 1.5 m	MS33736	183, 244
RD SET intermediate cable for Masimo [®] SET Mcable, 3.6 m	MS33737	183, 244
Recording paper for Infinity [®] R 50, 50 mm, approx. 30 m, 10 rolls	4711201	173, 417
Red cap, 7 mm	8401645	39, 51, 392
Restraining strap, disposable, 10 pcs.	MU06011	363, 368, 370
Reusable Breathing circuit (A), 900MR784	MP17129	28
Reusable Breathing circuit (N), 900MR782	MP17128	28, 392
Rubber plug, for immersion tube of water vacuum gauge	2M85139	317
S Sample gas return line for Fabius[®], with filter	M32692	25, 172
Sample line, disposable, 10 pcs.	8290286	25, 172
Sealing rings for spiral hoses, 20 pcs.	M07152	95
Seattle PAP plus	MP02260	376
Secretion control glass	M07582	300
Secretion jar 3 L Medela	MP03230	305
Secretion jar set with clamp	MP00646	295
Secretion jar set, with 700 ml secretion jars (2 pcs.) and 1.5 m suction hose, maximum capacity up to the overflow safeguard 2 × 1000 ml, scale in steps of 10 ml up to max. 700 ml	2M85056	279, 294, 429
Secretion jar set, with rail holder and secretion jar 700 ml, for continuous drainage units	2M85163	315
Secretion jar VacuSmart [®] Gel, PE (Polyethylene), disposable, 700 ml, 25 pcs.	MX23025	297, 430
Secretion jar VacuSmart [®] , PE (Polyethylene), disposable, 700 ml, 25 pcs.	MX23023	297, 430
Secretion jar, PSU (Polysulfone), reusable, 300 ml, autoclavable up to 134°C	2M86207	300, 431
Secretion jar, PSU (Polysulfone), reusable, 700 ml, autoclavable up to 134°C	2M85594	300, 431
Set of 5 red caps, 7 mm	8402953	44, 393
Set VacuSmart [®] Gel and Tube, 25 pcs.	MP03211	297, 430
Shut-off valve for VarioVac [®] , for temporary opening and closing of suction flow	MP00570	300
Silicone hose 6 × 2.5, 2 m	5001330	299
Silicone mask, reusable, size 0, small neonate	MP01910	77, 361, 404
Silicone mask, reusable, size 00, premature baby	MP02900	77, 361, 404
Silicone mask, reusable, size 1, large neonate	MP01911	77, 361, 404
Silicone mask, reusable, size 2, small child	MP01912	77
Silicone mask, reusable, size 3, large child	MP01913	77
Silicone mask, reusable, size 4, small adult	MP01914	77

INDEX BY PRODUCT NAME

Silicone mask, reusable, size 5, large adult	MP01915	77
Simple cylinder trolley, for one gas cylinders up to a volume of 11 L	MP00038	278, 280
Skin temperature probe Premi-Probe® 1, disposable, 10 pcs.	MU06933	363, 366, 369
Skin temperature probe Premi-Probe® 1, reusable	MU06951	363, 366, 369
Skin temperature probe Premi-Probe® 3, disposable, 10 pcs.	MU11228	360
Skin temperature probe Premi-Probe® 3, reusable	MU11225	360
Skin temperature probe Premi-Probe® 5, disposable, 10 pcs.	MU12551	355
Skin temperature probe Premi-Probe® 5, reusable	MU12533	355
Skin temperature probe Premi-Probe® for Isolette® 8000 plus, disposable, white, peripheral, 10 pcs.	MU26042	353
Skin temperature probe Premi-Probe® for Isolette® 8000 plus, disposable, yellow, central, 10 pcs.	MU26041	353
Skin temperature probe ThermoTrace, disposable, white, peripheral, 5 pcs.	MX11001	347, 349, 351, 357
Skin temperature probe ThermoTrace, disposable, yellow, central, 5 pcs.	MX11000	347, 349, 351, 357
Skin temperature probe, reusable, with 7-pin connector, 1.5 m	4329822	140, 192, 231, 415
Skin temperature probe, reusable, with 7-pin connector, 3 m	5204669	140, 192, 231, 415
Skin temperature probe, single-patient use, adult, 7-pin connector, 1.6 m	MP00995	139, 191, 230
Skin temperature probe, single-patient use, adult, 7-pin connector, 3 m	MP00996	139, 191, 230
Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 1.6 m	MP00997	139, 191, 230, 415
Skin temperature probe, single-patient use, child/neonatal, 7-pin connector, 3 m	MP00998	139, 191, 230, 415
Soda lime Drägersorb® 800+, for CO ₂ absorption, canister, 5 l, 2 pcs.	MX00001	99
Soda lime Drägersorb® Free, for CO ₂ absorption, canister, 5 l, 2 pcs.	MX50050	99
Soft connector N, disposable, 20 pcs.	MP03826	378, 380
SpO ₂ adhesive sensor Dräger, cloth, disposable, baby, 24 pcs.	MS16447	123, 411
SpO ₂ adhesive sensor Dräger, cloth, disposable, neonate, 24 pcs.	MS16446	123, 411
SpO ₂ adhesive sensor Dräger, foam, disposable, adult, 24 pcs.	MS16445	122, 214
SpO ₂ adhesive sensor Dräger, foam, disposable, child, 24 pcs.	MS16444	122, 214
SpO ₂ adhesive sensor Dräger, vinyl, disposable, adult, 24 pcs.	MS16449	122, 214
SpO ₂ adhesive sensor Dräger, vinyl, disposable, child, 24 pcs.	MS16448	122, 214
SpO ₂ adhesive sensor Masimo® LNCS Axt, disposable, adult, 20 pcs.	MP00790	126, 181, 215, 242
SpO ₂ adhesive sensor Masimo® LNCS Inf, disposable, baby, 20 pcs.	MP00791	126, 181, 215, 242
SpO ₂ adhesive sensor Masimo® LNCS Neo, disposable, neonate, 20 pcs.	MP00792	126, 181, 242, 365, 413
SpO ₂ adhesive sensor Masimo® LNCS NeoPt, disposable, premature baby, 20 pcs.	MP00794	126, 181, 242, 365, 413
SpO ₂ adhesive sensor Masimo® LNCS Pdt, disposable, child, 20 pcs.	MP00793	126, 181, 215, 242
SpO ₂ adhesive sensor Nellcor® Oximax® Max-A, disposable, adult, cable 46 cm, 24 pcs.	MX50065	136, 190, 218
SpO ₂ adhesive sensor Nellcor® Oximax® Max-AL, disposable, adult, cable 91 cm, 24 pcs.	MX50071	136, 190, 218
SpO ₂ adhesive sensor Nellcor® Oximax® Max-Fast, disposable, application on forehead, 24 pcs.	MX50070	137, 226, 239
SpO ₂ adhesive sensor Nellcor® Oximax® Max-I, disposable, child, 24 pcs.	MX50067	136, 190, 411
SpO ₂ adhesive sensor Nellcor® Oximax® Max-N, disposable, neonate/adult, 24 pcs.	MX50068	136, 190, 218, 411
SpO ₂ adhesive sensor Nellcor® Oximax® Max-P, disposable, child, 24 pcs.	MX50066	136, 190, 218
SpO ₂ adhesive sensor Nellcor® Oximax® Max-R, disposable, application on nose, 24 pcs.	MX50069	137
SpO ₂ adhesive square Masimo®, 12 pcs.	7496875	132, 187, 217
SpO ₂ adhesive wrap Masimo® Multisite-Wrap, for LNOP YI SpO ₂ sensor, 100 pcs.	7496917	132, 187, 217, 414
SpO ₂ adhesive wrap Masimo® Petite-Wrap, for LNOP YI SpO ₂ sensor, 100 pcs.	7496909	132, 187, 217, 414
SpO ₂ adhesive wrap Masimo® Posey-Wrap, for NeoPt SpO ₂ sensor, 12 pcs.	7496933	414
SpO ₂ adhesive wrap Masimo®, for NeoPt SpO ₂ sensor, 100 pcs.	7496925	414
SpO ₂ extension cable, 2 m, reusable	2606487	225, 238
SpO ₂ finger sensor Dräger, reusable	MS13235	122, 214
SpO ₂ finger sensor Masimo® LNCS DCI, reusable, adult	MP00796	125, 180, 215, 241

INDEX BY PRODUCT NAME

SpO ₂ finger sensor Masimo® LNCS DCIP, reusable, child	MP00795	125, 180, 215, 241
SpO ₂ finger sensor Nellcor® DURA DS 100A	7262764	135, 190, 218, 226, 239
SpO ₂ finger sensor, adult, 1 m, reusable	2606484	225, 238
SpO ₂ finger sensor, adult, 2.5 m, reusable	2606483	225, 238
SpO ₂ foam wrap Masimo® Foam-Wrap, neonate, LNOP YI SpO ₂ sensor, 12 pcs.	7496891	132, 187, 217, 414
SpO ₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 1 m	3368433	121, 124, 134, 410
SpO ₂ intermediate cable Dräger/Nellcor®, for MultiMed® 5/6 and NeoMed pods, 2 m	3375834	121, 124, 134, 410
SpO ₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 1.2 m	MS18683	121, 124, 134, 213, 218, 410
SpO ₂ intermediate cable Dräger/Nellcor®, for MultiMed® Plus/Plus OR, 3 m	MS17330	121, 124, 134, 213, 218, 410
SpO ₂ intermediate cable Masimo® LNCS, 1.2 m	MP00800	179, 240
SpO ₂ intermediate cable Masimo® LNCS, 1.2 m	MS24303	124, 215, 412
SpO ₂ intermediate cable Masimo® LNCS, 3 m	MP00797	179, 240
SpO ₂ intermediate cable Masimo® LNCS, 3 m	MS17522	124, 215, 412
SpO ₂ intermediate cable Masimo® LNCS, 4.3 m	MP03194	179, 240
SpO ₂ intermediate cable Masimo® RD SET, 1.2 m (3.9 ft)	MS34436	124, 216
SpO ₂ intermediate cable Masimo® RD SET, 3 m (9.8 ft)	MS34437	124, 216
SpO ₂ intermediate cable Masimo® Red LNCS, 1.2 m	MP02991	179
SpO ₂ intermediate cable Masimo® Red LNCS, 3 m	MP02992	179
SpO ₂ intermediate cable Masimo® Red LNCS, 4.3 m	MP03196	179
SpO ₂ Masimo® sensor LNCS-DBI, reusable, adult	MP02993	125, 180, 215, 241
SpO ₂ Masimo® sensor LNCS-TC-I Ear, reusable, pediatric/adult	MP00788	125, 180, 215, 241
SpO ₂ Masimo® sensor LNCS-TF-I, reusable, pediatric/adult	MP00799	125, 180, 215, 241
SpO ₂ Masimo® sensor LNCS-YI Multisite, reusable, neonate/pediatric/adult	MP00789	125, 180, 215, 241, 413
SpO ₂ multisite sensor Nellcor® DURA-Y™, reusable	MX01004	135, 226, 239
SpO ₂ Nellcor® FLEXMAX Adt	MS34677	137, 190
SpO ₂ Nellcor® FLEXMAX Pdt	MS34678	137, 190
SpO ₂ sensor, adult, 0.5 m, disposable, 25 pcs.	2606208	225, 238
SpO ₂ sensor, infant, 0.5 m, disposable, 25 pcs.	2606210	225, 238
SpO ₂ sensor, neonate, 0.5 m, disposable, 25 pcs.	2606211	225, 238
SpO ₂ sensor, pediatric, 0.5 m, disposable, 25 pcs.	2606209	225, 238
SpO ₂ silicone Soft-tip sensor, adult, 1 m, reusable	2606485	225, 238
SpO ₂ silicone Soft-tip sensor, pediatric, 1 m, reusable	2606486	225, 238
Suction bottle with filter and tube, 800 ml (27.05 fl. oz.), 20 pcs.	MU10918	350, 354, 356, 362
Suction bottle with holder for Resuscitaire, disposable, 40 ml	MU10916	362
Suction bottle, disposable, 40 ml, 10 pcs.	MU10917	362
Suction canister holder for GCX	3702751	350
Suction hose holder, suitable for various diameters	2M85239	296, 306, 339, 432
Suction hose, silicone, reusable, 1.5 m, for medical suction	M25780	299, 431
Suction set, disposable, VacuSmart®, 700 ml and VacuSmart® Tube, 25 pcs. each	MP00061	297, 430
Suction tube VacuSmart® Tube, PE (Polyethylene), disposable, 2 m, w/ finger tip, 25 pcs.	MX23026	297, 430
Switchover valve Medela	MP03232	305
Swivel table, reusable, for poles with Ø 38 mm	2M21186	352, 359
T T-adapter neonatal for NO therapy, 10 pcs.	8422215	23, 388
T-Adapter, adult, 10 pcs.	MP01085	94
T-Adapter, pediatric, 10 pcs.	MP01086	94, 409
T-connector for Aerogen Pro, 12 mm, with plug, 5 pcs.	MP01037	93, 408
T-connector, disposable, Ø 15 mm, 25 pcs.	MU10902	362
T-piece adapter for MX 300 and MX 300-i, reusable	MP01445	88
T-piece adapter for MX 300/MX 300i and Resuscitaire® breathing circuit, disposable, 20 pcs.	MP01880	361

INDEX BY PRODUCT NAME

T-piece AIR NIST	M36056	335
T-piece connector, reusable	8411044	51
T-piece set with plug for Aerogen Pro, N, 15 mm (0.6 in), w/ adapter 10 mm (0.4 in), 5 pcs.	MP01032	93, 408
T-piece with plug for Aerogen Pro, A, ISO 22 mm, 5 pcs.	MP01031	93
Tape for Eyebrow Sensor, 96 pcs.	MS34209	208
TAT-5000S Probe Caps, single-use, 1,000 pcs., Vista 120 SC	2600948	250
TAT-5000S Sheats, single-use, latexfree, 250 pcs., Vista 120 SC	2600949	250
TAT-5000S Temporal Scanner, Vista 120 SC	2600946	250
Tcore® adapter, reusable	MP00999	141, 193
Tcore® extension cable	MP02201	141, 193
Tcore® mount for M540	MS33239	193
Tcore® sensor, single patient use	MP00989	141, 193
Telemetry transmitter pouch, disposable, 100 pcs.	MS22905	219
TEMP Adapter 2-pin	MP00944	256
TEMP Adapter Jack	MP00943	256
Temperature probe adapter 900 MR 532, for F&P humidifiers	MP01007	31, 50
Temperature probe skin, adult, 3 m, reusable	2601197	229
Temperature probe skin, neonate/infant, 3 m, reusable	2601200	229
Temperature probe, rectal/oral, adult, 3 m, reusable	2601198	229
Temperature probe, rectal/oral, neonatal/Infant, 3 m, reusable	2601199	229
Temperature probe/flow sensor for F&P MR 850, for all breathing circuits, 1.1 m	8414989	30, 391
Temperature probe/flow sensor for F&P MR 850, for all breathing circuits, 1.5 m	8414966	30, 391
Temperature probe/flow sensor for F&P MR 850, for HFV reusable breathing circuits, 0.7 m	8418507	30, 391
Temperature sensor mount, reusable, N	8403056	39, 51, 392
Test filter	6870281	169, 205
Test gas cylinder 5 Vol.% CO ₂ , 95 Vol. % N ₂	6850435	169, 205
Test lung SelfTestLung, latex-free, reusable	MP20020	53
Test lung SelfTestLung, reusable, max. 1 L, not made with natural rubber latex	MP02400	53
Test Lung, reusable, max. 1.5 L, not made with natural rubber latex	8403201	53
ToFscan® Kit configurator	MS32318	208
Transition Set, disposable, 20 pcs.	MX17018	348, 352
Trolley	5790148	376
Trolley, for OR suction unit, with protection jar, w/o 4 L secretion jars	2M85441	305
Tubing grommet for Caleo®, reusable, angle	2M50385	351
Tubing grommet for Caleo®, reusable, flap	2M50412	351
Tubing grommet for Isolette® C2000, 8000 and 8000 plus, reusable	MU12609	354, 356
Tubing grommet for TI 500 and TI 500 Globe-Trotter, reusable	MU04366	363, 366, 369
Tubing grommet, reusable	2M20434	358
Turn knob, compl.	MP00587	435, 341
Two-sided pump ball for air management	MP01590	72, 78, 405
U Universal holder flex 15 for wall rails 25 × 8 to 35 × 10 mm with bottle screw	MP00637	302
Universal holder for wall rails 25 × 8 to 35 × 10 mm with M5 drill hole 15 mm distance	MP00635	302
Universal holder for wall rails 25 × 8 to 35 × 10 mm with M5 drill hole 22 mm distance	MP00636	302
Universal holder profile 30 × 5 for wall rails 25 × 8 to 35 × 10 mm	MP00634	302, 306, 377
Universal holder with DIN rail	G93140	29, 390
Universal holder with groove 25 × 5 for wall rails 25 × 8 to 35 × 10 mm, for compact suction system and secretion jar set	2M86757	302, 306

INDEX BY PRODUCT NAME

Universal pole mount	7485621	112, 147, 148
V VAC CS hose 1.5 m (4.9 ft) DISSB EN HAND	M36299	332
VAC CS hose 1.5 m (4.9 ft) NF EN NF	M35431	322
VAC CS hose 1.5 m (4.9 ft) NIST BS straight	M35404	324
VAC CS hose 1.5 m (4.9 ft) NIST EN AS	M36826	328
VAC CS hose 1.5 m (4.9 ft) NIST EN BS	M35395	323
VAC CS hose 1.5 m (4.9 ft) NIST EN Carba	M35477	326
VAC CS hose 1.5 m (4.9 ft) NIST EN DIN	M35498	319
VAC CS hose 1.5 m (4.9 ft) NIST EN NF	M35422	321
VAC CS hose 1.5 m (4.9 ft) NIST EN SS	M35386	325
VAC CS hose 1.5 m (4.9 ft) NIST EN UNI	M35413	327
VAC CS hose 1.5 m (4.9 ft) NIST NEUTRAL DIN	M36013	320
VAC CS hose 1.5 m (4.9 ft) NIST NEUTRAL NF	M35455	321
VAC CS hose 3 m (9.8 ft) DISSB EN HAND	M36300	332
VAC CS hose 3 m (9.8 ft) NF EN NF	M35432	322
VAC CS hose 3 m (9.8 ft) NIST BS straight	M35405	324
VAC CS hose 3 m (9.8 ft) NIST EN AS	M36827	328
VAC CS hose 3 m (9.8 ft) NIST EN BS	M35396	323
VAC CS hose 3 m (9.8 ft) NIST EN Carba	M35478	326
VAC CS hose 3 m (9.8 ft) NIST EN DIN	M35499	319
VAC CS hose 3 m (9.8 ft) NIST EN NF	M35423	321
VAC CS hose 3 m (9.8 ft) NIST EN SS	M35387	325
VAC CS hose 3 m (9.8 ft) NIST EN UNI	M35414	327
VAC CS hose 3 m (9.8 ft) NIST NEUTRAL DIN	M36014	320
VAC CS hose 3 m (9.8 ft) NIST NEUTRAL NF	M35456	321
VAC CS hose 5 m (16.4 ft) DISS JIS KAWA	M35820	329
VAC CS hose 5 m (16.4 ft) DISSB EN HAND	M36301	332
VAC CS hose 5 m (16.4 ft) NF EN NF	M35433	322
VAC CS hose 5 m (16.4 ft) NIST BS straight	M35406	324
VAC CS hose 5 m (16.4 ft) NIST EN AS	M36828	328
VAC CS hose 5 m (16.4 ft) NIST EN BS	M35397	323
VAC CS hose 5 m (16.4 ft) NIST EN Carba	M35479	326
VAC CS hose 5 m (16.4 ft) NIST EN DIN	M36000	319
VAC CS hose 5 m (16.4 ft) NIST EN NF	M35424	321
VAC CS hose 5 m (16.4 ft) NIST EN SS	M35388	325
VAC CS hose 5 m (16.4 ft) NIST EN UNI	M35415	327
VAC CS hose 5 m (16.4 ft) NIST NEUTRAL DIN	M36015	320
VAC CS hose 5 m (16.4 ft) NIST NEUTRAL NF	M35457	321
VAC CS hose DISS Hand 10L	4199590	330
VAC CS hose DISS Hand 15L	4199589	330
VAC CS hose DISS N/DISS N 10L	4199602	330
VAC CS hose DISS N/DISS N 15L	4199601	330
VAC CS hose DISS N/DISS N 20L	4199600	330
Vacuum aspirator DISS Canada, version for Atlan® A300/A350 and Perseus® A500	MP00687	287
Vacuum aspirator DISS, version for Atlan® A300/A350 and Perseus® A500	MP00686	287
Vacuum aspirator NIST, version for Atlan® A300/A350 and Perseus® A500	MP00689	287
Vacuum regulator dual VarioVac® B/OR Rail, for bronchial and OR suction, gas connection NIST VAC	MP00516	285, 304
Vacuum regulator VarioVac® B AGA, for bronchial suction, gas connection VAC, Swedish standard	MP00504	284
Vacuum regulator VarioVac® B BS, for bronchial suction, gas connection VAC	MP00502	284

INDEX BY PRODUCT NAME

Vacuum regulator VarioVac® B CARBA, for bronchial suction, gas connection VAC	MP00514	284
Vacuum regulator VarioVac® B compact, for bronchial suction, gas connection NIST VAC, including secretion jar set and integrated vacuum regulator	MP00506	292
Vacuum regulator VarioVac® B DIN, for bronchial suction, gas connection VAC	MP00501	284
Vacuum regulator VarioVac® B NF, for bronchial suction, gas connection VAC	MP00503	284
Vacuum regulator VarioVac® B/OR Rail DISS, for bronchial and OR suction, gas connection DISS VAC	MP00631	285
Vacuum regulator VarioVac® B/OR Rail, for bronchial and OR suction, gas connection NIST VAC	MP00505	285, 304
Vacuum regulator VarioVac® D AGA, for continuous and wound drainage, gas connection VAC	MP00510	308
Vacuum regulator VarioVac® D BS, for continuous and wound drainage, gas connection VAC	MP00508	308
Vacuum regulator VarioVac® D CARBA, for continuous and wound drainage, gas connection VAC	MP00515	308
Vacuum regulator VarioVac® D DIN, for continuous and wound drainage, gas connection VAC	MP00507	308
Vacuum regulator VarioVac® D NF, for continuous and wound drainage, gas connection VAC	MP00509	308
Vacuum regulator VarioVac® D Rail, for continuous and wound drainage, gas connection VAC	MP00511	308, 423
Vacuum regulator VarioVac® MRI B/OR Rail, for bronchial and OR suction, gas connection NIST VAC, for connection to all standard wall rails	MP03222	286
Vacuum regulator VarioVac® MRI P Rail, for bronchial suction, gas connection NIST VAC, for connection to all standard wall rails	MP03223	286, 424
Vacuum regulator VarioVac® P compact, for bronchial suction, gas connection NIST VAC, including secretion jar set and integrated vacuum regulator	MP00661	292, 425
Vacuum regulator VarioVac® P Rail DISS, for bronchial suction, gas connection DISS VAC, for connection to all standard wall rails	MP00632	286
Vacuum regulator VarioVac® P Rail, for bronchial suction, gas connection NIST VAC, for connection to all standard wall rails	MP00660	286, 424
Vacuum regulator VarioVac® T DIN, for thoracic drainage, gas connection VAC, for connecting to the central vacuum supply	MP00664	309
Vacuum regulator VarioVac® T Rail WV, for thoracic drainage, gas connection VAC, water vacuum gauge included, for connection to all standard wall rails, double-scaled immersion pipe (kPa and cmH ₂ O)	MP00512	310, 315
Vacuum regulator VarioVac® T Rail, for thoracic drainage, gas connection VAC, for connection to all standard wall rails	MP00663	309
Vacuum suction DISS Canada, for bronchial suction, gas connection VAC	MK03320	287
Vacuum suction DISS, for bronchial suction, gas connection VAC	MK03140	287
Vacuum suction NIST, for bronchial suction, gas connection VAC	MK01422	287
Ventilation hose holder, reusable, 40 cm	8411075	348, 352
Ventilation valve for Oxylog® 1000, reusable	2M86800	45
Ventilation valve for Oxylog® 2000/3000, reusable	8412001	46
VentStar Autobreath N for Resuscitaire® and TN300, N, disposable, 25 pcs.	MU10841	350, 361
VentStar® Aquapor dual heated, Adult and Paediatric, 10 pcs.	MP07130	27
VentStar® Aquapor, dual heated, N plus, 10 pcs.	MP07150	27, 388
VentStar® Aquapor, heated, Adult and Paediatric, 10 pcs.	MP07131	27
VentStar® Aquapor, heated, N plus, 10 pcs.	MP07151	27, 388
VentStar® Breathing Bag 0.5 L, 25 pcs.	MP02710	15, 386
VentStar® Breathing Bag 1 L, 25 pcs.	MP02711	15
VentStar® Breathing Bag 2 L, 25 pcs.	MP02712	15
VentStar® Breathing Bag 3 L, 25 pcs.	MP02713	15

INDEX BY PRODUCT NAME

VentStar® Helix, dual heated, 10 pcs.	MP02606	24
VentStar® Helix, heated, 10 pcs.	MP02607	24
VentStar® Resus Neo with pressure measurement line, disposable, N, 1 m, 25 pcs., not made with natural rubber latex	MP00311	350, 361, 386
VentStar® Resus Neo, disposable, N, 1 m, 25 pcs., not made with natural rubber latex	MP00310	350, 361, 386
Vista 120 printing paper, 4 rolls	2606231	237
Vista 120 rechargeable Lithium-Ion battery	2606229	237
W Wall holder, with fixing screws, for compact suction unit and secretion jar set	M26859	296, 432
Warming mattress TransWarmer®, disposable, self warming for approx. 2 h, 6 pcs.	MP01414	373
Water bag holder 900 MR 290	MP01008	29, 390
Water container with screw thread + gaskets, for VarioVac® T Rail WV (MP00512), VarioAir® T Rail WV (MP00544) and water vacuum gauge (MP00518)	MP00565	318
Water container, for continuous drainage system 2 – 4, w/o screw thread	M20668	317
Water tank and water transfer hose, for Caleo®, reusable	2M50040	352
Water trap for Evita® 2 dura, Evita® 4, Evita® XL, reusable	8413125	85
Water trap Infinity® ID WaterLock® 2, 12 pcs.	6872020	25, 172
Water trap with outlet, reusable N/P	8404760	39, 392
Water trap, reusable, A	8404985	40, 43, 51
Water trap, reusable, N/P	8409627	39, 51, 392
Water vacuum gauge, for thoracic drainage, combination w/ VarioVac® D and VarioAir® D, for connection to all standard wall rails, double-scaled immersion pipe (kPa and cmH2O)	MP00518	315
Watertrap WaterLock® 2, 12 pcs.	6872130	25, 172
Y Y-cable CO ₂ and analog/sync	MS22259	168, 205
Y-cable, reusable, 7-pin connector, for 2 temperature probes	5592154	139, 191
Y-piece 90° with Luer Lock, reusable, A	M33278	42, 43, 50
Y-piece for F&P, reusable, A	8412888	50
Y-piece straight, reusable, N/P	M27077	50, 396
Y-piece with sensor mount, reusable, A	8405435	40, 50
Y-piece, polyamide, reusable, straight, A	M25650	41, 50
Y-piece, reusable, N, 90°, with suction opening	8403075	39, 44, 50, 392, 393, 396

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

REGION EUROPE

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION MIDDLE EAST, AFRICA

Drägerwerk AG & Co. KGaA
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA / PACIFIC

Draeger Singapore Pte. Ltd.
61 Science Park Road
The Galen #04-01
Singapore 117525
Tel +65 6872 9288
Fax +65 6259 0398

REGION CENTRAL AND SOUTH AMERICA

Dräger Indústria e Comércio Ltda.
Al. Pucurui - 51 - Tamboré
06460-100 - Barueri - São Paulo
Tel. +55 (11) 4689-4900
relacionamento@draeger.com

CANADA

Draeger Medical Canada Inc.
2425 Skymark Avenue, Unit 1
Mississauga, Ontario, L4W 4Y6
Tel +1 905 212 6600
Toll-free +1 866 343 2273
Fax +1 905 212 6601
Canada.support@draeger.com

Locate your Regional Sales
Representative at:
www.draeger.com/contact

